

Teamaze

QN-MUS-1

Octopus Music Massager



Instruction

Thank you so much for choosing us,
please read the instruction carefully before using.

Product Introduction

As a small-town youth, the fast pace of first-tier cities and the pressure of life often make me tired and anxious, so I always listen to music to relax. One day I was working with headphones and a colleague put a massager on my head, an ASMR-like experience hit me instantly.

At the same time, inspiration came to mind, "If you listen to the electric sound, you can massage according to the rhythm, wouldn't it be great?" Based on curiosity is the theory that drives human progress, this product was born.

Product Parameters

Product Name	Octopus Music Massager
Minimum Voltage	3.45V
Maximum Voltage	4.2V
Charging Current	500ma
Product Material	ABS+ silicone
Battery Capacity	1100mAh
Charging Time	2.5h
Operating Temperature	10 °C ~50 °C
Endurance Time	18h

*Tips


Thank you for purchasing our Teamaze Octopus Music Massager, please read the instructions carefully before using, pls keep it , especially pay attention to the contents of safety matters.

Product list

Name	quantity
host	1
USB charging cable	1
Product instruction	1
Product warranty card	1
Product certificat	1

The main structure of the product

The main structure and components of the QN-MUS-1 octopus music massager are shown in the following figure:

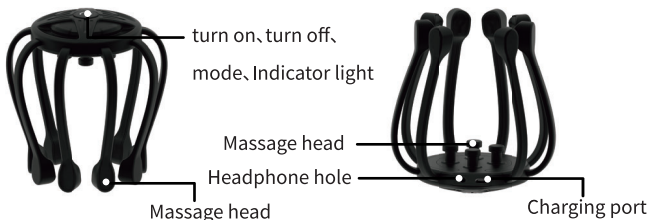

Operating instructions

(Functional operations)

1. Tap the  key to turn on the system and press and hold the 3-second  key to turn it off.

2. Default power-on Bluetooth mode.

3. Tap the  key to start switching modes.

① Bluetooth mode→② Sleep mode→③ Casual mode→④ Sport mode


4. Bluetooth function with headset use safely)

(Bluetooth function)

According to the product voice prompt, open the phone - settings - open Bluetooth, search for "light energy" click to connect


As shown in the figure

It must be connected under Bluetooth mode

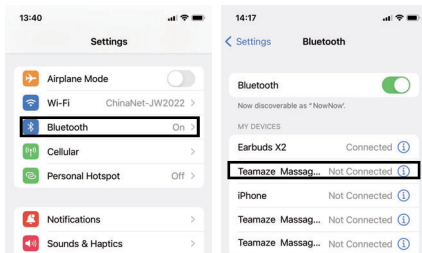

(Safe Use of Headphones) The massager playing volume has been reduced, the maximum volume of the massager is equivalent to 70% of the maximum volume of the mobile phone, and long-term wearing does not damage the ear.

Apple device users need to use in the - settings - sound and touch - click on the headphone safety - turn off to reduce the high volume, otherwise when plugging in the headphones to use the massager, the volume is insufficient will cause the massager to fail to capture the song frequency, so that the vibration frequency and force cannot be automatically adjusted, as shown in the figure:

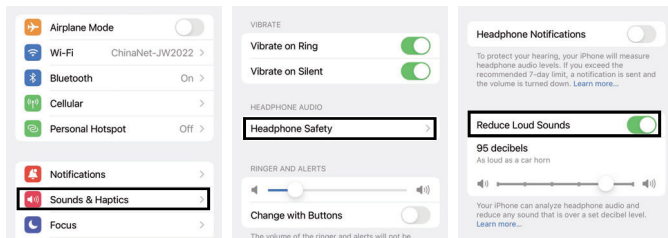

Notes and descriptions

- 1.Power on automatic detection of power: if the power is too low, the display light will display red and blue at the same time, and then automatically shut down, please charge immediately, do not use when charging.
- 2.This product uses lithium ion battery. If the device is not used for more than three months, it needs to be recharged once to avoid failure to start the device due to discharge of lithium battery after a long time no use.
- 3.In the case of normal use, the life of the rechargeable battery is not less than 500 times of charging and discharging, and the life of lithium batteries is generally about 3 years.
- 4.All modes are fixed functions, and the timing time and function are not adjustable.

safety matters

- 1.Do not use this product in places with high heat, flammable, electromagnetic radiation, or humidity.
- 2.Do not use this product in the presence of thunder or lightning.
- 3.Do not disassemble, repair or modify this product without permission, which may cause accidents or malfunctions.
- 4.This product is not waterproof, it is strictly forbidden to enter the water inside the massager, and it is forbidden to use it when the massager enters the water.
- 5.Do not use on infants and young children, people with head trauma, or people who have recently had head surgery.
- 6.Patients with deep stroke, hypertension and hypotension should use with caution, please consult a professional physician, and if there is discomfort during use, please discontinue use immediately.
- 7.Use with caution in patients with slow mobility, physical and mental discomfort or intellectual disability, please use under the supervision of the accompanying staff, and if there is discomfort during use, please discontinue use immediately.
- 8.Please avoid using the massager when drinking alcohol or being unwell.
- 9.This product is a personal health product and is not a substitute for pharmaceuticals and medical devices.
- 10.This product is a household electronic health product and is not intended for diagnostic treatment and other commercial uses.

WARN!  **Not allowed in the following groups**

infants and young children, people with head trauma, people who have recently had head surgery.

Troubleshooting

If the fault occurs during operation, please make the following checks. If the fault cannot be removed or the product still cannot be operated ideally, please contact the after-sales service center.

The basic fault	Possible triggers	Basic exclusion method
Unable to boot	Is the battery insufficient	If yes, please charge the massager
The vibration intensity is too weak	Whether the power is weak	If yes, please charge the massager
The overhead tremor is weak	The four massage heads on the top of the head do not touch the skin well	Please let the massage head fit against the scalp

Maintenance

Take moisture-proof, rat-proof, anti-pollution measures; Cleaning of the product should be carried out in a shutdown state; the massager should be placed in a safe place after use, and heavy objects should not be placed on the massager; When there is a fault should ask a professional repair.

Storage and transportation

- 1.Ambient temperature : -20° C - + 40°C.
- 2.Atmospheric pressure : 860Hpa~ 1060hpas.
- 3.Relative humidity : No more than 80%, no corrosive gas and good ventilation room.
- 4.Cannot be close to the floor, walls and roof.
- 5.Indoor ventilation should be good, avoid strong sunlight and corrosive gas invasion.
- 6.Transportation: the product is equipped with simple anti-vibration facilities; Suitable for air, rail, road and ship transportation, should avoid rain, snow splashing, inversion and collision.
- 7.When the storage period is more than 3 months, the machine should be checked out of the box and recharged before use.
- 8.This product is harmless to human body and can be used safely for a long time.

Product certificate


Light Energy
Octopus Music Massager

NO. _____
model:QN-MUS-1

The date of the inspection

Inspector

Shenzhen Tianyin Juyue Information Technology Co. LTD

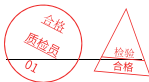

Product warranty card

Product name : Octopus Music Massager

Product model : QN-MUS-1

Customer name :

Telephone :

Date of purchase :

Promise

1.Within one year from the date of sale of the product, in the event of a malfunction in the case of normal use; after testing non-human destruction, The company will provide free maintenance, parts replacement and other services according to the fault situation.) (The free maintenance promised by our company only covers the cost of working hours and materials, excluding the cost of express delivery.

2.This warranty card is the company's certificate to provide after-sales service to customers, this card needs to be carefully filled in, the dealer stamped after the valid.

3.In any of the following cases, free warranty service is not implemented and maintenance is required for a fee.

- (1) Beyond the warranty period.
- (2) Failure to use according to the requirements of the product instruction manual, improper maintenance or storage caused damage.
- (3) Self-disassembly and maintenance, use under special circumstances, etc.
- (4) Failure or damage of the massager caused by force majeure, man-made damage.
- (5) Replacement of wearing parts and random accessories.

This card is issued with the goods, in order to ensure that you fully enjoy the company's free warranty service rights, please keep this card properly, lost and will not be replaced.

Card number

Customer name	Telephone
Address	
Date of purchase	



Teamaze

**Counterfeiting of patented products
will be banned.**

Utility models : 202222074904.1

Product name : Octopus Music Massager

Product model : QN-MUS-1

Company name : Shenzhen Tianyin Juyue Information Technology Co. LTD

Company address : Wanfuyuan Industrial Park, Dongsheng Town, Zhongshan City, Guangdong Province

After-sales hotline : 177-4859-2546

*The final interpretation of this specification belongs to Shenzhen Tianyin Juyue Information Technology Co., Ltd., if there are any product changes and design changes, the actual product shall prevail without prior notice.

FCC Warning:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 0cm between the radiator and your body.