



H1发箍形头颈部按摩器

H1 Head Massager

使用说明书
Instruction Manual



按摩
Mini Size

便携
Light Weight

蓝牙APP
App Control

图片仅供参考，产品以实物为准！

Dear Hi5 Customer,

Thank you for choosing the Hi5 Hairband Head Massager! H1 Head Massager is a mini and light weight massager with hairband-shape. Its flexibility and elasticity make the product fit the head contour like a headband. It can also be used for neck massage.


This product has 11 independent vibration massage points, the elasticity of each massage head can adapt to different head contours.

Head massage can promote the blood circulation of the head, improve the ability of the brain to absorb oxygen, and improve the function of the cerebral cortex. Stimulating hair follicles can strengthen hair roots and prevent hair from falling off, which is beneficial to hair growth. Through massage, the acupuncture points on the head can be effectively stimulated to refresh the mind, eliminate fatigue, tension and anxiety, and make the brain energetic.

Please read this manual carefully before use, and use this product correctly.

Hi5 Team

Product Configuration

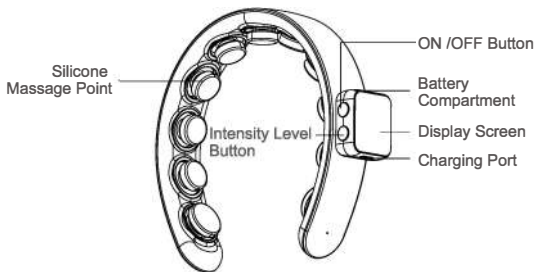

Product Parameter

Product Name: Head Massager

Model: H1-Deluxe

Rated Input Power: 3.5W

Rated Charging Voltage: 5V

Battery Capacity: 700mAh

Smart Timer: 15 Minutes

Jiangsu Mingteng Medical Equipment Co., Ltd

Executive Standard: GB4706.1-2005 GB4706.10-2008

Product Features

- 11 Independent vibration massage points
- Three intensity level adjustment
- 5 Working Modes: -Awake mode,
-Soothing mode,
-Sleep mode,
-Active mode
-Neck massage mode
- Mobile APP Control

Wearing Method :



Middle of Head



Neck



Forehead



Thigh



Back Head



Top of Head

Instructions for Use

- **ON/OFF:** Press and hold the [ON/OFF] button for 1.5s to start up with default awake mode vibration and medium intensity level.
- **Massage Mode Adjustment:** Short press the [ON/OFF] button to switch the massage mode.
- **Intensity Level Adjustment:** Short press [Intensity Level] button to adjust the intensity level.
- **Shutdown:** Press and hold the [ON/OFF] button for 1.5s to shutdown equipment when its working.
- **Charge:** Please use the 5V1A or 5V2A charger and the original type C charging cable for charging. It is recommended to separate the battery compartment from the equipment when charging, and charge the battery compartment separately until it is fully charged.
- **Battery Compartment:** The battery compartment can send instructions to equipment and supply power to equipment. If the battery compartment is removed during massage, the massage will be suspended, and the massage will resume after installation.
- **Monitoring Function of Battery Compartment:** After the battery compartment is separated from the headband for 120 seconds during use, it will automatically shut down.
- **Display Screen:** The display screen will display: current massage mode, current massage intensity, remaining battery value, Bluetooth connection status.

App Installation

iPhones:

Step1, Enter the Apple's App store and search for "Hi5 Health" to download and install.

Step2, Or, open a browser with scanning function to scan the QR code to download and install.



Android Phones:

Step1, Enter the Andriod's App store and search for "Hi5 Health" to download and install.

Step2, Or, open a browser with scanning function to scan the QR code to download and install.



Reminding: Since the App is constantly being upgraded, if you need the latest version, please uninstall it and reinstall it.

Bluetooth Connection

Connect the Device

1. Power on the device.
2. Turn on the Bluetooth function of the mobile phone.
3. Open the "Hi5 Health" APP in mobile phone.
4. Enter the H1 Headband massager operation page.
5. Click the Bluetooth icon at the bottom of the operation page to connect with the H1 headband device;
6. The bluetooth name is : H1 Hairband Head Massager


App Control

Please follow the FIG below to operate the equipment

Note1: When using the APP to control the device as well as the buttons on the battery compartment to control the device, which may cause conflict.

Note 2: When device is more than a certain distance from the mobile phone, the App control may fail.

Bluetooth Connection

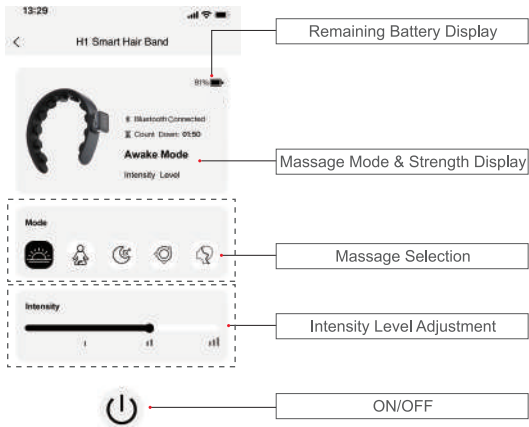

Disclaimer

Warnings & Precautions

a) Prohibition

1. Patients with head diseases and injuries are prohibited to use;
2. Patients with osteoporosis, vertebral fractures are prohibited to use;
3. Patients with hypertension caused by diabetes are prohibited to use;
4. Patients with a history of epilepsy are prohibited to use;
5. Patients with moderate and severe cardiovascular diseases and severe heart diseases are prohibited to use;
6. Patients with head skin disturbance or abnormal head skin are prohibited to use;
7. Patients who have had craniotomy on the head and have intracranial hemangioma are prohibited to use
8. Patients in the rehabilitation stage after head surgery are prohibited to use;
9. Pregnant women are prohibited to use;
10. Patients whose surgical wounds have not healed or have not recovered are prohibited to use;
11. Patients who are being treated and feel physically abnormal are prohibited to use;
12. Patients with hypertension, arteriosclerosis or cerebral thrombosis are prohibited to use;

Disclaimer

13. Patients with high fever and malignant tumors are prohibited to use;
14. Personnel, infants and children who cannot take care of themselves are prohibited to use.

b) Safety Precaution

1. Please turn off the power supply when not using this product;
2. For your safety, please power off the product when its not in use;
3. Keep away from fire, corrosive liquids and various unfavorable environments;
4. Do not wash the product. If the product is wet, it must be completely dried before use;
5. Clean the surface with a soft cloth or neutral detergent and dry it naturally;
6. Do not use the product during charging;
7. This product contains lithium batteries. To avoid danger, do not assemble, disassemble, repair by yourself, nor throw it into fire;
8. Do not use the product continuously for a long time, otherwise it may cause your massage area to a fatigue or discomfort state;
9. Please stop using the massager immediately when you feel uncomfortable.
10. The massager is designed as a head massager and is not for diagnostic, therapeutic and other commercial purposes.

Disclaimer

c) Service Environment

1. Do not place the product in the bathroom or use it in a humid place;
2. Do not put the product into water or place it in a high-temperature or humid environment such as balcony or bathroom for a long time;
3. Do not place the product in an electric heater, stove or place exposed to direct sunlight;
4. Do not use it when lightning.

d) Power Plug and Power Cord

1. Do not use a power supply other than the rated voltage for charging. Please use the original power cord of the product to avoid fire or damage to the product.
2. When pulling out the power cord, hold the plug of the power cord instead of the power cord.
3. If the power cord is damaged, do not replace the unqualified cord by yourself to avoid danger.

Warranty

If the product fails due to non-human factors within one year from the date of purchase, we will provide services including repair, replacement or necessary parts for free. However, this warranty does not apply for the below conditions:

- a) Damage caused by accident, misuse, improper handling or transportation.
- b) Unauthorized repair.
- c) Damage caused by improper use.
- d) Damage that exceeded the original amount of the product.
- e) Damage caused by improper storage conditions;
- f) Evidence of purchase date cannot be provided.

Trouble Shooting

If you encounter problems when using this product, please refer to the following solutions. If this product is still in a faulty state or cannot be used, please unplug the power cord and contact us for inspection and maintenance.

Trouble Shooting

Fault	Possible Causes	Methods
Massager does not work	Low battery	Recharge
	Automatic shutdown	Reboot
	Battery compartment not installed	Install battery compartment
Intensity got weak	Low battery	Recharge
	Battery compartment is not in good contact with the device	Install battery compartment
Display screen is abnormal	Circuit failure	Contact sales personnel for after service
Missing battery compartment	/	Contact sales personnel for after service

- Users are not allowed to disassemble or repair this product by themselves. Damage to the product caused by disassembly or repair without authorization of the company will invalidate the product warranty, and the company is not responsible for any consequences caused by this situation.

FCC Warning

15.19 Labeling requirements.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

15.21 Information to user.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

15.105 Information to the user.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are

designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception,

which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC RF Radiation Exposure Statement:

- 1.This Transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.
- 2.This equipment complies with RF radiation exposure limits set forth for an uncontrolled environment.

This equipment should be installed and operated with minimum distance 5mm between the radiator and your body.