



QUICK START GUIDE — SM

Accordion Fold

Product: Quick Start Guide
Material: Matte Paper — 80 lb Text
Printing: Double-sided
Style: Accordion Fold
Details: Rounded Corners — 0.125 Inches
Color: CMYK
Dimensions: 5.125 (W) x 3.75 (H) Inches

—

Part Number: 05-040-05
Description: Quick Start Guide, SM
Revision: C



C0 / M32 / Y100 / K5

FRONT

What's in the box

- 1x SM grip Device
- 1x USB-C Cable Charger
- 1x Quick Start Guide

Device diagram

1. Finger Plungers (Index, Middle, Ring, Small)
2. Thumb Pinch/Strain Gauge
3. Status Light
4. USB-C Charging Port

Download the GRIP mobile app

Download on the App Store

GET IT ON Google Play

BACK

Set up with your phone

1. Download and install the grip mobile app:

- Apple App Store for iPhones
- Google Play Store for Android phones

2. Account Setup:

- If you already have a grip account, log in to your account.
- If you do not have a grip account, tap the "Get Started" button to begin the account setup process. A series of steps will assist you throughout the process.

3. Pair grip device:

- Enable Bluetooth on your smartphone and open/sign-in on the Umajin app.
- With your SM grip device, press anyone of the figure plungers down to turn on the device, this will automatically put the device in pairing/advertising mode.
- In the Umajin app, once your device is connected, you will be prompted to start your baseline testing and begin your care plan.

How to use

Charging the grip Device

- Charging: Alternating pulses between battery indicator color and blue for pairing.
- Charged: Pulses battery indicator color then alternates with blue LED for pairing.

Sleep Mode

- The device will go to sleep mode ONLY when you disconnect it from the app.
- Wake up or turn on the device by pressing any figure plunger.

Product features

- Handheld device for strengthening hand muscles
- Pairs to smartphone for added functionality
- Offers interactive hand therapy exercises with real-time feedback
- Track and share progress over time with a hand therapist
- Works for both right and left hands
- Comes in two different sizes "SM and LM"

Technical specifications

Bluetooth 5.1

Battery Life up to 14 days

Charging Type USB-C

Compatible iPhone and Android

Product Dimensions 105.0 x 65.25 x 36.5 mm (SM)

Device Description
Smart hand therapy and exercise device

Intended Use
To increase flexion and the extension range of motion of the joints of the fingers and hand.

Indications for Use
For use in the home or clinical setting.

Indications for Use

- Open wounds on hand or fingers
- Allergy to materials that comprise the GRIP device
- Spastic contracture of hand
- Immediate fracture
- Inflammatory arthritis with diffuse deformity

Facebook @OrthoRPM

Twitter @ortho_rpm

LinkedIn @orthorpm

Manufacturer

Date of manufacture

Catalog number

Serial number

Consult instructions for use

Protect from heat and radioactive sources

Recycle the lithium-ion battery

Do not dispose of in the trash



QUICK START GUIDE

Accordion Fold

- Dieline
- - - Valley Fold
- - - Mountain Fold

FRONT

1

2

3

4

BACK

5

6


7

8

[5.125 (W) Inches]

[3.75 (H) Inches]

1

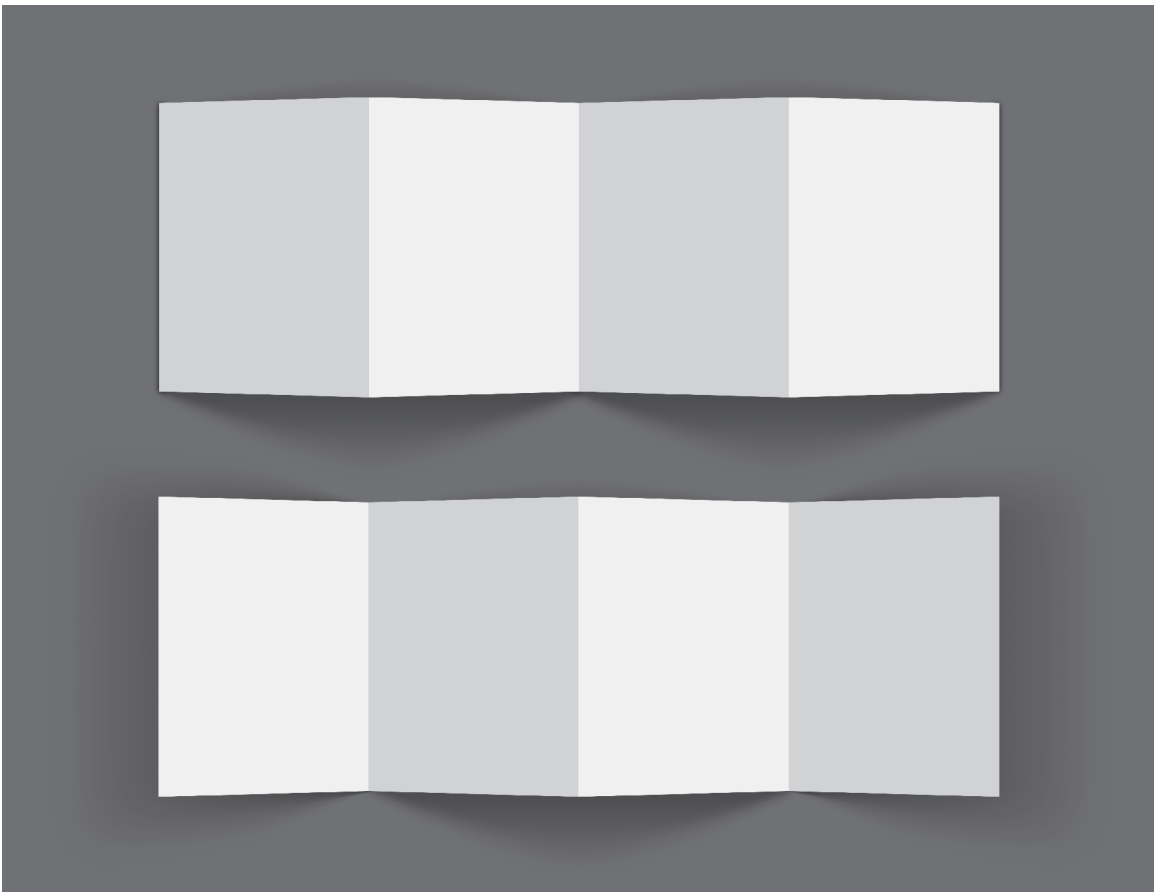

Please create the actual die-line with all the details that are shown from the reference die-line on the right.

1. Rounded corners — 0.125 Inches