

Basic Product Information

Products


Intelligent head eye massager adopts microcomputer chip control technology, through the air pressure/hot compress/vibration effect on the head and eye area, it helps the head and eye muscles to relax, and brings you an all-round comfortable experience.  
This product is suitable for all kinds of people who use computers, cell phones and other electronic devices for a long time.

Product Parameters

[Product Features] Head massage  
[Massage Mode] Five modes adjustable  
[Rated Voltage] DC 5V  
[Packing List] Host \*1 USB cable \*1 Instruction manual \*1

[Battery Capacity] 1000mAh  
[Working Time] 15 minutes  
[Rated Current] 1A

Product mix

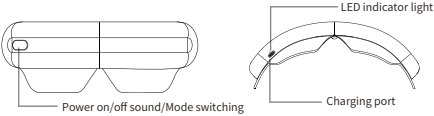

(Note: The figure above is a schematic diagram, product appearance shall be subject to the actual)

Preparation and charging before use

Pre-use Preparation

Aim the head massager at your forehead, put on the machine, adjust the elastic band to the most comfortable size, and adjust the machine up and down to the most comfortable position.

Please remove glasses, contact lenses, hair clips, hair accessories, etc. before wearing to avoid discomfort during use. Do not use this product with wet hands. Do not fall or step on this product.

Charging

The USB plug of the USB charging cable is inserted into the USB socket of the power adapter (5V-1A), then the power adapter (5V-1A) into the power outlet, and finally insert the plug of the USB charging cable into the charging port of the product.

Warm reminder

During charging, the indicator light flashes, and the indicator light is always on when fully charged.  
During charging, the product can not be turned on and work, after charging is completed, please disconnect the power adapter.  
If the product is charging while working, the product will stop working automatically.

Usage

Turn the product on or off

Press the " " key for 2 seconds to turn on the machine. After the machine is turned on, the indicator light is lit and the machine starts to work, default is automatic mode. The working time is 15 minutes, playing music.

While the product is working, keep pressing the " " button for 2 seconds to turn off the machine or the machine will turn off automatically after 15 minutes of use.

Switching Mode

Press the " " key for less than one second to switch modes. Operate once to switch to the second mode; operate twice to switch to the third mode; and so on.

Music control

Press the " " key twice in a row for less than - seconds to turn the music off or on.

Mode introduction

Automatic mode = Air pressure massage + hot compress massage + vibration massage + music  
Decompression mode = Air pressure massage + hot compress massage + music  
Vitality mode = Air pressure massage + vibratory massage + music  
Comfort mode = Hot compress massage + music  
Sleep mode = Air pressure massage + music  
(Note: Because the product version is constantly updated, the actual mode may differ from the manual)

Connecting Bluetooth

After powering on the machine, the machine Bluetooth is turned on by default. When you need to link Bluetooth, please open your cell phone's Bluetooth and search for Bluetooth (default Bluetooth name is:Head Eye Massager). Play cell phone music after successful connection. After disconnecting the Bluetooth, the machine plays the built-in meditation music.

Transportation and storage

Storage environment

Room temperature, no corrosive gases, well ventilated, dry room.

Transport

This product has simple shock-proof facilities, suitable for air, railway, road and ship transport, should avoid rain, exposure, inversion and collision and other bad conditions.

Cleaning and storage

Please use a dry soft cloth to clean the product (please turn off the machine when cleaning, do not use benzene or diluent and other corrosive cleaning agents, otherwise it will cause fading, surface aging, structural damage and other problems.  
Please put the product in a well-ventilated, dry room, and cover it with a protective cover when it is not used for a long time. When the storage time exceeds 3 months, please open the box for inspection and charge for more than half an hour.

Caveat

- Do not use when there are symptoms of eye disease. People who have performed eye surgery, suffered from glaucoma, cataract, retinal detachment and other major eye diseases are prohibited from using.
- Do not wash this product directly with water. Please wipe it with cotton soaked in alcohol or water.
- This product has a hot surface, please pay attention to the use of people who are not sensitive to heat. Minors, the disabled, the elderly, etc., should be used under adult supervision.
- For sensitive skin, atopic dermatitis, cosmetic materials placed on the site and such special groups, please use this product under the guidance of a doctor.
- Do not disassemble this product by yourself. This product contains batteries, do not put dry high temperature (above 50° C) or into the fire, otherwise there will be explosion risk.
- This product is a personal health care massage product, can not replace medicine or medical equipment.
- This product is not a toy and should not be used by children alone unless parents or guardians have set it up in advance, or have given children adequate guidance on its use.
- The use of this product will block the line of sight, do not use in walking, driving and other behaviors.
- After using this product, the eyes need to return from the relaxed state to the tense state, and there will be short-term blurred vision, which is a normal phenomenon. Please close your eyes for 1-3 minutes, and your vision will gradually become clear.
- The storage and use temperature of this product is -10°C ~ 50°C. It is prohibited to store and use in extreme environment (low dry temperature negative 10°C or more than 50°C).
- This product should be charged at least once every 3 months when it is not used for a long time.
- Due to product hardware and software upgrade, product functions, please refer to the actual.

Common problems and solutions

Problems encountered	Possible cause	Treatment method
Automatic shutdown unable to turn on	Low battery	Please charge
The charging indicator light is off	The USB port is improperly connected to the power supply	Plug in the USB charging cable/power cable

If the above processing is carried out and there is an abnormal phenomenon, please turn off the power and contact the dealer you purchased or the company.

Please do not disassemble or repair the product without authorization, and the company will not be responsible for any product damage, personal injury or other losses.

After-sales and maintenance

This product enjoys free maintenance within one year from the date of purchase. Free maintenance only includes the cost of working hours and materials, not including it includes the transportation fee and does not provide on-site service. This product comes with accessories such as charging cable and other vulnerable products, not within the scope of three guarantees.

Free warranty

- Can provide "warranty card"
- Can provide the purchase invoice or purchase certificate, and non-human causes of failure or damage.
- In the warranty period (has purchased the invoice or the date of issuance of the purchase of the machine, no purchase invoice or voucher, to the body bar code factory date).

Non-free warranty content (material fee, labor fee)

- Machine failure or damage caused by human causes.
- More than three guarantee period or no "warranty card".

Warranty card

Thank you for purchasing this product, we will be in accordance with the "People's Republic of China Consumer Rights and interests Protection Law", sincerely to provide you with services, for your rights and interests, please fill in the following details and true.

Customer name		Birth time	
Correspondence address			
Purchase model	Contact number		
Buying shop	Date of purchase		
Local service telephone	Clerk's signature		

Fault/Problem description:

Maintenance record

Maintenance date	Fault content and troubleshooting	Maintenance unit	Maintenance personnel signature

Note: This warranty card is used when you need after-sales service after purchasing this product, please keep it properly. When you need to repair, please show this warranty card and the proof of purchase.

Statement

Due to product technology innovation or product improvement, some parameters may be changed without prior notice. Please refer to the actual product received.  
The legend shown in this manual is based on the appearance of the standard model and is for demonstration purposes only. The actual appearance is subject to the purchased product.

The right to interpret the relevant regulations of this manual belongs to the manufacturer

Product certificate


Product Name : Massager

Production date :

Inspector :

Inspection conclusion : Qualified

Version number: V2.0



MASSAGER

(Air pressure + Vibration + Hot compress + Bluetooth + Built-in music + Voice broadcast)

INSTRUCTION MANUAL

Please read this manual carefully before using the product  
This manual contains service instructions, please keep it in a safe place.

1

2

3

4

5

6

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception,

which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction