

尺寸：120*70g
材质：128g双铜

Enjoy life and health

Smart Head Massager
User's Manual




Introduction

Dear users:
Hello!
Welcome to use our products, in order for you to more fully experience our products, please be sure to carefully read and save this instruction.


- 01 -

Matters Needing Attention



Danger Warning

- Do not use this product when there are medical electronic devices such as pacemakers, artificial heart and lung, and electrocardiogram devices in the body;
- When patients use high-frequency surgical equipment on their bodies, do not use this product at the same time;
- It is prohibited to use this product in flammable and explosive places such as gas stations, gas stations, high-dust workshops;
- Do not use this product while driving, bathing or showering;
- Do not use this product in humid, high temperature, high electromagnetic environment;
- Do not use or clean this product while charging.




Attention

Consult your doctor before use if you have any of the following conditions:

- Pregnant women, sick or in poor health;
- Patients undergoing treatment with sensory abnormalities
- Patients with high fever and infectious diseases;
- Patients with hypertension, heart disease, cerebral thrombosis, arteriosclerosis or malignant tumors;
- Skin perception disorder or abnormal skin sensitivity;
- Postoperative rehabilitation stage;
- People who are uncertain about their physical health.

Matters Needing Attention



Attention

- Do not use this product when the head, eyes, neck and other parts are injured, redness, inflammation, congestion or other discomfort;
- This product has a heating function. For users who are prone to heat allergies or are not sensitive to heat sources, please use the heating function with caution;
- If you feel uncomfortable during or after use, please stop using this product and consult a doctor in time;
- This product is a personal care product, not a substitute for professional treatment or health care, and should not be used for diagnostic or other commercial purposes;
- Do not use with other medical devices;
- Do not use it in places with heavy humidity such as bathrooms.

Do not allow liquids such as water to come into contact with the product;

Do not use on the same part for more than 30 minutes, otherwise it may cause discomfort;


Do not use the heating function for more than 30 minutes, otherwise the skin may be scalded by low temperature.

Even if the temperature is relatively low (38℃~42℃), if you contact the same part of the skin for a long time, you will still have burns without feeling heat or pain.

- 02 -

- 03 -

Matters Needing Attention




Attention

- This product contains a battery, do not disassemble, modify, water the product or throw it into the fire;
- In case of abnormal situation that can not work, please immediately remove and cut off the power supply, check and eliminate the basic fault or consult customer service
- Do not charge, use or place in high temperature places such as near the stove or under the hot sun;
- The electric control box of this product is not waterproof, please do not wash and pay attention to waterproof;
- This product can only be decomposed at the time of disuse, if you need to repair, please consult customer service in time;
- Do not leave the product in front of unattended children;
- Do not allow children or persons with limited ability to use this product without adult supervision;
- This product cannot replace any medical treatment;
- This product is not intended for any sick or incapacitated person.

Packing Instructions

Please check whether the cloth cover is damaged before use. If found damaged, do not use. If you find a broken leather case, do not repair it yourself. Please contact the store or dealer.

The list of contents is as follows:
01. Product
02. Type-C charging cable
03. Instruction Manual




Attention

- Place this product in a dry location, away from prolonged exposure to sunlight.
- Carefully inspect the product for any damage prior to use
- Cannot be turned on during charging. The product works for 15 minutes by default and turns off automatically.

- 04 -

- 05 -

Matters Needing Attention

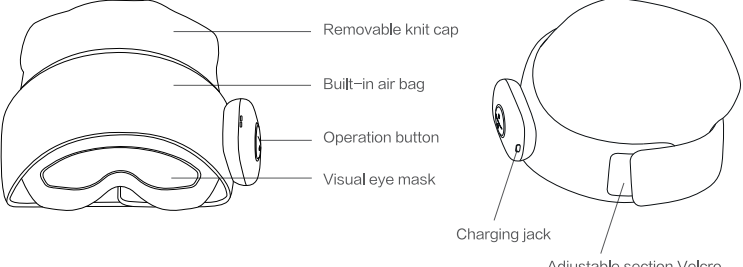

Removable knit cap

Built-in air bag

Operation button


Visual eye mask

Charging jack

Adjustable section Velcro


Preparation Before Use

Ensure the product is powered on. If the battery level is too low, the machine will either turn off or fail to start. Please charge it before use.

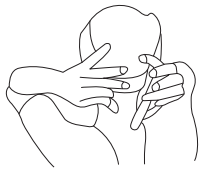


Wear a massager

Remove the massager from the box and place it on your head in a comfortable position.



Align the position and adjust the wear

The eyes and the bridge of the nose are aligned, and the visual area is flush with the line of sight of the eyes; Adjust the rear elastic strap to personal comfort Place, wear complete.

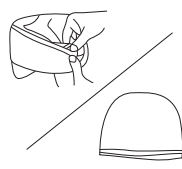


Remove the cap body

Find the zipper head of the messenger cap body and pull the leather case along the zipper to remove the cap body,

- 06 -

- 07 -


Functional operation indication



Power key:

Long press the power button for 1.5 seconds to switch the machine on and off, short press the switch mode, and cycle between [Relieve mode], [Vitality mode], [Sleeping mode] and [Pause mode] in turn;

Double click: ON/OFF the light.



Heating key:

By default, the hot compress [low temperature] is started at startup. Short press the heating key to adjust the temperature gear, and switch between [high temperature], [temperature off], and [low temperature] in turn. The voice will remind you accordingly. The change of hot compress gear does not affect the massage mode. Press and hold Heating key for 1.5 seconds: Turn on/off voice announcement.Temperature resistance (38±3℃) voice reminder (low temperature), the temperature of high-grade (42±3℃) voice reminder (high temperature).

Mode introduction:


① Relieve mode: Relieve Air Pressure + hot compress (Default Low Temperature)

② Vitality mode: Vitality Air Pressure + hot compress (Default Low Temperature)

③ Sleeping mode: Gentle Air Pressure + hot compress (Default Low Temperature)

④ Pause mode: hot compress (Default Low Temperature)


Support Bluetooth, Bluetooth name:RLXPROO, connected Bluetooth devices can play music and audio sources in any mode .

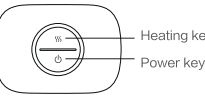

Heating key

Power key

Charging

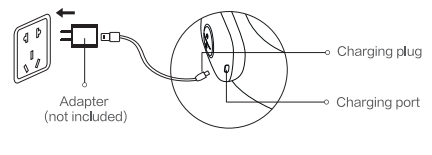

Charging plug

Charging port

Adapter (not included)

1. Use the charging cable delivered with the product to connect to the USB port for charging.

2. The power button indicator is steady white and blinks red when the battery is low.

3. During the charging process, the power button light blinks red. After the charging is complete, the power button light keeps on green.

4. During the charging process, the product does not support startup work, and the product starts the self-protection program and stops working.

Multiple charging modes

Computer USB interface/Mobile power supply /5V ≡ 1A power adapter

Note: The default working time is 15 minutes. If you want to extend the working time, you need to start the host again.

- 08 -

- 09 -

Product Information

Product name: Smart head massager

Product model: FB-02

Rated input: 5V ≡

Rated power: 5W

Battery capacity: 1500mAh

Charging interface: Type-C

Product weight: about 315g

Product size: L670mm*W150mm*H50mm

Troubleshooting and Maintenance

1. If a fault occurs during operation, perform the following check. If the fault cannot be removed or the product still cannot operate normally, it should be repaired by professional and technical personnel, do not disassemble the product by yourself!

2. The basic fault of the product still cannot be solved, please consult customer service staff.

Basic Failure	Possible causes	Basic exclusion method
Unable to boot	1. Low battery power 2. It's charging 3. Power is not turned on	1. Charge before use 2. Do not use this product while charging 3. Make sure the massager is on
No charging	Undercharge voltage	Replace the charging method or charging device
Products to stop working at runtime	Security timing function start, 15 minutes products automatically stop. It's not a malfunction.	Just press the on button again, and it'll work again.

- 10 -

- 11 -

FCC Warning

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

NOTE 1: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

NOTE 2: Any changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.