


High-end Homecare Device

Lymfit.Plus

MANUAL

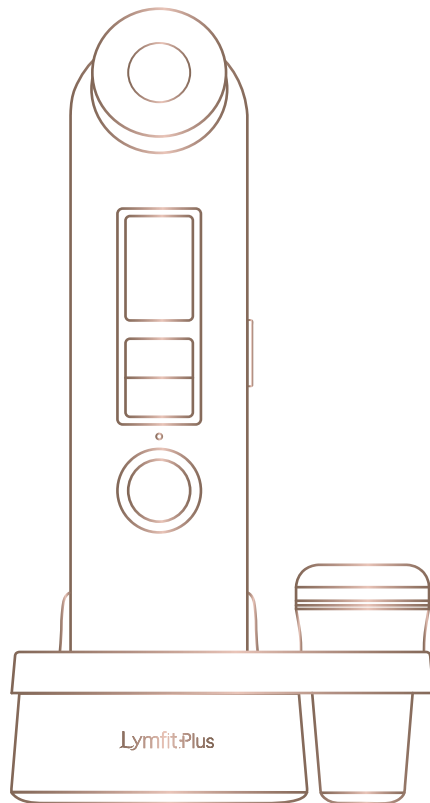

Lymfit.Plus

A high-end homecare device
containing the performance of aesthetic equipment

[Information in English →](#)

1) Before Use	25
Precautions / Product Care	
2) Lymfit.Plus Introduction	26
Product Introduction, Specifications	
Product Components, Name / Head Attachment and Detachment Method / Head Holder Assemble Method / Charging Method	
3) SONO Usage Method	30
Product Main Function / Product Operation Method	
Location of Lymph Nodes / Must Read Before Use, Mode	
Lymph Care Method / Ampoule Absorption Method	
4) HIFU Usage Method	36
Product Main Function / Product Operation Method	
Area Specific Shots / Must Read Before Use	
Lifting Care Method / Handling Usage Tip	
5) Warranty	42
Warranty Regulations / Product Warranty	
6) Manufacturer Information	44

✳ For your safety, please follow these guidelines.

Please read this manual carefully before using the product, and keep it for future reference.

✳ CAUTION


1. This product is not a medical equipment.
2. Remove the film from the head before using the device.
3. Please keep out of reach of children.
4. Not intended for use by individuals with reduced physical-sensory, or mental capabilities(including children) or those lacking experience and knowledge, unless supervised or given instructions on device use by a responsible person.
5. Avoid use on children-pregnant-breastfeeding women.
6. Remove metal accessories such as necklaces-earrings etc., before use.
7. Avoid use on areas with skin diseases such as wounds-infections-boils.
8. It is recommended to consult with a medical professional when using with dental treatments such as implants-nerve treatments-orthodontics.
9. It is recommended to consult with a medical professional before undergoing any cosmetic procedures.
10. Avoid use on people who require health precautions, such as having pacemaker-device implants-history of severe illness-outpatient treatment.
11. This product is a skin care device intended for use on the face only. Avoid using it on any area other than the face.
12. Avoid use on pupils(eyes), around the eyes and where the thyroid gland passes.
13. During the procedure, be sure to apply a sufficient amount of gel and use the device in close contact with the skin to ensure there is no space between the skin and the device.
14. Use carefully, as prolonged use on one area may cause skin irritation.
15. Set the appropriate power based on the skin sensitivity.
16. We are not responsible for any side effects that may occur as a result of not following our recommendations.

✳ Product Care

1. If a detachable power supply is provided, only use the device with the power supply provided with it.
2. Do not modify or disassemble-repair the product. If the user modifies-disassembles-repairs the product, free A/S service cannot be provided within the warranty period, therefore be sure to use the designated A/S center.
3. Do not wash or immerse the product in water. Be careful as it may cause electric shock-malfunction.
4. Disinfect the SONO head with alcohol after use. When disinfecting, use a soft cloth to avoid scratches or damage, as rough cloths or excessive wiping may cause surface damage.
5. Disinfect the HIFU head with alcohol after use. When disinfecting, use a soft cloth.(Alcohol prohibited for black film) The black film may harden upon contact with alcohol, causing energy output issues and film damage.
6. Use wet tissues or a soft cloth to clean the main body, as using alcohol may cause discoloration-contamination.
7. For detailed product care-A/S services, follow the manufacturer's warranty regulations.

Product Introduction

Lymfit.Plus is a high-end homecare device that delivers the performance of professional aesthetic equipment. It combines high-density multi-ultrasound (SONO) for epidermal lymph care with high-intensity focused ultrasound (HIFU), making it effective for lymph circulation-lifting-anti-aging treatments.

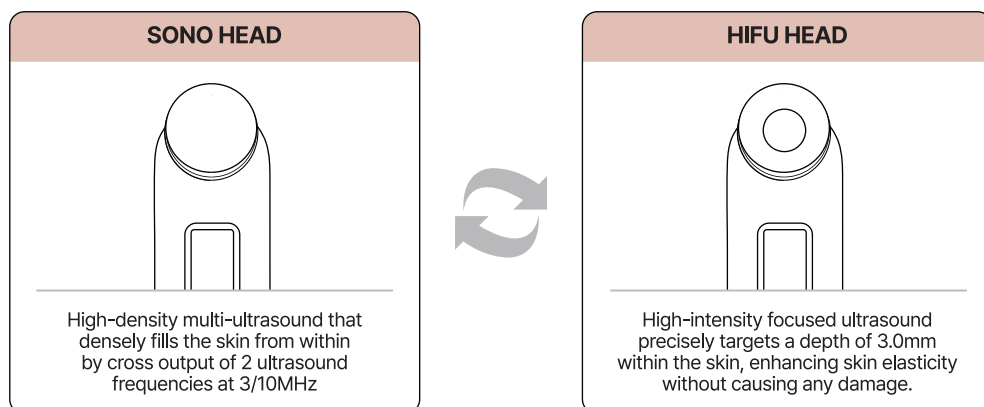

Product Specifications


This product's specifications may be subject to change without separate notice for performance improvement.

Model Name	JHS200B	Product Name	Skin care device
Product Size	W48 x H198 x D57 mm	Product Weight	About 274g
SONO Weight	About 47g	HIFU Weight	About 28g
Input	DC 5V2A	Battery	2400mAh
Charging Cable	C-TYPE	LCD	1.14inch (135 x 240 px)
Charging Time	About 180 minutes	Usage Time	About 150 minutes
Manufacturer/Seller	JLU Co., Ltd.	Date of Manufacture	Separate Notation
A/S Center	+82 10-4505-7989	Customer Center	+82 51-331-7185


2) Lymfit.Plus Introduction

27


Product Components

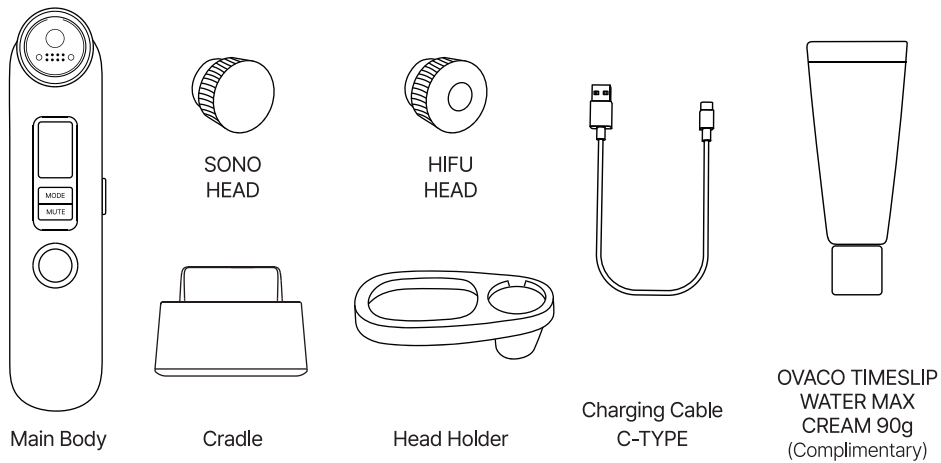


Product Name

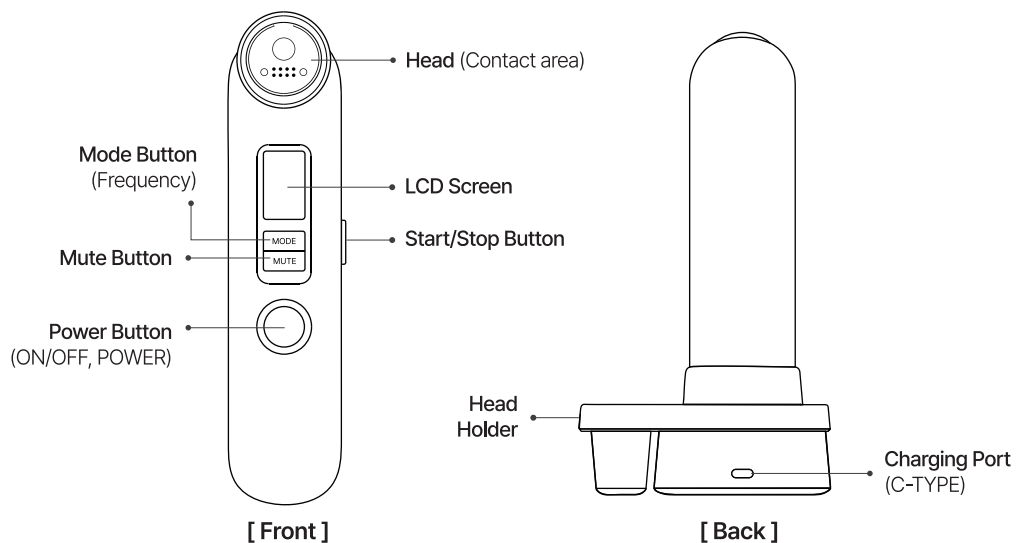


Head Attachment and Detachment Method

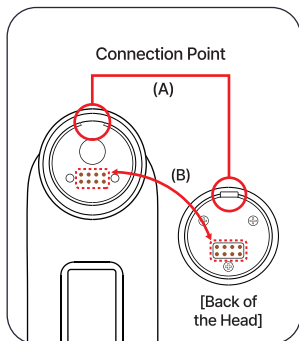

1. Magnets are inserted in the main body and the back of the head, allowing for easy attachment and detachment.

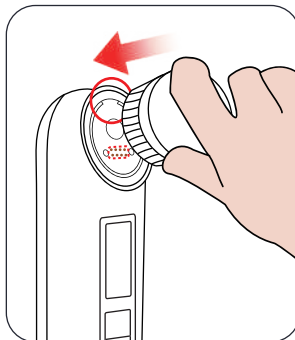


2. During the attachment, align the main body and head with the connection points (A) and (B) and secure them together.

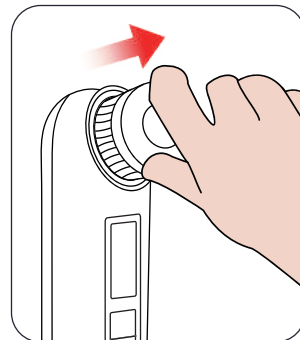


3. During the detachment, pull the head away to separate the magnets.


Head Holder Assemble Method

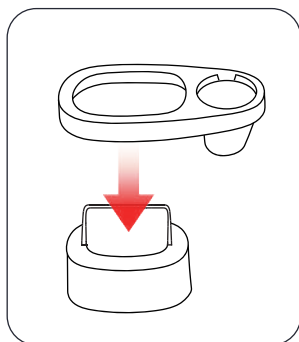

1. Assemble the head stand and cradle by aligning them according to their shape.

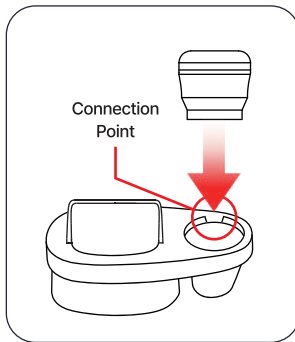

2. Store the head in its designated position.

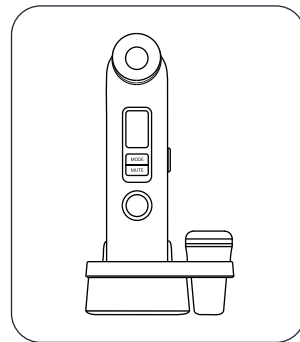


3. Place the main body and head as shown in the picture.

Charging Method

When connecting the charging cable to the cradle, push the connector in completely.

If not, charging may not proceed smoothly.

If the device does not power on, check the connection at the cradle charging port.

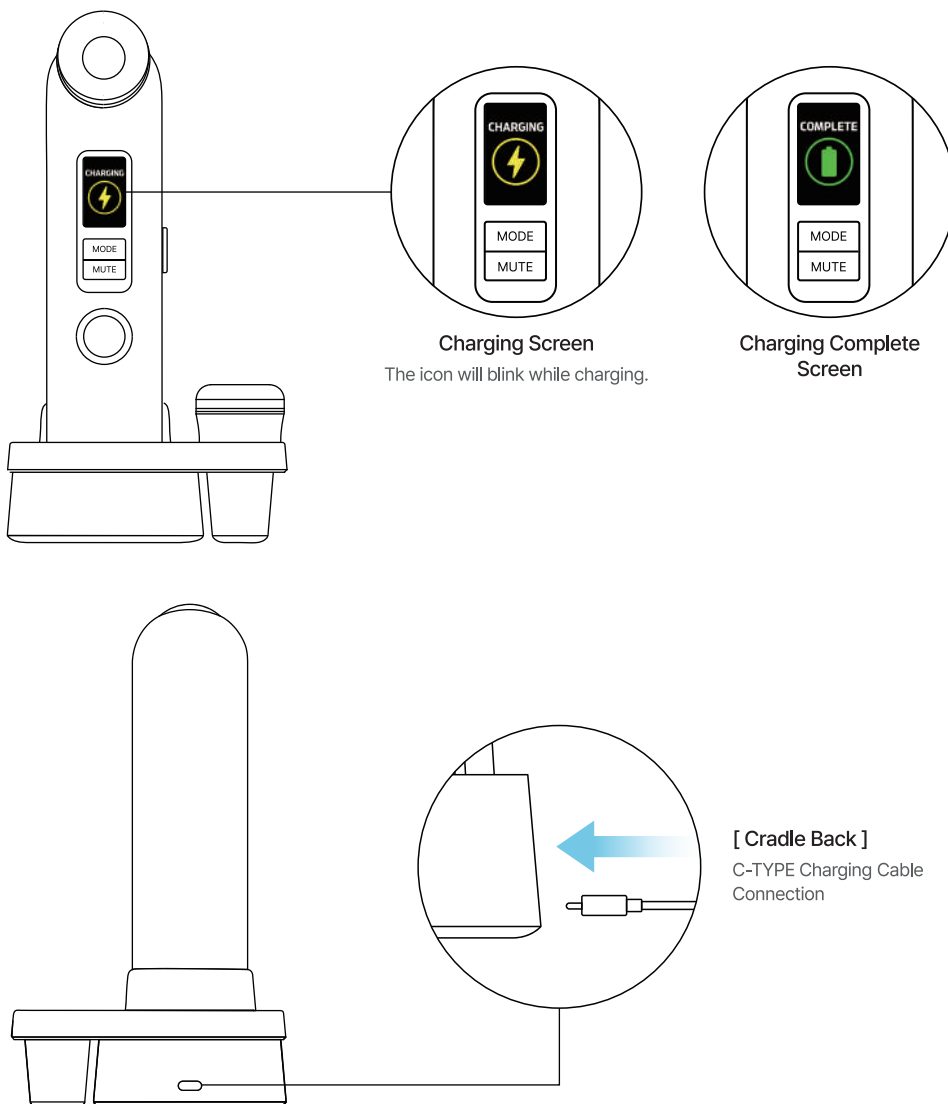


3) SONO Usage Method

30

Product Main Function

■ High-density multi-ultrasound with 2 depths (3/10MHz)

By cross output of ultrasound at 2 frequencies(wavelengths) of 3/10 MHz, energy penetrates from the epidermis to the SMAS layer, helping to improve skin condition.

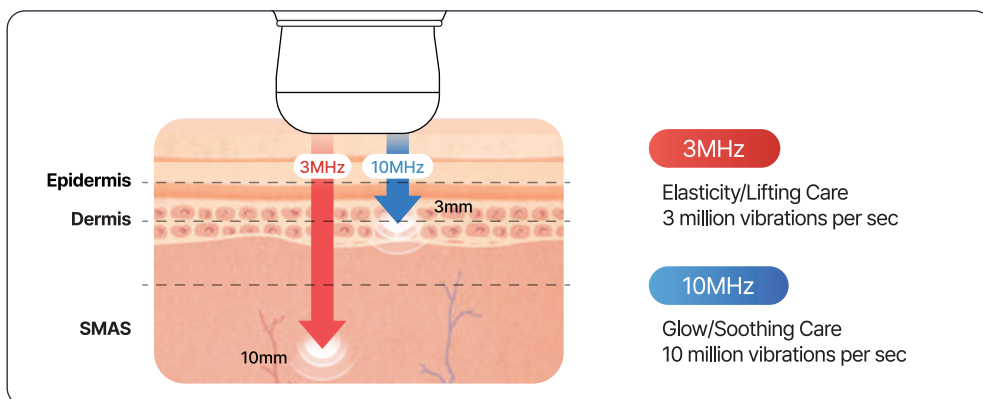


■ Gyro Sensor that Detects Movements

The device detects movement to prevent malfunction and overheating.


If no movement is detected, the output automatically stops.

When the gyro sensor is activated and movement is detected again, the output resumes automatically.

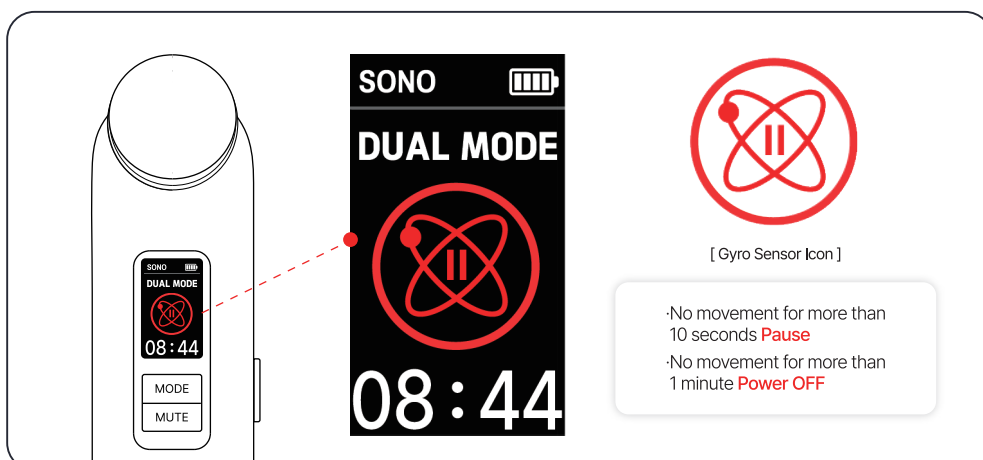

Product Operation Method

- Power ON** Press the power button for more than 3 seconds and the device will turn on with the screen.
- Power OFF** When it turns on, press the power button for more than 3 seconds to turn it off.
- Mode Settings** Set the mode by briefly pressing the MODE button. (DUAL MODE / LIFTING / MOIST)
- Power Settings** Set the power by briefly pressing the power button. (Level 1~5)
- Silent Settings** Press the MUTE button briefly to enable mute mode.
- Screen Information**

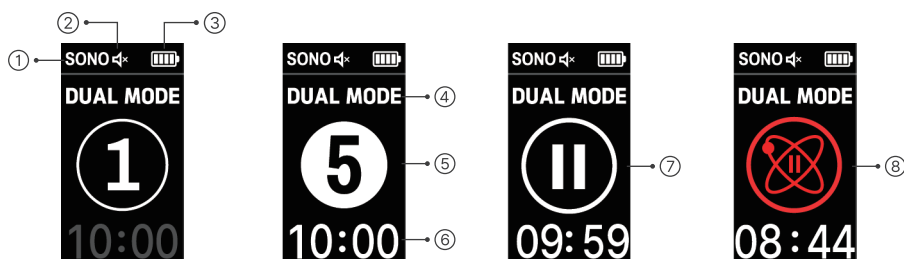

- ① **SONO** The currently attached head is displayed.
- ② **Mute** The icon will be displayed when the device is in mute mode.
- ③ **Battery** The battery status is displayed.
- ④ **Mode** The currently set mode is displayed.
- ⑤ **Energy Level** The currently set power is displayed.
- ⑥ **Treatment Time** The device is set to a default of 10 minutes, and during operation the remaining time will be displayed.
- ⑦ **Pause** It is possible to pause during use with the start/stop button.
- ⑧ **Gyro Sensor** If no movement is detected for more than 10 seconds during operation, the device will enter a pause state where energy output automatically stops.

3) SONO Usage Method

32

*Reference for lymph nodes locations



The Importance of Lymph Care

The lymph system is an important part of our body, playing a crucial role in blood circulation, fluid exchange, and the immune system. When the circulation of lymph is blocked, it can lead to side effects such as swelling, fatigue, weakened immunity, worsening skin condition, and weight gain. Therefore, it is important to first relieve blocked lymph before starting any skin care.

*Must Read Before Use

Please make sure to carefully read the precautions in this manual before use.

 Prohibited Areas of Use Pupils(eyes), around the eyes	 Prohibited Areas of Use The area where the thyroid gland passes	 Usage Cycle Can be used daily, recommended for 10 minutes a day
--	--	---

| Mode (Select the frequency)

It can be used in 3 modes.

SONO DUAL MODE 5 09:59 DUAL MODE 3/10MHz 3-10mm penetration - Lymph circulation / skin regeneration - 2 ultrasound frequencies for the optimal combination - Cross output (recommended use)	SONO LIFTING 5 09:59 LIFTING 3MHz 10mm penetration - 3 million vibrations per second - elasticity / lifting / anti-aging / troubles / swelling care	SONO MOIST 5 09:59 MOIST 10MHz 3mm penetration - 10 million vibrations per second - glow / soothing / moisturizing / pores / redness / brightening care
---	--	---

| Be careful of damage to the SONO head

Disinfect the SONO head with alcohol after use. Use a soft cloth for cleaning.

*Using a rough cloth or wiping the surface too hard may cause scratches or damage. Product damage due to user negligence is not eligible for exchange or refund.

3) SONO Usage Method

34

Lymph Care Method

■ Recommended Usage (Based on full face)

DUAL MODE (3/10MHz)

Level 3~5

10~15 minutes

Can be used daily



STEP 01

Apply Gel

Apply a sufficient amount of gel to the area to be treated.




STEP 02

Select Mode-Power

Select the mode(frequency) and power before starting to use the device.

*Refer to the recommended usage.

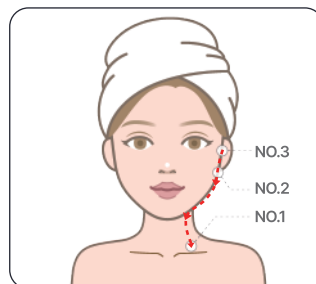

STEP 03

Lymph Nodes Treatment


Rub each lymph node and surrounding area in the order of 1>2>3 for 30 seconds each. Then, handle by sweeping in the direction of 3>2>1.



STEP 04

Full Face Care


Rub the entire face starting from the chin area then cheeks>eyes>forehead for about 6 minutes.

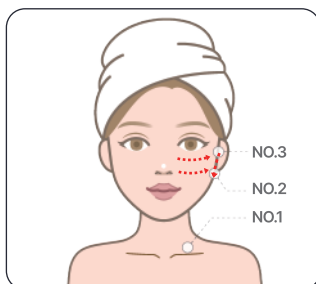

STEP 05

Lymph Nodes Treatment

Handle by sweeping from the center of the face towards lymph nodes 3 and 2 for about 1 minute.

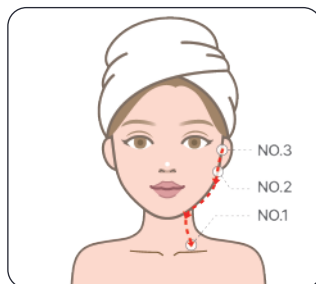

STEP 06

Lymph Nodes Treatment

Handle by sweeping towards lymph nodes 3>2>1 for about 1 minute.

3) SONO Usage Method

35

Ampoule Absorption Method


■ Recommended Usage (Based on full face)

DUAL MODE (3/10MHz)

Level 3~5

About 5 minutes


Can be used daily

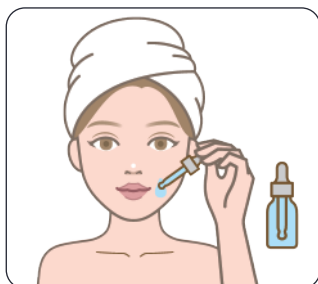

STEP 01

Apply the Ampoule

Apply a liquid type of the ampoule to the skin.

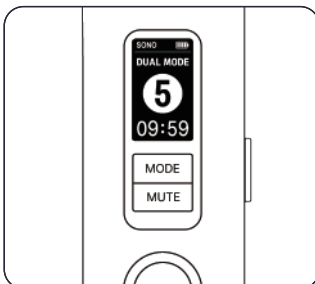

STEP 02

Select Mode-Power

Select the mode and power before starting to use the device.
*Refer to the recommended usage.




STEP 03

Start Treatment

Rub the entire face and let the ampoule absorb as much as you want.

Usage TIP!

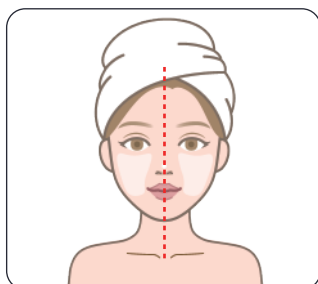

Section Care

Divide the face into half to prevent the gel or ampoule from drying out during treatment.



Use Special Gel

For synergistic effects with ultrasound penetration (lymph care), try using it together with Ovaco Timeslip Water Max Cream.




Use of Ampoule-Serum-Mask Pack

Please use a liquid type formula, not a cream. Please use a sheet mask, not a jelly type mask.

Product Main Function

■ High-Intensity Focused Ultrasound (3.0mm)

High-intensity focused ultrasound energy precisely targets a depth of 3.0mm within the skin, helping to stimulate skin transformation without damaging the surface.


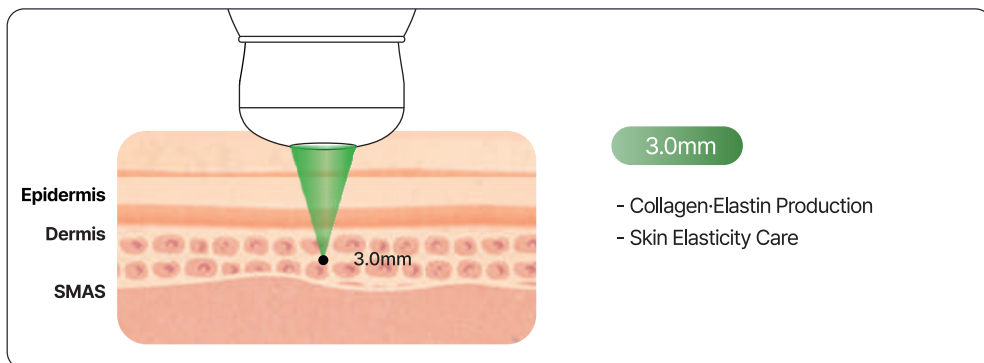


■ Gyro Sensor that Detects Movements

The device detects movement to prevent malfunction and overheating.

If no movement is detected, the output automatically stops.

When the gyro sensor is activated and movement is detected again, the output resumes automatically.

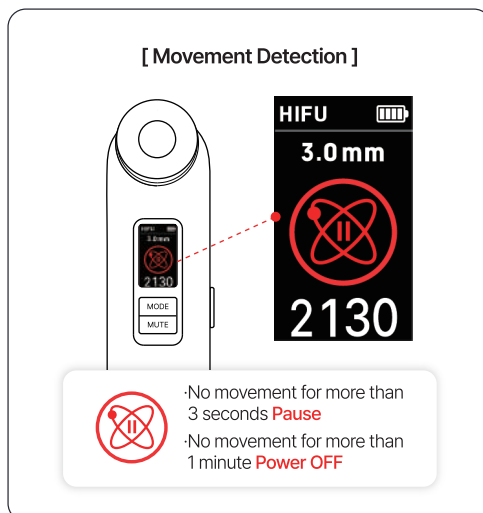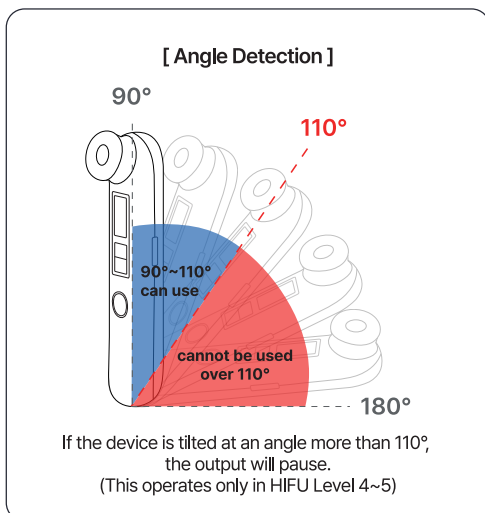

Product Operation Method

Power ON

Press the power button for more than 3 seconds and the device will turn on with the screen.

Power OFF

When it turns on, press the power button for more than 3 seconds to turn it off.


Power Settings

Set the power by briefly pressing the power button. (Level 1~5)

Silent Settings

Press the MUTE button briefly to enable mute mode.

Screen Information

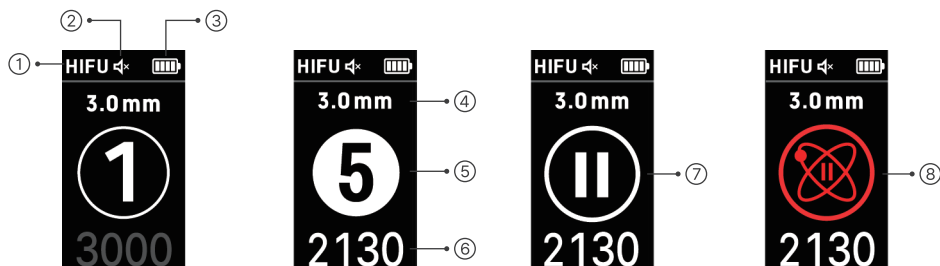

- | | |
|-----------------------|---|
| ① HIFU | The currently attached head is displayed. |
| ② Mute | The icon will be displayed when the device is in mute mode. |
| ③ Battery | The battery status is displayed. |
| ④ Mode | The currently set mode is displayed. |
| ⑤ Energy Level | The currently set power is displayed.
Levels 1~3 : 3Hz (3 shots per second) / Levels 4~5 : 5Hz (5 shots per second) |
| ⑥ Shot | The default set to 3,000 shots, and during operation the remaining shots will be displayed. |
| ⑦ Pause | It is possible to pause during use with the start/stop button. |
| ⑧ Gyro Sensor | If no movement is detected for more than 3 seconds during operation, or if the device tilts back more than 110° (Level 4~5), the device will enter a pause state where energy output automatically stops. |

4) HIFU Usage Method

38

***Area Specific Shots**
(Half of face 1,500 shots)



Forehead/Eye Corner
350 shots

Cheeks/Chin 800 shots



Below Jaw/Neck
350 shots

4) HIFU Usage Method


39

*Must Read Before Use

Please make sure to carefully read the precautions in this manual before use.

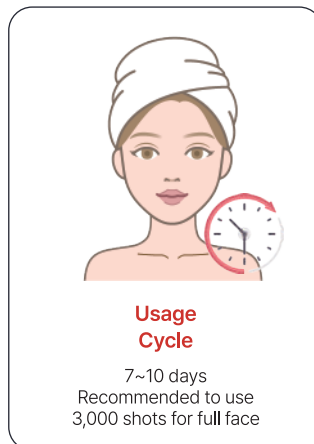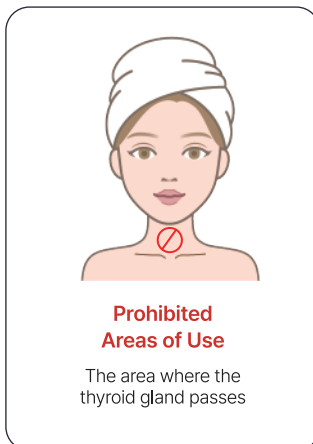

| HIFU Head Care Method

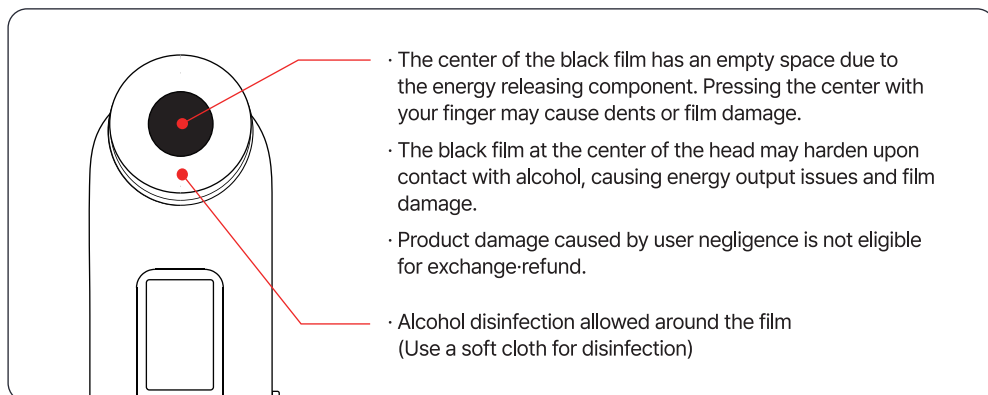

| HIFU Cartridge 300,000 Shots Output Guaranteed

- Recommended usage : 144,000 shots per year / 288,000 shots over 2 years
- Output of 300,000 shots guaranteed within 2 years from the date of manufacture for products that have completed product registration
(Free A/S service is not available in cases of element damage, product breakage, scratches, coolant leakage, or discoloration.)

4) HIFU Usage Method

40

Lifting Care Method

■ Recommended Usage (Based on full face)

3.0 Mode

Level 3~5

3,000 shots


7~10 day cycle



STEP 01

Apply Gel


Apply a sufficient amount of gel to the area to be treated.

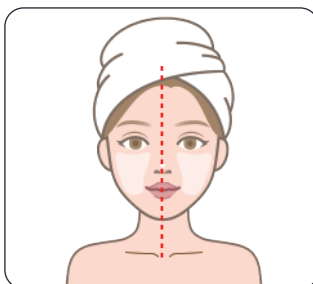

STEP 02

Section Care

To prevent the gel from drying, try dividing your face in half for treatment.

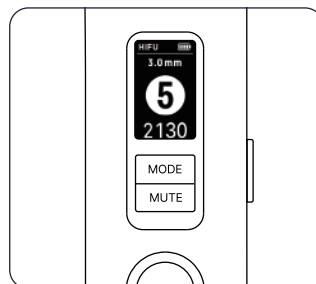


STEP 03

Select Power

Select the power before starting to use the device.


*Refer to the recommended usage.

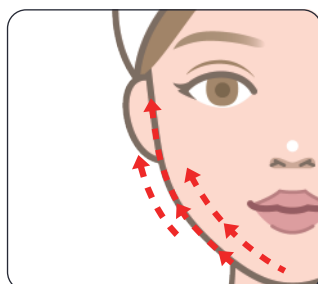

STEP 04

Lifting Care


There is no fixed order for the treatment areas. For effective contour lifting, handle the device as if you are lifting the skin.

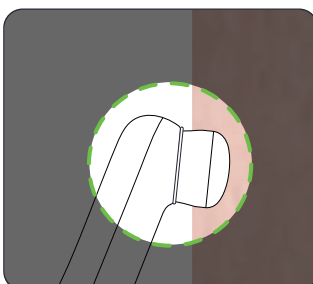

PROPER USE

Head Contact O

The head must be in full contact with the skin without any space.

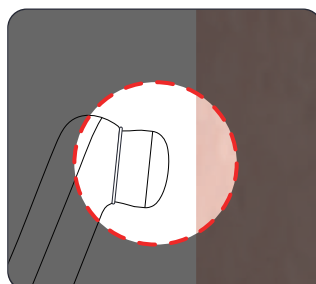

IMPROPER USE

Head Contact X


If there is space between the head and the skin, you may feel discomfort, and energy may not penetrate effectively.

4) HIFU Usage Method

41

Handling Usage Tip

■ Rubbing Handling

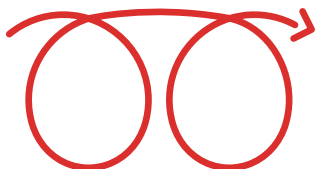

This method involves moving in a circular motion. Start from the center of the face and move outward, lifting from the bottom to the top.


It's a suitable method for first-time users who are not familiar with using HIFU.

■ Linear Handling




This method involves moving in a straight motion. Start from the center of the face and move outward, lifting from the bottom to the top.

It's a suitable method for users who are familiar with using HIFU for homecare.

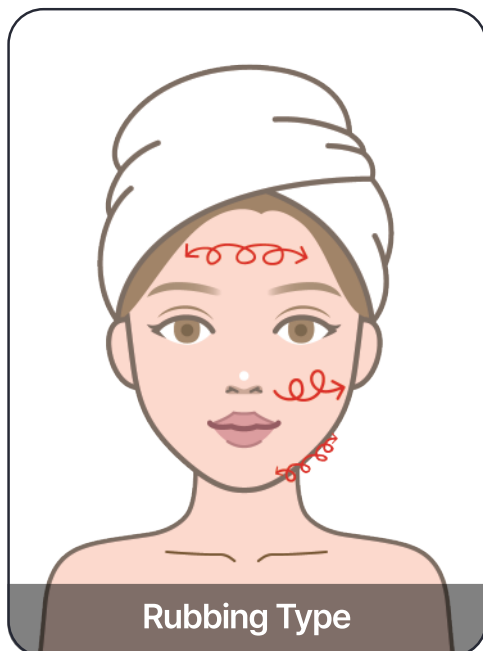

Rubbing Type

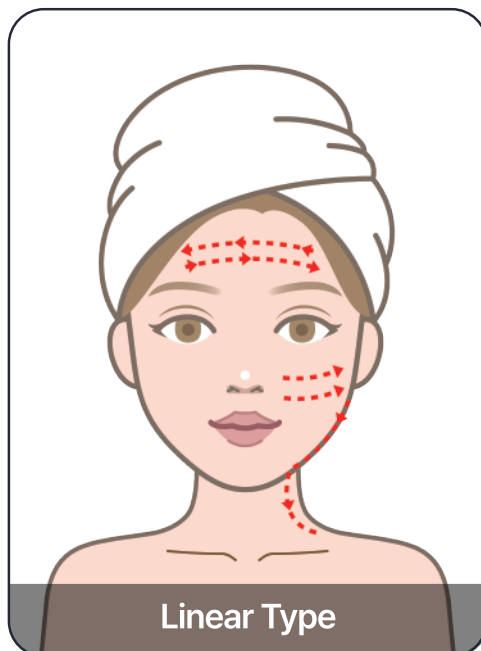

Linear Type

Warranty Regulations

※ Note

The product warranty period is 1 year from the date of manufacture.

These warranty regulations apply only to the scope of the warranty period and conditions provided by our company. For A/S service, please contact the place of purchase or the A/S center.

※ Within the warranty period

Free maintenance within 1 year from the day of production (product registration is required)

During the warranty period, if a problem occurs because of the production process of the equipment, free A/S service can be provided.

※ After the warranty period

Paid A/S service can be provided for products that are not discontinued after the warranty period, and A/S service may not be available for discontinued products.

※ Exceptions

In the following cases, paid A/S service is provided even within the warranty period.

- 1) Lost or stolen
- 2) Damage caused by user negligence
- 3) User arbitrarily modified, disassembled, or repaired by unauthorized personnel
- 4) Improper storage condition
- 5) Damage caused by natural disasters (fire, salt damage, flood damage, etc.)

Product Warranty

Product Warranty

Product Name : Lymfit.Plus

Thank you for purchasing JLU products.

This product is a genuine product manufactured by the headquarters and shipped after passing strict quality control and inspection processes.

Free service within the warranty period is available if you complete the product registration through the homepage.

(S/N : Separate notation on the product)



Please access the QR code and proceed with product registration.

A/S Center : +82 10-4505-7989

JLU Co., Ltd.

410~412 office, Busan Knowledge Industry Center, 111, Hyoyeol-ro, Buk-gu, Busan, Republic of Korea

 +82 51-331-7185

 +82 10-3752-8240 (EN line)

 +82 10-4505-7989

 jlukorea@naver.com

 jlu_korea_official

 JLU Co.,Ltd

 jlu_korea_official

 blog.naver.com/jlukorea

FCC Compliance Statement

This device complies with part 15 of the FCC rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Interference Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to correct the interference by one of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC Caution

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this equipment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

FCC Radiation Exposure Statement

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 5mm between the radiator and your body. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.





Lymfit.Plus