



Etee Toothbrush Instruction Manual

Patented semiconductor chip microcurrent and the synergistic effect of ultrasonic vibration

• PRECAUTIONS

- ① Use it after familiarizing yourself with the operation procedure and safety instructions. Be careful not to injure yourself due to carelessness.
- ② Use of this product in a way other than the manual is prohibited. The manufacturer and seller are not responsible for accidents such as damage and electric shock caused by arbitrary disassembly or change of product parts.
- ③ Do not arbitrarily repair this product and contact the place of purchase and receive A/S. Product warranty due to random repairs will not be applied.
- ④ If the power cord or plug is damaged, contact after-sales service and take actions.
- ⑤ Do not put the product in a microwave oven, high heat, or fire.
- ⑥ Do not charge, use or leave at high temperatures.

• WARNING

- ① Do not use this product for purposes other than cleaning teeth, gums and tongue.
- ② Store this product in a safe place so that it cannot be dropped or submerged in water.
- ③ Do not disconnect the power plug with wet hands.
- ④ Before connecting the power plug, check that the voltage of the charging part matches the operating voltage before use.
- ⑤ Do not damage, damage or modify the power cord. Do not forcefully pull, bend, or twist the power cord.
- ⑥ Keep this product (body, charger) away from fire.
- ⑦ Do not store in high temperature and high humidity and avoid direct sunlight.
- ⑧ Be careful not to drop the charging cradle into water.
- ⑨ If you do not use it for a long time, turn off the power and store it.

• Declaration of Conformity

SMLNS Co., Ltd. declares that this device complies with the essential requirements and relevant provisions of Directive 2014/53/EU.
Go to <https://www.smlns.com/Doc> to view the Declaration of Conformity.

Note : This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Modifications not expressly approved by the manufacturer could void your authority to operate the equipment under FCC rules.

FCC Part 15.19

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Part 15.21

Any changes or modifications (including the antennas) to this device that are not expressly approved by the manufacturer may void the user's authority to operate the equipment.

FCC RF Radiation Exposure Statement: This equipment complies with FCC RF Radiation exposure limits set forth for an uncontrolled environment. This device and its antenna must not be co-located or operating in conjunction with any other antenna or transmitter. This equipment should be installed and operated with a minimum distance of 0 cm between the radiator and your body.

• PRECAUTIONS FOR USE

- ① Use this product only with the dedicated toothbrush head and charger provided together.
- ② For patients with severe periodontal disease, use it after consulting with a specialist.
- ③ If other symptoms such as tooth shaking in the mouth or continuous bleeding of the gums occur during use, avoid using it and consult with a specialist before using it.
- ④ This product may not be suitable for children to use alone. Children and those who are physically uncomfortable using it by themselves must use it under the guidance and supervision of a guardian.
- ⑤ This product cannot be used as a toy for children and infants.
- ⑥ The toothbrush is a personal hygiene tool, so do not share it with others.
- ⑦ Regularly wipe the charger head with a clean, dry cloth. Otherwise, the charging effect may be reduced.



• READ BEFORE USE

► Power on

It works by briefly pressing the power button once.

► Transition mode


After operation, immediately press the power button briefly to switch modes.

► Smart Memory

If you use one mode for more than 5 seconds, the system automatically remembers the mode and saves it as the initial mode in the system.


► After 5 seconds of operation, press and hold the power button to turn off the power, or after 2 minutes of operation, the power turns off automatically.

• PRODUCT DESCRIPTION & COMPOSITION

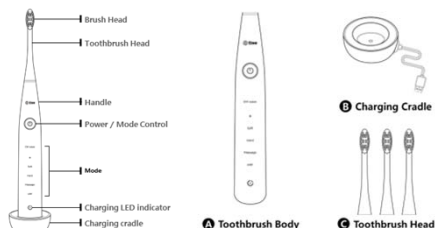

• HOW TO CHANGE MODE

After the power is turned on, press the switch again to select the desired mode. You can change to another mode by pressing the button even while in use.

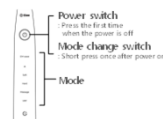

► Mode detailed function

EM wave	A microcurrent from the electrodes manages plaque and inflammation where a toothbrush cannot reach.
EM wave + Soft	It is a soft mode function that enables delicate cleaning with microcurrent.
EM wave + Hard	It is a hard mode function that enables powerful cleaning with microcurrent.
EM wave + Massage	It is a mode function that massages the gums with microcurrent.
EM wave + User	It is a function that allows the user to directly adjust the intensity of sonic vibration along with the microcurrent. * Usage TIP: Press the button for more than 3 seconds to select the intensity you want

He 4 sections of the oral cavity can be brushed evenly.

※ Do not rub it on your teeth! If you hold the bristles between the gums and teeth at an angle without scrubbing, the microcurrent removes plaque to the places where the toothbrush cannot reach, and the non-irritating sound wave vibration effectively removes large debris.

• Recommended brush head replacement cycle to keep your mouth clean and healthy, dentists recommend replacing your toothbrush head every 2-3 months.

• CHARGE

- ① When using an electric toothbrush for the first time, if there is remaining battery, charge it after using up the battery. From the next use, it can be charged whenever necessary.
- ② If the remaining battery charge is less than 10%, the charging indicator turns red and a charging alarm appears.
- ③ You can charge it by placing it on the charging cradle.
- ④ The rechargeable battery indicator light blinks during charging and turns off when fully charged. Battery charge time is 12 hours

- ⑤ While charging, the electric toothbrush does not work even if you press the button.
- ⑥ If you place the electric toothbrush on the charging cradle while it is working, it will stop working and switch to the charging state.

• CHARGER PRECAUTIONS

Please use only the specified adapter. In case of damage to the charger, please contact the A/S center.

• Precautions when disposing of batteries

- ① When disposing of this product, remove the battery through the bottom and dispose of it at a designated disposal site.
- ② Before removing the battery, the power must be turned off, and the remaining power of the battery must be exhausted before removing it.
- ③ After removing the battery, do not connect the product to the power line.

• SELF-DIAGNOSIS

Issue	Cause & Solution
Weak vibration	It's very likely that the bristles were pressed against the teeth. · Solution: Gently place the bristles on the teeth.
· Work too short · Not working	Most likely the battery is low or not charging properly. · Solution: Please charge enough
Other issue	Please contact the A/S center

• Precautions for A/S

- ① Warranty period: 1 year from the date of purchase
- ② Do not separate the product arbitrarily. AS processing is not possible
- ③ If you want free repair or exchange of the product within the warranty period, you must present data related to purchase details. (Related documents such as purchase receipts, etc.)
- ④ Purchase the toothbrush head for a fee)
- ⑤ Intentional or negligent damage, etc., and if purchase confirmation is not possible, free repair is difficult.
- ⑥ It may take about 2 weeks as it is carried out after confirming the collection of the product at A/S.

• RF function information

Function	TX Frequency	Power
WPC	190 kHz, 210 kHz (± 5kHz)	5 W