


zoesmiien

USER MANUAL

Eye Massager

FJ256B

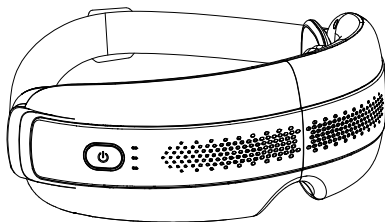

*

Thank you for purchasing this product.

Before use, please carefully read this user manual and pay special attention to safety precautions to ensure correct operation and use.

Please keep this applicable manual properly for future reference.

CATALOGUE

Safety precautions	1
Component names and functions	2
Product Structure	2
Control panel	3
Technical parameters	4
Instructions for use	5
Charging	6
Maintenance and Storage	7
Product component recycling instructions	8
Troubleshooting	9






➤➤➤ SAFETY PRECAUTIONS

zoesmlen

- ❶ Children or individuals with disabilities need special supervision to use this device;
- ❷ Individuals with eye diseases or injuries, as well as those who have undergone eye surgery, are prohibited from using it;
- ❸ Those who have medical electronic devices (such as pacemakers) installed in their bodies, malignant tumors, and heart diseases Patients with diseases, severe hypertension, and other physical discomforts are prohibited from using it;
- ❹ If you experience or feel any abnormalities in your body during use, stop using it immediately and consult a doctor;
- ❺ If the instrument or charging cable is damaged, do not use it;
- ❻ Before charging, please ensure that the charger and product charging socket are dry;
- ❼ Do not use the product or charger for charging when it may be submerged in water;
- ❽ Do not use other therapeutic devices at the same time as using the unit;
- ❾ Do not use for therapeutic purposes;
- ❿ Do not drag the product on the ground, or pull the power cord to move the product to avoid damage products;
- ⓫ Do not disassemble the product and accessories without permission, or disassemble, repair and modify the product to avoid the possibility of machine failure;
- ⓬ Do not use a power supply other than the rated voltage to avoid causing fire or damaging the product;
- ⓭ When unplugging the power cord, hold the plug of the cord and not the cord;
- ⓮ Do not damage, excessively bend, pull, twist, or knot the power cord.

>>>> COMPONENT NAMES AND FUNCTIONS


- Product Structure

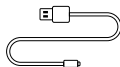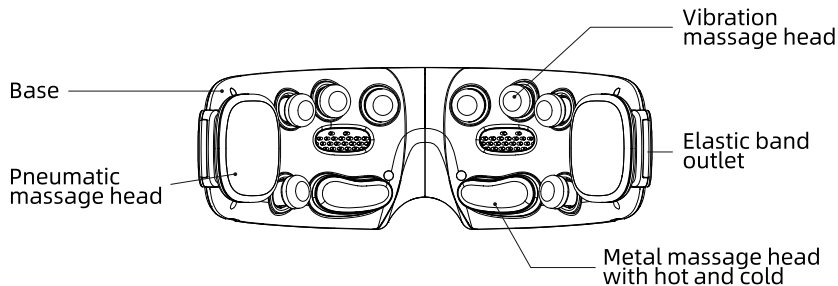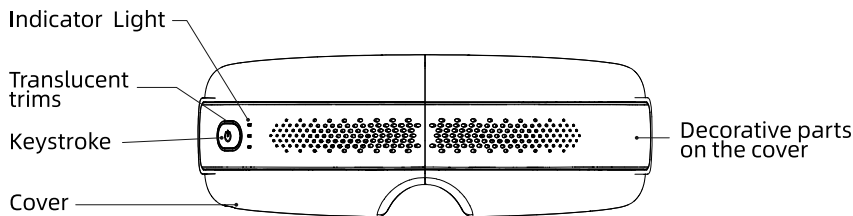

Type-c charging cable

COMPONENT NAMES AND FUNCTIONS

- Control panel

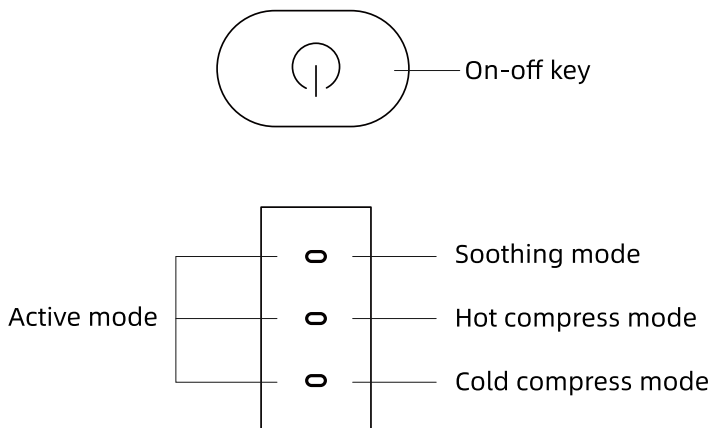

► Note

The default working time of soothing mode, hot compress mode, and active mode is 15 minutes, and the default working time of the cold compress mode is 6 minutes, if you want to continue using it, press the on/off mode button again.

»»»» TECHNICAL PARAMETERS

zoesmlen

Product name: Eye massager

Product model: FJ256B

Rated voltage: 3.7V 

Rated power: 5W

Rated frequency: 50/60Hz

Safety standards: III

- ▶ ① Press and hold the key product to turn on the power-on, and the voice broadcasts "welcome to the smart portable eye massager".

When the product is turned on, short press the on/off key to switch the mode of the product, and cycle in turn: soothing mode/hot compress mode/cold compress mode/active mode.

② Double-click the switch to turn on silent mode, and the voice prompt "silent mode is enabled" (local music and Bluetooth music are turned off, the voice broadcast is normal, and the music can be turned on when Bluetooth music is turned on manually). Double-click again to turn off the silent mode, and the voice prompt "silent mode is off".

③ On/off key: Short press: Mode switching, in order:

Mode 1: Soothing mode (vibration + air pressure), after switching to this mode, it will run automatically for 15 minutes.

Mode 2: Hot compress mode (vibration + air pressure + hot compress), after switching to this mode, it will run automatically for 15 minutes.

Mode 3: Cold compress mode (vibration + air pressure + cold compress), after switching to this mode, it will run automatically for another 6 minutes.

Mode 4: Active mode (vibration + air pressure + hot and cold alternation), after switching to this mode, it will run automatically for another 15 minutes.

- ▶ Bluetooth: After synchronously opening the product and turning on the Bluetooth of the phone, find EYE-MASSAGER Bluetooth and click on the connection to broadcasts "bluetooth connection successful". Third-party music can be played on the phone; After Bluetooth connection once, open the product and mobile phone Bluetooth synchronously again, Bluetooth will be connected automatically.

- ▶ Shutdown: Press and hold the power button for 2 seconds, the voice broadcasts "thanks for using bye", and the power button indicator light off.

- ▶ Ambient lighting: green, blue, purple, red; Mode indicator: white ;

① The product is not connected to Bluetooth when powered on, and the ambient light is blue;

② After the product is connected to Bluetooth, the ambient light is always purple;

③ The first indicator of mode 1 is on, the second indicator of mode 2 is on, the third indicator of mode 3 is on, and the three lights of mode 4 are all on.




- ▶ Please wear the product on your eyes before turning it on, and those who wear contact lenses, hair accessories and other items should be removed before use.

CHARGING

zoesmilien

- ▶ Rotate the product to the angle shown below to expose the Type-C charging port, and then plug the charging plug into the charging port of the product.
- ▶ Plug the charging plug into a 100-240V~, 50/60Hz power outlet.
- ▶ The ambient light is red when the product is charging, and the ambient light is always green when fully charged.
- ▶ It takes about 3.5 hours to fully charge the battery, and the charging ambient light is always green when fully charged.
- ▶ Low battery reminder: When the battery level is lower than 20%, the ambient light will flash red slowly, and the voice will be broadcast "the battery is low please charge it".

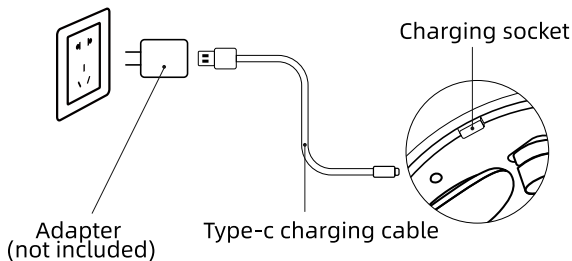

NOTE

- ▶ During charging, the product does not work.
- ▶ Please match the 5V1A adapter by yourself.
- ▶ Do not squeeze the power cord to avoid damage.

»»»» MAINTENANCE AND STORAGE

Maintenance:

If the product is dirty, wipe it with a dry, soft cloth.

(If the product is particularly dirty, wipe it with a cotton or cotton cloth soaked in alcohol).

Note:

It is forbidden to clean this product with water or corrosive liquid, and professional technicians should be called in when the product fails. Do not disassemble the product by yourself.

Storage:

After use, please put the product in a safe place and put it in the box with a bag.

Please do not place heavy objects on the product, and protect it from moisture, rodents, and stains.

Note:

Do not expose this product to high temperatures or water.


Do not give this product to children for play.

>>>> PRODUCT COMPONENT RECYCLING INSTRUCTIONS

zoesmlen

- Recycling indicative instructions

In order to protect the environment, it is your responsibility as a consumer to dispose of this product or its parts when they are discarded. It is separated from household waste and sent to a qualified recycling site, which is divided into national standards. Relevant provisions for classification, dismantling, recycling and other treatments.

PART NAME	Hazardous substance					
	lead Pb	mercury Hg	cadmium Cd	hexavalent chrome Cr(VI)	polybrominated biphenyl PBB	polybrominated biphenyl ethers PBDE
Hardware	×	○	○	○	○	○
Circuit boards	×	○	○	○	○	○
Plastic parts	○	○	○	○	○	○
Wire	○	○	○	○	○	○
<p>○ It means that the content of the hazardous substance in all homogeneous materials of the part is below the limit requirements specified in GB/T 26572.</p> <p>×</p> <p>It means that the content of the hazardous substance in at least one of the homogeneous materials of the part exceeds the limit requirements specified in GB/T 26572.</p>						
<div>  <p>Environment-friendly use period mark</p> </div> <p>The "environment-friendly use period" of this product is 10 years, and the company formulates the period according to the design characteristics of the product. The "Environment-friendly use period" is only valid when the product is used under normal conditions as described in the instruction manual of this product.</p>						

► NOTE

- Other parts contained in this product and their homogeneous materials are not listed in the above table and do not contain hazardous substances.
- All parts and homogeneous materials containing hazardous substances in the above table comply with the strict requirements of China RoHS on the restriction of the use of hazardous substances, so please feel free to use them.

>>>> TROUBLESHOOTING

Problems	Possible causes	Solutions
Won't turn on	On/Off is not on	Turn on/off eg long press to open
	The battery runs out	charge
	The time setting is off	Turn on/off again
Unable to charge	The plug is not plugged in	Plug it in
	The adapter is not applicable	Select the 5V1A adapter
The indicator light is dim or flashing	The battery is low	Please charge
The air pressure is too low		

► NOTE

If there are other abnormal phenomena, please contact the dealer for maintenance, please do not disassemble by yourself.

*
Product specifications, appearance, design, etc., will be improved without prior notice.

▶▶▶▶ FCC WARNING

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

NOTE: This device and its antenna(s) must not be co-located or operation in conjunction with any other antenna or transmitter.

▶ RF exposure

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.