

CERAGEM MASSAGE CHAIR USER MANUAL

CGM-EMCG 2403

Thank you for purchasing the CERAGEM PAUSE massage chair.
Please read this instruction manual carefully and use it correctly before using it.
Please always keep it where you can see it after reading it.

Product warranty is included in the instruction manual.
This product is a home massage chair,
so please do not use it for business or industrial purposes.
The headquarters may improve the design and specifications
of this product without further notice of circumstances.

CERAGEM LIFESTYLE MASSAGE CHAIR USER MANUAL

PAUSE M7 USER MANUAL

CONTENTS

1. safety precautions
2. name of each part
3. Remote control description

The following safety precautions are used to ensure that the product is safe and accurate
It is intended to prevent unexpected risks or injuries by using it.

When installing

Never renovate, disassemble, or repair except for service technicians
→ It can cause fire and electric shock.
Use it on a flat surface
→ It can cause a failure or accident.

Power related

If it doesn't move or if something is wrong, immediately remove the power plug from the outlet Remove and ask the customer center for inspection and repair
→ There is a risk of fire or electric shock due to short circuits,

Please guide the child not to touch the power plug.
→ It can cause electric shock or injury.

Avoid undue bending of the power cord or being crushed by heavy objects
→ It is a risk of fire and electric shock.

Do not use multiple devices at the same time when using multiple outlets
→ There is a risk of fire due to abnormal heating of the outlet.

When removing the plug, hold the plug and disconnect it from the outlet.
→ Disconnecting the power line poses a risk of fire and electric shock.

Do not touch the power plug with wet hands,
Disconnect the power plug from the outlet when cleaning or maintaining the
→ There is a risk of electric shock.

Be sure to use a grounded outlet when using a multi-tap.
→ Using a non-ground outlet can cause electric shock and injury.

In the event of a power outage, immediately turn off the power switch and unplug the power plug from the outlet.
→ It can cause electric shock and fire.

Make sure that the power cord wires do not get stuck between the products during operation of the product
→ Care must be taken when using because short circuits may occur due to damage to the sheath.

1. Name of each part of massage chair



1. Head cushion

2. Bluetooth speakers

3. Touch remote control

4. Shoulder massage unit

5. Back pad


6. Arm massage unit

7. Back massage unit

8. Calf massage unit

9. Foot massage unit

2. Built-in Quick Controller

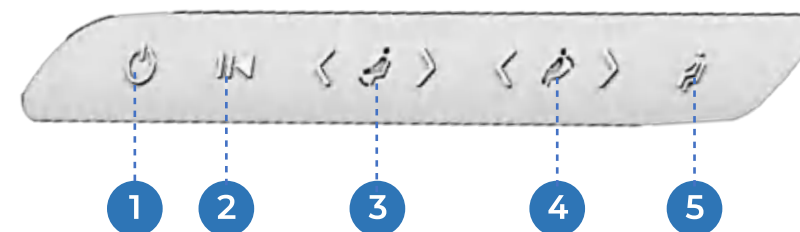

1. POWER On/Off

- you can turn the massage chair on and off by pressing the button for about 3 seconds.

2. PAUSE

- Press the button to temporarily stop or re-activate all functions of the massage chair.

3. Detailed leg angle adjustment

- < > You can use the (left and right) key to adjust the foot angle of the massage chair


3. 4. Detailed body angle adjustment

- < > You can use the (left and right) key to adjust the body angle of the massage chair

5. Quick adjustment of body angle


- You can quickly adjust the angle of the body. (1 weightless / 2 weightless / 3 weightless)

1. Start screen

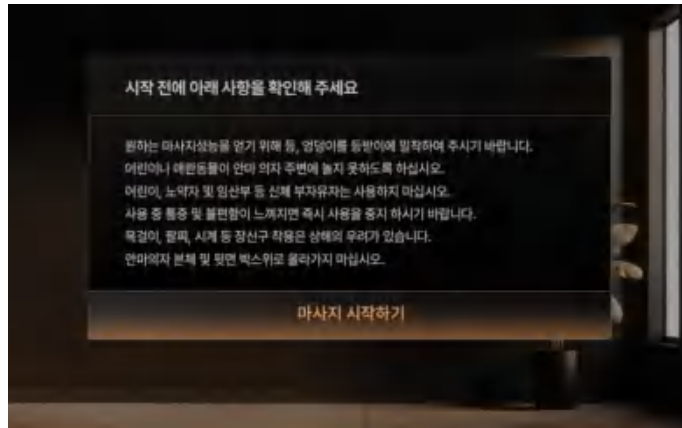


Press and hold the power button of the built-in quick controller for 3 seconds or the power at the top of the remote control. Press the button to turn on the power.

2. Automatic body scan

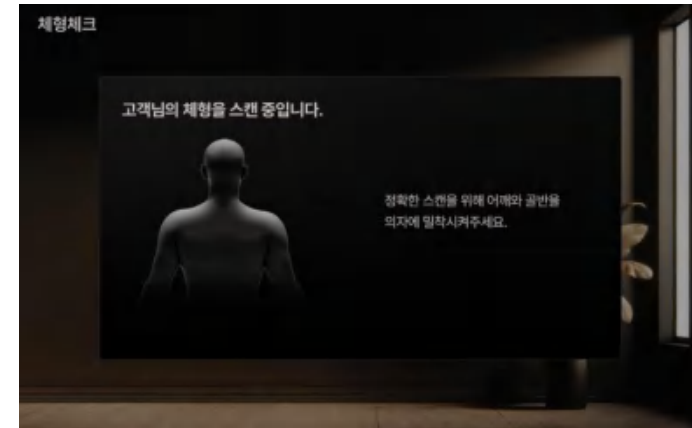

The body shape scan starts automatically after the power is turned on.

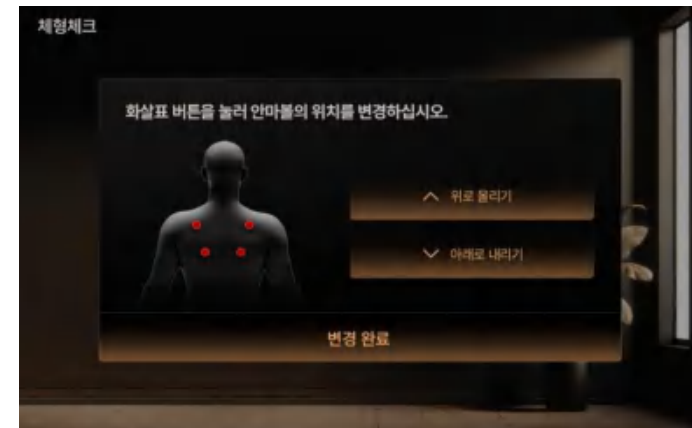

After the automatic body type scan, you can manually adjust the position of the user's shoulder.

3. Home screen



1. Automatic mode

If you choose automatic mode, you can choose 24 massage courses.

2. Manual mode

If you select manual mode, you can set it to receive a massage on the desired area using the desired massage method.

3. Air settings

If you select Air Settings, you can select the preferred Air Massage area and strength You can choose.


4. Device settings

You can set up a thermal / language / screen brightness / Bluetooth connection.

5. Quick button

You can activate the previous screen / main screen / massage chair pause.

4. Automatic mode

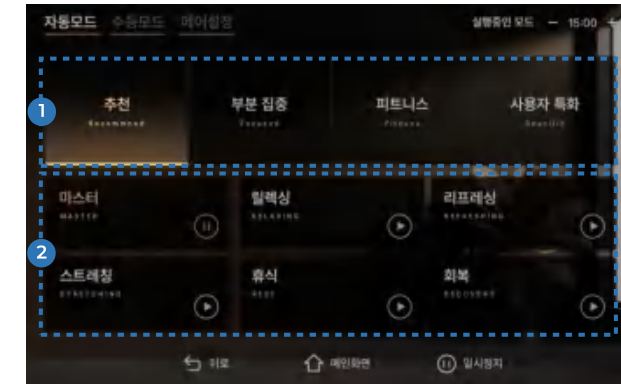

1. Message Theme: Choose four types of message themes (Recommendation / Partial Focus / Fitness / User Specialization)
2. Message Mode: Choose 24 types of message modes to match your message theme



1. Massage Strength: Adjustable massage module upper/lower strength (step 5)
2. Calf Roller: You can adjust the strength of the calf roller. (Step 3)
3. Foot Roller: You can adjust the strength of the foot roller. (Step 3)
4. Weightless Angle: You can adjust 3 weightless angles
5. Massage chair angle: You can adjust the upper/lower body angle and calf length of the massage chair.
6. Thermal function: You can control the heat of the massage module and the body of the massage chair. (Step 3)

5. Manual mode



1. Massage module: You can select one of the upper/lower massage modules.
2. Massage Method: You can choose a total of 6 massage methods.
3. Massage Strength: You can adjust the massage intensity of the massage module (step 5)
4. Massage Speed: You can adjust the massage speed of the massage module (step 3)
5. Calf Roller: You can adjust the strength of the calf roller. (Step 3)
6. Foot Roller: You can adjust the strength of the foot roller. (Step 3)
7. Massage section: You can select the massage section of the massage module.
8. Massage chair angle: You can adjust the upper/lower body angle and calf length of the massage chair.
9. Massage Width: You can adjust the massage width of the massage module.
10. Massage Module Location: You can manually adjust the massage module to the desired location.
11. Thermal function: You can control the heat of the massage module and the body of the massage chair. (Step 3)

6. Device Setting



Thermal function control: You can adjust the thermal function ON/OFF and intensity

Language Settings: You can choose between Korean and English.

Brightness/Sound: You can adjust the button sound, screen brightness, and mood lighting functions.

External device: Bluetooth function can be turned on/off.

FCC WARNING

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Note: The Grantee is not responsible for any changes or modifications not expressly approved by the party responsible for compliance. Such modifications could void the user's authority to operate the equipment.

The device has been evaluated to meet general RF exposure requirement.

This equipment complies with FCC's RF radiation exposure limits set forth for an uncontrolled environment. This device and its antenna(s) must not be co-located or in conjunction with any other antenna or transmitter.