



Product List

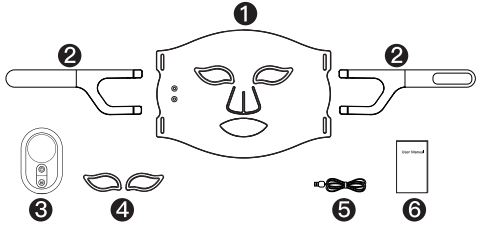


1. LED Light Therapy Mask x1
2. Head Strap x2 (Left & Right)
3. Remote Control x2
4. Blocking Eyeshield x2 (Left & Right)
5. USB Type-C Charging Cable x1
6. Usage Manual x1

Product Structure

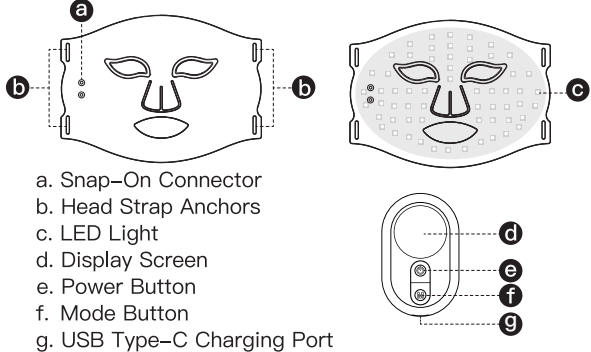

a. Snap-On Connector
b. Head Strap Anchors
c. LED Light
d. Display Screen
e. Power Button
f. Mode Button
g. USB Type-C Charging Port

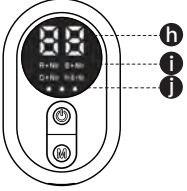

h. Treatment Time
i. Mode
j. Intensity Level

Usage Instructions

1. Power On:
a) The device turns on automatically when magnetically connected.
b) Manual startup: Press and hold the power button for 1.5 seconds to turn on.
c) After powering on, it automatically enters Red Light + Infrared mode with a default timer of 10 minutes.

2. Operation:
a) After startup, briefly press the power button to switch between 4 modes:(R+Nir, B+Nir, O+Nir, R+B+Nir)
b) Briefly press the mode button to adjust intensity (3 levels). Press and hold the mode button to set duration (cycles through 10/15/20/25/30 minutes).

3. Power Off:
a) During operation, press and hold the power button for 1.5 seconds to turn off.

4. Standby & Auto-shutdown:
If the controller is detached during use, the screen will blink and enter standby mode. The device will automatically power off if not reconnected within 1 minute.

5. Timer Reset:
The timer resets after each adjustment.


Mode Settings

Mode	Mode 1	Mode 2	Mode 3	Mode 4
Function	Red + Infrared	Blue + Infrared	Orange + Infrared	Red + Blue + Infrared
Visible Color	Red	Blue	Orange	Purple
Wavelength	630nm+850nm	460nm+850nm	590nm+850nm	630nm+460nm+850nm
Irradiance	19.0 mW/cm²	29.0 mW/cm²	17.5 mW/cm²	35.0 mW/cm²
Time of use	10-20min	10-20min	10-20min	5-15min


The Key to Phototherapy Effectiveness:
Your dose = Irradiance × Time
Different skin goals require tailored dosages. Your skin's tolerance also matters. Our products are designed to be both safe and effective.

Usage Steps


Attentions
1. Please avoid eye exposure to LED Light during use.
2. The frequency we recommend: 3-5 times per week
Use consistently for 4 week as cycle.




1. Clean and dry your face.

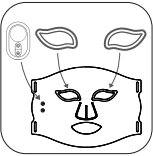


2. Attach eye cover and controller to face mask.




3. Attach straps. Adjust mask for comfort.
Power on and choose mode.



4. After use, remove mask. Resume normal skincare. Keep your skin moist.



LenGen



封底

Product Specifications

Product Name	Wireless LED Light Therapy Mask
Model	SKB-2418
Input	DC 5V 2A
Battery Capacity	1200mAh
Number of LED Beads	272
Wavelengths	Red, Blue, Orange, Near-Infrared
Operating Modes	4
Charging Time	≤2.5h
Product Dimensions	265*195*4 mm
Weight	145g

Power Instructions

1. Charge using a DC 5V/2A adapter and the included USB Type-C cable (adapter not included).
2. Low battery alert: Mode indicator flashes before shutdown.
3. Remove the controller for charging. Device cannot power on while charging.
4. Indicator flashes during charging, stays lit when fully charged.
5. Fully charge before first use or after prolonged storage.
6. Unplug after charging to extend battery life.
7. Turn off when not in use.
8. Recharge every 2 months to maintain battery health.
9. Avoid charging or using in high temperatures.
10. Dispose of used batteries according to local regulations.

Important Notes

1. For personal facial care only. Do not use for purposes not specified in this manual.
2. Stop immediately if skin redness, itching, or discomfort occurs.
3. Who Should NOT Use:
a) Minors under 18.
b) Pregnant, nursing, or women trying to conceive.
c) Who with photosensitivity or taking those medications.
d) Those with eye/genetic skin conditions (eg:lupus, albinism) or facial wounds.
* Consult a doctor if uncertain about suitability.
4. Do Not disassemble, repair, or modify the device. Tampering may void warranty and risk fire, electric shock, or injury.

Warning

1. Keep away from children/pets.
2. Avoid heat sources,open flames, moisture,or direct sunlight.
3. Do not drop or impact the device.
4. Avoid folding the mask to prevent LED damage.
5. Avoid extended direct LED exposure to prevent eye strain or damage.
6. Attach and securely fit the light-blocking eye shield to the mask. Ensure no LED light is visible during use.
7. Handle plugs with dry hands. Insert fully to avoid shock.
8. If damaged or malfunctioning, stop use immediately.Contact customer support.

FCC Warning

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) this device may not cause harmful interference, and
(2) this device must accept any interference received,including interference that may cause undesired operation.
Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
NOTE: This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is noguaranteee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the device off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
-Reorient or relocate the receiving antenna.
-Increase the separation between the device and receiver.
-Connect the device into an outlet on a circuit different from that to which the receiver is device.
-Consult the dealer or an experienced radio/TV technician for help.

Cleaning & Maintenance

1. Please Clean this product when it is powered off.
2. Do not use alcohol, gasoline, detergents, hot water.
3. Wipe gently with a damp cloth. Air-dry in a ventilated area before storage.
4. Do not soak in water.

Troubleshooting


Issue	Cause	Solution
Failed to start	Held power button too briefly	Press and hold power button for 2 second
	Weak magnetic connection (controller & mask)	Remove and reattach controller firmly
	Low battery	Recharge controller
Unchargeable	Power source off	Check if outlet/power is on
	Loose charging cable	Fully insert cable into controller port
Dim LED lights	Low battery	Recharge controller

These are common issues and solutions. If problems persist, contact our support team. Unauthorized disassembly by non-professionals is prohibited.

Instructions

LenGen

Wireless LED Light Therapy Mask

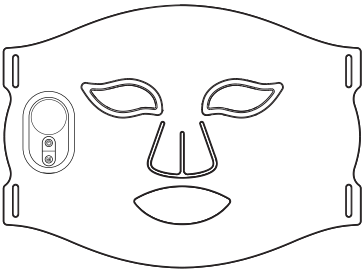

SKB-2418
Read instructions before use and retain for reference

封底

封面

5

6

7