

Eye Massager

User Manual




Power
display



Cooling



Hot

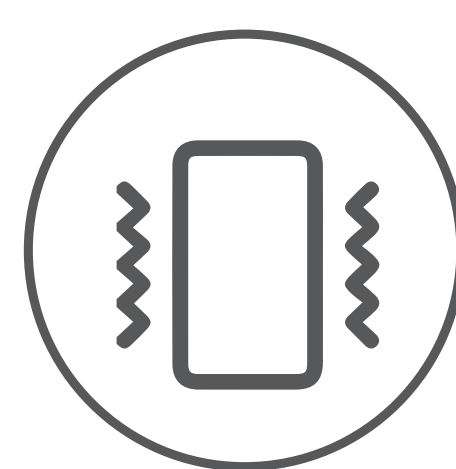

Air
compression



Bluetooth



Music

Thank you for choosing our product. Please read this manual before using it. Please keep it properly for future reference.

For the appearance of the product, please refer to the real object.

Must read before use


Unsuitable For

1. For people nearsighted with more than 8.0 diopters, please use the product under the guidance of a doctor;
2. Patients with severe glaucoma and cataract;
3. Children under 7 years old are suggested not to use the eye massager;
4. People with stroke sequelae are forbidden to use it;
5. Patients with injuries around the eyes and patients who have just gone through surgeries (like double eyelid surgery) are forbidden to use it during the rehabilitation period;
6. Those who have undergone major eye surgeries such as glaucoma, cataract and retinal detachment are forbidden to use it; For people who have done the myopic surgery for one year, please consult your doctor for advice;
7. People with severe inflammation of eyes are forbidden to use it.

Safety Instructions

1. It is recommended to appropriately adjust the product with suitable tightness when using it (the tightness can be adjusted slowly based on the user's experience in use). When it is in use, the user should lean against the sofa or chair for better fitting, more uniform stress, and better effect and comfort!
2. It is recommended to use it 1-2 times a day for 10-20 minutes each time. Remember to close your eyes and rest for 3-5 minutes after use, and your eyes will be more comfortable in this way.
3. The product can be used after charging.
4. Product cleaning: Wipe the mask and cover of the product with a wet towel or a cotton ball dipped in 75% medical alcohol.
5. Keep away from fire source, corrosive liquids and various kinds of adverse environments

Product Description

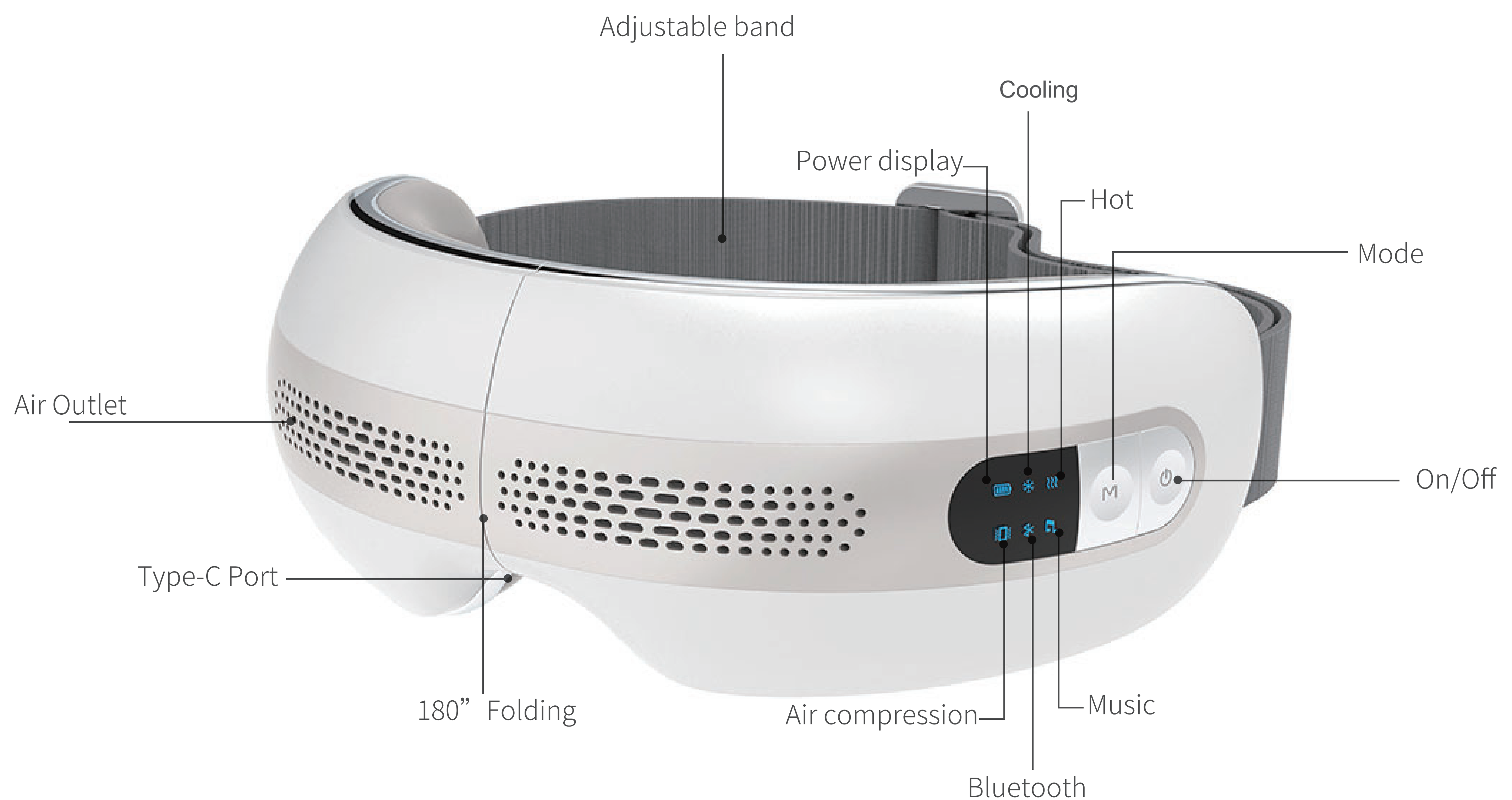

Product Functions

1. Four modes
2. Voice
3. Built-in music
4. Bluetooth connection
5. One-key function operation
6. Type-C charging interface
7. Battery capacity:2500mAh
8. Rated power:3.5W
9. Rated voltage:3.7V

Operation Instructions

1. Press the button  for 3 seconds to turn on and off.
2. Short press the button  to switch between different modes, Dynamic mode, Warm mode, Cooling mode, Automatic mode.
3. Double-click the button  to pause/resume music playback.
4. Short press the button  to adjust the air pressure intensity, turn on and off the air pressure. The air pressure cannot be turned off in Dynamic mode.
5. When the Bluetooth  lights up, it means that the mobile phone has been successfully connected to the product Bluetooth and can play mobile music; When the music  lights up, it means that the product's own music is being played and the mobile phone Bluetooth is not connected.

Function Description

- 1.Dual channel output speaker 8R0.8W
- 2.Charging current 900mA.
- 3.Two hot and cold tablets.
- 4.Two 10KR thermistor.
- 5.Horn, air pump, air valve using 1.0MM terminals, battery, hot and cold pieces, temperature resistor using 1.25MM terminals.
- 6.When the battery power is 3.3V, the battery quantity detection and alarm are carried out.
- 7.Charging indicator light charging red light, fully charged blue light.
- 8.The default working time is 10 minutes.
- 9.Warm compress temperature: 39 degrees.
- 10.Cold compress mode temperature is 18 degrees.

Mode	Air compression	Hot compress	Cold compress	Music
Dynamic mode	√			√
Warm mode	√	√		√
Cooling mode	√		√	√
Automatic mode	√	√	√	√

Voice Prompt

- 1.Boot: Welcome to use the smart eye massager.
- 2.Shut down: Thank you for using.
- 3.Working mode: Dynamic mode, Warm mode, Cooling mode, Automatic mode.
- 4.Air pressure: Air pressure soft, air pressure strong, air pressure off.
- 5.Low power alarm: Low battery power.
- 6.Bluetooth Connected: Successfully connected Bluetooth.
- 7.Bluetooth Disconnect: Bluetooth has been disconnected.

FAQ

Status	Potential causes / treatment methods
Unable to start the machine	Please confirm whether the battery is low; Please confirm whether the product switch has been turned on;
Unable to charge or no charging indicator	Please check whether the type-C cable is connected properly;
The machine suddenly stopped when it is working	Is the product low on power? Or is the massage time up?



The product is designed with simple shockproof facilities. During transportation, please avoid violent shaking, throwing, falling and exposure to rain, snow and sun.

Warning Note

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction