




FAR-INFRARED Hair Growth Device

HL800 User manual

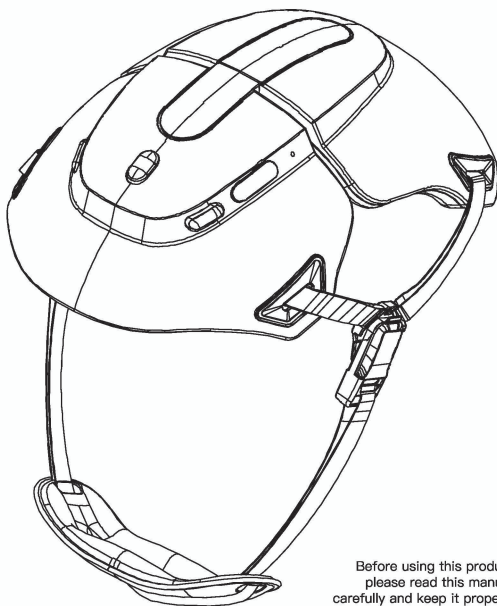

Before using this product,
please read this manual
carefully and keep it properly.



Notice: Please read and follow the instructions in the Instruction Manual carefully before using this product. The Company will not be responsible for any consequences caused by incorrect use or non-compliance with the instructions.

Warning:

1. Please avoid using the instrument after drinking or sleeping.
2. The instrument cannot be used as a head protection device such as a safety helmet.
3. The instrument should be avoided from falling to the ground and severe impact.
4. The instrument should avoid strong sunlight and gas corrosion.
5. The instrument should be placed away from children and pets.
6. When the instrument is not used, please turn off the power supply in time to avoid long-term work.
7. The instrument is prohibited to be used in places with anesthetic gas, oxygen, hydrogen and other flammable gases or other flammable chemicals, otherwise there is a risk of explosion or fire.
8. Do not drop or insert objects into the vent or opening of the instrument, and attention should be paid to prevent insects or small animals from entering the device.
9. When the instrument is not used for a long time, it should be stored in the required environment, and charged at least once a month to maintain the performance of the battery.
10. The scalp temperature will increase when the instrument is used, and it will not burn the skin under normal circumstances. If you feel uncomfortable due to high temperature, please immediately stop using it.

※ Special Reminder:

Using pure graphene heating products, long-term heating will lead to local pore expansion. If the cold air is removed immediately, it will aggravate the symptoms before maintenance. Therefore, after use, please disconnect the power supply and wait for the temperature to drop before removing the instrument.

Use With Caution According To Doctor's Advice:

1. Heat allergy, not suitable for hot compress people.
2. People without self-care ability should be used under the care of family members.
3. People who are pregnant or people who are forbidden by doctors.
4. People who are sensitive to temperature or have delayed reaction to temperature, or are receiving treatment, please consult a doctor before using this product.
5. Patients with acute inflammation, suppurative inflammation and bleeding tendency are not suitable for use.
6. Patients who have not recovered from brain surgery, please consult a doctor before using this product.

Packing List

Equipment * 1 The manual * 1 QC pass card * 1

Adapter * 1 TYPE-C quick charging cable * 1

Cleaning cloth * 1 Small column guide comb * 1

※The product packaging include power adapter.The power adapter applicable to this product is: all types of TYPE-C 5V==1A /12V==1.67A power adapter.

Specification

Material: silicone + ABS + PC

Size: 220mm*180mm*110mm

Weight: 430g

Working voltage: 3.7V==

Working power: 5W

Battery capacity: 3000mAh

Charging volatage:DC5V/PD12V

Fully charged battery: about 60–90 minutes for fast charging,
about 300–360 minutes for slow charging.

Product Key Position Diagram

A: Forced shutdown button


D: Type–C charging port

B: Massage adjustment button

E: Heating adjustment button

C: Charging indicator light

F: Power on/off button

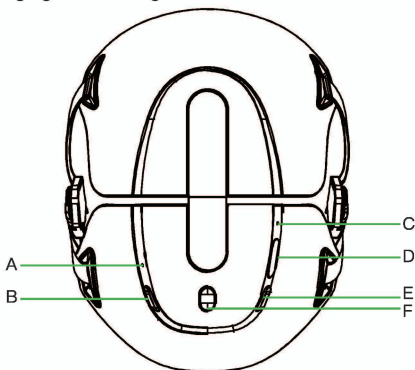

Indicator Light Description

- White light of power on/off button always on: the device is in the power on/off state.
- Heating button green light: Low temperature setting.
- Heating button blue light: Medium temperature setting.
- Heating button red light: High temperature setting.
- Heating button red light flashing: Heating system malfunction.
- Massage button green light: Low intensity setting.
- Massage button blue light: Medium intensity setting.
- Massage button red light: High intensity setting.
- Massage button red light flashing: Massage system malfunction.
- Charging white light blinking: Device is charging.
- Charging white light steady: Device charging complete.

User Manual

I. Steps for Wearing

1. Preparation: Prepare the penetrator and remove the left/ right magnetic buckle. Press the button on top of the magnetic buckle and move it outward to remove it (Refer to Figure 1).
2. Place the penetrator on the head, then fasten the lettered magnetic buckles on the headband together (they automatically attach when placed on the same horizontal plane). (Refer to Figure 2).
3. According to the condition of ones head , the instrument can be manually stretched back and forth to adjust to match the head shape.((Refer to Figure 3).
4. Depending on the user's comfort, adjust the tightness of the headband as needed (the protruding part of the headband can be adjusted by contraction).

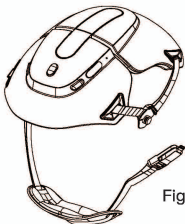


Figure 1

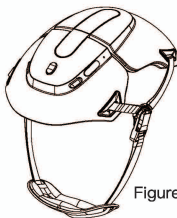


Figure 2

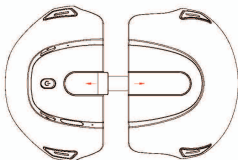


Figure 3

II. Power On and Usage

1. Press and hold the power button for 3 seconds to turn on the device. Upon startup, the massage and heating functions will activate simultaneously, the default green light low level work(Refer to Figure 1).
2. Use the massage adjustment button to adjust the intensity to your comfort (Refer to Figure 2).
3. Use the heating adjustment button to adjust the temperature to your comfort (Refer to Figure 3).
4. Massage and heating gear short press in (green) low – (blue light) middle – (red light) high – close, cycle.
5. Turn off the massage and heating function, the power indicator is on and enters the standby mode, and automatically shuts down after 5 minutes.
6. In standby state (turn off heating and massage), press any key to turn on the massage and heating function at the same time, and automatically restore the default low-grade heating and massage after starting the work.
7. The product does not support the function of charging while discharging, standby or working state access charger, charging priority, heating or massage automatically off.
8. The equipment starts to work normally, work 1 hour automatically shut down.

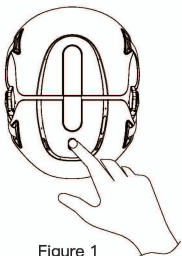


Figure 1

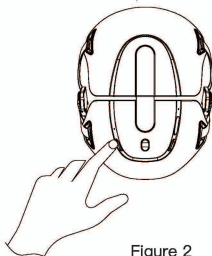


Figure 2

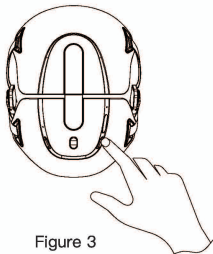


Figure 3

III. Device Charging

1. Remove the silicone plug from the charging port of the device (refer to Figure 1).
2. Prepare a Type-C charging cable (refer to Figure 2).
3. Insert the Type-C charging cable into the Type-C port. The indicator light(White) will blink, indicating charging is in progress. When the indicator light(White) stays on continuously, it means charging is complete (refer to Figure 3).

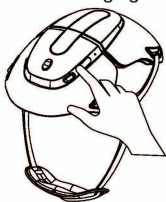


Figure 1

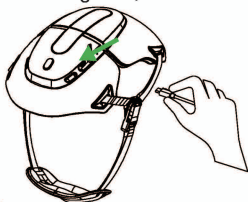


Figure 2

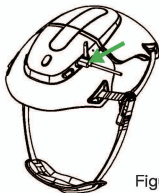


Figure 3