


UREVO Care Wireless Recovery Massager

User Guide



US
Model. UCRM002

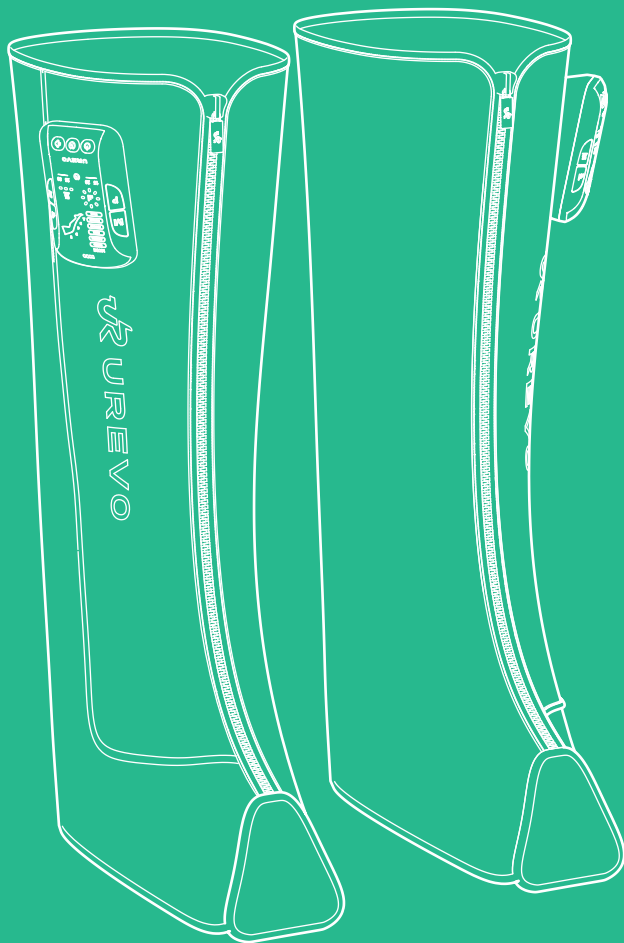

REGISTER NOW TO ACTIVATE 24-MONTH WARRANTY!

You'll receive a 24-month warranty and lifetime customer service after activation. This warranty applies only to eligible purchases from UREVO authorized distributors. For more details, visit urevo.com/warranty.



[WWW.UREVO.COM](https://www.urevo.com)

Customer Service

Email:

- **US:** services@urevo.com
- **EU:** service_eu@urevo.com (Response within 24 hours)

Tel: (1) 844-998-2473

Hours: Mon-Fri, 7 AM - 6 PM PST

Sat-Sun, 7 AM - 11 AM, 12 AM - 4 PM PST



Scan QR code to download
and install UrevoFitness APP!

IMPORTANT: Read all instructions carefully before using this product.

The specifications of this product may vary from this photo, subject to change without notice.



Scan now for the User Manual!!

For more instructions, please scan the QC Code.

IMPORTANT SAFETY INSTRUCTIONS

Read all instructions before using (this appliance).

DANGER – To reduce the risk of electric shock:

1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

1) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.

2) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

3) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.

4) Keep the cord away from heated surfaces.

5) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.

6) Never drop or insert any object into any opening.

7) Do not use outdoors.

8) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

9) **WARNING:** Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

10) **WARNING:** Do not stand on or in appliance. Use only while seated.

11) This appliance contains batteries that are non-replaceable.

12) This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

SAVE THESE INSTRUCTIONS

	Note: There is a possibility of personal injury and damage when used incorrectly		Handle fragile items with care
	Please read the instructions carefully before using this product		Moisture-proof and rainproof product
	Do not throw rubbish carelessly. Please dispose of it according to local regulations		



Manufacturer: ZHILE US CO., LTD.
Address: 22635 NE MARKETPLACE DR STE 100, REDMOND, WA,
98053-5898, UNITED STATES
Version:1.0



IMPORTANT SAFETY INSTRUCTIONS

- 1.Charging with adapters that do not meet the requirements of this product is prohibited.
- 2.Non-professional repair technicians are prohibited from decomposing and repairing this product.
- 3.If the power cord is damaged,in order to avoid danger,it must be replaced by the manufactur-
er or its maintenance department and similar full-time personnel.
- 4.Do not use in wet areas or when the product has penetrated into water or other liquids.
- 5.After use,please place the massager in a dry place, away from humidity and direct sunlight.
- 6.Do not place in fusible,explosive,chemical solution,damp place.
- 7.Do not touch the product with sharp objects to prevent damage to the product parts.
- 8.If you feel uncomfortable or abnormal during use,please stop using immediately.
- 9.Do not use in high temperature environment.
- 10.Please keep the product clean and use a dry towel during maintenance.
- 11.Before using this product,it is recommended to consult a doctor or use under the guidance
of a professional:
 - High-risk pregnant women and children should not use.
 - Mental patient
 - Suffer from skin infection,skin swelling,scalding burns,etc
 - High-risk patients who have a serious illness and are not suitable for massage
 - The injured have weak or damaging bones
 - People with cardiac pacemakers
 - The body is fitted with other metal devices,such as prosthetics
 - People with serious heart disease or abnormal blood pressure

SAVE THESE INSTRUCTIONS

Specification

Product Name	UREVO Care Wireless Recovery Massager
Model No.	UCRM002
Pressure level	8 levels
Noise	≤60dB
Charging port	Type C
Battery	5000mAh/7.4V
Charging current	5V ⚡ 2A
Charging time	<6H
Working time	➤4H
Product Weight	1.75 KG
Input voltage	5V

FCC STATEMENT:

Operation is subject to the following two conditions:
(1) This device may not cause harmful interference, and
(2) This device must accept any interference received, including interference that may cause undesired operation.
(3) Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
Note:This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try the interference by one or more of the following measures:
to correct

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with RSS247 of Industry Canada. Cetappareil se conforme à RSS247 de Canada d'Industrie. This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device. appareils radio exempts de licence. Son fonctionne-mentest sujet aux deux conditions suivantes: (1) le dispositif ne doit pas produire de brouillage-
préjudiciable, et (2) cedispositifdoit accepter tout brouillagerecu, y compris un brouillage susceptible de provoquer un fonctionnementindésirable.
Radiation Exposure Statement: The product comply with the Canada portable RF exposure limit set forth for an uncontrolled environment and are safe for intended operation as described in this manual.
The further RF exposure reduction can be achieved if the product can be kept as far as possi-
ble from the user body or set the device to lower output power if such function is available.
Déclarationd'exposition aux radiations : Le produitestconforme aux limitesd'exposition pour les
appareils portables RF pour les Etats-Unis et le Canada établies pour un environnement non
contrôlé.
Le produitest sûr pour un fonctionnementtelquedécritdanscemanuel. La réduction aux exposi-
tions RF peutêtreaugmentéesil'appareilpeutêtreconservéaussi loin que possible du corps de
l'utilisateurouque le dispositifestrégulé sur la puissance de sortie la plus faiblesiunetellefonction-
estdisponible.
This Class B digital apparatus complies with Canadian ICES-003.
Cetappareilnumérique de la classe B estconforme a la norme NMB-003 du Canada.