CrocBird

Fig.1

Fig. 2

Fig.3.1

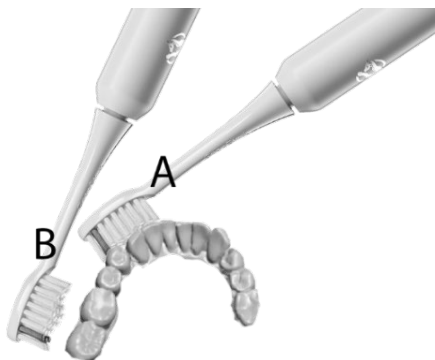
Fig.3.2

Fig.3.3

Fig. 4

Fig.5

Fig.6.1

Fig.6.2

Fig.7

Fig.8.1

Fig.8.2

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Developed in Canada

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180 Northfield Dr. Suit 4, Waterloo, N2L 0C7
Ontario, Canada

Customer service\Service client:

Email: David@crocbird.com,

Tel: 1-877-249-7122 ext. 1001

Website: www.crocbird.com

Special reminder:

- ❖ **This product features a switchless design. The brush head will automatically vibrate upon contact with teeth or any other objects**
- ❖ **The power ramps up to its maximum after operating for a total of 15 minutes**

1. Product Definition (Fig.1)

- 1.1 **Crocbird Overview:** Crocbird is an advanced electronic toothbrush features a proprietary cutting-edge linear motion technology to deliver a superior cleaning experience. It utilizes a translational brushing movement where the brush head sweeps side-to-side. This action more effectively cleans the surfaces of teeth and the gumline while minimizing irritation to both teeth and gums.
 - 1.2 **Switchless design:** There is no switch button on the handle surface. Crocbird activates or deactivates seamlessly as bristles contact or leave the teeth.
 - 1.3 **Detachable Brush Head:** The Crocbird toothbrush features a detachable brush head. Replacement is indicated by a color reminder on three specific tufts of bristles.
 - 1.4 **Ergonomic Design:** A slim handle allows both adults and kids grip and operate it effortlessly.
 - 1.5 **Travel Case with Charger:** The Crocbird comes with a travel case that includes a built-in wireless battery charger and a standard Type-C USB port allowing for connecting to a safety certificated 5V DC power supply, such as a smartphone charger.
- ### **2. How Crocbird Works (Fig.2)**

- 2.1 **Motor Functionality:** Crocbird is powered by a flat motor that provides linear reciprocal motion.
- 2.2 **Automatic Activation:** The motor starts and stops automatically, initiated by the brush head. It activates when the brush head contacts teeth or any other physical objects and deactivates upon detachment.

3. Specifications

- 3.1 **Innovative Motor Technology:** The flat motor is designed based on Lorentz force law principles, ensuring maximum energy efficiency. This unique technology is proprietary to Crocbird Innovation.
- 3.2 **Battery Life:** The normal model lasts 2.5 to 3 months on a single charge, assuming 2 minutes of daily usage. The kid's edition provides double the usage period per charge.
- 3.3 **Motor Power:** The normal model operates at about 1.35 W, while the kid's edition operates at about 0.7 W.
- 3.4 **Battery Specifications:** It uses a rechargeable Lithium-ion battery with an 800mAh capacity, operating between 3.3–4.2V.
- 3.5 **Charger input:** The average power input to the wireless charger is about 1 W.

4. Main Features of Crocbird

- 4.1 **Deep Cleaning:** Help with a thorough and efficient tooth cleaning experience.
- 4.2 **Whitening Effect:** help whiten teeth and maintain their natural whiteness with daily use, making it especially beneficial for coffee and tea drinkers or smokers, showing visible results before and after brushing.
- 4.3 **Smart Control:** Activates or deactivates seamlessly as bristles contact or leave teeth.

4.4 Gentle on Gums: Ideal for sensitive teeth and gums; suitable for kids and the elderly.

4.5 Ergonomic Design: Features a slim body for effortless grip and operation.

5. How to Use Crocbird

5.1 Attaching the Brush Head (Fig.3.1)

- **Occasions for Installation:** When removed from the travel case or when installing a new brush head.
- **Installation Steps:** Align the brush head with the metal shaft, double-check that the bristles face the side with the logo, and push until the gap is about 1.5mm (see Fig. 1).

5.2 Applying Toothpaste

- Wet the bristles, if necessary, before applying toothpaste.
- Quickly and carefully apply toothpaste to the bristles to prevent triggering vibrations that could shake it off.

5.3 Starting and Stopping Brushing

- The toothbrush starts automatically when the brush head is placed in your mouth, contacting the teeth (Fig. 3.2-A).
- It is programmed to operate for two minutes, the recommended brushing duration to ensure optimal dental care.
- The toothbrush pauses its brushing action if the bristles lose contact with the teeth during the two-minute cycle (Fig. 3.2-B).
- It resumes brushing when the bristles contact the teeth again within the two-minute frame.
- To initiate a new brushing session, the user must release the handle and grip it again, resetting the process for a fresh cycle.

5.4 Brushing Technique

- Manually control the position of the brush head, moving it slowly and repeatedly across and along the teeth.

- Thanks to its linear motion feature, the bristles can glide across the gingival line and reach the gum areas without causing stinging or discomfort.

5.5 **Cleaning the Brush Head (Fig.3.3):** Rinse with water, scrubbing bristles with fingers if necessary. Vibration during rinsing is normal.

5.6 **Storing Crocbird:** Place on a countertop without other objects contacting the brush head.

6. Brush Head and Replacement

6.1 **Types of Brush Heads:** Available in normal and kid's editions, both fit the same handle.

6.2 **Time for Replacement:** Replace the brush head every 3 months or when the color on the reminder tufts fades significantly. The reminder consists of 3 different colored tufts of bristles near the tip (Fig.4).

6.3 **Detaching Brush Head:** Pull out gently from the handle.

6.4 **Attaching New Head:** Reference section 5.1 for guidance.

7. The Handle

7.1 **Handle Variations:** Handles differ by power level, including a kid's edition marked by a distinctive label (Fig.5).

7.2 **Model Information:** The product model and power level information are located on the handle surface, denoted as Pxx, where xx ranges from 00 to 99.

8. Battery Recharging

8.1 **Low Battery Indicators and Recharge Requirement:** Recharge when performance is intermittent and the battery indicator at the bottom of the handle blinks. Ensure the battery is charged before a normal brushing (Fig.6.1).

8.2 **Preparations for Charging**

- Remove the brush head.

- Place the handle in the travel case properly and ensure it can close completely.

- Connect to Power:** Connect the travel case using the provided cable or a standard USB Type-C cable to a certified 5V DC power supply, such as a smartphone battery charger (Fig.6.2).

- Charging Duration:** The charging can last up to 10 hours, depending on the battery's depletion level.

- Charging complete:** The battery charging is complete once the charging indicator on the travel case stops blinking

9. Travel Preparation (Fig.7)

9.1 **Packing for Travel:** Store the Crocbird and brush head in the travel case; ensure it's closed securely and carried in your luggage without pressure.

9.2 **Noise Check:** Ensure no noise from the package; the metal shaft must not contact objects.

9.3 **Noise Cessation:** Abnormal noise stops within four minutes after the handle is released from holding.

10. Disposal

10.1 **Battery Disposal:** Follow local regulations regarding Lithium-ion battery disposal. If necessary, remove the battery from the handle.

10.2 **Removing the Battery:**

- Caution:** Attempt this procedure only if you have a basic understanding of mechanics and electronics. Prior to beginning the removal process, ensure you are equipped with safety gear to protect your eyes, face, and hands. Refer to Fig. 8.1 for the battery's location inside the handle via a break-out view.

- Tools preparation:**

- Have a hammer and a flat-tip screwdriver ready for use.

- Ensure you have a table with a sturdy corner to provide adequate support during the procedure.

10.3 suggested steps for removing battery

- Place the handle on a corner of the table, ensuring that the long axis of its elliptic cross-section is vertical.
- Use the hammer to strike the end portion of the handle from the top. Ensure the point of impact is within 15mm from the end section (Fig.8.2).
- Observe for a gap forming between the housing and the plastic end part. This indicates that the impact was sufficient.
- Insert the screwdriver into the gap and pry the plastic end part away. Once loosened, the housing can be easily removed.
- Forcefully pull the battery away from the PCB to detach it.

10.4 Important Warning:

- Throughout the process, be cautious not to damage the battery body to prevent fire or explosion risk.
- Ensure each step is followed accurately to maintain safety and effectiveness.

11. Important Safeguards

11.1 **Child Safety:** Make sure children are supervised to prevent them from using the appliance as a toy.

11.2 **Operational Care:** Do not have children replace the brush head by themselves.

12. Dangers and Cautions

12.1 Pay special attention to any instructions highlighted in bold font within the manual.

12.2 Ensure no object touches the brush head or the metal shaft when the Crocbird is not in use.

12.3 When traveling, ensure the Crocbird is properly enclosed within the travel case.

- 12.4 To recharge the battery, make sure the brush head is removed from the handle and the handle is properly positioned inside the travel case so that the case can be closed completely.
- 12.5 Always place and store the travel case at a safe distance from water.
- 12.6 Do not use a wet travel case for charging until it is completely dry.

13. Warnings

- 13.1 Do not clean the appliance, whether in whole or in part, in a dishwasher.
- 13.2 Any charger with a standard USB Type-C adapter for a smartphone can be used to connect the travel case for Crocbird battery charging.
- 13.3 Additionally, only UL or CE certified 5V DC power supplies with a standard USB Type-A adapter, either as a stand-alone product or a wall outlet, can be used for Crocbird battery charging.
- 13.4 If you lose the provided charging wire, any other charging or data transfer wire with standard USB Type-A to Type-C adapters can be used as a substitute.
- 13.5 The product is not designed to be repaired by service personnel in the event of failure or component defect; if this occurs, the unit should be disposed of.
- 13.6 Discontinue use immediately if the product appears damaged in any way, including the brush head, handle, or travel case.
- 13.7 This product contains no user-serviceable parts. Refer to the "Warranty and Support" section if the product no longer functions properly.

13.8 FCC Warning

This device complies with Part 18 of the FCC Rules.

**Refer to FCC Part 18.213(c) and KDB
784748 Annex A section A.5:**

- (a) Max. Electric Field Emissions: 6.19V/m
- (b) Max. Magnetic Field Emissions: 0.74A/m
- (c) The maintenance of the system is the responsibility of the manufacturer.
- (d) To prevent interference, the product should be as far away from devices operating at the same frequency as possible.

RF exposure considerations

This equipment complies with the FCC RF radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20cm between the radiator and any part of your body.

14. Medical Warnings

- 14.1 Consult your dental professional before using this product if you have had oral or gum surgery within the past 2 months.
- 14.2 Contact your dental professional if excessive bleeding occurs after using this product, or if bleeding continues after 1 week of use.
- 14.3 Discontinue use of this product and contact a physician or dentist if you experience discomfort or pain.
- 14.4 Consult your physician if you have any other medical concerns.

15. Warranty and Support

15.1 Warranty and Support

Crocbird Innovation warrants its products for two years from the date of purchase. Defects due to manufacturing and transportation will be replaced at Crocbird Innovation's expense, provided that valid proof of purchase within the qualifying period is presented. Replacements can be completed through your online or offline purchase platform, stores, or other agencies. If you encounter any difficulties with the replacement process any other issues, please

visit www.crocbird.com or send an email to David at david@crocbird.com, or call the toll-free number 1-877-249-7122 ext. 1001.

15.2 Warranty Exclusions

- Brush heads.
- Damage caused by improper use and neglect, including parts replacement and repair.
- Normal wear and tear, including scratches, abrasions, discoloration, or fading.

15.3 Limitation of Remedies

The remedies provided under this warranty are exclusive and in lieu of all other remedies and warranties, whether expressed or implied. Crocbird Innovation will not be liable for any incidental or consequential damages arising from the use, misuse, or inability to use the product, even if Crocbird Innovation has been advised of the possibility of such damages. In no event shall Crocbird Innovation's liability exceed the purchase price of the product. This warranty gives you specific legal rights, and you may also have other rights which vary from state to state or country to country. Some jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusions may not apply to you.

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