

CERAGEM MASSAGE CHAIR USER MANUAL

CGM-EMCG 2404

Thank you for purchasing the CERAGEM PAUSE massage chair.
Please read this instruction manual carefully and use it correctly before using it.
Please always keep it where you can see it after reading it.

Product warranty is included in the instruction manual.
This product is a home massage chair,
so please do not use it for business or industrial purposes.
The headquarters may improve the design and specifications
of this product without further notice of circumstances.

CERAGEM LIFESTYLE MASSAGE CHAIR USER MANUAL

PAUSE M8 USER MANUAL

CONTENTS

1. safety precautions
2. name of each part
3. Remote control description

The following safety precautions are used to ensure that the product is safe and accurate
It is intended to prevent unexpected risks or injuries by using it.

When installing

Never renovate, disassemble, or repair except for service technicians
→ It can cause fire and electric shock.

Use it on a flat surface
→ It can cause a failure or accident.

Power related

If it doesn't move or if something is wrong, immediately remove the power plug from the outlet Remove and ask the customer center for inspection and repair
→ There is a risk of fire or electric shock due to short circuits,

Please guide the child not to touch the power plug.
→ It can cause electric shock or injury.

Avoid undue bending of the power cord or being crushed by heavy objects
→ It is a risk of fire and electric shock.

Do not use multiple devices at the same time when using multiple outlets
→ There is a risk of fire due to abnormal heating of the outlet.

When removing the plug, hold the plug and disconnect it from the outlet.
→ Disconnecting the power line poses a risk of fire and electric shock.

Do not touch the power plug with wet hands,
Disconnect the power plug from the outlet when cleaning or maintaining the product
→ There is a risk of electric shock.

Be sure to use a grounded outlet when using a multi-tap.
→ Using a non-ground outlet can cause electric shock and injury.

In the event of a power outage, immediately turn off the power switch and unplug the power plug from the outlet.
→ It can cause electric shock and fire.

Make sure that the power cord wires do not get stuck between the products during operation of the product
→ Care must be taken when using because short circuits may occur due to damage to the sheath.

1. Name of each part of massage chair



1. Head cushion

2. Bluetooth speakers

3. remote control

4. Shoulder massage unit

5. Back pad


6. Arm massage unit

7. Back massage unit

8. Carf massage unit

9. Foot massage unit

2. Built-in Quick Controller

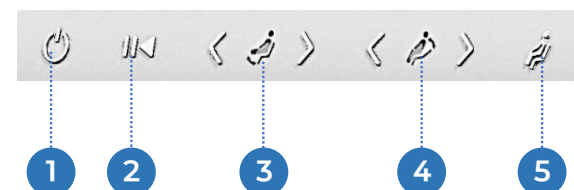

1. POWER On/Off

- you can turn the massage chair on and off by pressing the button for about 3 seconds.

2. PAUSE

- Press the button to temporarily stop or re-activate all functions of the massage chair.

3. Detailed leg angle adjustment


- < > You can use the (left and right) key to adjust the foot angle of the massage chair

3. 4. Detailed body angle adjustment

- < > You can use the (left and right) key to adjust the body angle of the massage chair

5. Quick adjustment of body angle

- You can quickly adjust the angle of the body. (1 weightless / 2 weightless / 3 weightless)

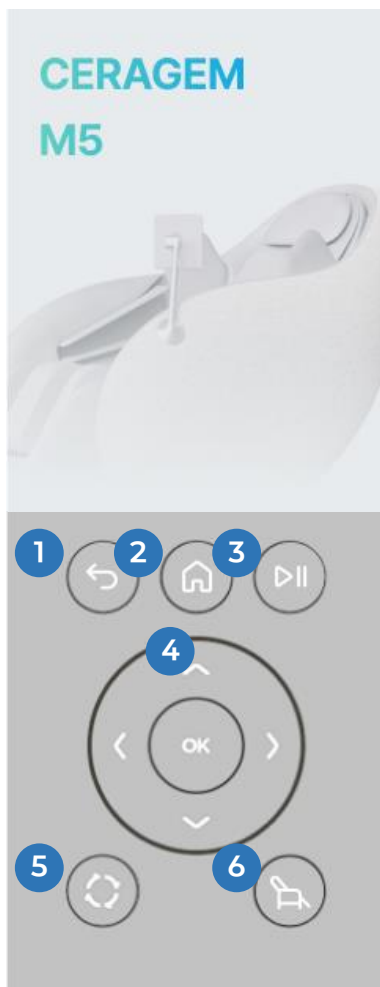

(1) Back button: Return to the previous screen.

(2) Home button: Enter the home screen where you can operate automatic mode, manual mode, setting, recliner.


(3) Play/Pause button: Play or pause the massage.

(4) Move up, down, left, and right buttons: Enter the screen you want to select or complete the settings.

(5) Auto mode button: Enter the menu in auto mode with 18 massage modes.

(6) Recliner button: Enter the recliner menu where you can adjust the angle of the massage chair.

1. Home Screen

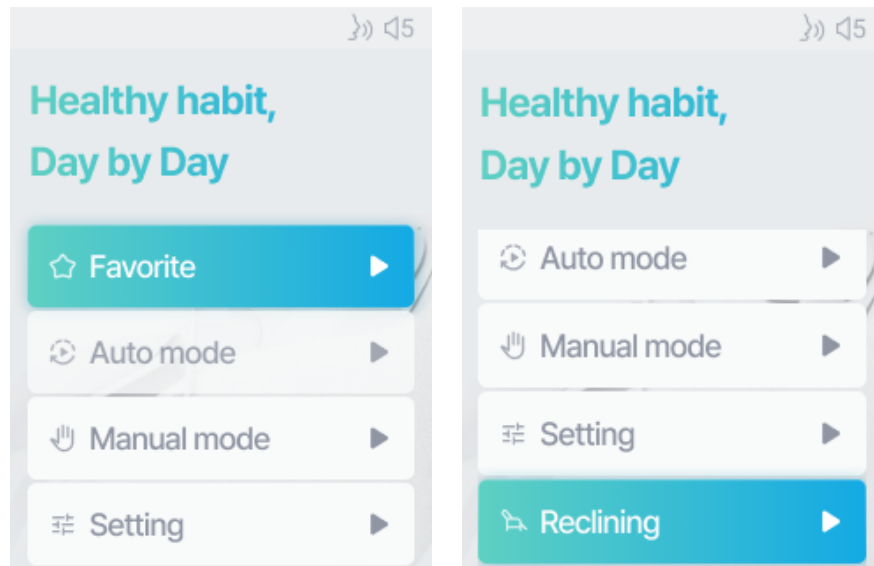


(1) You can control the home screen by turning on the massage chair or pressing the home key.

(2) On the home screen, you can go to the Automatic Mode, Manual Mode, Setting, and Reclining menus.

2. Auto mode


1) Auto mode select

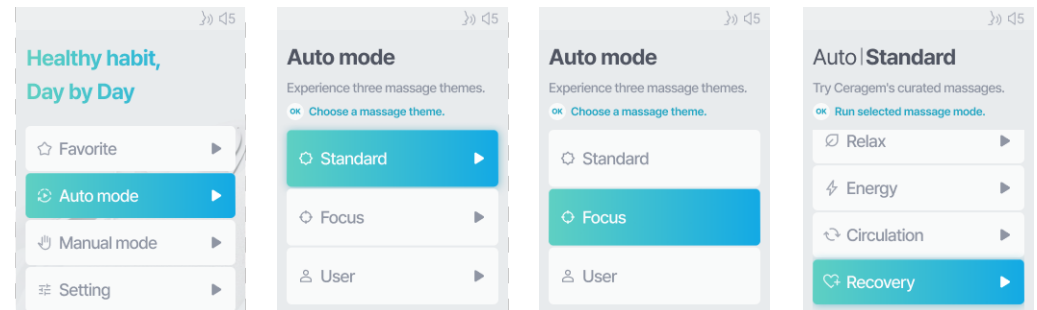

(1) You can select the automatic mode by pressing the up/down button on the home screen.

(2) In auto mode, you can choose 18 automatic modes that fit the theme.


2) Message setting

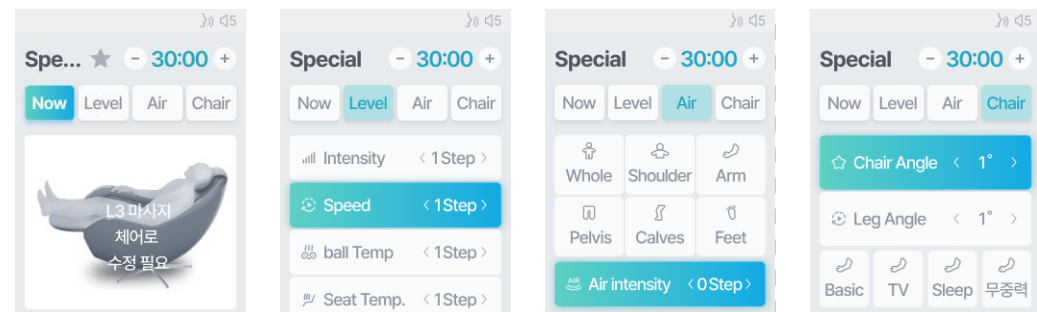

(1) In automatic mode, you can check the current status, detailed message settings, air message settings, and recliner settings.


(2) You can adjust the message time through the + - button at the top of the screen.

3. Manual mode

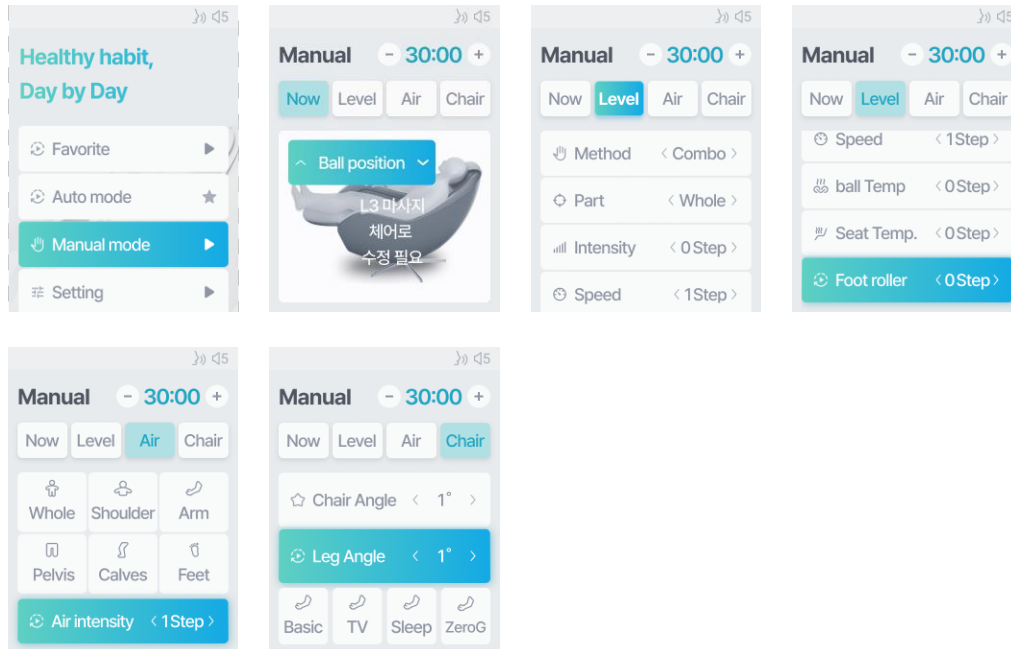

- (1) You can enter the manual mode screen through the up/down button on the home screen.
- (2) The Now menu on the manual mode screen allows you to adjust the status of the massage and the position of the massage module.
- (3) The Level menu on the manual mode screen, you can adjust the massage intensity, massage width, massage method, massage speed, and leg massage settings.
- (4) The Chair menu on the manual mode screen allows you to adjust the angle of the massage chair.


4. Favorite

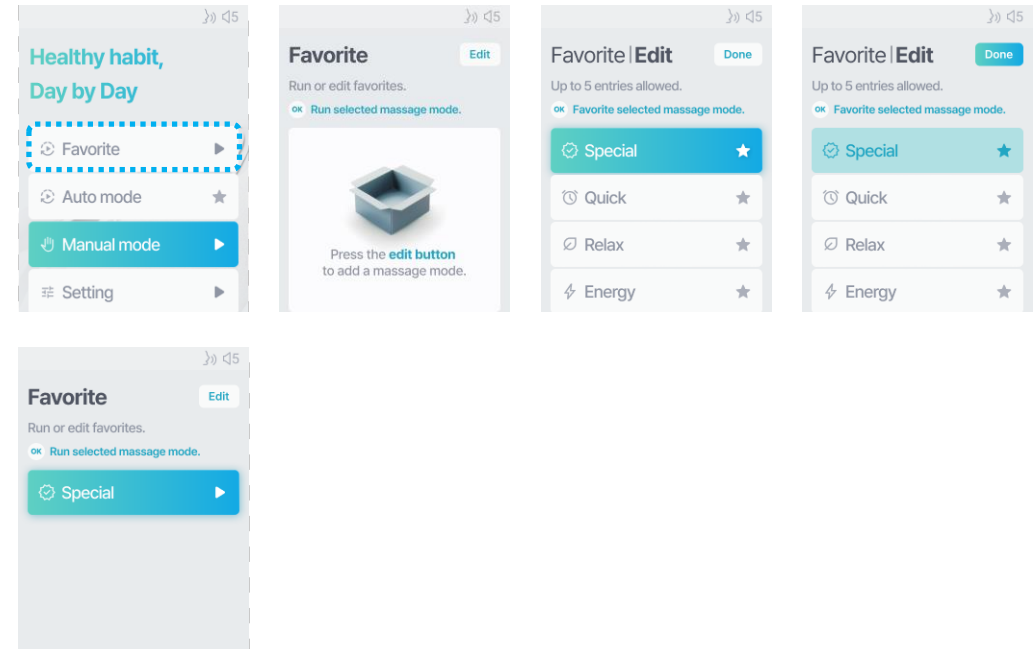

- (1) On the home screen, you can enter the favorite menu through up/down keys.
- (2) Use the Edit button at the top of the Favorite menu to save your preferred message mode..

3. setting

1) Sound

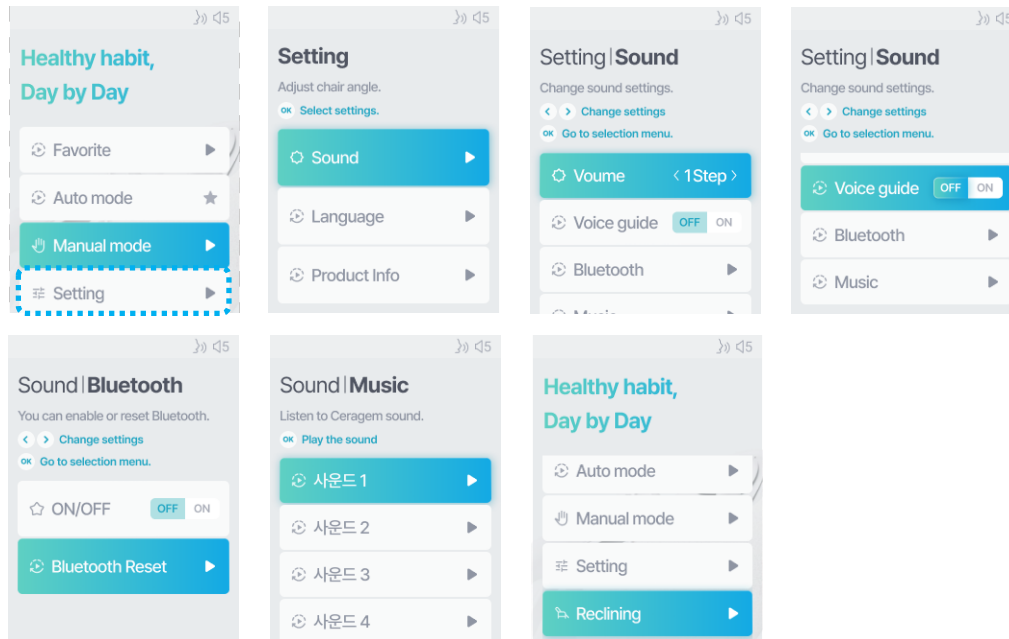


(1) From the home screen, you can enter the Setting menu through the up/down button.

(2) From the Sound menu on the Setting menu, you can adjust the sound volume, select voice guidance ON/OFF, set Bluetooth, and select healing sound.

4. Setting


2) Language

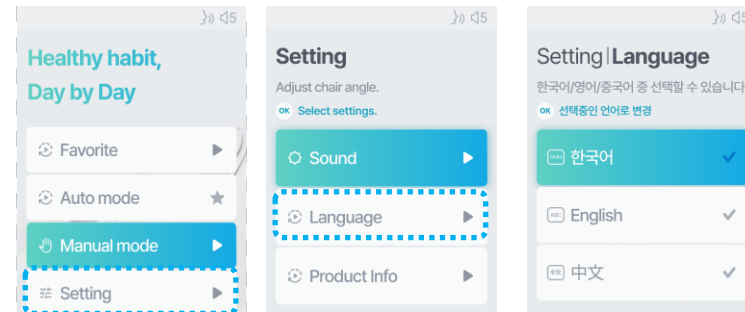

(1) On the home screen, you can enter the Setting menu through up/down buttons.

(2) You can select a language from the Setting menu. (Korean/English/Chinese)


3) Product info

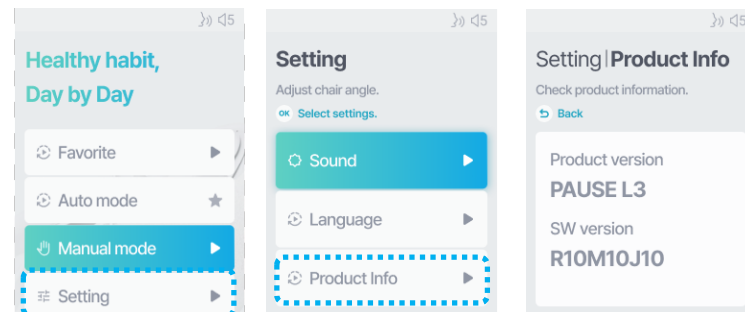

(1) On the home screen, you can enter the Setting menu through up/down buttons.

(2) On the Setting screen, you can view information about the product from the Product info menu.

5. Reclining

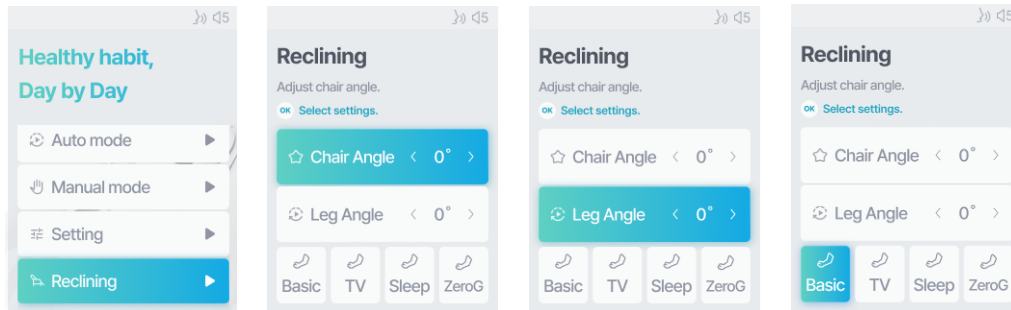

(1) You can enter the recliner menu by pressing the up and down button on the home screen.

(2) In the recliner menu, you can adjust the body angle and leg angle of the massage chair.

04. Product Specifications

제품명	Massage Chair
모델명	CGM EMCG-2404
안마의자 크기	
순중량	
사용 전압	120V, 50/60Hz
소비전력	280W
전자파 적합 인증번호	
전기 안전 인증번호	
고객 센터	1577-5570

05. FCC WARNING

CGM EMCG-2404

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Note: The Grantee is not responsible for any changes or modifications not expressly approved by the party responsible for compliance. Such modifications could void the user's authority to operate the equipment.

The device has been evaluated to meet general RF exposure requirement. This equipment complies with FCC's RF radiation exposure limits set forth for an uncontrolled environment. This device and its antenna(s) must not be