


材料表			
零件名称	物料编码	规格描述	备注


记号: MARK:	日期: DATE:	ECN编号: ECN NO:	说明: DESCRIPTION:	登 记 者: BOOKER:
	2024.04.02		受控版本	韦诗捷
⚠	2024.09.13		文字内容改为英文	韦诗捷

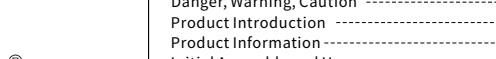

## Smart Capsult Toothbrush

### Dimmi-Q3

**User Manual**



Scan the QR code to install  
the Dimmi Life App

<b>Contents</b>	
Danger, Warning, Caution .....	3
Product Introduction .....	4
Product Information .....	5
Initial Assembly and Usage .....	6
Assembly Precautions .....	7
Travel Mode .....	7
How to Connect the Smart Capsule Toothbrush to Your Phone .....	8
Getting Started with the Dimmi-Q3 Smart Capsule Toothbrush .....	8
Usage Precautions .....	9
Charging .....	10
Low Battery Alert .....	10
Charging Precautions .....	11
Common Issues and Troubleshooting .....	11
Battery Removal Method .....	12
Battery Disposal .....	13
Warranty and Service .....	13
Warranty Card .....	14

### Dear Customer

Thank you for choosing the Dimmi-Q3 Smart Capsule Toothbrush. Please read this user manual carefully before using the product and keep it for future reference.

### Danger

- Keep the charger away from water. Do not immerse the charger in water or any other liquids.
- Ensure the charger is completely dry before connecting it to a power outlet.

### Warning


- Do not remove the brush head during the automatic toothpaste working states to avoid splashing.
- Do not disassemble the smart toothpaste during the automatic toothpaste working states to prevent leakage or damage to the device.
- Keep this product away from children and do not use it as a toy.
- Do not disassemble the device or charging base. If the charger base, power cord, or other components are damaged, replace them with original parts to avoid danger.
- The product contains a built-in rechargeable battery. Do not expose it to fire, high temperatures, or charge it near high temperatures.
- Keep the power cord away from high temperature surfaces.
- Do not hit or drop the product. If the product is damaged due to impact or falling, stop using it.
- This product is intended for cleaning teeth and gums only. Do not use it for other purposes.

### Caution

- Only use accessories specifically designed for the Dimmi-Q3 Smart Capsule Toothbrush. Dimmi is not responsible for any consequences resulting from the use of non-Dimmi certified products.
- It is normal for the toothbrush body to feel slightly warm during charging. If it feels hot, disconnect the power immediately and contact Dimmi Customer Service.
- This product is intended for cleaning teeth and gums only. Do not use it for other purposes.

### Product Introduction

Dimmi-Q3 Smart Capsule Toothbrush Accessories  
Toothbrush Handle x1 Motor Shaft Protection Cover x1  
Brush Heads x2 Brush Head Covers x2 Charging Base x1 Capsule Holder x1

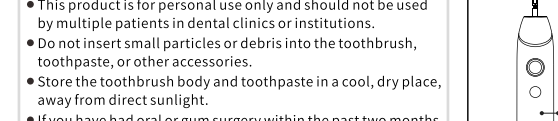

not use it for other purposes. If you feel discomfort or pain, stop using it immediately and consult a doctor.

- This product is for personal use only and should not be used by multiple patients in dental clinics or institutions.
- Do not insert small particles or debris into the toothbrush, toothpaste, or other accessories.
- Store the toothbrush body and toothpaste in a cool, dry place, away from direct sunlight.
- If you have had oral or gum surgery within the past two months, please consult your dentist before using this product.
- If you have braces, dental restorations, or implants, use this product under the guidance of your dentist.
- Initial use of this product may cause itching or minor bleeding due to the gums not being accustomed to this high-efficiency cleaning. This typically subsides within 1 to 2 weeks. If it persists beyond 2 weeks, stop using the toothbrush and consult your dentist.
- It is recommended to start with the lowest vibration intensity when using for the first time. Gradually use other intensities as you become accustomed.
- Replace the Dimmi-Q3 Smart Capsule Toothbrush head once daily if the bristles become worn, deformed, or bent due to normal use. For optimal oral hygiene, replace the toothbrush head every three months.
- The recommended usage time is 2-3 minutes per session, not exceeding 3 minutes.
- Do not place any metal objects on the wireless charging base.
- Disconnect the charging base from the power source once the toothbrush is fully charged.
- Replace the Dimmi smart toothpaste promptly when it runs out. Do not refill the empty bottle with other brands of toothpaste to avoid damaging the device.

### Initial Installation and Use

- Remove the motor shaft protector.
- Take out the brush head compatible with the Dimmi-Q3 Smart Capsule Toothbrush. You can choose between the Whitening or Comfort brush head based on your preference. Align the front of the brush head with the front of the handle and insert it vertically onto the motor shaft at the top of the handle, ensuring it is securely attached.
- Rotate the bottom capsule compartment to counter clockwise as indicated in the diagram to remove it.
- Choose a Dimmi smart toothpaste, unscrew the bottle cap, and remove the silicone seal from the nozzle.
- Fill the capsule compartment with the Dimmi-compatible toothpaste.
- Align the T-shaped icon on the capsule compartment with the T-shaped icon on the handle. Insert the capsule compartment upward in the direction indicated by the arrow on the metal ring. Then rotate it to the right to lock it in place. The smart toothpaste is now installed.

### Installation Precautions

- It is recommended to keep the motor shaft protector in a safe place. After removing the brush head, please cover the motor shaft with the motor shaft protector to ensure the interior remains clean and hygienic.
- When installing the brush head, ensure there is a slight gap between the brush head and brush to ensure the proper vibrating.
- During the installation of the toothpaste, avoid shaking the toothpaste capsule vigorously to prevent leakage and ensure a smooth usage experience.

### Travel Mode

- Quickly press the power button twice to enter Travel Mode. The toothbrush will enter deep sleep, and both the power button and the toothpaste button will be disabled. It is recommended to activate Travel Mode during transport, such as when checking your toothbrush in luggage, to prevent accidental activation of the power or toothpaste buttons.
- To exit Travel Mode, quickly press the power button again. Placing the toothbrush on the charging base will also automatically exit Travel Mode.

[illegible]

## Wk02英文说明书 材质与规格描述

材质:100G双面胶纸

尺寸:115mmX62mm

工艺:单色印刷 风琴折页(7折)

### 工艺和材质要求：



1.色相准确、印迹牢固、套印准确、各色套印不露杂色；

2.制件表面不允许有明显的异色、划痕、模糊等缺陷；

3.套印误差 $\leq 0.1\text{mm}$ ;

4.啤切走位 $\leq 0.2\text{mm}$ ;

5.材料需符合欧盟ROHS要求,以及符合PAH、REACH、销售地食品卫生等相关认证要求。
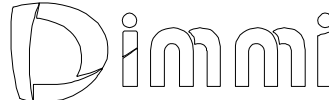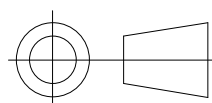
<div></div> <div>深圳迪米生活科技有限公司 SHENZHEN DIMMI LIFE TECHNOLOGY CO.,LTD</div>																																		
<div>除非特别说明，默认单位为mm，公差如下： <i>UNLESS OTHERWISE SPECIFIED : THE UNIT IS MM, THE TORLERENCE IS AS BELOW</i></div> <table><tr><td>尺寸</td><td>公差</td><td>尺寸</td><td>公差</td><td rowspan="4">未注尺寸公差</td></tr><tr><td>3以下</td><td>±0.08</td><td>50~80</td><td>±0.25</td></tr><tr><td>3~10</td><td>±0.10</td><td>80~120</td><td>±0.30</td></tr><tr><td>10~30</td><td>±0.15</td><td>120~160</td><td>±0.35</td></tr><tr><td>30~50</td><td>±0.20</td><td>160~200</td><td>±0.40</td><td>单位 (mm)</td></tr></table> <div>本文档的所有权归深圳迪米生活科技有限公司所有，未经允许任何人不得复制。 <i>THIS DOCUMENT IS THE PROPERTY OF SHENZHEN DIMMI LIFE TECHNOLOGY CO.,LTD ANYONE CAN NOT COPPY UNLESS ALLOWED.</i></div>					尺寸	公差	尺寸	公差	未注尺寸公差	3以下	±0.08	50~80	±0.25	3~10	±0.10	80~120	±0.30	10~30	±0.15	120~160	±0.35	30~50	±0.20	160~200	±0.40	单位 (mm)	物料编码		P05.01WK02.001		产品型号		WK02	
					尺寸	公差	尺寸	公差		未注尺寸公差																								
					3以下	±0.08	50~80	±0.25																										
					3~10	±0.10	80~120	±0.30																										
					10~30	±0.15	120~160	±0.35																										
30~50	±0.20	160~200	±0.40	单位 (mm)																														
零件名称		说明书		产品名称		智能胶囊牙刷																												
材 质		100G双面胶纸		比 例		1:1																												
版 本		V1.3		幅 面	A3	共 1 页	第 1 页																											
制 定		确 认		审 核																														
韦诗捷																																		

**FCC Warning:**

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

**Caution:** Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 0cm between the radiator and your body.