


Eye Massager

User Manual

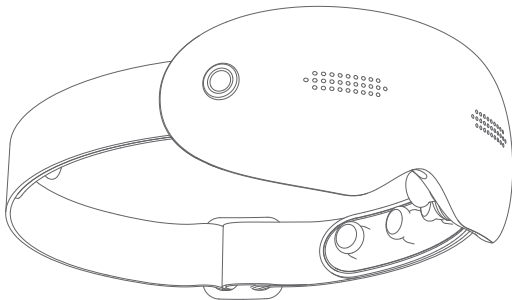

*Before using this product, please read it carefully and keep it properly.

Thank you for using the product provided by us, please read this manual carefully and correctly understand and carefully read the safety cautions before using this product. The company reserves the right to improve the product and user manual. If there is any change without information, the appearance of the product is subject to the actual product. This manual has been carefully checked, but it does not rule out that there may be typographical errors or misunderstandings in the content, please understand.

This manual provides users with safe and effective operation methods. Please save it properly after reading it for future reference.

Contents

Safety Cautions.....	01~04
Production Specification	04
Name and Functions of Machine Parts.....	05
Preparation before Usage.....	06
Operation Instructions.....	07~08
How to Use Bluetooth.....	08
Clean and Maintenance of Product.....	09~10
Problems and Methods.....	10~11
Warnings for Product Discarding.....	12
Recycle Guidance.....	13

Safety Cautions

Warnings

01. When using it for the first time, or after the product has been put on hold for a long time, the battery power is depleted, please charge it according to the instruction manual to ensure that it can work normally.
02. During the use, if your body appears or feels abnormal, please stop using it immediately and consult a doctor.
03. If the product is damaged, do not use it; if the charger/charging cable is damaged, please replace it before charging the product.
04. Before charging, make sure the charger and the charging jack of the product are dry.
05. Do not charge the product or charger when it may be immersed in water.
06. Do not use other massage equipment while using this unit.
07. This product is for relief only, not for therapeutic purposes.
08. Do not use other massage equipment while using this unit.
09. Do not drag the product on the ground, or pull the power cord to move the product, so as not to damage the product.
10. When use it, please tilt your head slightly to make it easier to wear the eye protector.
11. This product contains a lithium battery, please do not disassemble, modify, immerse the product in water or throw it into fire, otherwise it may cause damage or explosion of the product.

Safety Cautions

12. Do not disassemble the internal components of the product without permission to avoid machine failure.
13. When an abnormal situation occurs that cannot operate, please cut off the power and remove it immediately, otherwise it may cause high temperature or short circuit, resulting in an accident. At the same time, please contact customer service for after-sales consultation.
14. Do not use a power supply other than the rated voltage to avoid fire or damage to the product.
15. Do not damage, excessively bend, pull, twist, or knot the power cord.
16. Do not forcefully bend/reversely bend or twist the product (otherwise the product will be damaged).
17. When the USB charging power cord is damaged, please do not disassemble and replace it without permission. Please contact our company or qualified professionals for replacement and repair.
18. When the product works for a long time, the surface will heat up, which is a normal phenomenon.
19. When disposing of this product, please deal with it in accordance with the relevant environmental protection regulations of the city, otherwise it may cause environmental pollution.

Prohibition

01. For myopia of more than 800 degrees, please use it under the guidance of a doctor.

Safety Cautions

- 02. Pregnant women are forbidden to use.
- 03. Not recommended for children under 7 years old.
- 04. The use of stroke sequelae is prohibited.
- 05. Those with injuries around the eyes and those who have just undergone surgery (double eyelid) are not allowed to use it during the recovery period.
- 06. Those with eye disease symptoms cannot use it, and it is forbidden for patients with severe glaucoma and cataracts. It is forbidden for patients who have had eye surgery, glaucoma, cataract, retinal detachment and other major eye diseases; please consult a doctor for advice after one year of myopia surgery.
- 07. Those with severe inflammation of the eyes are forbidden to use.
- 08. Children or persons with disabilities require special supervision to use this product.
- 09. Do not use this product while intoxicated or unwell.
- 10. People with a body temperature of 38 degrees or higher (fever period).
- 11. Those who are not sensitive to heat and those who are prone to scalds are forbidden to use.
- 12. People with trauma around the eyes and people who need to be quiet should not use it.
- 13. Those who use medical electronic devices such as pacemakers should not use this product.

Safety Cautions & Product Specification


14. Patients with malignant tumors, heart disease, severe hypertension, etc., please use it with caution (need to be monitored by a normal person).
15. Patients who cannot express their wishes freely should use it after consulting a doctor.
16. Do not use in places with high humidity, such as bathrooms, or use the product while taking a shower, as this may damage the product.
17. Do not use it for people who do not have the ability to express consciousness, otherwise it may cause an accident or cause physical discomfort.

※ Product Specification

Product name	Eye Massager	Charging interface	Type-C
Rated input voltage	5V $\overline{=}$	Lithium battery	800mAh
Product size	178×80×72mm	Rated input power	5W
Massage duration	10 minutes	Massage modes	5 kinds
Body material	High-quality PC + high-quality plastic + skin-friendly silicone		

Name and Function of Machine Parts

※ Name of Machine Parts

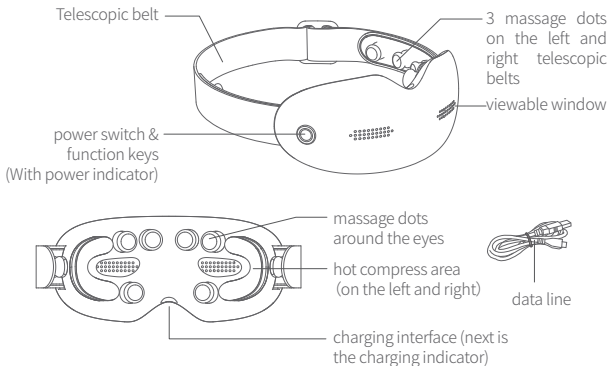

※ Product Functions

5 massage modes

Hot compress

Included soothing music


Voice broadcast

Bluetooth music

Dot massage technique

Preparation before Usage

- Please check the product carefully before use.
- Before using it for the first time, please charge the product for about 2-3 hours to fully charge it.
- When charge it, connect the port of the adapter to the charging port of the handle and plug the adapter into an AC outlet. (As shown below)

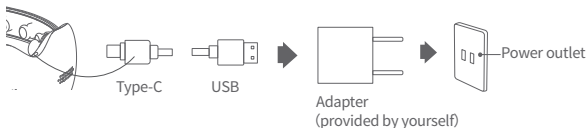

- Be sure to fully charge it before using it for the first time.
- While charging, the product cannot be turned on and used.
- When charging, the red light will be on; when fully charged, the light will turn green .
- When the battery is low, the red light will flash, and at the same time the voice broadcasts “low battery, please charge”, then please charge as soon as possible.
- The default working time of each gear is 10 minutes. After running for 10 minutes, the product will automatically shut down. If you want to continue the massage, please restart.
- Please remove glasses, contact lenses, hair accessories, etc. before wearing the product to avoid discomfort during use.
- Put the massager on your eyes and adjust the elastic strap to the most comfortable position.

Operation Instructions

※ Operation Instructions

◎ On/Off/Function Keys

- 1 - Press and hold the button on the main unit, the main unit power on/off, after power on, the power indicator is always on blue light, and after hearing the voice broadcast “Welcome to the head and eye massager, comprehensive mode”, the product starts to work, and the default Bluetooth is turned on and comprehensive mode.
- 2 - Click to switch the massage mode cyclically, there are 5 modes in total (comprehensive mode, smart mode, eye protection mode, head mode, sleep mode), there will be a corresponding voice broadcast for each mode switched.
- 3 - Double-click, when not connected to Bluetooth, pause/start the background music; when connected to Bluetooth, pause/start the mobile phone music.

Massage Modes

1. Comprehensive mode: eye plus head vibration massage mode 1+hot compress; play music 1.
2. Smart mode: eye plus head vibration massage mode 2; play music 2.
3. Eye protection mode: eye vibration massage+hot compress; play music 3.
4. Head mode: head vibration massage; play music 4.

Operation Instructions & How to Connect Bluetooth

5. Sleep mode: ye plus head vibration massage mode 3; play music 5.

Shutting down the product

Press and hold the button on the main unit, after hearing the voice broadcast “Thank you for using, goodbye!”, the main unit shuts down.

● Notes

1. The default working time of each gear is 10 minutes (Switching the mode midway will restart the timer by 10 minutes). After running for 10 minutes, the product will automatically shut down. If you want to continue the massage, please restart.
2. During sleep mode: After automatic shutdown, “Thank you for using, goodbye” will not be broadcast.

※ How to Connect Bluetooth

1. Connect Bluetooth to mobile phone: Turn on the Bluetooth function of your mobile phone, search for Bluetooth devices on the mobile phone, and select AS-EM005 to complete the connection according to the mobile phone prompts (after the first successful connection to Bluetooth, subsequent use of the product will automatically connect as long as the mobile phone's Bluetooth is turned on mobile phone bluetooth).

How to Connect Bluetooth & Maintenance

2. After the connection is successful, turn on the music on the mobile phone and start the music journey (cannot answer the call talk).
3. Disconnect the Bluetooth connection: Turn off the Bluetooth of the mobile phone, and the eye massager will be directly disconnected from the Bluetooth of the mobile phone.

※ Clean and Maintenance of Product

- Before cleaning, make sure that the power switch is turned off and the power plug is unplugged from the outlet.
- Please clean the surface of the product with a damp paper towel.
- Please clean the dirty place with a soft cloth dipped in water or neutral solvent, and then dry it with a dry cloth.

Notes

- Do not place the product in a place where there is direct sunlight or high temperature, so as not to cause the product to fade.
- Do not plug or unplug the power plug with wet hands.
- Never rinse the product with water, or get water or other liquids inside the product.
- Do not use corrosive cleaning agents such as benzene or thinner to clean the product, otherwise it may cause the product to discolor and corrode, and crack the product structure.
- Do not disassemble the product when it is wet, just let the product dry naturally .

Maintenance & Problems and Methods

- Make sure the power switch is turned off and the power plug has been unplugged from the outlet.
- Please keep the product clean and away from places with high humidity or high temperature and environments with volatile gases.
- If you do not use this product for a long time, please cover the product with a dust cap (please buy it yourself) to prevent dust from falling into it.

Special Warnings

1. Make sure to check for the following symptoms with year-round use:

There is burnt smell,sometimes power on and sometimes power off after touching the wire,the power cord is hot,there the power cord is damaged and other exceptions.

2. Turn off the power switch, unplug the power plug, and contact the dealer or our company where you purchased this product.
3. Do not disassemble and repair by yourself.

※ Problems and Methods

- If there is still abnormal phenomenon after the above treatment, please contact the dealer and request for maintenance.
- Please do not disassemble and repair the product without authorization, otherwise, the company is not responsible for any damage to the product, personal injury or any other loss.

Problems and Methods

Phenomenon	Possible problems	Methods
Automatic shutdown after startup	low battery	Please charge
No working	low battery	Use after charge it
	Working time is too long,motor overheating protection	Allow the motor to cool down naturally before continuing to use it
	Power switch not turned on	Turn on the power switch
Stop during use	Default working time is up	Restart it
	low battery	Use after charge it
When the product is working, there is an odor or smoke inside	Product malfunction	Please stop using, turn off the power and send for repair
Battery life becomes short and cannot be charged	Aging of battery	Please replace the battery
The charging light does not light up	Poor contact of Type-C socket/poor power contact	Unplug and then plug the Type-C charging cable/power cable
Hot compress temperature is not high	Battery is low	Please charge

Warnings for Product Discarding

Please support the recycling of rechargeable batteries (lithium-ion batteries).The lithium-ion battery used in this product is a valuable resource that can be recycled and regenerated. Please take out the lithium-ion battery when it is used up and discarded, and hand it over to a professional organization that can recycle and reprocess rechargeable batteries.

Dangerous

The battery used in this product is a professional rechargeable battery and cannot be used in other products. Do not charge the removed rechargeable battery.

- Do not throw in fire to overheat/do not pierce with nails, gravity shock, or dismantle/do not connect the positive and negative poles of the battery directly with metal/do not store it with necklaces, hairpins, etc/do not charge, use or place it in high temperature places such as the stove or the hot sun.

The above conditions can cause heat, fire, rupture and other phenomena.

Dangerous

- In case of accidental ingestion, please seek medical attention immediately, when the battery fluid leaks, do not touch the fluid with bare hands.
- When battery fluid gets into your eyes, there is a risk of blindness. Do not rub your eyes directly with your hands. Wash your eyes with clean water immediately and seek medical attention. When battery fluid gets on your body or clothes, it may cause skin irritation or injury. Please Rinse with water, if you still feel unwell, seek medical attention immediately

Recycle Guidance

In order to protect the environment, after this product or its parts are scrapped, you, as a consumer, are responsible for separating it from domestic waste and sending it to a qualified recycling site, and it will sort, dismantle, recycle and reuse in accordance with relevant national regulations and so on.

Name and Content Table of Hazardous Substances in Products

Name of Machine Parts	Hazardous Substances				
	Pb	Hg	Cr (VI)	PBB	PBDE
Outer shell	○	○	○	○	○
Circuit board assembly	×	○	○	○	○
Metallic parts	○	○	○	○	○
Rubber parts	○	○	○	○	○
Connecting line	○	○	○	○	○

* Circuit board assemblies include components made up of printed circuit boards such as resistors, capacitors, integrated circuits, connectors, etc.

○/×: Indicates that the content of the hazardous substance in all homogeneous materials of the part is below the limit requirement specified in GB/T26572 / exceeds the limit requirement specified in GB/T26572.



This electronic and electrical product is used under the conditions of use described in the manual, and the harmful substances or elements contained in it will not leak, and will not cause pollution to the environment or cause serious damage to people and property. The term is 10 years.

FCC Statement

This device complies with Part 15 of the FCC rules. Operation is subject to the following two conditions: 1) this device may not cause harmful interference, and 2) this device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC RF Radiation Exposure Statement

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.