

ibreoTM See 7

SIMPLY RELAX

Table of Contents

Cautions	03,04
Necessary knowledge	05
Box Contents	06
Main Features	07
Product Overview	08,09
Charging the Massager	10
How to Wear the Massager	11
How to Use the Massager	12,13
Bluetooth Connection	14
Troubleshooting	15
Maintenance & Care	16,17
Specifications	17
Disposal & Recycling	17
FCC Regulatory Compliance	18,19
Warranty	20, 21, 22

Cautions

To ensure the correct and safe use of this product, please read this User Manual carefully before use. iBreo is not responsible for any consequences resulting from improper use or failure to follow the instructions provided.

- This product is a household electronic health device and is not a substitute for medical equipment. It is not intended for diagnosis, treatment, or any commercial use.
- This product obstructs vision. Do not use it while walking or playing.
- Remove glasses or any objects that may interfere with the massage before use.
- Do not use this product if you have undergone surgery in the massage area or have skin injuries, inflammation, swelling, bruising, or other discomfort.
- Do not use this product if you have heart disease, a history of heart disease, acute illness, skin disease, malignant tumors, blood disorders with a tendency to bleed, osteoporosis, or soft tissue conditions.
- Individuals with mobility impairments, mental or physical discomfort, or cognitive disabilities should not use this product unless advised by a physician.

- It is not recommended to use this massager during pregnancy or menstruation. If necessary, consult a doctor before use.
- Do not use this product if you have a pacemaker, fractures, or metal implants in your body.
- Do not use this product for entertainment purposes. Store it in a place out of children's reach.
- Do not scratch or puncture this product with sharp objects.
- Children under 8 years old are prohibited from using this product. Children aged 8–12 must use it under the guidance and supervision of a parent or guardian.
- If any abnormal reactions occur during use, stop using the product immediately.
- This device has a heated surface. Users who are insensitive to heat should use it with caution.
- Use a charging device that matches the product's input voltage and follow electrical safety precautions. It is recommended not to mix and use different charging devices.
- The power supply used must not exceed an output power of 15W under normal or abnormal operating conditions.
- It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- This appliance contains batteries that are non-replaceable.

Necessary knowledge

※ Usage Frequency

- It is recommended to use the device twice daily, each session should be limited to 10 minutes.
- Allow an appropriate interval before using the device again.
- Use the device in an environment with a temperature range of 15°C to 30°C (59°F to 86°F) for optimal performance.


※ Precautions Before Use

- Do not operate the product with wet hands.
- Remove contact lenses, glasses, and false eyelashes before use to avoid discomfort during the massage.
- If you experience pain or a burning sensation in the massage area, stop using the device immediately and consult a doctor.


※ Product Cleaning

- To clean the exterior, use a slightly damp cloth or cotton to wipe the surface.
- Do not rinse the product with water.


Box Contents

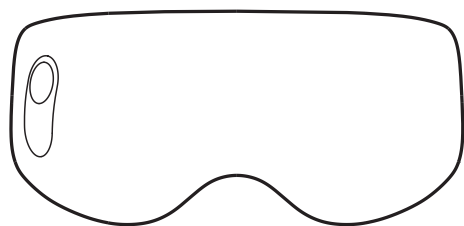


Eye Massager

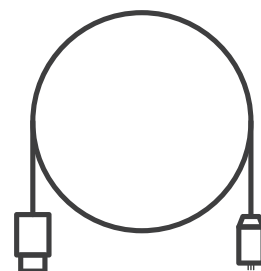

Charging Cable

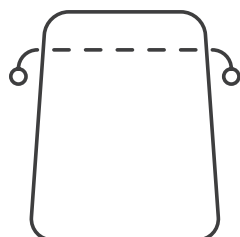


Pouch

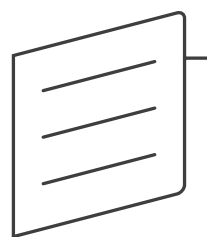

User Manual
& QSG

Main Features

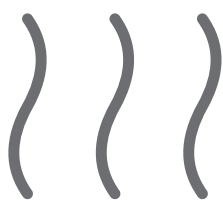

Heating
Comfort



Cooling
Relief



Vibration
Massage



Bluetooth
Music Play


Product Overview

Power Button
(On/Off, Mode Switch,
Voice Prompt Switch)

LED
Indicator

Multifunction Button
(Heating / Cooling /
Vibration Control)

Product Overview

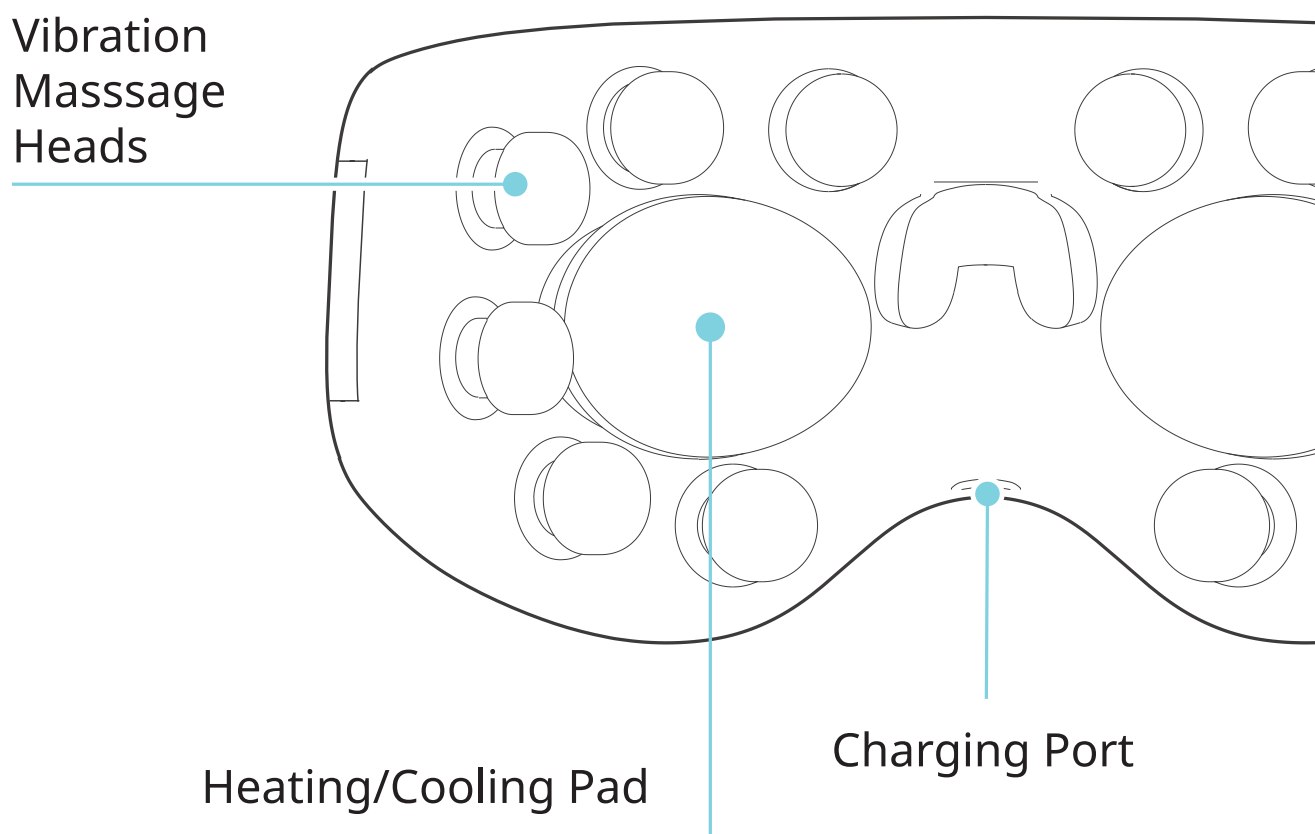

Charging the Massager


※ Charging adapter is not included.

※ Charging Time

- The device takes approximately 4 hours to fully charge.
- It is normal for the device and adapter to slightly heat up during charging.

※ Battery Indicator

- Low battery: The red indicator light blinks.
- Charging: The red indicator light breaths.
- Fully charged: The LED indicator light turns green and stays on.

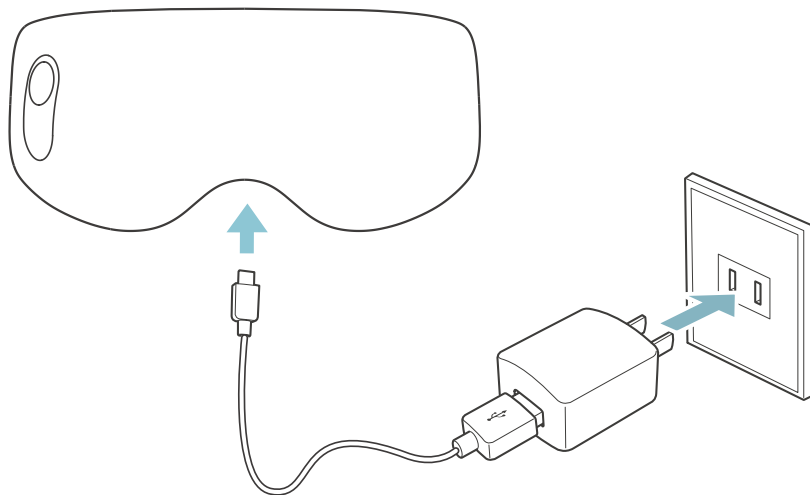

How to Wear the Massager

1. Remove the silicon covers from the heating/cooling pads.
2. Loosen the adjustable head strap.
3. Place the device on your head and adjust for a comfortable fit.
4. Power On: Press and hold the Power Button for 2 seconds until you hear a voice prompt, indicating the device has started.

※ Note:

- Do not wear the device too tightly.
- Ensure you have removed contact lenses, glasses, and false eyelashes before use.



How to Use the Massager

※ Power On

- Press and hold the Power Button for 2 seconds to turn on the device.

※ Change Massage Mode

- Press the Power Button briefly to switch between:
Eye Care Mode → Eye Massage Mode → Sleep Mode

※ Change Vibration Level

- Quickly single-click the Multifunction Button to cycle through:
Low → Medium → High → Turn Off Vibration

※ Change Heating/Cooling Level

- Quickly double-click the Multi-Function Button to cycle through:
Low Heat → High Heat → Low Cool → High Cool →
Turn Off Temperature

※ **Enable/Disable Voice Prompts**

- Quickly double-click the Power Button to toggle voice prompts ON or OFF.

※ **Power Off**

- Auto Power Off: The device will automatically shut down after 10 minutes of use.
- Manual Power Off: Press and hold the Power Button for 2 seconds to turn off the device.

※ **Note:**

- The cooling fan may continue running for one minute after the massager is powered Off. This is an intentional feature designed to prevent overheating.

Bluetooth Connection

- The device can stream music from your mobile device via Bluetooth.
- When the massage device is powered ON, it will enter pairing mode and await connection.
- Go to your mobile device's Bluetooth settings, locate "ibreo_see7", and connect.
- Once connected, the device will automatically reconnect to your mobile device each time it is powered ON.
- You can play, pause, and adjust the volume using your mobile device.

※ Note:

- The speaker in the eye massager is designed solely for relaxation during massage sessions. Its volume is moderate, and the sound quality is not comparable to a professional music player. Therefore, we recommend using it for soothing audio rather than playing popular music.

Troubleshooting

If the device malfunctions during use, check the following:

Issue	Possible Cause	Solution
Device does not turn on	Out of battery	Charge the device
Weak massage intensity	Low battery	Charge the device
Device does not heat	No heat in certain mode	Switch to other massage mode
Device does not cool	Temperature protection	Turn off device for at least 30 min

※ **Note:** If the issue persists, disconnect the power and seek professional repair. Do not attempt to disassemble the device yourself.

Maintenance & Care

- Keep the device dry, free from pests, and uncontaminated.
- Always turn off the device before cleaning.
- Store the device in a safe place after use. Do not place heavy objects on the device.
- If the device malfunctions, seek professional repair services. Do not attempt self-repair.


※ Storage:

- If the device is stored for more than 3 months, it should be inspected and recharged before use.
- Do not place the device directly on the floor, against walls, or against the ceiling.
- Store in a well-ventilated indoor space, avoiding direct sunlight and exposure to corrosive gases.

※ Transportation:

- The product packaging includes basic shock protection and is suitable for air, rail, road, and sea transport.
- Avoid exposure to rain, snow, overturning, or impacts during transportation.

Specifications

Item Name:	Eye Massager
Model No.:	See 7
Power Supply:	3.7V  2000mAh rechargeable lithium battery
Rated Power:	5W
Item Weight:	290g
Dimensions:	175mm (L) x 74mm (W) x 83mm (H)

Disposal & Recycling



- At the end of the product life cycle, do not throw this product into the normal household garbage, but bring it to a collection point for the recycling of electronic equipment and battery.
- Waste Electrical and Electronic Equipment may have potentially harmful effects on the environment.
- Improper disposal may lead to the accumulation of harmful toxins in air, water and soil, which is harmful to human health. You have the obligation to dispose of the device correctly. Consult your municipal authority for information about disposal.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Warning Statement

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the

equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Warranty

※ One-Year Limited Hardware Warranty

Your Massager ("Product") is covered by a One-Year Limited Hardware Warranty ("Warranty") against defects in materials and workmanship under normal use. This Warranty applies only to original purchasers and is limited to residents of the country/region of sale. It provides specific legal rights, which may vary by state.

※ Warranty Period

- The Warranty is valid for one (1) year from the date of purchase.
- The Warranty does not extend if the Product is repaired or replaced.
- Any changes to this Warranty will not apply retroactively.

※ What's Covered

If a hardware defect arises within the Warranty Period and a valid claim is made, iBreo will, at its discretion:

1. Repair the Product using new or refurbished parts that match new performance and reliability standards.

2. Replace the Product with a new or refurbished item of equivalent function and value.

Note: This Warranty allows only one replacement per purchase. If the Product is discontinued, an equivalent replacement may not be available.

※ Limitations & Exclusions

This Warranty does not cover damages caused by:

- Improper handling: Transport, storage, accidents, abuse, or external causes beyond our control.
- Misuse: Failure to follow product instructions or perform preventive maintenance.
- Unauthorized modifications or repairs.
- Normal wear and tear.

Important: Do not disassemble the Product. Unauthorized disassembly will void the Warranty. Repairs must be performed only by iBreo or an authorized service provider.

※ Legal Disclaimer

- The remedies under this Warranty are your sole and exclusive solutions.

- Our liability will not exceed the original purchase price of the Product.
- We are not responsible for any incidental, consequential, special, or punitive damages.
- Some states/countries do not allow limitations on implied warranties or exclusions of incidental damages, so these limitations may not apply to you.

All implied warranties, including but not limited to merchantability and fitness for a particular purpose, are limited to the duration of this Warranty.

ibreoTM

Copyright© 2025. iBreo. All Rights Reserved.