

尺寸: 110*145mm

Product Information

Product Description

This product is a multi-polar RF home beauty instrument (hereinafter referred to as RF instrument), using RADIO FREQUENCY (hereinafter referred to as RF).ELECTRICAL MUSCLE STIMULATION (hereinafter referred to as EMS), plasma hardening three major technologies, designed to improve a modern woman's Some facial skin problems.

Product Specifications

- Product Name: JS-02 frequency conversion RF instrument
- Product Model: JS-02
- Product Dimensions:181mm*40mm*40mm
- Product weight: 200g
- Rated input: 5V/1A
- Battery Capacity: 2000mAh
- Production country: Shenzhen, China


Product List

Radiofrequency instrument×1
Charging cable × 1
Packaging bags × 1
Instruction manual×1
User manual × 1

Instructions for use

Product diagram

1. Electrode piece (9 pieces)
- 2.LED light (internal)
- 3.Temperature sensor
- 4.Mode indicator
- FACE (Facial Mode)
- EYE (Eye Mode)
5. On/Off Mode key
6. Strength adjustment key
7. Strength adjustment key/power indicator
- 8.Type-C charging port

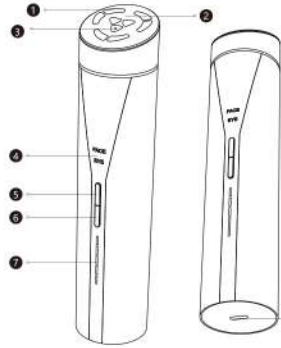

Operation Introduction

1. Charging

- Mode 1: Insert the charging cable Type-C port completely into the charging port of the host.
- Mode 2: Plug the charging cable Type-C connector into the charging dock charging port and place the host on the charging dock.
- When the intensity/power indicator 5 lights remain on, that is, the power is full.
- **The machine cannot be used in the charging state.
- ***The product should be in a flat state when charging to prevent tipping.

2. Switching on and off

- After power on, if the instrument does not have any operation within 5 minutes and does not touch the skin, it will turn off automatically, if you need to use it, please turn it on again.
- In the power on state, press and hold the switch/mode button for 2S to turn off the device directly.

3. Mode switching

- In the power-on state, press the On/Off/Mode button briefly to switch between [FACE (facial Mode)] and [EYE (eye mode)].

4. Strength selection

- The default intensity of the first boot is 1 gear (1 light); after the boot, short press the gear key to increase the gear; the highest gear is 5 gears (5 lights), the lowest grade is 1 grade (1 light).

5. Low power, lack of power alert

- When the power is on or in standby mode, the corresponding gear indicator will keep flashing when the power is low, so please charge it as soon as possible to avoid affecting the normal use.Please charge it as soon as possible to avoid affecting the normal use.
- When the power is off, when the switch/mode key cannot be turned on and the 1st gear indicator flashes 5 times and then goes off, it means that the power is insufficient, please charge it timely.Please recharge the battery in time to avoid affecting the normal use.

Nursing advice

- FACE (Face Care): Recommended to be used once every 2 days.
- EYE (eye care): Recommended to use once every 2 days.
- FACE(Face Care) mode uses the following areas of the face: left cheek, right cheek, and forehead.
- EYE (eye care)mode uses the eye area: eye bags, eye corners.

Facial Treatment Steps

Step 1: Makeup removal and cleaning

- Remove makeup and cleanse your face of oil and dirt to ensure your skin is in a clean state.

Step 2: Apply the gel

- Take out the gel, squeeze out an appropriate amount of gel and apply it to the center of the RF electrode piece or apply it to the area of the face (when you feel the machine pushing on the face during use, or feel the skin dry, please replenish the gel).When you feel the resistance of the machine on the face increase, or feel the skin dry, please replenish the gel in time).

Step 3: Facial Treatment

FACE (face mode) the whole face beauty process about 5-10 minutes, EYE (eye mode) the whole face beauty process about 5 minutes.

- Start the RF device and select the desired facial beauty mode.
- Press the RF device against the skin of the face and start working after a short vibration (refer to the face beauty technique illustration and operation instruction).
- If you want to continue to use another mode, short press the mode switch button and repeat the above steps, you can not use it on your face for more than 30 minutes per day.

Step 4: Cleaning and hydration

- After the facial, use a slightly damp soft cloth to wipe the gel left on the RF device thoroughly, dry it with a dry cloth and put itPut back in the box or storage bag.

Step 5: Machine maintenance

- After finishing, please use a slightly damp soft cloth to wipe the gel left on the RF meter thoroughly, dry it with a dry cloth and put it back into the box or storage bag. Put it back in the box or storage bag.

Nursing techniques (illustration and introduction)



Forehead manipulation techniques

- Left and right cheeks: press the machine firmly against the skin and glide it slowly over the skin in the inside-out direction, repeat several times.
- Forehead: press the machine against the skin and glide it slowly over the skin in the direction from bottom to top.
- Eyes: press the machine firmly against the skin and glide it slowly over the skin in the direction from the eye bags to the temples.

**During the use of the instrument, please keep the RF instrument always in a slow moving state, do not stay still.

Safety Precautions

Warm Tips:

Before use, be sure to read all [safety precautions] carefully and use this product in strict accordance with the instructions in this manual, failure to follow instructions for use may result in adverse effects,

Operation Note

- The recommended ambient temperature for this product is 10° C-30° C. Please keep the charging environment in the temperature range of 10-30° C.
- Be sure to charge the product first and make sure the buttons and indicators are normal when using it for the first time or after a long period of time.
- For the first time or before your skin gets used to the facial, start using it at intensity 1.
- If skin abnormalities occur, please stop using immediately and contact after-sales service or consult a professional doctor.
- When using, please pay attention to the radio frequency electrode piece to fit the skin, please be sure to slide on the skin to use, the recommended sliding speed: 2cm/s, if you feel too hot, you can slightly accelerate the sliding speed.
- Do not stick to the same position for a long time, otherwise it may cause burns, please slide on t be skin to use.
- Do not look directly at the LED light for a long time, otherwise it may cause eye injury.
- Do not apply the RF electrode to the eyelids and around the eyes, as this may cause eye damage.
- It is normal to feel different temperature when using RF because of the difference of individual skin perception.
- When the RF device is connected to the power supply for charging, it cannot be turned on or work, and all buttons are invalid.
- Before discarding the product, the battery must be removed from the product and disposed of safely. Take out the battery step is first make the product does not work to disconnect w ith the aid of a tool, pull the bottom cover off and pull out the internal battery. Batteries should be recycled in a specified manner.

Dangerous Matters

- This product is suitable for ages 18 years and older; please do not use this product on minors; please place this product out of reach of children ;keep this product out of the reach of children.
- Do not use for pet care as it may cause injury to your pet.
- Do not use the product for purposes other than those specified in the instructions, as this may result in personal injury or mechanical failure.
- Do not use for physically or mentally handicapped persons, as this may result in accidents or bodily injury.
- Do not use the machine in a malfunctioning state, as this may cause adverse effects on the body and pose a danger.
- Do not use near flammable substances (such as alcohol, thinner, gasoline, nail polish remover or styling spray), as this may result in fire/electric shock or bodily injury.
- When charging, please do not put in the fire, heating device, high temperature, use or place, otherwise it may cause the instrument to overheat, fire, burst.
- Do not use near water, in the bathroom or in the shower, as this may result in electric shock, short circuit and fire.
- Please avoid dropping or colliding with the main unit, otherwise it may cause malfunction or injury due to the fall of the main unit.Do not use the device if it has been damaged by dropping or collision.
- Do not let the electrode contact with metal, otherwise it may cause malfunction of the instrument and personal injury.
- Do not use solvents with corrosive properties (alcohol, thinner, gasoline, nail polish remover, etc.) to wipe the main unit, as this may damage the surface of the unit.
- Please turn off the power when idle and cleaning, otherwise it may cause electric shock, short circuit and fire.
- Do not disassemble or modify, otherwise it may cause fire, electric shock or physical injury, if you need to repair, please contact the aftermarket.
- Please unplug the charging cable immediately if any abnormal situation occurs during charging. If there is abnormal smoke, noise or odor when running, please stop using.

Unsuitable people

Do not use this product if you have any of the following conditions, otherwise it may cause accidents, skin and body problems.

- During menstruation, pregnancy or breastfeeding;
- during development or puberty;
- Metallic implants in the area of use;
- You are or have been suffering from the following diseases,such as cancer (especially skin cancer, precancerous nevus), acute illness, tuberculosis, febrile illness, infectious disease infectious diseases, blood diseases, visceral diseases, skin diseases, allergic dermatitis);
- The use of this product may cause discomfort or adverse effects on people with allergies, people with skin diseases and atopic dermatitis, people with sensitive skin, etc. discomfort or adverse effects, if the use of adverse effects, then the use of this appliance should be stopped immediately;
- Impaired immune system due to immunosuppressive diseases, such as AIDS or taking immunosuppressive drugs;
- Positive permanent fillers have been placed within the scope of use, or temporary subcutaneous fillers or botox have been implanted within the scope of use, please consult first the opinion of a physician;
- Any disease within the scope of use due to thermal irritation, such as recurrent rash virus;
- Physical irritation within the scope of use, etc. that can cause symptoms of blemishes;
- Allergic, sensitive or fragile skin;
- The body is in a state of excessive fatigue, sickness, after drinking or taking sleeping pills;
- Dilated capillaries caused by long-term use of steroids, hormones and other drugs, as well as liver dysfunction;
- People who use these medical devices: artificial lungs, cardiac pacemakers, internal defibrillators and other assistive devices implanted in the body, or electrocardiographs, etc.portable devices.

Not applicable to body parts

Do not use this product for the following body parts or if the following conditions exist

- Nodules, eyeballs, eyelids, inside the mouth, pubic area, mucous membranes, etc.;
- Severe acne, sores on the application site;
- Purulent and inflamed skin, previous scarring or abnormal wound healing;
- Undergoing surgery, such as laser peels, chemical peels or wounds healing within 3 months of the application area;
- Inflammation or redness of the skin due to sunburn;
- Itchy or pathological pigmentation on the skin;
- Floating moles, viral blemishes on the skin;
- Abnormalities of the skin such as rash, inflammation, itchy cancer;
- Areas with damaged skin integrity (e.g. with burns, blisters, incisional wounds, active skin disease or infection);

Note: The above list of safety precautions is not exhaustive. If you are unsure whether you can use this product, please consult your doctor or contact our company before making a decision.

Product Warranty Terms

Mainland China

- The company provides 1 year warranty service for customers who purchase this product in mainland China, the warranty period starts from the date of receipt (if there is no evidence to the contrary.The warranty period is calculated from the date of receipt of the goods (in the absence of evidence to the contrary, the date of receipt of the goods shall prevail for electronic transactions, and the date of physical delivery shall prevail for physical transactions).

- The company follows the "Electronic Commerce Law of the People's Republic of China", "Product Quality Law of the People's Republic of China", "Consumer Rights Protection Law of the People's Republic of China Consumer Rights Protection Law" and other relevant laws and regulations. For offline physical transactions, if quality problems occur within 7 days (inclusive) of the purchased products, you can contact the seller for return or exchange with the invoice or similar valid documents: within 15 days (inclusive) of the quality problems, you can contact the seller for exchange with the invoice or similar valid documents.If the quality problem occurs within 15 days (inclusive), the seller can be contacted for exchange with the invoice or similar valid documents. For online electronic transactions, the customer can return the goods within 7 days from the date of receipt without reason; the returned goods should be The customer shall bear the freight cost of the 7-day return, except for the return of goods caused by product quality problems.

- If you purchase this product in a country or region other than mainland China, the product warranty period, return policy, maintenance and other related matters Please contact your local authorized dealer for specific matters.

- Please ask for the invoice and other proof of purchase when you purchase the product and keep it properly. During the warranty period, if the product warranty period cannot be determined due to the inability to produce the invoice or valid purchase receipt, the production date on the product body production label will be used as the starting time of the warranty period.

Please note that free warranty service will not be available for the following cases during the warranty period.

- Improper use caused by man-made damage, such as access to improper power, the use of improper accessories, not in accordance with the instructions for use, wrong use or due to damage caused by negligence, etc. Repair or modification not authorized or approved by the Company; other damage caused by force majeure (such as natural disasters, voltage abnormalities, etc.)Damage caused by natural disasters, voltage abnormalities, etc.

- The product aging, wear and tear caused by normal use, but does not affect the normal use of the product.

- Free service does not include consumable parts and accessories, such as charging cables, etc: free service does not include shipping costs, etc.

Troubleshooting and Maintenance

Common fault problem table

Question	Possible reasons	Solutions
Unable to start the device	Battery depletion	Charge is complete before use
	Battery life exhaustion	Contact the manufacturer for repair or replacement
Automatic device shutdown	No operation for more than 5 minutes	Rebooting
	Battery depletion	Charge is complete before use
	Radiofrequency instrument electrode piece is not tightly attached to the skin	Keep the electrode of the RF instrument close to the skin
	Unused gel	Use of gel
	Insufficient amount of gel	Replenish the amount of gel on the face
Feeling excessive heat on the face	Too high gear	Adjusting the intensity of the instrument
	Instrument stays localized and does not move	Slightly accelerated instrument movement
	Insufficient amount of gel	Replenish the gel on the face
Face not hot enough	Too low gear	Boost intensity gears
	Instrument sliding speed is too fast	Slightly reduce the speed of RF meter movement
Current tingling sensation	Eye care mode being used	The normal mode is dominated by current sense
	The device is not pressed against the skin	Please hold the device close to the skin
Redness of the skin after use	The gear is too high	Downshift
	RF warming generated by the phenomenon	Short-term recovery is normal, if not recovered or other uncomfortable reactions, please seek medical attention promptly

JS-02 Variable Frequency RF Instrument - Instruction Manual



Note: Please use the product in strict accordance with this instruction manual. If you have any questions about this manual or the use of this product, please consult your doctor or contact our company before use.