


Yoga massager instruction manual

Thank you for purchasing our products to support your health needs

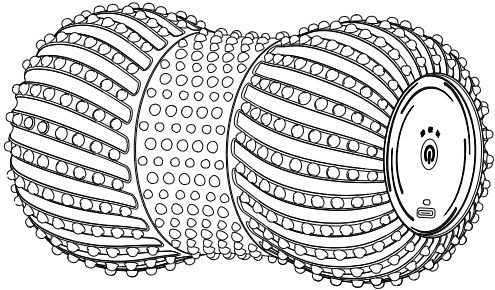

Please read this user manual carefully before use, so as to use the product correctly. Please carry the user manual with you for easy reference.

Critical safety note

When using this product, follow all instructions and basic precautions which should be followed at all times including the following:

Danger

1. To reduce the risk of electric shock: After use and before cleaning, be sure to unplug the product from the electrical outlet immediately.

warn

To reduce the risk of burns, fire, electric shock or personal injury:

- 1) This product should not be left unattended when plugged in. Unplug the socket when not in use and before using the product.
- 2) Do not operate under blankets or pillows. Overheating can cause fire, electric shock or personal injury
- 3) Do not carry this product through the power cord or use the power cord as a power cord as a handle.
- 4) Do not operate this product in the case of blocked pores. Keep the outlet free of fluff, hair, etc.
- 5) Do not drop or insert any object into any opening.
- 6) When using in outdoor humid environment, be sure to take precautions.
- 7) Do not operate in areas where aerosol (spray) products are used or where oxygen is available.
- 8) To disconnect, turn off all controls and remove the plug from the socket.
- 9) When children, persons with disabilities or disabilities use this device, or use it near them, they must be closely supervised to ensure that they do not misuse this product.

10) After prolonged use, the surface of the product may be hot. People who are sensitive to heat must be careful when using this product

11) Keep the power cord away from heated surfaces as it can cause damage, power shortages, fire and personal injury.Do not use in areas of the body with poor circulation.

13) This device is not intended for use by persons (including children) with physical, sensory or intellectual impairment or lack of experience and knowledge, unless They are supervised or instructed to use the device by the person responsible for its safety.

Be careful when using this product

1. Before each use or charge, carefully check the power unit, cable, and device for damage. Do not try to use Defective equipment.
2. Do not use electrical appliances or power cords if they show signs of damage, are in poor working condition, are dropped or are damp. To preventIn case of danger, be sure to take the equipment to the repair center.
3. Before cleaning, ensure that the switch is OFF and the power plug is removed from the socket. Please use dry cloth on the product Do the cleaning.
4. Keep the product clean and away from wet environment.
5. Please do not use benzene or diluents and other corrosive cleaners to clean the product: it may fade and cause damage.

You should consult your doctor if you have any of the following health problems:

1. Sensitive skin, circulatory diseases, blood clots or inflamed leg veins, deep vein thrombosis, open hands or feet Sexual wounds, thyroid problems, seizures, low platelet counts, or other blood problems that may make you more prone to bruising Injuries and bleeding.
2. Pregnant women should take precautions because some pressure points may cause contractions.

- ⚠ This device is not intended for commercial or medical use. If you have health problems, consult your doctor before using a smart massager
- ⚠ Don't wash! Don't dry clean! Can be wiped with a soft dry cloth!

Product parameter


Product name: Yoga massager	Product size:17*10.5*10CM
Range of use: hips/waist/limbs	Massage gear: 3 modes and 6 gears
Color: Black/blue/gray/orange	Enter: 5V1A
Material: Plastic. Hardware	Battery capacity: 3000 mah
Horn power: 4Q5W	Bluetooth name: Yoga massager

Product function

1. Press and hold the key for 2S to turn on/off. Start the M1 mode vibration,English prompt Bluetooth is turned on by default, charging position flashing blue light reminder.
2. Short press the open key to adjust the vibration gear. Double click on the switch to turn off Bluetooth.
- 3 Automatically turn off the motor 30 minutes after starting. Bluetooth doesnot affect, such as pressing the button can wake up the motor.
- 4 When the battery is low, the three indicators blink at the same time, indicating that charging is required. The red light is displayed when charging, and the green light is displayed when fully charged.

3 modes and 6 gears	flicker	No flicker
First gear M1(connected vibration)	Mode 1 The speed is slow	Mode 2 Fast
Second gear M2(vibrate 2 seconds, stop 1 second)	Mode 1 The speed is slow	Mode 2 Fast
Third gear M3(shake 1 second stop 1 second)	Mode 1 The speed is slow	Mode 2 Fast

Same function (two models)

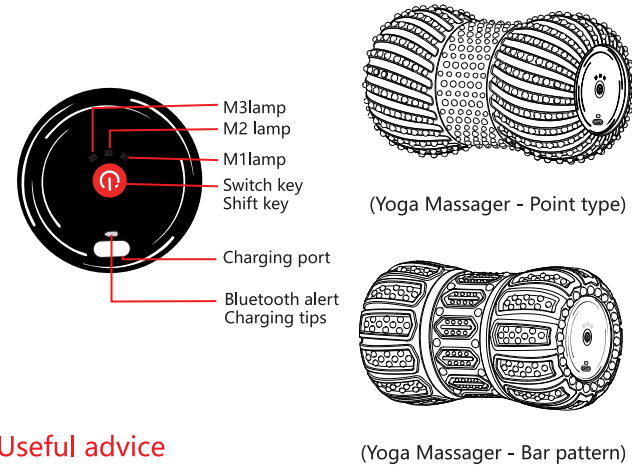

Useful advice

- In some applications, a non-slip surface such as a yoga mat is recommended for better traction.Topical (for in vitro) essential oils, such as lavender, peppermint, rosemary, and cedar oils, will help to slide the knob To the area being treated.

After-sales service

Name of consignee	
Delivery address	
Maintenance reason	
Product name	
Contact number	
Date of purchase	

1. If the product has quality problems, please contact us on working days, we will solve it for you as soon as possible
2. If the product can not be used due to quality problems, we support return and refund and do not accept any surface mail delivery.
- 3.The company reserves the right of final interpretation of after-sales service terms.

symptom	Lead to	solution
Machine not working	The machine is not charged or the charging plug is not properly inserted	Please insert the charging plug correctly and safely
	The Power On button is not pressed correctly	Press the on button for 1.5 seconds to start
	After 30 minutes of massage, the machine will automatically Off	For further massage time, please again Press the open key to use
Hot or unusable	You've had several 30-minute massages	If charging, unplug the charging plug and let set Ready cooling
Certificate of conformity		
Production date		
Quality check		

FCC Warning Statement: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

FCC Radiation Exposure Statement

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.