

HOMEDICS®

UNIVERSAL HANDLE

USER INSTRUCTIONS MANUAL

MODEL NO: HT100

FCCID: NS4HT100

PREFACE

The universal handle is a creative design in a series of HOMEDICS® dental product family.

Compared to other electric toothbrush already on the market, this product is more attractive and easy to be used.

Thank you for choosing HOMEDICS® UNIVERSAL HANDLE. Please fill in and mail the warranty card enclosed with this manual to ensure registration of your universal handle warranty.

FCCID: NS4HT100

IMPORANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING

1. **DANGER** - To reduce the risk of electric shock:

- **DO NOT** handle the plug with wet hand
- **DO NOT** place or store the charger base where it can be pulled into a tub sink or toilet.
- **DO NOT** reach for a charger base that has fallen into water. Unplug immediately.
- **DO NOT** submerge the charger base in water or any other liquid.
- This charger base is designed for use in 110 volt, 60 Hz outlet.

2. **WARNING** - To reduce the risk of burns, electric shock, fire or injury

- **DO NOT** excessively bend or pull the cord. Unplug the cord by grasping the plug-not the cord.
- Keep the cord away from heated surfaces
- Use this product only for its intended use as described in this manual. **DO NOT** use the product for any other purpose.
- **DO NOT** force over-sized tooth brush handle into the product.
- Remember to apply reasonable inserting force to prevent any damage to the mechanism, and to be able to pull your toothbrush out as well.
- **NEVER** stick your toothbrush in or out of the universal handle while the product is on.
- **NEVER** operate the product if it is not working properly, if it has been dropped or damaged, or it is immersed in water.
- **NEVER** use while sleeping or drowsy. .
- Close supervision is necessary when this product is used by, on or near children or disabled persons.
- Avoid using a brush with crushed or bent bristles.
- Contact your dentists prior to use if you have had oral or periodontal surgery in the previous twelve months.

SAVE THESE INSTRUCTIONS


FCCID: NS4HT100

PRODUCT FEATURES

1. The HOMEDICS® UNIVERSAL HANDLE is compatible for most *standard* types of toothbrushes known to us on the market.
2. The movement of the toothbrush is linear, which means that the user brushes his tooth in an horizontal movement back and forth, the universal handle will vibrate in the same direction.
3. Change your personal toothbrush into an electric toothbrush at once.
4. Rechargeable batteries.
5. Replaceable toothbrush, especially suitable in family use.

PRODUCT SPECIFICATION




- Weight:
 - (a) Handle: 119 gram (0.26 lb.)
 - (b) Charger: 146 gram (0.32 lb.)
- Volts: 110
- Watts: 3w
- Frequency: 60 hz.
- Dimension:
 - (a) Handle: 7inch long x 1.5 inch width.
 - (b) Charger: 4 inch long x 3 inch width x 2.5 inch height.
 - (c) when handle placed into charger, the gross height is 7.5 inch.

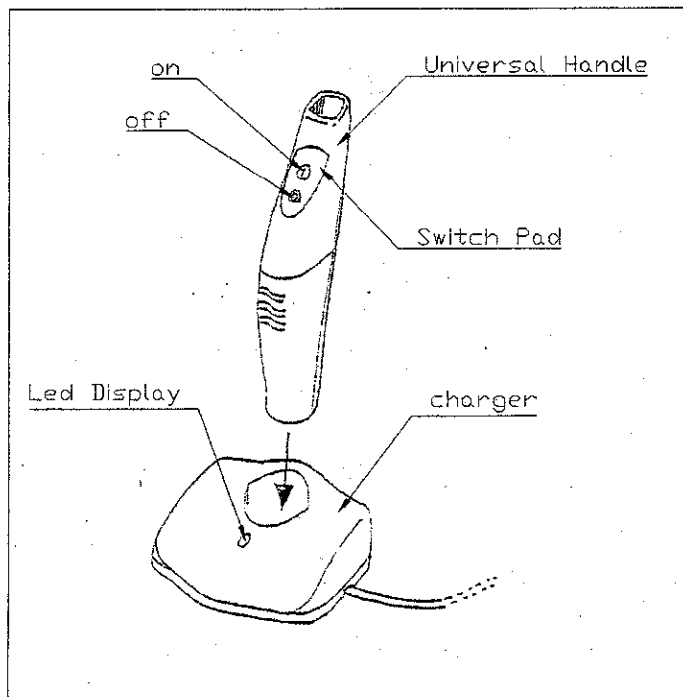


FCCID: NS4HT100

HOW TO USE

1. Choose a suitable size personal toothbrush. Over-sized or irregular one is inadequate.
2. With moderate force, insert the toothbrush with bristles at the same side of the switch pad, so that about $\frac{2}{3}$ to $\frac{1}{2}$ of its stem will penetrate the universal handle.
3. To assure the toothbrush is working properly. Press the 'I' symbol on the switch pad to start for a few seconds and then press the 'O' symbol to terminate.
4. Apply toothpaste to the bristle. Any toothpaste or gel may be used.
5. Turn on the universal handle. Your toothbrush is vibrating, ready to brush your teeth.
6. Wash the bristle of your toothbrush under water stream, so that the rinse will not pour into the opening of the universal handle.
7. Dry the handle. Use a towel to dry the universal handle. Replace the dry handle in the base for charging and safe storage.

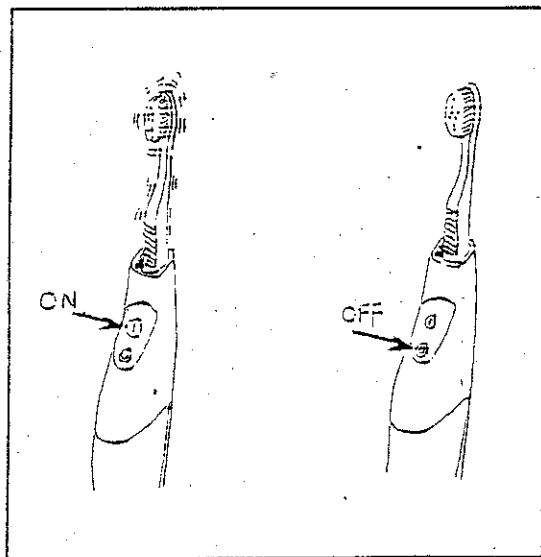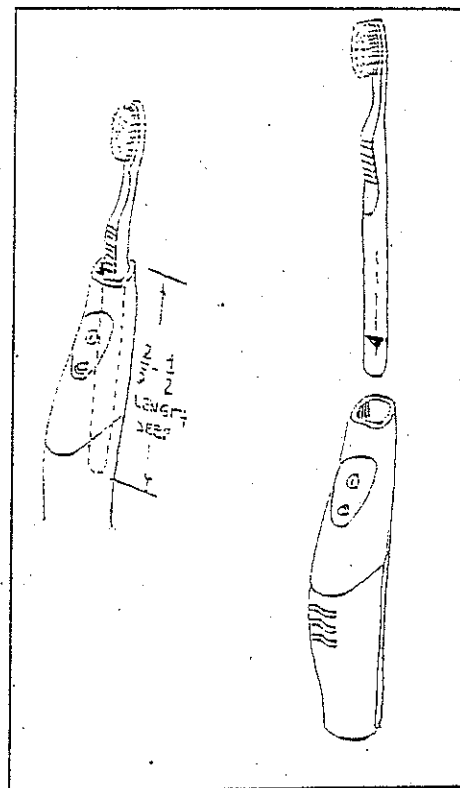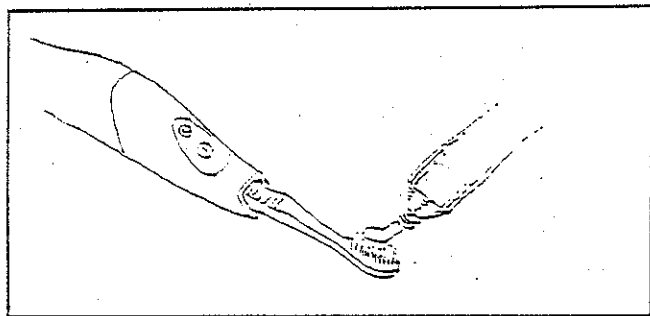

FCC ID: NS4HT100

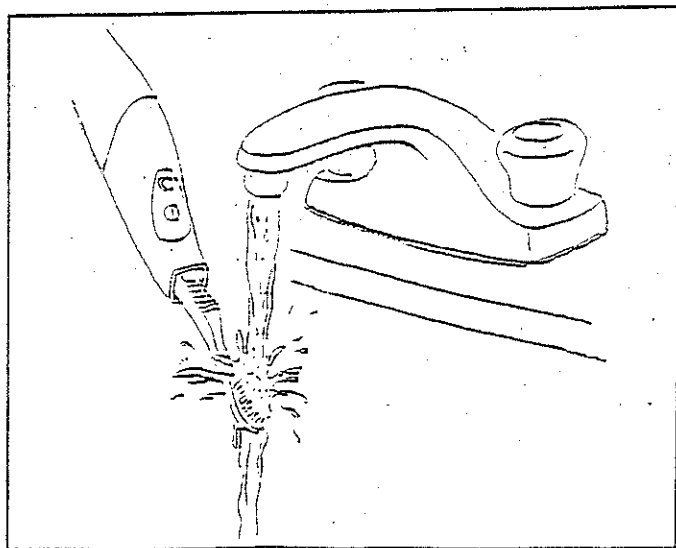

Charging

To ensure that the charger will function correctly, it is important that you follow the procedure below you first use the appliance.

1. **Plug the charger base into a suitable mains outlet** and turn it on. Place the HOMEDICS® UNIVERSAL HANDLE in the base and leave to charge for 24 hours.
2. **The charger led will be lighted up in red when the handle is placed inside the charger in right position.** Light out when the universal handle disengages the base.
3. **The light may blacken or dark red if the universal handle is not placed in right position**, this may cause inadequate charging.
4. **The best voltage to indicate the light function is 110 v.** It is more bright if the voltage ascend to 120v, there is nothing in properly charging.
5. **Unplug the base and use the universal handle as normal.** Do not recharge it until the batteries are completely discharged.
6. **When the batteries are completely discharged**., plug in the charger base again and recharges for 24 hours before use.

Follow points 5 & 6 every six months to help preserve the batteries.

NICKEL- CADMIUM BATTERIES REMOVAL

To help protect our environment, the nickel-cadmium batteries within the HOMEDICS® UNIVERSAL HANDLE should be removed at the end of its life, for safe disposal.

To remove the batteries

1. Split the universal handle open by hitting it with a hammer along the cosmetic line. The handle should then be cracked into top and bottom halves. Discard the bottom half.
2. Once split open, completely submerge the top half of the handle in water for a minimum of 10 minutes to completely discharge the batteries.
3. Take the universal handle out of the water and pull the battery chassis out and disengages the battery.
4. Once the battery is removed , cover any exposed metal with tape.

For environmental reasons, do not discard the batteries with your regular refuse. Contact your local authority for more details on recycling and safe disposal facilities in your area.

MEDICAL SAFEGUARDS

Please note the following instructions for your continued safety.

1. Contact your dentist before using the HOMEDICS® UNIVERSAL HANDLE if you have had oral or periodontal.
2. Contact your dentist if excessive bleeding occurs or if bleeding continuous after the first week of use.

FCCID: NS4HT100

FCC Information to the User

This device complies with Part 18 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference that may cause undesired operation.

Emissions Note: this equipment has been tested and found to comply with the limits pursuant to part 18 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is needed.
- Consult the dealer or an experienced radio/TV technician for help.

FCCID: NS4HT100