Important Precautions

Please read before use of product.

Beam toothbrushes are safe when used within the safety guidelines. Beam toothbrushes are an electrical appliance and should be handled cautiously especially around children.

Dangers

To lessen the risk of injury or fire:

1. Do not use any power source other than the provided AC/DC adapter.
2. Do not store the charger near bathtubs, toilets, and sinks.
3. Do not attempt to grab or use a charger if it has fallen into the water. Unplug from the outlet immediately.
4. Do not use the brush while bathing.
5. Do not attempt to clean the charger while plugged into an electrical source.
6. Do not attempt to plug the charger into any other voltage other than specifically stated on the product.
7. Do not forcefully plug the charger into the outlet.
8. Do not use a charger if the plug does not easily fit into the outlet.
9. Do not attempt to take apart the charger.

Caution

To lessen the risk of injury or fire:

1. If the brush or charger appears to be damaged in any way, please discontinue use of the product. As further use could result in injury.
2. If the brush is damaged and can no longer be used the battery must be removed from the product before it is discarded. Do not attempt to remove the battery if the brush is not ready to be discarded. Please refer to the battery disassembly below for more information on how to remove the battery.
3. Once the battery has been removed please discard safely.
4. The Beam toothbrush has been specifically designed to clean teeth, gums, and tongue. Do not attempt to use this product in any other way than as previously described.
5. Do not store brush or charger near a heat source.
6. Beam toothbrush is NOT dishwasher safe. Do not attempt to wash this product in the dishwasher.
7. Please refrain from storing or using this product outdoors.
8. Supervision is required for children or people with reduced physical, sensory, or mental capabilities to safely use this product. Without these specific supervision requirements this product could lead to injury.
9. The Beam toothbrush contains a battery that is non-replaceable.

Medical Caution

If you have received any injuries from this product please consult your dental professional.

How to remove and recycle the battery

Please note: This process will fully disable the product and will not be able to be undone. The steps to remove the battery from the brush are listed below.

Steps to remove battery

1. Before you start the removal process please make sure that the brush has no charge left and the battery is completely drained. You can accomplish this by turning the brush on and letting it cycle. Repeat until the battery is drained.
2. Once you have made sure that the battery has been fully drained and the brush has no power you will need to remove the battery from the brush. For this you will need a screwdriver and channel locks.
3. Take the channel locks and place it 1/4 inch away from the bottom of the brush approximately where the beam logo is on the front and proceed to apply pressure with the channel locks. Rotate brush 90° and apply pressure, repeat until cap can be removed.
4. Once the cap is off, then proceed to take the screw driver and pull back on 2 locking tabs inside the bottom of the brush.
5. Take off the brush head and proceed to turn the cap upside down and press it against a table to be able to remove the carrier from the handle. The carrier houses the motor, circuit board, and the battery.
6. Once the carrier is free please push the battery from the back and pop the battery out of the carrier. Twist the battery wires until the wires come loose from the circuit board, watch to make sure the positive and negative wires do not come in contact with one another. Once the battery is free it needs to be properly recycled please read below to learn the proper way to recycle the battery.

How to properly recycle the battery

These batteries can not be thrown away since they contain harmful materials that can be a hazard to the environment. To recycle the battery it will need to be disposed of at your local waste management facility. For more information please ask your local waste management facility for information on how to dispose of the battery.

Setup Guide

Beam Brush Pearl

beam®

Welcome to your new brush

Unlock the full power of the Beam Brush

Included in box:

- 1 Beam Brush handle
- 2 Removable brush heads
- 1 Charging base
- 1 AC/DC charging adaptor

Compatible devices:

- iOS 10+
- Android 5.0+

Using your new Beam Brush

Charge

1. Connect your brush to the charger. The button light pulses white while charging. When the lights turn solid white, your brush is fully charged.
2. When your battery is low, you will receive an audible low battery tone along with a red light indicator that lasts for 10 seconds.
3. We recommend charging your battery for a full 16 hours before the first use. A fully charged battery should provide at least 2 weeks of regular use (2 mins, 2x daily).


Note: Keep the charger in a dry area, and make sure the brush is clean and dry before charging.

Connect to the App

1. Download the Beam on the Apple App Store or Google Play. Sign into your Beam account to get started.
2. Place your brush on the charger. Navigate to the Brush tab and click **“Set up your Beam Brush.”** When the screen indicates it is searching for your brush, press the button on the brush for 3 seconds until the light flashes blue to enter provisioning mode.

Note: You can only pair 1 brush at a time.

3. When the app has discovered your brush, *Beam Brush 3.0*, select it and follow the in-app instructions to finish set-up.
4. Once your brush is connected to the app, you can customize your brushing experience, access your digital insurance card, view claim information and more!

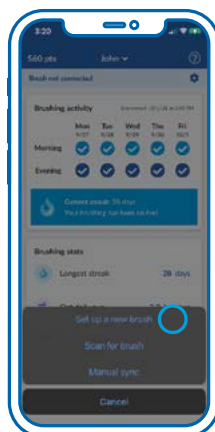Customize your brush

Picking your brush mode

When turned on, your brush defaults to the standard Clean brushing mode. To switch brushing modes, navigate to the **“Brush Settings”** tab in your Beam app and select **“Brush Modes.”**

Brush Modes


Clean: our standard brushing mode

Whitening: made to combat surface stains (*note: this mode runs for 2 minutes and 30 seconds*)

Sensitive: a softer mode for sensitive gums

Kids: a gentle mode for children’s teeth

Note: Intensity can be adjusted for each mode. Setting changes may take up to 15 seconds after the brush is turned on to update.

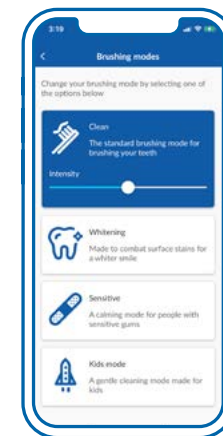Picking brush colors

Pick out your favorite color for your brush light ring in the Brush light settings. You can also adjust the brightness to get your perfect choice.

Brush Indicator Lights

Brushing: solid white

Low Battery: fade to red followed by solid red

Charging 0%-99%: slow white pulse

Fully charged: solid white

Pairing: pulse blue

Pairing Failure: solid red

Pairing Successful: solid white

Firmware Update: blink white and blue

A healthier, happier smile

Taking care of your Beam Brush

Keep it clean. After brushing, rinse the brush head and handle separately, and wipe dry. Never place the Beam brush in the dishwasher or submerge in water. Clean with a damp cloth.

Keep it fresh. Dentists recommend replacing your brush head at least every 90 days with normal use. Using the Beam app, you can track the life of your brush head, simply navigate to the Brush tab. You can order new brush heads via our app or shop.beam.dental.

Brushing tips

1. Brush for two minutes twice a day
2. Brush each quadrant (upper left/right, lower left/right) for 30 seconds.
3. Brush along the gum line
4. Angry gums? Try our Sensitive brush mode to help clean gentle areas.
5. Don't forget to brush the roof of your mouth and tongue.
6. No need to press hard - let the brush do the work!

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help.beam.dental | (800) 648-1179 | help@beam.dental

Beam Brush — Instructions for Use

SAVE THESE INSTRUCTIONS. PLEASE READ ALL INSTRUCTIONS BEFORE USE.

FCC Statement

This device is equipped with Transmitter Module FCC ID: O7P-4343.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: **(1)** this device may not cause harmful interference, and **(2)** this device must accept any interference received, including interference that may cause undesired operation.

The Beam Brush's wireless operation is safe and complies to RF Exposure requirements. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Any unauthorized modification to the devices included in this Beam product could void the user's authority to operate the product.

This device with transmitter module has been tested to SAR and complies with FCC exposure requirements for portable devices. SAR testing has been done at a distance of 0mm from the body.

Industry Canada Statement

This device contains Transmitter Module IC: 10147A-4343.

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This device complies with Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions: **(1)** this device may not cause interference, and **(2)** this device must accept any interference, including interference that may cause undesired operation of the device.

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The Beam Brush's wireless operation is safe and complies to RF Exposure requirements. This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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Under Industry Canada regulations, this radio transmitter may only operate using an antenna of a type and maximum (or lesser) gain approved for the transmitter by Industry Canada. To reduce potential radio interference to other users, the antenna type and its gain should be so chosen that the equivalent isotropically radiated power (e.i.r.p.) is not more than that necessary for successful communication.

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This device with transmitter module has been tested to SAR and complies with ISED exposure requirements for portable devices. SAR testing has been done at a distance of 0mm from the body.

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