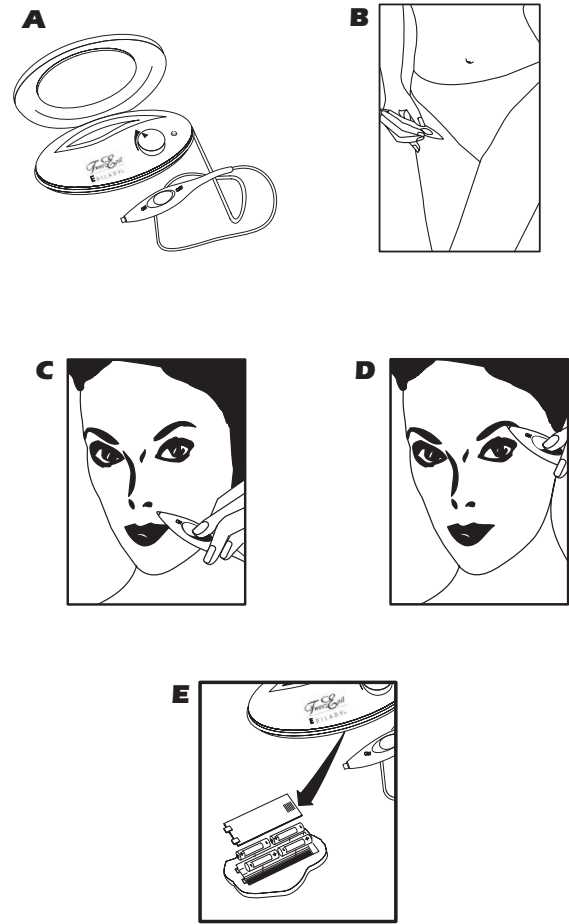
67405001

EPILADY®

Instructions for use

TweezEpil

English
Français
Deutsch
Español
Italiano
Nederlands

ENGLISH

TweezEpil is an electric tweezer especially designed for painless and long-lasting removal of hair from the root (Picture A).

The hair's root is destroyed by the transmission of radio waves through the hair to the root, dissolving it by electrolysis. After about 30 seconds the hair can be easily removed. The tweezer can remove hairs from all parts of the body. TweezEpil is battery operated, uses four AA alkaline batteries and does not need any electric outlet. Its strength selector monitors the intensity of the radio waves. The tweezer's locking device relieves you from keeping the operating button pressed throughout each hair treatment (Picture B). After 30 seconds a buzzing tone will be heard. A gentle pull will then remove the hair painlessly.

Please read the following instructions carefully.

Safety Instructions

- Use the device only for its intended purpose and in accordance with the user instructions. Any damage caused by disregarding these instructions shall be at the user's sole responsibility.
- DO NOT use the tweezer to remove hair from nostrils, ears or eyelashes.
- DO NOT use the device in the bath or shower.
- Pregnant women, patients after heart transplants or with pacers, and patients under medical supervision should consult their physician before using the device.
- If the device becomes damaged, take it to one of the service stations. **Do not attempt to repair it yourself.**
- Prevent children from playing with or handling the tweezer or any of its parts.

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ENGLISH

Guarantee

Mepro Epilady Ltd. makes every effort to ensure that its products are manufactured without any defect. However, in the rare event of defects due to faulty workmanship or materials, we guarantee to repair or replace the product, or any part thereof, free of charge.

- This guarantee is granted to the original purchaser only and is valid for one year from the date of purchase indicated below by the sales agent.
- This guarantee does not apply if the product is altered or repaired or otherwise tampered with by unauthorized persons and/or is used other than in accordance with the manufacturer's instructions and/or is subjected to other than reasonable use of a product of its the nature.
- This guarantee in no way affects your rights under statutory law.

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