



SIMPLY RELAX

Breo iSee M

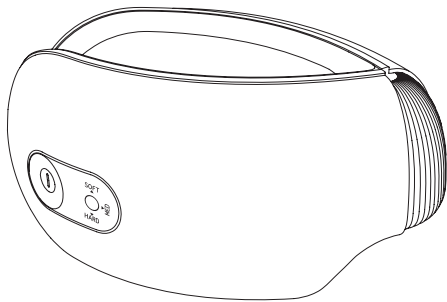

breo®

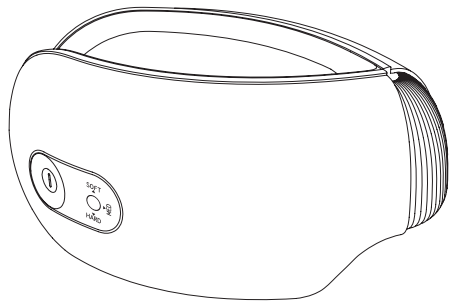

MENU

Caution	4
Product Configuration	5
Product Characteristics	6
Product Functions	7
Troubleshooting	8
Protection	9
Product Information	10



Thank you for your purchasing Breo!

Please read this user's manual carefully prior to use and, keep it properly for future reference.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio/TV technician for help.
- changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.


Caution

Note: To ensure safety in operation, please read the Use Instruction Manual carefully prior to operation. The Company should not be liable for the consequences of improper use or non-compliance with the instructions for the product.

1. Do not use the product when driving a vehicle.
2. Prior to operation, remove the jewelry or articles that may affect the massage from the body part to receive massage.
3. Do not use the product in a humid environment such as a bathroom.
4. Please keep the product away from water source, source of ignition and corrosive environment.
5. Observe safety regulations when you use the product in a specific environment.
6. Do not use the massager on body parts subjected to surgical operation.
7. Do not use the massager on body parts suffering from skin injury, swelling, inflammation or blood stasis.
8. Those with heart disease or a history of heart disease, acute disease, skin disease, malignant tumor or physical abnormality should not use the product.
9. Those with blood disease showing bleeding tendency, osteoporosis, soft tissue disease and apoplexy should not use the product.
10. Those with leaden feet, physical and mental illness, or mental retardation should not use the product, and should follow the doctor's advice if they need to use it.
11. If you have any discomfort during use, stop using the product.
12. Pregnant women and menstrual women should not use the massager.
13. Children under 6 years of age are prohibited from using the product, and children aged 6-12 can use the product only when they are subject to guidance and care from their parents or supervisors.
14. To clean the appearance, use alcohol and cotton instead of water.
15. Do not scratch or stab the product using a sharp object.
16. This product is a domestic electronic health-care product. The product is not designed to work as a medical device, and it does not apply in diagnosis, treatment and other commercial uses.
17. Never fail to charge the product using charging equipment provided by the Company. Pay attention to the electrical safety.
18. Once any exception occurs during operation, power off the product immediately. Please do not disassemble, reassemble and repair the product without consent.

Product Configuration

Massager, charging cable, storage bag, and instruction manual

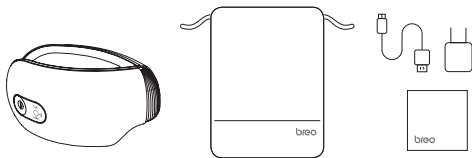

Product Characteristics



kneading-pressing acupoints



hot compress with
constant temperature



intelligent gas pressure



meditation guidance

Product Functions

1. **Tragen:** Der Benutzer stellt die Länge des elastischen Bandes entsprechend der Größe des Kopfumfangs auf die bevorzugte Position des Benutzers ein;
2. **Einschalten:** Im ausgeschalteten Zustand drücken Sie lange auf die Taste, das Produkt gab einen "Tropf" Ton aus;
3. **Modusauswahl:** Drücken Sie kurz den Modusschalter, das Produkt kann zwischen MED, HARD, SOFT und benutzerdefinierten Modus umgeschaltet werden;
4. **Musik-Player:** Nachdem das Bluetooth mit dem Mobiltelefon verbunden wurde, können Musik in der APP oder lokale Musik zur Wiedergabe an das Produkt gesendet werden;
5. **Ausschalten:** Drücken Sie lange die Taste, und das Produkt wird ausgeschaltet, nachdem es "tropft";
6. **Strommenge:** Wenn die Stromversorgung des Produkts niedrig ist, beginnt das Batteriesymbol zu blinken;
7. **Aufladen:** Während des Ladevorgangs blinkt das Batteriesymbol. Wenn das Batteriesymbol aufleuchtet, zeigt es an, dass es voll ist;

Troubleshooting

If a fault occurs during operation, do check it. In case the fault cannot be eliminated, or the product still fails to work satisfactorily, please turn off the massager and ask a professional technician to repair it. Please do not disassemble the product without consent!

Basic fault	Possible cause	Basic troubleshooting method
Failure to start up	Check whether the battery level is low	If yes, please charge in time
Low massage strength	Check whether the battery level is low	If yes, please charge in time

Protection

Storage

To use the machine when a storage period of over 3 months expires, never fail to perform an open-case inspection and recharge it. Never allow the machine to abut against the floor, wall, and roof when you place it. The room should be well ventilated, without intense sunlight exposure and corrosive gas erosion.

Transportation

The product package with simple anti-vibration functions applies in air, railway, highway and ship transportation, and should be protected from rain, snow, splash, inversion and collision.

Power Check

Check the power supply. If the supply voltage exceeds the specified range, the machine should not be turned on.

Maintenance

Damp-proof, anti-mouse and anti-pollution measures should be taken. The machine shall be shut down before it is cleaned. Place the product in a safe place after use, and do not put any heavy object on it. The machine shall be repaired by professional technician in case any fault occurs. The non-professional are not allowed to disassemble and reassemble the product without consent.

Abandonment and Recovery

This product should not be discarded along with household waste. Please follow local environmental laws and guidelines when you decide to dispose the product.



Product Information

Product name: Breo iSee M

Model: see m1

Voltage: 5V---

Power: 5W

Battery: 1150mAh

Net weight: 275g

Dimensions: 210mm × 75mm × 110mm

