breo.iDream5

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USER MANUAL

Thank you for purchasing Breo

Please read this manual carefully before using and keep it properly for future reference

Contents

Product function	1
Inside the package	2
Product Deployment	3
Operation	4-5
Precautions	6
Troubleshooting	7

Product function

Point Massage

Heating

Air Pressure


Music

1

Inside the package


Your massager includes main unit, power adapter, earphone, user manual

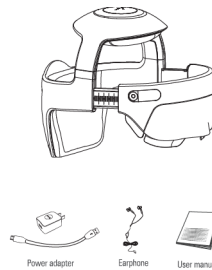2

Product Deployment

Main unit

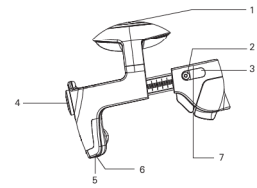1. Message head height adjust bottom
2. On/off /mode switch
3. OLED display
4. Size adjust Knob
5. Earphone jack
6. DC power interface
7. Disassemble button for eye part

3

Operation

Preparation before use

1. Check if there are missing parts of the package.
2. Make sure the main unit with enough power.
3. Make sure there is no damage of the adaptor.

Charge the unit

1. Open the cover of ⑥ DC power interface. Insert the adaptor.
2. When charging, the battery status on ③ OLED display will flashing when fully charged, it will stop flashing.

Attention

1. Make sure the main unit is power off when charging. Charging time is about 4 hours.
2. Do not use the product when charging. The main unit doesn't work while charging.
3. Download APP for customizing massage.

How to wear

1. Press the button in the middle of ④ Size adjust Knob and wear the main unit.
2. Turn ④ Size adjust Knob to adjust the size to fit your head and eyes.

4

Operation

3. The top side size could be adjusted after machine on.
Short press ① Message head height adjust bottom for shorter the size(maximum 3 times), and long press for higher size.
4. The eye part could be disassembled before you wear the main unit by ⑦ Disassemble button for eye part.

Operation& Mode Selection

Short press ② On/off /mode switch, MED mode will on, Short press ② On/off /mode switch to switch the massage mode.

Operation& Mode Selection

MED	HARD
Air pressure	Air pressure
Heating	Heating
Music	Kneading
15 min	15min

Shutdown

Long press ② On/off /mode switch.

5

Precautions

1. Do not use the head massager when driving.
2. Do not use the massager in places too hot and wet.
3. Check before using the massager: No hairpin, hair accessory on your head and no wear contact lenses.
4. Wear the massager before switch on to avoid overpressure damage to the air bag.
5. The massager is suitable for adults of 56cm-64cm head circumference.
6. For those who suffer from heavy inflammation, blood stasis or brain surgery, do not use the massager.
7. For those who suffer from heavy slow action, poor physical and mental function or disturbance of intelligence, use under the supervision of companion, and stop using immediately if there is any discomfort.
8. For those who suffer from apoplexy, high blood pressure and low blood pressure, please use with caution and consult professional physician, and stop using immediately if there is any discomfort.
9. Use the main blood OCSU adaptor, and use safely.
10. The product is not suitable for diagnosis, treatment and other commercial purposes.
11. Use sterile cotton with neutral detergent for cleaning, rather than water.
12. The massager has heating function, those who are insensitive to temperature shall pay attention.

6

Troubleshooting

If there is any problem when you use the device, please follow the instructions below to check. Do not disassemble or repair the device yourself.

Faults	Possible reasons	Solutions
Unable to boot	Low battery	If yes, please charge the main unit
		If not, check the connection of charging cable
Air pressure or vibration is over low	Low battery	Please charge the machine

Technical parameters

DC	5V/1A
Power	5W
Use environment	Temperature: 5°C-40 °C, Relative Humidity: ≤ 80%, Atmospheric Pressure: 860hpa-1060hpa
Dimension	L254mm*W210 mm*H42mm
Net weight	1075g

Storage

1. If the storage period exceeds 3 months, please check and recharge it.
2. Do not store close to floor, wall or roof.
3. Keep in good ventilation condition and away from direct sunshine and corrosive gases.

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27

· This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

· **The statements should be displayed in the user manual :**
changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Product name: breo iDream 5S

Model No.: idream 5s

Voltage : 5V --- 1A

Power : 5W

Net Weight : 960g

Dimensions : L254×W210×H242mm

FCC ID: PXUIDREAM5S