

Congratulations...

on purchasing the Crest SpinBrush Pro Recharge – its dual moving head and ergonomic design have been developed with Crest Dental Professionals to:

- Deliver effective cleaning
- Whiten stained teeth*
- Fight gingivitis, an early form of gum disease*

With proper brushing, SpinBrush Pro Recharge is designed to help you keep your smile healthy-looking and beautiful for life.

* With regular brushing

The Crest SpinBrush Pro Recharge has been type tested and found to comply with the limits for Part 18 of the FCC rules, which are designed to provide reasonable protection against such interference. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try and correct the interference by one or more of the following measures:

- 1) Re-orientate the receiving antenna.
- 2) Relocate the microwave with respect to the receiver.
- 3) Plug the microwave into a different outlet so that the receiver and microwave are on different branch circuits.

If necessary the user should consult the dealer or an experienced radio/television technician for additional suggestions.


Note: Modifications should only be carried out by the manufacturer or authorized representative to ensure continuing conformance.

Recharging Base

- The base of your SpinBrush Pro Recharge comes with a brush head holder that can conveniently hold 2 extra replaceable heads (not included).

- To attach the brush head holder, simply place it onto the bottom of the recharging base at a 45 degree angle and then turn it to lock in place.

- You can also mount your recharging base to the wall using 2 screws (not included).

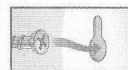

IMPORTANT SAFEGUARDS

Crest SpinBrush Pro Recharge has been evaluated by Underwriters Laboratories (UL) and conforms to UL standards for quality electrical products.

READ ALL INSTRUCTIONS BEFORE USING

DANGER:

To reduce the risk of electrocution:

- Do not store the charging unit in a position where it can fall into a tub or sink.
- Do not place or drop charging unit into water or other liquid.
- Do not touch charging unit that has fallen into water. Unplug unit immediately.
- Do not use while bathing or showering.

WARNING:

To reduce the risk of damage to unit, burns, electrocution, fire, or injury to persons:

- Close supervision is necessary when this product is used around, on or by children under 12.
- Use as instructed, and do not use attachments not recommended by the manufacturer.
- Do not operate this unit if the cord has been damaged, the plug has been damaged, or if the charger has been dropped into water. Return the product to a service center for examination & repair.
- Keep cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

SAVE THESE INSTRUCTIONS

It is recommended that the handle be stored in the charging unit to maintain full power and maximize battery life.

Battery removal / disposal instructions:

The nickel-metal hydride batteries inside your Crest SpinBrush Pro Recharge will eventually wear down and no longer recharge. They cannot be replaced but can easily be removed for disposal or recycling. When the brush needs to be disposed of, simply remove the battery by following the instructions below.

1. Discharge the battery completely by turning on the brush outside of the recharging base and letting it run until it stops.
2. Unscrew the bottom battery cap of the brush using the bottom of the brush-head holder as a tool.
3. Remove the battery cap and pull out the batteries from the brush.
4. Dispose of batteries according to local environmental regulations, or at an appropriate recycling collection site.

Caution: Opening the handle will destroy the appliance.

Note: You will not be able to place new batteries in the handle.

Customer Service Information

Thank you for purchasing a Crest SpinBrush Pro Recharge brush. If you are in the US or Canada and have any questions, please call: 1-800-561-0752 or visit:

www.spinbrush.com

Total Satisfaction Guarantee

We guarantee that you will be totally satisfied with the oral care you receive from your SpinBrush Pro Recharge or we will refund your money.

If you are not fully satisfied with your Crest SpinBrush Pro Recharge, simply call 1-800-266-2209 within 90 days of purchase with the original UPC and receipt.

Getting Started

Charging instructions:

Crest SpinBrush Pro Recharge is electrically safe when used as directed, and can be used in the bathroom.

1. Plug Charging unit into an electrical outlet.
2. Place the brush handle in the charging unit.
3. The green charging light (6) will pulse when the brush is charging and remain on when the brush is fully charged.
4. Charge your brush for approximately 16 hours before your first use.
5. You're ready to brush – just like you normally would with your manual brush.

Using Your Crest SpinBrush Pro Recharge:

- You can use any toothpaste with your Crest SpinBrush Pro Recharge


- To avoid splattering, place the brush head in your mouth before turning on the brush

Turning the Brush On and Off

- Turn on the brush by gently pressing the on/off button (4). When you take the brush out of the base, it may feel warm to the touch.

- You can turn the brush off manually by pressing the on/off button at any time.

- The brush will automatically turn off when placed in the charging base (7).

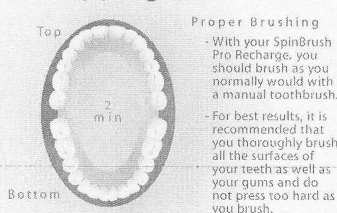

Proper Brushing

- With your SpinBrush Pro Recharge, you should brush as you normally would with a manual toothbrush.

- For best results, it is recommended that you thoroughly brush all the surfaces of your teeth as well as your gums and do not press too hard as you brush.

Brushing-Timer

- Your SpinBrush Pro Recharge is equipped with a brushing timer to help you brush for a recommended 2-minutes.
- The brush will pause momentarily when 2 minutes of brushing is complete.

Your Replaceable Dual Action Brush Head

- Your SpinBrush Pro Recharge features a dual moving head designed to deliver great cleaning and care.

- It features colored bands that allow you to customize your brush head. This is useful for those who share a rechargeable brush.

- Change brush head every 3-4 months due to normal wear.

- For optimum performance, replace your brush head with the Crest SpinBrush Pro head specially designed for your Pro Recharge-Brush. These are available at retail stores.

CLEANING AND MAINTENANCE

- After brushing, turn on the brush and rinse the brush head for several seconds under running water.

- Periodically, clean the brush handle and charging base with a damp cloth.


- To extend the life of your brush head, periodically remove the head and rinse it out by running water through the neck of the brush head.

- Do not submerge the toothbrush handle or the charger into water for cleaning.

- Store the product in dry cool place when it is not in use.

- Do not wrap the cord around the charger unit.


- This product has no user serviceable parts.

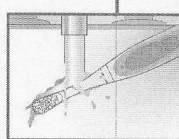

DESCRIPTION

of Your Crest SpinBrush Pro Recharge:

- 1 Easy-replace brush head
- 2 Powerful dual moving head
- 3 Custom color band
- 4 On/off button
- 5 Ergonomic, soft grip handle
- 6 Charge indicator
- 7 Compact charging base
- 8 Convenient brush head holder

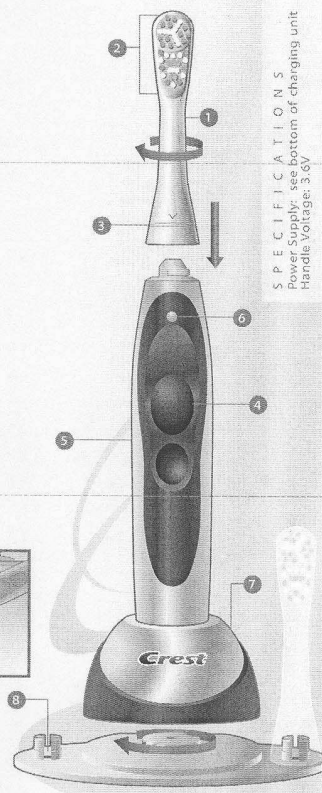

SPECIFICATIONS
Power Supply: AC 120V
Handle Voltage: 5.0V

00000000

FINAL ART 03.07.04

Interbrand

111 East High Street, Cincinnati, Ohio 45202

PROJECT INFORMATION		STUDIO INFORMATION	
Product: Crest Rechargeable SpinBrush (U.S.) – Franklin Description/Size: SpinBrush (English/French) User Manual Rev	Job No: 202-2000	File No: 1201-114010 (User Manual Rev) 348524	Studio Contact: Tom Aberg
Project #: P-050204	Blueprint #: 001	Brand Code: 001	
Ref/PMs #: 94327081	UPC #: 001	PLU/Item #: 001	
New IPMS #: 001	Last Revision Date: 03/07/04		
PROCTOR & GAMBLE		TECHNICAL INFORMATION	
PRODUCTION MANAGER	INITIALS	Application/Version:	Initials/Revision: 0.0
LEGAL		Layers/Pages:	000001
REGULATORY		Placed Images:	Interbrand Breakdown Front
TRANSLATION		Fonts Used:	Helvetica, Helvetica Cond., Helvetica Neue, Times, Helvetica 277
Signatures of approvers are on file with the originals kept in the P&G Design Department.		Separator: Kim Berg	Release Date:
		Approved For Release	/ /
COMPUTER EQUIVALENTS			
PMS 100	PMS 105	PMS 150	PMS 155
PMS 160	PMS 165	PMS 170	PMS 175
PMS 180	PMS 185	PMS 190	PMS 195
PMS 200	PMS 205	PMS 210	PMS 215
PMS 220	PMS 225	PMS 230	PMS 235
PMS 240	PMS 245	PMS 250	PMS 255
PMS 260	PMS 265	PMS 270	PMS 275
PMS 280	PMS 285	PMS 290	PMS 295
PMS 300	PMS 305	PMS 310	PMS 315
PMS 320	PMS 325	PMS 330	PMS 335
PMS 340	PMS 345	PMS 350	PMS 355
PMS 360	PMS 365	PMS 370	PMS 375
PMS 380	PMS 385	PMS 390	PMS 395
PMS 400	PMS 405	PMS 410	PMS 415
PMS 420	PMS 425	PMS 430	PMS 435
PMS 440	PMS 445	PMS 450	PMS 455
PMS 460	PMS 465	PMS 470	PMS 475
PMS 480	PMS 485	PMS 490	PMS 495
PMS 500	PMS 505	PMS 510	PMS 515
PMS 520	PMS 525	PMS 530	PMS 535
PMS 540	PMS 545	PMS 550	PMS 555
PMS 560	PMS 565	PMS 570	PMS 575
PMS 580	PMS 585	PMS 590	PMS 595
PMS 600	PMS 605	PMS 610	PMS 615
PMS 620	PMS 625	PMS 630	PMS 635
PMS 640	PMS 645	PMS 650	PMS 655
PMS 660	PMS 665	PMS 670	PMS 675
PMS 680	PMS 685	PMS 690	PMS 695
PMS 700	PMS 705	PMS 710	PMS 715
PMS 720	PMS 725	PMS 730	PMS 735
PMS 740	PMS 745	PMS 750	PMS 755
PMS 760	PMS 765	PMS 770	PMS 775
PMS 780	PMS 785	PMS 790	PMS 795
PMS 800	PMS 805	PMS 810	PMS 815
PMS 820	PMS 825	PMS 830	PMS 835
PMS 840	PMS 845	PMS 850	PMS 855
PMS 860	PMS 865	PMS 870	PMS 875
PMS 880	PMS 885	PMS 890	PMS 895
PMS 900	PMS 905	PMS 910	PMS 915
PMS 920	PMS 925	PMS 930	PMS 935
PMS 940	PMS 945	PMS 950	PMS 955
PMS 960	PMS 965	PMS 970	PMS 975
PMS 980	PMS 985	PMS 990	PMS 995
PMS 1000	PMS 1005	PMS 1010	PMS 1015
PMS 1020	PMS 1025	PMS 1030	PMS 1035
PMS 1040	PMS 1045	PMS 1050	PMS 1055
PMS 1060	PMS 1065	PMS 1070	PMS 1075
PMS 1080	PMS 1085	PMS 1090	PMS 1095
PMS 1100	PMS 1105	PMS 1110	PMS 1115
PMS 1120	PMS 1125	PMS 1130	PMS 1135
PMS 1140	PMS 1145	PMS 1150	PMS 1155
PMS 1160	PMS 1165	PMS 1170	PMS 1175
PMS 1180	PMS 1185	PMS 1190	PMS 1195
PMS 1200	PMS 1205	PMS 1210	PMS 1215
PMS 1220	PMS 1225	PMS 1230	PMS 1235
PMS 1240	PMS 1245	PMS 1250	PMS 1255
PMS 1260	PMS 1265	PMS 1270	PMS 1275
PMS 1280	PMS 1285	PMS 1290	PMS 1295
PMS 1300	PMS 1305	PMS 1310	PMS 1315
PMS 1320	PMS 1325	PMS 1330	PMS 1335
PMS 1340	PMS 1345	PMS 1350	PMS 1355
PMS 1360	PMS 1365	PMS 1370	PMS 1375
PMS 1380	PMS 1385	PMS 1390	PMS 1395
PMS 1400	PMS 1405	PMS 1410	PMS 1415
PMS 1420	PMS 1425	PMS 1430	PMS 1435
PMS 1440	PMS 1445	PMS 1450	PMS 1455
PMS 1460	PMS 1465	PMS 1470	PMS 1475
PMS 1480	PMS 1485	PMS 1490	PMS 1495
PMS 1500	PMS 1505	PMS 1510	PMS 1515
PMS 1520	PMS 1525	PMS 1530	PMS 1535
PMS 1540	PMS 1545	PMS 1550	PMS 1555
PMS 1560	PMS 1565	PMS 1570	PMS 1575
PMS 1580	PMS 1585	PMS 1590	PMS 1595
PMS 1600	PMS 1605	PMS 1610	PMS 1615
PMS 1620	PMS 1625	PMS 1630	PMS 1635
PMS 1640	PMS 1645	PMS 1650	PMS 1655
PMS 1660	PMS 1665	PMS 1670	PMS 1675
PMS 1680	PMS 1685	PMS 1690	PMS 1695
PMS 1700	PMS 1705	PMS 1710	PMS 1715
PMS 1720	PMS 1725	PMS 1730	PMS 1735
PMS 1740	PMS 1745	PMS 1750	PMS 1755
PMS 1760	PMS 1765	PMS 1770	PMS 1775
PMS 1780	PMS 1785	PMS 1790	PMS 1795
PMS 1800	PMS 1805	PMS 1810	PMS 1815
PMS 1820	PMS 1825	PMS 1830	PMS 1835
PMS 1840	PMS 1845	PMS 1850	PMS 1855
PMS 1860	PMS 1865	PMS 1870	PMS 1875
PMS 1880	PMS 1885	PMS 1890	PMS 1895
PMS 1900	PMS 1905	PMS 1910	PMS 1915
PMS 1920	PMS 1925	PMS 1930	PMS 1935
PMS 1940	PMS 1945	PMS 1950	PMS 1955
PMS 1960	PMS 1965	PMS 1970	PMS 1975
PMS 1980	PMS 1985	PMS 1990	PMS 1995
PMS 2000	PMS 2005	PMS 2010	PMS 2015
PMS 2020	PMS 2025	PMS 2030	PMS 2035
PMS 2040	PMS 2045	PMS 2050	PMS 2055
PMS 2060	PMS 2065	PMS 2070	PMS 2075
PMS 2080	PMS 2085	PMS 2090	PMS 2095
PMS 2100	PMS 2105	PMS 2110	PMS 2115
PMS 2120	PMS 2125	PMS 2130	PMS 2135
PMS 2140	PMS 2145	PMS 2150	PMS 2155
PMS 2160	PMS 2165	PMS 2170	PMS 2175
PMS 2180	PMS 2185	PMS 2190	PMS 2195
PMS 2200	PMS 2205	PMS 2210	PMS 2215
PMS 2220	PMS 2225	PMS 2230	PMS 2235
PMS 2240	PMS 2245	PMS 2250	PMS 2255
PMS 2260	PMS 2265	PMS 2270	PMS 2275
PMS 2280	PMS 2285	PMS 2290	PMS 2295
PMS 2300	PMS 2305	PMS 2310	PMS 2315
PMS 2320	PMS 2325	PMS 2330	PMS 2335
PMS 2340	PMS 2345	PMS 2350	PMS 2355
PMS 2360	PMS 2365	PMS 2370	PMS 2375
PMS 2380	PMS 2385	PMS 2390	PMS 2395
PMS 2400	PMS 2405	PMS 2410	PMS 2415
PMS 2420	PMS 2425	PMS 2430	PMS 2435
PMS 2440	PMS 2445	PMS 2450	PMS 2455
PMS 2460	PMS 2465	PMS 2470	PMS 2475
PMS 2480	PMS 2485	PMS 2490	PMS 2495
PMS 2500	PMS 2505	PMS 2510	PMS 2515
PMS 2520	PMS 2525	PMS 2530	PMS 2535
PMS 2540	PMS 2545	PMS 2550	PMS 2555
PMS 2560	PMS 2565	PMS 2570	PMS 2575
PMS 2580	PMS 2585	PMS 2590	PMS 2595
PMS 2600	PMS 2605	PMS 2610	PMS 2615
PMS 2620	PMS 2625	PMS 2630	PMS 2635
PMS 2640	PMS 2645	PMS 2650	PMS 2655
PMS 2660	PMS 2665	PMS 2670	PMS 2675
PMS 2680	PMS 2685	PMS 2690	PMS 2695
PMS 2700	PMS 2705	PMS 2710	PMS 2715
PMS 2720	PMS 2725	PMS 2730	PMS 2735
PMS 2740	PMS 2745	PMS 2750	PMS