Congratulations...

on purchasing the Crest SpinBrush Pro Recharge – its dual moving head and ergonomic design have been developed with Crest Dental Professionals to:

- Deliver effective cleaning
- Whiten stained teeth\*
- Fight gingivitis, an early form of gum disease\*

With proper brushing, SpinBrush Pro Recharge is designed to help you keep your smile healthy-looking and beautiful for life.

\* With regular brushing

IMPORTANT SAFEGUARDS

Crest SpinBrush Pro Recharge has been evaluated by Underwriters Laboratories (UL) and conforms to UL standards for quality electrical products.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

To reduce the risk of electrocution:

- Do not store the charging unit in a position where it can fall into a tub or sink.
- Do not place or drop charging unit into water or other liquid.
- Do not touch charging unit that has fallen into water. Unplug unit immediately.
- Do not use while bathing or showering.

WARNING:

To reduce the risk of damage to unit; burns; electrocution, fire, or injury to persons:

- Close supervision is necessary when this product is used around, on or by Children under 12 or invalids.
- Use this product only for its intended use as instructed in this Manual, and do not use attachments not recommended by the manufacturer.
- Do not operate this unit if the cord has been damaged, the plug has been damaged, or if the charger has been dropped into water. Return the product to a service center for examination & repair.
- Keep cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

SAVE THESE INSTRUCTIONS Information

It is recommended that the handle be stored in the charging unit to maintain full power and maximize battery life.

Battery removal / disposal instructions:

The nickel-metal hydride batteries inside your Crest SpinBrush Pro Recharge will eventually wear down and no longer recharge. When the brush needs to be disposed of, simply remove the battery by following the instructions below.

1. Discharge the battery completely by turning on the brush outside of the recharging base and letting it run until it stops.
2. Hold brush handle in right hand, bristles facing down, grab battery cap with left hand, apply downward force until cap snaps.
3. Remove the battery cap and pull out the batteries from the brush.
4. Dispose of batteries according to local environmental regulations, or at an appropriate recycling collection site.

Caution: Opening the handle will destroy the appliance. Note: You will not be able to place new batteries in the handle.

Customer Service Information

Thank you for purchasing a Crest SpinBrush Pro Recharge brush. If you are in the US or Canada and have any questions, please call: 1-800-561-0752 or visit:

www.spinbrush.com

Total Satisfaction Guarantee

We guarantee that you will be totally satisfied with the oral care you receive from your SpinBrush Pro Recharge or we will refund your money.

If you are not fully satisfied with your Crest SpinBrush Pro Recharge, simply call 1-800-266-2209 within 90 days of purchase with the original UPC and receipt.

This product has no user serviceable parts. The Crest SpinBrush Pro Recharge has been type tested and found to comply with the limits for Part 18 of the FCC rules. There is no guarantee that interference will not occur in a particular installation. If this equipment does cause interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try and correct the interference by one or more of the following measures:


- 1) Re-orientate the receiving antenna.
- 2) Relocate the brush with respect to the receiver.
- 3) Plug the brush into a different outlet so that the receiver and brush are on different branch circuits.

If necessary the user should consult the dealer or an experienced radio/ television technician for additional suggestions.

Note: Modifications should only be carried out by the manufacturer or authorized representative to ensure continuing conformance.

DESCRIPTION of Your Crest SpinBrush Pro Recharge:

- 1 Easy-replace brush head
- 2 Powerful dual moving head
- 3 Custom color band
- 4 On/off button
- 5 Ergonomic, soft grip handle
- 6 Charge indicator
- 7 Compact charging base

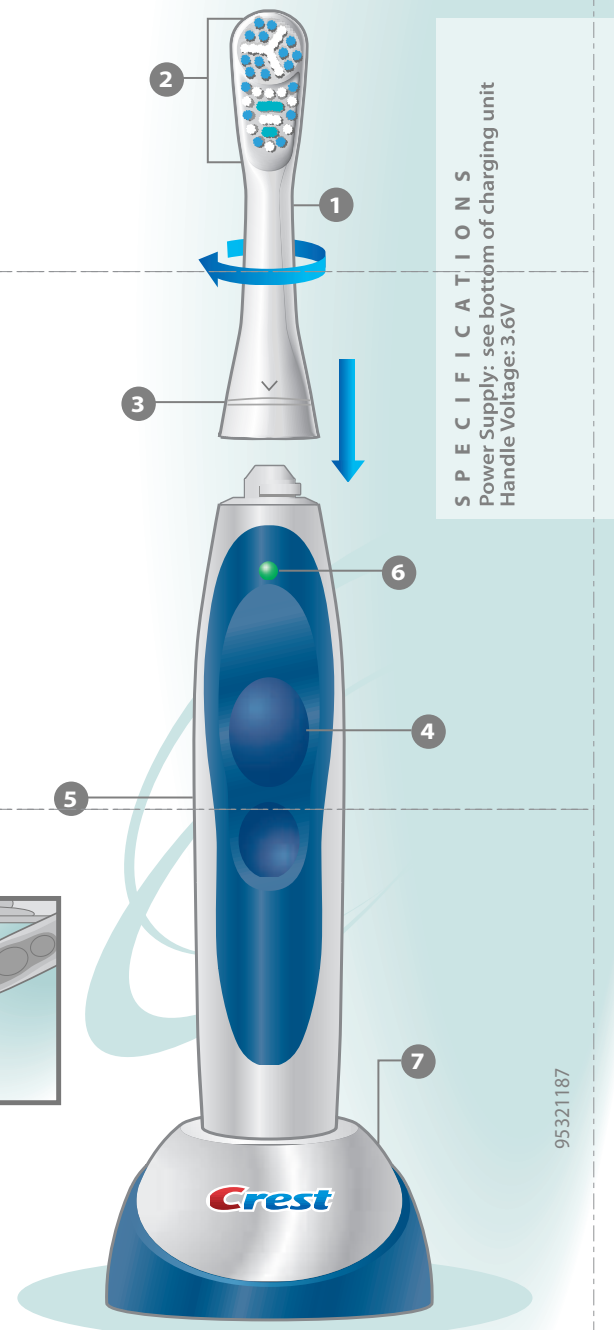Power Supply: see bottom of charging unit  
Handle Voltage: 3.6V

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