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2024_12_25

activate your warranty and meet your product

Scan the QR code **within 7 days** of receiving your product to activate your warranty and experience all the benefits!

body composition glo smart scale
with ambient illumination

INSTRUCTION MANUAL AND WARRANTY INFORMATION

SCL-BBC200, SCL-BBC250, SCL-BBC350, SCL-XBC200 | L-05257, Rev. 1

IMPORTANT SAFETY INSTRUCTIONS:

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ AND SAVE THESE INSTRUCTIONS.


- DO NOT use if you have implanted medical devices such as a cardiac pacemaker or artificial lung.
- DO NOT use while wearing devices such as a wearable electrocardiograph.
- DO NOT use if you are pregnant.
- Keep batteries away from children.
- DO NOT use on slippery floor.
- This product is intended for use as a body composition analyzer for personal fitness and health monitoring only.
- This product should not be used for treating or diagnosing any medical conditions.
- Always consult your physician or healthcare professional prior to making any changes to your diet or exercise regimen.
- DO NOT reach for a unit that has fallen into water.
- DO NOT use while bathing or while in a shower.
- DO NOT place or store appliance where it can fall or be pulled into a tub or sink.
- DO NOT place in or drop into water or other liquid.
- DO NOT use outdoors. This product is intended for household use only.
- DO NOT strike or drop the scale.
- DO NOT stand on the edge of the scale.
- DO NOT jump on the scale.
- DO NOT step onto a wet scale. Wet glass is slippery. Always ensure the scale surface is dry and feet are dry before use.
- DO not place in direct sunlight, near heating equipment, air outlet of air conditioning.
- Always place the scale on a hard, dry, flat surface before measurement.
- Check the scale before use. DO NOT use the scale if the glass is cracked or broken. If you see a crack, contact Customer Support.
- DO NOT overload the scale. Maximum weight capacity is 400lbs (180kg or 28st:7lb).
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- The measurement data may not be accurate if using the scale immediately after intense exercises.
- Replacement of a battery with an incorrect type that can defeat a safeguard.
- Do not disposal of a battery into fire or a hot oven, or mechanically crushing or cutting of a battery, that can result in an explosion;
- Do not leaving a battery in an extremely high temperature surrounding environment that can result in an explosion or the leakage of flammable liquid or gas.
- Do not expose the battery to extremely low air pressure that may result in an explosion or the leakage of flammable liquid or gas.

Operating Environment	Temperature: 59°F - 104°F (15°C - 40°C) Humidity: 30 - 80% RH Atmospheric Pressure: 970hPa - 1050hPa
Storage Environment	Temperature: -31°F - 140°F (-35°C - 60°C) Humidity: 10 - 90% RH Atmospheric Pressure: 700hPa - 1060hPa


SYMBOL	EXPLANATION
	This symbol means “Manufactured by” and provides the name and address of the legal manufacturer.
	Date of manufacture.
	This symbol refers to the batch code of the device.
	Type B Applied Part, B indicates that it is used for the body. This device is connected to the main voltage.
	This symbol indicates to read the manual before use
	This symbol indicates the product may not be discarded as household waste.
	This symbol indicates the product is designed for indoor use only.

GET TO KNOW YOUR GLO DIGITAL SCALE WITH AMBIENT ILLUMINATION

SCL-BBC200, SCL-BBC250 AND SCL-BBC350

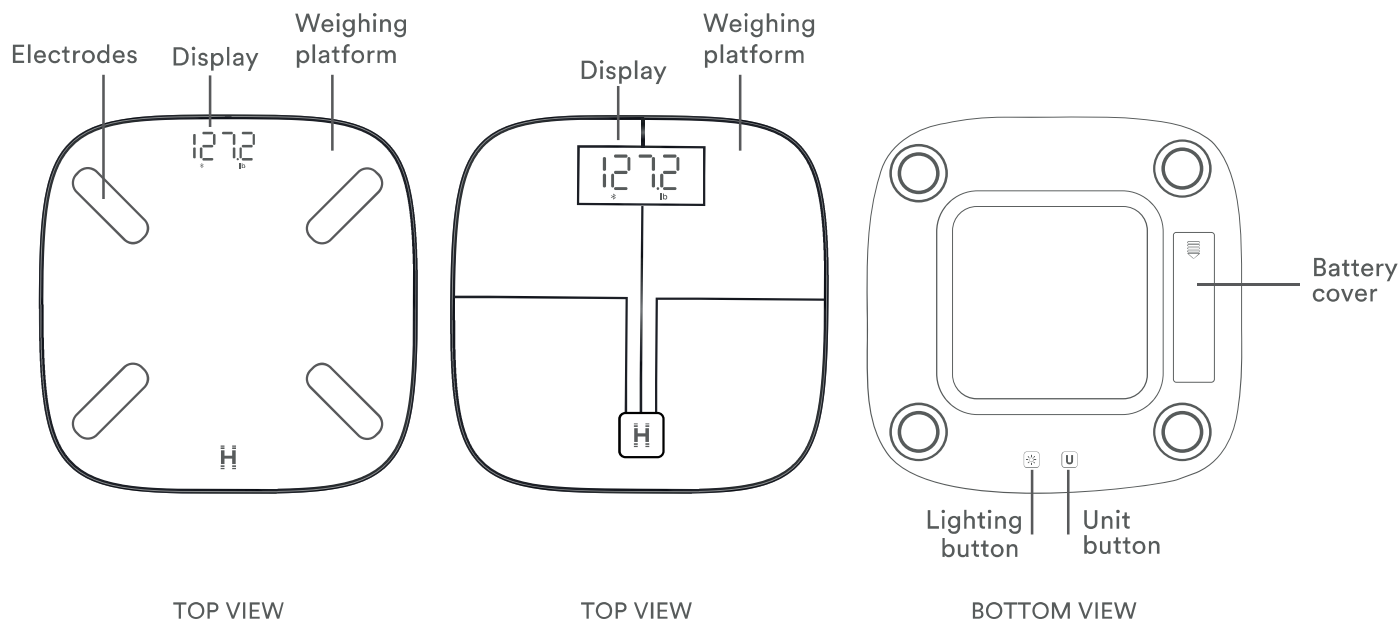FOR SCL-XBC200

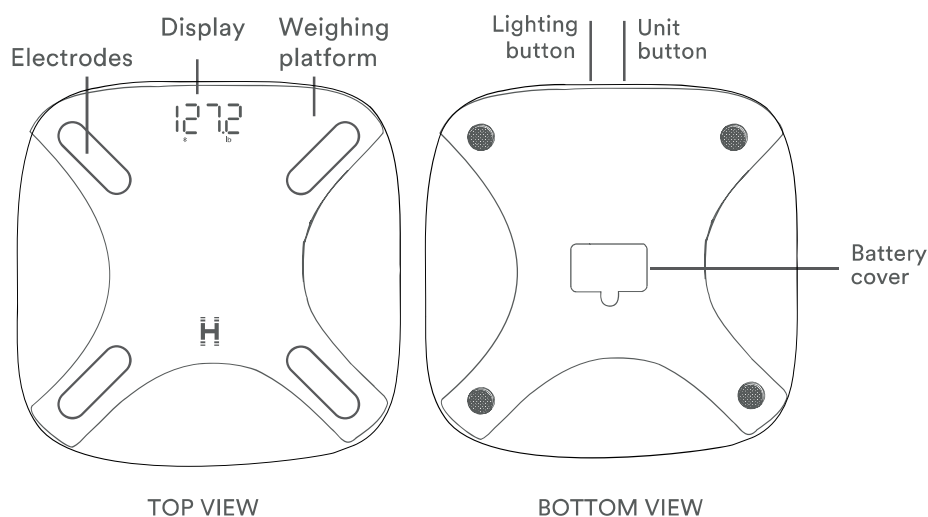BEFORE TURNING THE UNIT ON

This scale is shipped with 4 AAA batteries for the SCL-BBC200, SCL-BBC250 and SCL-BBC350 or 3 AAA batteries for the SCL-XBC200. Please install the batteries inside the battery compartment to power on your scale before use. After power on, you can take your first measurement.

BACK PANEL BUTTONS

LIGHTING BUTTON

Press the lighting button to choose your lighting brightness: high, low, or off

UNIT BUTTON


Press to choose your desired units of measure: lb, kg, or st:lb

DISPLAY SCREENS


Upon measurement, weight, body mass index (BMI) and body fat percentage (%) will display on the scale’s display screen for 4 seconds each. Download the Homedics Scales App for all 13 essential biometrics.

NOTE: If you have not connected your scale to the Homedics Scales App prior to measurement and input your height, date of birth, and sex, the scale will only display weight. BMI and body fat % will only show on the scale screen after connecting to the Homedics Scales App and inputting your height, date of birth, and sex.


Weight

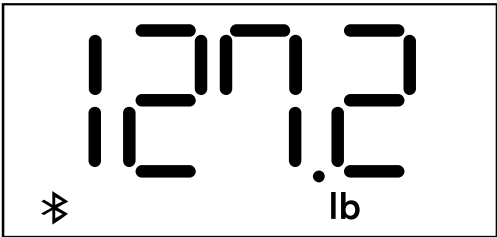Body Mass Index

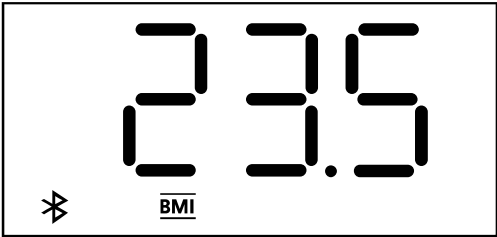Body Fat Percentage

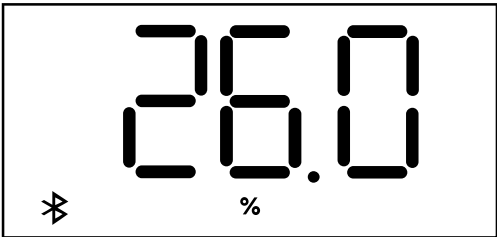
HOW TO CONNECT THE MOBILE APP

1. Download the Homedics Scales app in the App Store® or Google Play store or scan the QR code.
2. Open the Homedics Scales app. If you already have an account, log in. To create a new account, tap Sign Up.
3. To add a new device, tap “+ Pair the device”.
4. Follow the in-app instructions to add your scale.

HOW TO TAKE A MEASUREMENT

NOTE: To use athlete mode or baby mode, follow the prompts in the app BEFORE taking your first measurement.

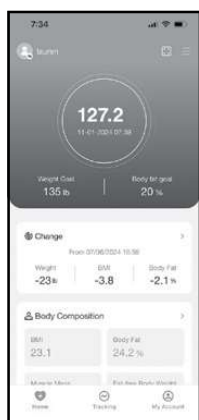
1. Open your Homedics Scales App. Then step on the scale.
2. Wait while the app syncs your measurement data.
3. Once the app shows all 13 measurements, you can step off the scale.
4. Tap your measurements to get more information.

AMBIENT ILLUMINATION

1. Install 4 AAA batteries for the SCL-BBC200, SCL-BBC250 and SCL-BBC350, or 3 AAA batteries for the SCL-XBC200 and place scale on a hard, flat surface.
2. Once you are within the proximity zone, the scale will turn on and stay on until you leave the proximity zone, or until motion is not detected.
3. If the LEDs don't turn on when you are close to the scale, click the lighting button located on the bottom of the scale to turn the night light on.
4. If you want to dim the motion sensing night light, click the lighting button once located on the bottom of the scale.
5. If you ever want to turn the motion sensing night light off, click the lighting button twice located on the bottom of the scale.

HOW TO USE

1. Place scale on a hard, flat, dry surface. Avoid carpet or other soft surfaces.
CAUTION: Placing the scale on soft surfaces such as shag carpet may cause tipping or inaccurate results.
2. Step on the scale with both feet. Use bare feet for best results.
3. Stand on the scale until the display flashes several times. This is your final measurement.

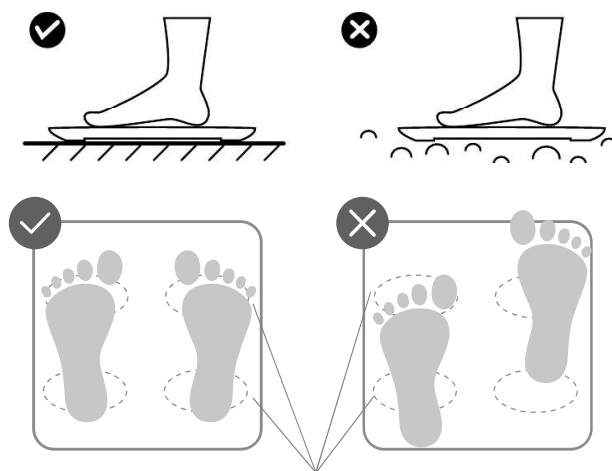MEASUREMENT TIPS

This product employs BIA (Bio-electrical Impedance Analysis) technology. To ensure accurate measurement results, follow these instructions:

- Use the scale at the same time each day.
- Measure under the same conditions as much as possible.
- Measure with clean, bare feet. Wearing socks or having soiled feet will impact the accuracy of your results.

Avoid measuring under the following conditions:


- During or immediately after strenuous exercise
- Immediately after showering, bathing, or using a sauna
- After heavy alcohol consumption
- During or immediately after eating large meals
- While experiencing a fever.

CLEANING THE DIGITAL SCALE


- Use a slightly damp cloth or disinfecting wipe to clean the surface of the scale, followed by a soft, dry cloth to wipe away any remaining moisture.
- DO NOT use abrasive cleaning agents to clean the scale.

BATTERY REPLACEMENT


FOR SCL-BBC200, SCL-BBC250 AND SCL-BBC350

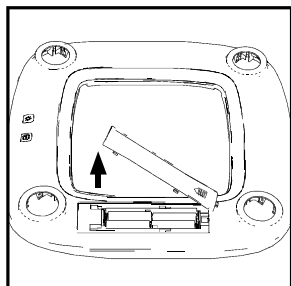
1. Remove the battery cover by sliding the cover in the direction of the open arrow and lifting the top edge.

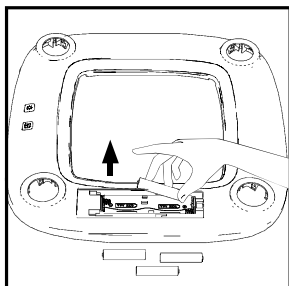2. Remove the batteries and dispose of per local regulations.

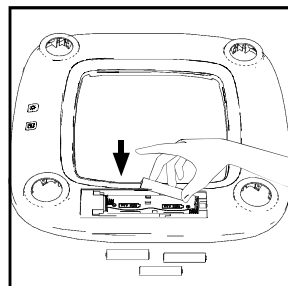
3. Install 4 new AAA batteries following the polarity as shown.

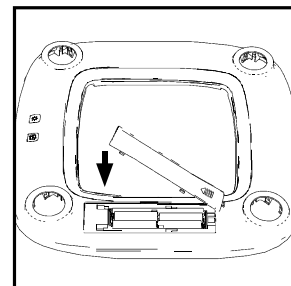
4. Replace the battery cover.


FOR SCL-XBC200

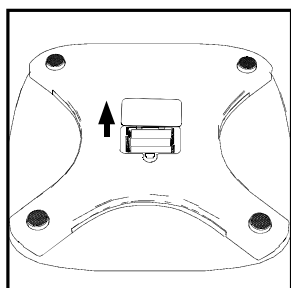
1. Remove the battery cover by pinching the tab and lifting the cover.

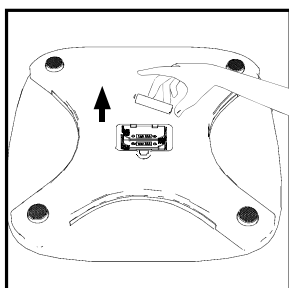2. Remove the batteries and dispose of per local regulations.

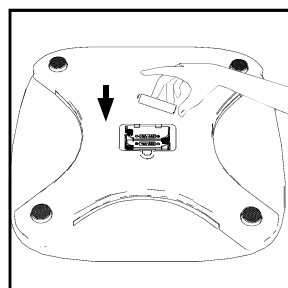3. Install 3 new AAA batteries following the polarity as shown.

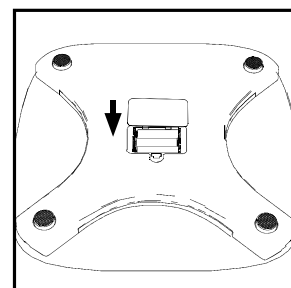4. Replace the battery cover.

BATTERY PRECAUTIONS

⚠ WARNING — KEEP BATTERIES OUT OF REACH OF CHILDREN.

1. Replace all batteries at the same time. Do not mix new and used batteries. Do not mix alkaline, standard, or rechargeable batteries. Use only the size and type of batteries specified.
2. Clean the battery contacts and the contacts of the device with a dry cloth prior to installing the batteries.
3. Ensure the batteries are installed correctly with regard to polarity (+/-).
4. Remove batteries from the device when the device is not used for an extended period.
5. Remove batteries promptly from the device when the batteries are depleted.
6. Keep batteries away from children.
7. Batteries are small objects and could be ingested. If swallowed, contact a physician at once.
8. Remove depleted batteries from the unit promptly. Dispose of used batteries safely.
9. Recycle or dispose of batteries properly in accordance with all local, state, province, and country regulations. Do not open batteries or dispose of in fire.

BIOMETRICS DESCRIPTIONS

WEIGHT

Body weight is a simple measure of total body components and is compared to a “standard” or desirable weight.

BODY FAT %

Body fat (BF) percentage refers to the ratio of fat (adipose tissue) in the body in comparison to overall body weight. While body mass index (BMI) is often used to estimate body fat, it’s not always an accurate method of doing so. According to the American Council on Exercise (ACE), a body fat range of 25% to 31% for women or 18% to 24% for men is typical for the average nonathlete. A BF below 14% for women or 6% for men may be considered dangerously low and can lead to health risks.

	Essential	Athletes	Fitness	Acceptable	High	Excess
MEN	2%-6%	6%-13%	13%-17%	17%-25%	25%-32%	≥32%
WOMEN	10%-14%	14%-21%	21%-25%	25%-32%	32%-38%	≥38%

BMI

A standardized ratio of weight to height, used as a general indicator of health.

The Body Mass Index (BMI) is a widely used health indicator. It can be roughly calculated by dividing your body weight (in kilograms) by your height (in meters) squared. If the resulting figure is less than 18.5, you are underweight. A number between 18.5 and 25 indicates a healthy weight. A number above 25 is overweight and a number above 30 is obese. Although BMI is a generally accepted health indicator, it is not the only one. For example, someone with a lot of muscle mass can have a high BMI without being unhealthy. By looking at the ratio of muscle to fat mass, among other things, you get a much better understanding of your health.

LOW	STANDARD	HIGH	EXCESS
< 18.5	18.5 - 25	25 - 30	≥ 30

SPECIFICATIONS

Model	SCL-BBC200, SCL-BBC250, SCL-BBC350, SCL-XBC200
Weight Units	lb / kg / st:lb
Weight Capacity	400lb / 180kg / 28st 7lb
Increments	lb: 0.1lb kg: 0.05kg for less than 100kg total weight, 0.1kg for 100kg or greater total weight st:lb: 0.1lb for less than 20st, 1lb for 20st or greater
Platform	Tempered Glass
Weight Display Auto-Off	10 Seconds
Dimensions	SCL-BBC200: 11.8×11.8×1.1in / 30×30×2.9cm SCL-BBC250: 11.8×11.8×1.1in / 30×30×2.9cm SCL-BBC350: 11.8×11.8×1.1in / 30×30×2.9cm SCL-XBC200: 11.8×11.8×1.1in / 30×30×2.8cm

TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
Scale won't turn on or no display	- Batteries are dead - Batteries aren't fully installed - Scale isn't on a hard, flat surface	- Use new batteries - Remove battery cover to check if all batteries are installed in the correct direction - Place your scale on a hard, flat surface
Lights won't turn on	- Batteries are dead - Batteries aren't fully installed - Lights are turned off	- Replace batteries with new batteries - Remove battery cover to check if all batteries are installed in the correct direction - Click the light button located on the bottom of the scale
The scale readings are unstable or inaccurate	Scale isn't on a hard, flat surface	Place your scale on a hard, flat surface
Display shows: 	Batteries are dead	Replace batteries with new batteries.
Display shows: 	Weight applied to the scale is greater than 400 lbs / 181 kg / 28 st:lb	Step off the scale. DO NOT apply weight greater than 400lbs / 180kg / 28st 7lb to the scale.

SUPPLIER'S DECLARATION OF CONFORMITY

Product Description: BODY COMPOSITION GLO SMART SCALE WITH AMBIENT ILLUMINATION

Model Number: SCL-BBC200, SCL-BBC250, SCL-BBC350, SCL-XBC200

Trade Name: Homedics

FCC Compliance Statement

Caution: The user is cautioned that changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s) and

Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause interference.

(2) This device must accept any interference, including interference that may cause undesired operation of the device.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two

U.S. Contact Information

Company: Homedics, LLC.

Address: 3000 N Pontiac Trail, Commerce Township, MI 48390

8:30am-7:00pm EST Monday-Friday 1-800-466-3342

Homedics is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This product has been tested and complies with the requirements for Federal Communication Commission, Part 15. Although this product is tested and complies with FCC, it may interfere with other devices. If this product is found to interfere with another device, separate the other device and this product. Conduct only the user maintenance found in this instruction manual. Other maintenance and servicing may cause harmful interference and can void the required FCC compliance.

CAN ICES-003(B)/NMB-003(B)

RF Exposure Statement

This equipment complies with FCC and ISED RF radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operated in conjunction with any other antenna or transmitter.

5-YEAR LIMITED WARRANTY

Homedics sells its products with the intent that they are free of defects in manufacture and workmanship for a period of 5 years from the date of original purchase, except as noted below. Homedics warrants that its products will be free of defects in material and workmanship under normal use and service. This warranty extends only to consumers and does not extend to Retailers.

To obtain warranty service on your Homedics product, contact a Consumer Relations representative for assistance. Please make sure to have the model number of the product available.

Homedics does not authorize anyone, including, but not limited to, Retailers, the subsequent consumer purchaser of the product from a Retailer or remote purchasers, to obligate Homedics in any way beyond the terms set forth herein. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; improper installation; unauthorized repairs or modifications; improper use of electrical/power supply; loss of power; dropped product; malfunction or damage of an operating part from failure to provide manufacturer's recommended maintenance; transportation damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Homedics.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adoption to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. THERE SHALL BE NO OTHER WARRANTIES EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS OR ANY OTHER OBLIGATION ON THE PART OF THE COMPANY WITH RESPECT TO PRODUCTS COVERED BY THIS WARRANTY. HOMEDICS SHALL HAVE NO LIABILITY FOR ANY INCIDENTAL, CONSEQUENTIAL, OR SPECIAL DAMAGES. IN NO EVENT SHALL THIS WARRANTY REQUIRE MORE THAN THE REPAIR OR REPLACEMENT OF ANY PART OR PARTS WHICH ARE FOUND TO BE DEFECTIVE WITHIN THE EFFECTIVE PERIOD OF THE WARRANTY. NO REFUNDS WILL BE GIVEN. IF REPLACEMENT PARTS FOR DEFECTIVE MATERIALS ARE NOT AVAILABLE, HOMEDICS RESERVES THE RIGHT TO MAKE PRODUCT SUBSTITUTIONS IN LIEU OF REPAIR OR REPLACEMENT.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on internet auction sites and/or sales of such products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate as to any products or parts thereof which are repaired, replaced, altered, or modified, without the prior express and written consent of Homedics.

This warranty provides you with specific legal rights. You may have additional rights which may vary from state to state and country to country. Because of individual state and country regulations, some of the above limitations and exclusions may not apply to you.

For more information regarding our product line in the USA, please visit www.homedics.com. For Canada, please visit www.homedics.ca.

FOR SERVICE IN THE USA

cservice@homedics.com

8:30am-7:00pm EST

Monday-Friday

1-800-466-3342

FOR SERVICE IN CANADA

cservice@homedics.ca

8:30am-5:00pm EST

Monday-Friday

1-888-225-7378

©2024 Homedics, LLC. All rights reserved. Homedics is a registered trademark of Homedics, LLC. The Bluetooth word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Homedics is under license. The App Store is a service mark of Apple, Inc. Android, Google Play, and the Google Play logo are trademarks of Google LLC.

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