

1. This massage product has a wireless remote control with 15, 20, 30 minutes timer. There are three modes, 9 intensity, 9 speeds and 9 sensitive adjustment. There are Zone-1, Zone-2 massage positions for selection.
2. This product has 4 parts: Wireless remote control, Middle control, Vibration units and Charger.
3. Before using, please connect the middle control with the 120VAC/60Hz power source, the vibration units and the charger. Make sure the product is be on the set on position.
4. Detail functions for Picture 4 is as below:

Name of the press button	functions
POWER/TIMER	<p>On/off set up</p> <p>After the power is on</p> <p>Press this button one time, the relative 15 minutes timer light will turn on. This means the product has been set “start on” and the it will turn off automatically after 15 minutes.</p> <p>Press the button the second time, The product is till on the “start on” position and the 20 minutes timer light(yellow) will turn no. This means it will turn off automatically after 20 minutes.</p> <p>Press the button the third time, the TIMER 30(Red light) will turn on, this means the product will turn off automatically after 30 minutes.</p> <p>Press the button the fourth time, The product will change from the “start on” position to the “power off” position. All lights will turn off. Press the button once again, it will repeat the former 4 situations.</p> <p>(In addition: 1.When the product change from the original Power on situation to the start on situation, the WAVE is default mode(Default intensity is middle, default speed is middle, default sensitive is middle).Below buttons are only valid when the product is set to the “start on” position.</p>
ZONE-1	To turn on/off for Massage position 1.


		Press the button one time and the second, the position will change from the on and off situation.
ZONE-2		To turn on/off for Massage position 1. Press the button one time and the second, the position will change from the on and off situation
WAVE/PULSE		To change the wave and pulse massage modes. Press the button one time and the second, massage position 1 and 2 will change between WAVE and PULSE.(When the modes changes, It will reserve the same intensity, speed and sensitive of each mode)
STEADY		To enter into continuous vibration massage position. When the product is set on, press the button to enter into continuous vibration massage position. (When in this steady situation, only intensity will change, the speed and the sensitive will not change.
INTENSITY	+	Increase the intensity. (Valid for WAVE、PULSE and STEADY) Press this button, Both zone-1 and zone-2 will increase one level. Maximum to be 9 levels. When it reaches the maximum intensity, there will not change even keep pressing this button.
	-	Decrease the intensity. (Valid WAVE、PULSE and STEADY) Press this button, Both zone-1 and zone-2 will decrease one level. Maximum to be 9 levels. When it reaches the minimum intensity, there will not change even keep pressing this button.
SPEED	+	Increase the speed. (Valid for WAVE and PULSE) Press this button one time, both zone-1 and zone-2 will increase the speed of their massage, wave or pulse by one level, maximum to 9 levels. When it reaches the maximum intensity, there will not change even keep pressing this button.
	-	Decrease the speed. (Valid for WAVE and PULSE) Press this button one time, both zone-1 and zone-2 will decrease the speed of their massage, wave or pulse by one level, maximum to 9 levels. When it reaches the minimum intensity, there will not change

		even keep pressing this button.
SENSATE	+	Increase the sensitive level (Valid for WAVE、PULSE) Press this button one time, massage zone-1 and zone-2 will decrease the trend of their massage, wave and pulse, that means the user must have higher sensitivity to feel the change of the wave or pulse's change. Maximum 9 levels. When reaches the highest level, there will be no changes even keep pressing this button.
	-	Decrease the sensitive level. (Valid for WAVE、PULSE) Press this button one time, massage zone-1 and zone-2 will increase the trend of their massage, wave and pulse, that means the user even have a lower sensitivity, they can still feel the change of the wave or pulse. Maximum 9 levels. When reaches the lowest level, there will be no changes even keep pressing this button.


5. When the timer is due, the product will quit all the working, turn to the “power on” situation. If you want to use again, you need to restart the product and press the relative button to fulfill the functions.
6. If the remote control get mal-function , Press the reset key on the back of the remote control to restart it (picture 2)
7. How to use the charger: (picture 3)
 - A. Put the remote control into the charger, make sure they are connect well.
 - B. When the red light is on, the remote control is being charged.
 - C. When the green light is on, it means the charge is done or the remote control is not connected well.

FCC STATEMENT


1. This device complies with Part 15 of the FCC Rules.
 Operation is subject to the following two conditions:
 - (1) This device may not cause harmful interference, and
 - (2) This device must accept any interference received, including interference that may cause undesired operation.
2. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

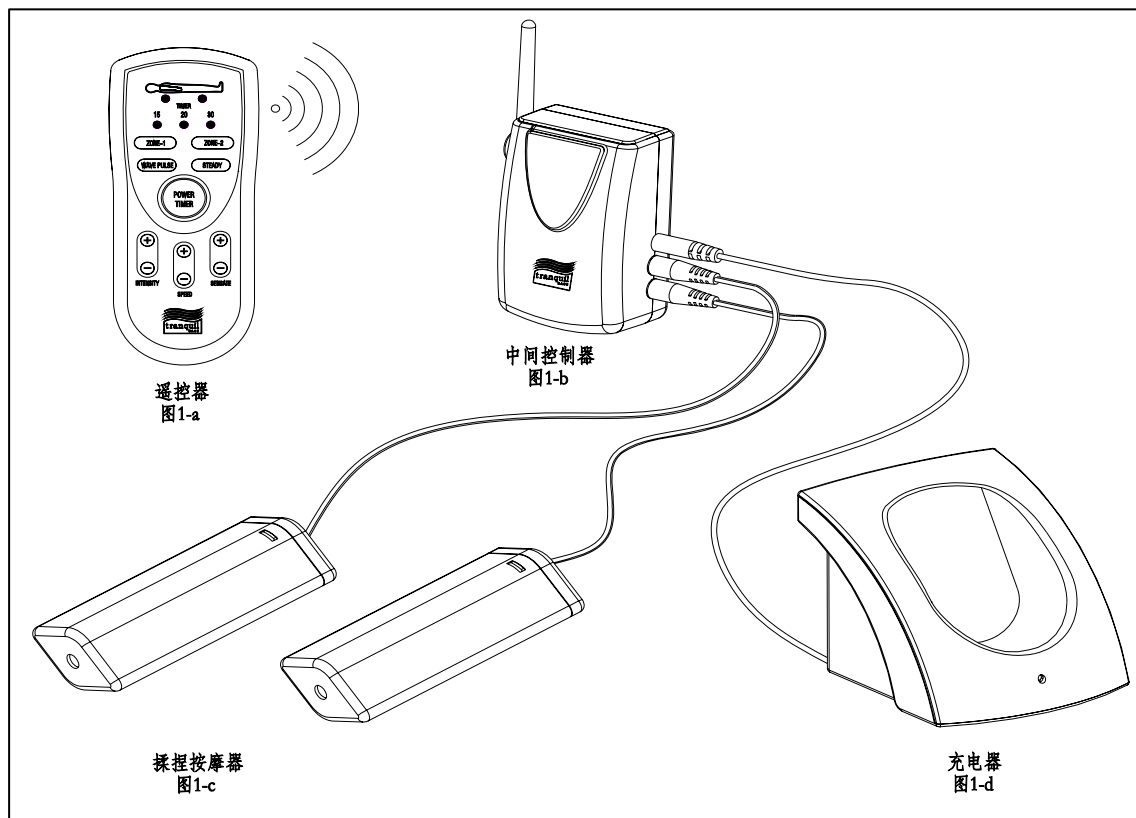


Picture 1

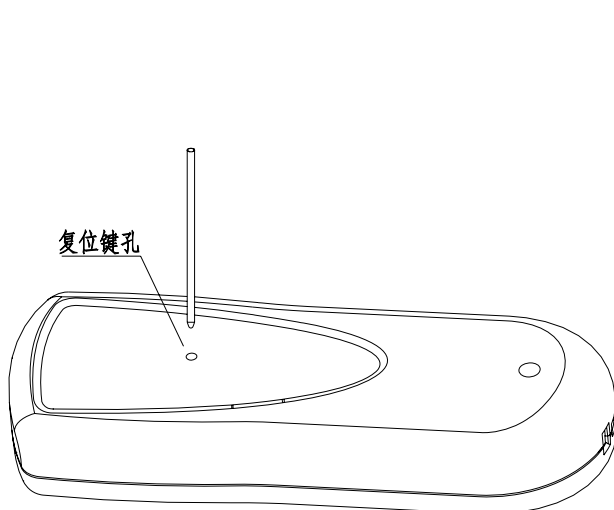

Picture 2

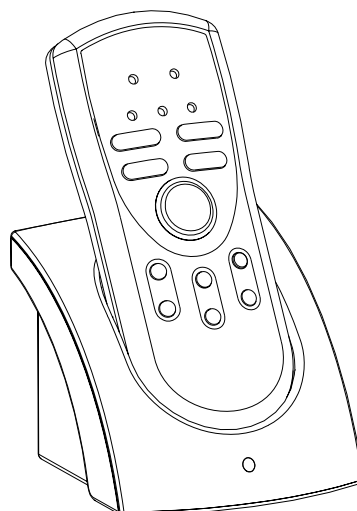

Picture 3

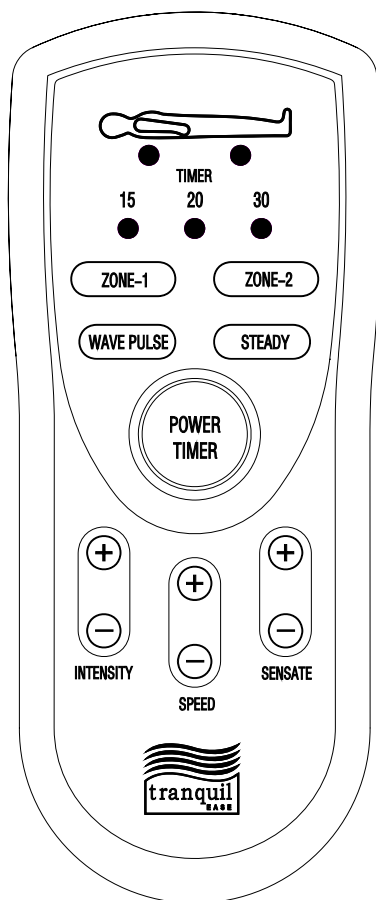

Picture 4