

LEDs Status	Red LED	Green LED	Yellow LED	Alarm
Event	Bat.	Alert	Pair	
Unit power on		3 sec Solid on		
Unit connected to Integrated Alarm			Blink once / sec	
Ready to Pair		Blink 8 times / sec		
Pairing successful		Blink 3 times / sec		
Pairing failed or time-out	Blink 3 times / sec			
Alarm (Unit on)		Solid on		Alarm / Melody1 / Melody 2 / Mute
Unit Activated				Beep 3 times
Press Reset				Beep 3 times
Incorrect Setup				Long Beep, turn off device and restart
Low battery warning	Blink once / sec			
Very low battery warning*	Solid on			

* No function



WARNING / CAUTION

Before using this device, read entire instruction sheet and save for future reference

- The device will not prevent falls.
- This device is not a substitute for visual monitoring of patient. Make routine inspections for proper use of the device. Test alarm daily to verify the alarm is operating properly. Failure to do so may result in injury or death to a person in your care.
- Use of alarm unit may not be suitable for "High Risk" patients. Other safety measures may still be necessary to ensure patient safety.
- The Silent (Mute) mode should only be used in conjunction with a nurse-call system. Do not use the device in silent mode for stand-alone operation, as no alarm tone will sound. When using the device in silent mode, be sure to reset the alarm each time it is triggered. Resetting only the nurse-call system will dramatically shorten the alarm's battery life.
- This product is warranted for 30 days from date of purchase.

The Gleen Medical Incorporated's Alarm is designed to be installed by the end user. As such, it is the entire responsibility of the buyer to insure that the system is properly installed and tested. Further, the alarm is not designed to replace visual monitoring of patient. Where permissible, it is the seller's intent to limit any liability for special, consequential or incidental economic damages to refund of the purchase price or replacement of then goods. Seller disclaims any liability for compensatory or other damages arising out of the use of the goods.

FCC Compliance information for IA/EA/PIR/TX ID:


This device complies with Part 15 of the FCC rules . Operation is subject to the following two conditions: (1) This device may not cause harmful interference; and (2) this device must accept any interference received, including interference that may cause undesired operation. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:


- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Contact the dealer or an experienced radio/TV technician for help.

GleenMed Line of Products

*Each Sold Separately

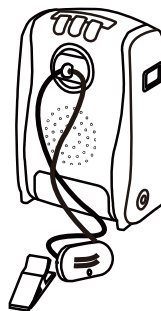


GLN0161
POM Wireless Emergency Alarm

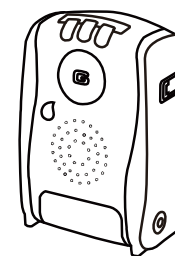


GLN0111
POM Wireless Integrated Alarm
GLN0110
POM Wired Alarm



GLN0171
POM Wireless Integrated PIR Light

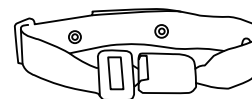

GLA5030
Seat Belt

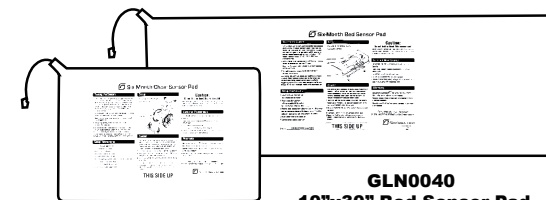

GLN0030
10"x15" Chair Sensor Pad

GLN0040
10"x30" Bed Sensor Pad

Follow us on



<http://facebook.com/gleenmedical>



<http://youtube.com/NETBcWhrpQo>



GLEENMED



For question or comments, please email:


INFO@GLEENMED.COM

Or Visit: www.gleenmed.com



Getting Start 1


Install 3 x AA Batteries into the Emergency Alarm (not included)
(An optional AC Adapter is available)

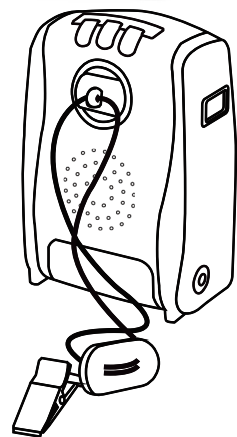

POM Wireless Emergency alarm

Quickstart Guide

GLN0161

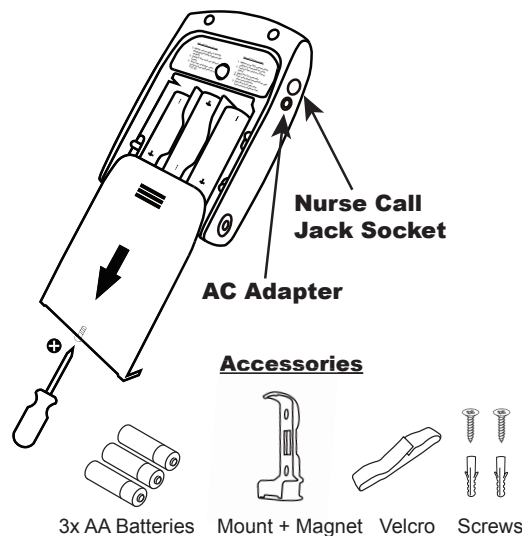


Setting 2

ON / OFF: Slide the delay switch to 0 second or 2 seconds

Volume: 85 dB (low) or 100 dB (High)


Sound: Melody 1, 2 or 3 or mute*

*Only use mute with nurse call jack

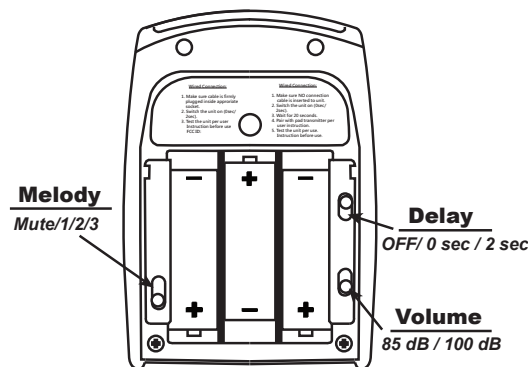

Test the system 3

Remove Magnet from Alarm. It is ready for use if Alarm Sound

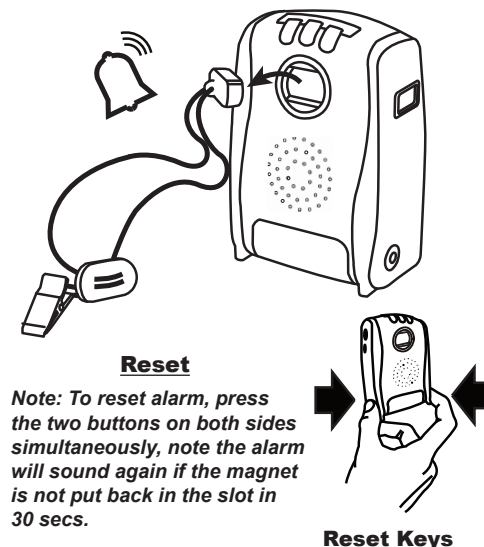


Mounting 4

- Magnet
- Screws
- Velcro




Bathroom / Toilet

If paired with GLN2001 Wireless Integrated Alarm (See section 5), Wireless Integrated Alarm can be mounted outside bedroom or bathroom to alert caregiver.

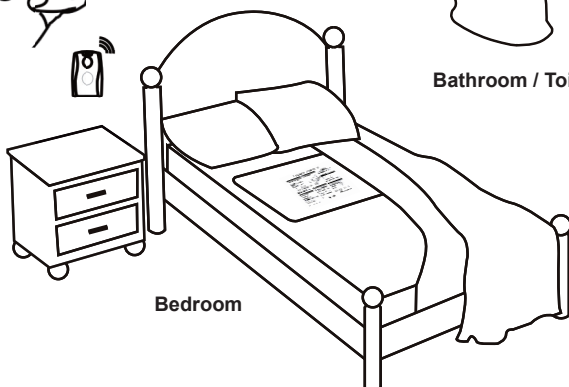


Bedroom

Wheel Chair
(Attach Alligator clip to patient's clothing around the shoulder / neck area)

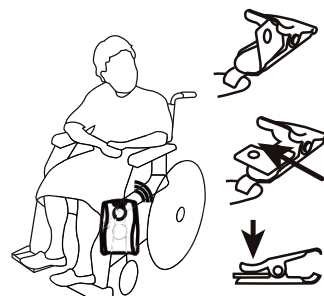

Pairing Option 5

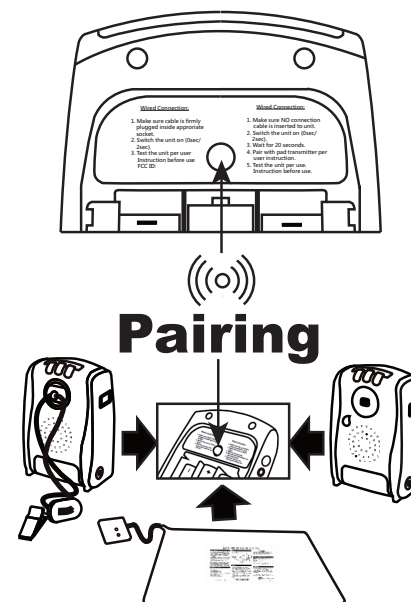

1. Turn the GLN0161 POM Wireless Emergency Alarm and GLN0111 POM Wireless Integrated Alarm ON.
2. When the green light flashes rapidly on both devices, press the Blue pairing button in the middle of the POM Wireless Integrated Alarm . (It is looking for other devices to pair with).
3. Press the pairing button on GLN0161 POM Wireless Emergency Alarm until the central light rapidly flashes green as well.
4. Make sure Integrated Alarm and Emergency Alarm are rapidly flashing green at the same time before.
5. When pairing is successful, the light on both devices should slowly flash green 3 times.
6. Test the System by removing magnet from Alarm, both alarms will sound now.
7. The green light of the POM Wireless Integrated Alarm will flash rapidly again to search for the next device to pair. If no pairing device is found or pairing is not successful after 15 sec, the red light will flash slowly 3 times. (max 4 devices can be paired with 1 POM Wireless Integrated Alarm)
8. If pairing is unsuccessful, red light will flash slowly 3 times. Turn all devices off and repeat from step 1 to step 7.