

Grid Pad 13 Manual

Document: Grid Pad 13 - Manual - UK - A
Device: GP13A
Language: English
Date of issue: <month> 2024

Copyright © 2024 Smartbox Assistive Technology


This guide provides the basics for getting started with your Grid Pad 13.

You can find a complete guide to all the features of your Grid Pad on the Smartbox Hub:

<Hub.thinksmartbox.com/GridPad13-2nd-gen>

Your Grid Pad and Grid

Your Grid Pad has many features that can be customised and controlled in Grid.

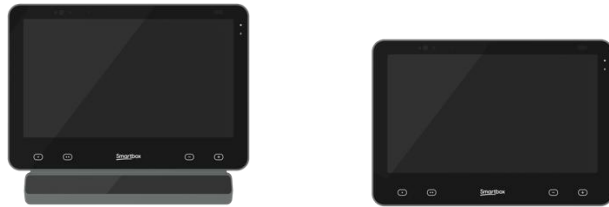

For more resources and guides on how to get started with Grid visit the Smartbox Hub.

<qr code>

- Get started with Grid
- Access methods and setup
- Troubleshooting guides
- Videos and webinars
- Resources for grid sets
- ...and much more

Power

To power on your Grid Pad press the power button on the side of the device for just under a second.



A white light on the front of the device will come on and the system will start up.

If the device is in sleep mode it will wake up.

Shutting down your Grid Pad

Close any apps you have open and make sure you have saved any work.

Tap the Start menu
Select the Power icon
Select Shut down

The Grid Pad will shut down and any pending Windows updates will be installed.

You can also Restart your device from here, which will apply any updates.

Force shut down

If you press and hold the power button for five seconds, your Grid Pad will shut down. Please note that this is only recommended in emergencies, and that any unsaved work will be lost.

Remote Power button

Your Grid Pad has a Remote Power Button included.

The Remote Power Button allows you to turn your Grid Pad on and off.

It can also be used to setup switches away from the device.

The switch port S3 also acts as a power button. This can be toggled in your Grid settings.

The Remote Power Button comes paired with your Grid Pad out of the box.

Charging

To charge your Grid Pad, plug the power lead into the charging port on the left side of the device.



When the Grid Pad is charging the charging indicator is illuminated white.

When the battery is full the charging indicator is illuminated green.

When your Grid Pad requires charging, the charging indicator will be illuminated orange.

Mounting

Adjustable stand

Your Grid Pad has an adjustable stand that can be folded out, to position your device on a desk or table comfortably. Pull the stand from the bottom and adjust the angle to suit you. When finished, simply push the stand back in.

Using a desk, floor, or wheelchair mount


Your Grid Pad is supplied with a mounting plate for both Rehadapt and Daessy, so that you can securely mount and position your device.

The mounting plate is attached with three screws (supplied with the mounting plate). Ensure all three screws are secure before mounting your Grid Pad to a stand or mounting system.



Smart Buttons

The front of your Grid Pad has five Smart Buttons that are accessible with eye gaze and touch.

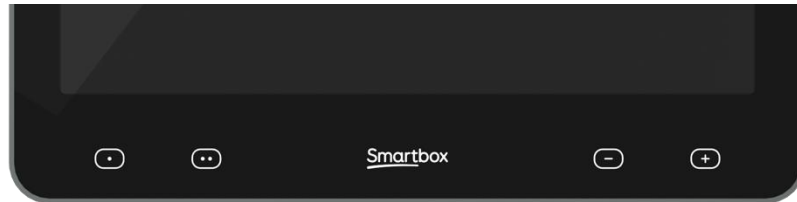

Buttons 1, 2 and the Smartbox button can help you control and navigate your Grid Pad device when using Grid, or they can be customised using Grid's set of powerful commands.

You can customise the colour of the smart button highlight as well as add audio feedback.

The two Smart volume buttons cannot be customised and offer eye gaze access when using Grid.

Connecting switches

Your Grid Pad has two switch ports which can be configured in Grid.



Switch ports S1 and S2 have an LED indicator above that illuminate when a switch is activated.

The Remote Power Button also has two additional switch ports, S3 and S4.

Connecting Eye gaze

If your Grid Pad is supplied with an eye gaze camera it will be connected and ready to use.

You can remove or swap the eye gaze camera.

- Place your Grid Pad face down on a flat surface, ideally a table so the camera is easy to access.
- Raise the stand.
- Loosen the two screws on the back of the camera unit.
- Pull the camera unit away from the Grid Pad.

Changing the stand

Note the additional stand is included when the EG camera is supplied separately.

- Place your Grid Pad on a flat surface with Raise the stand so it is vertical.
- Remove the screw from each side of the stand.
- Pull both sides of the stand outwards at the same time to unhook the stand.

Second screen

The second screen on the back of the Grid Pad can show your Grid messages and Symoji's to the people around you.



The second screen can be toggled on and off in the Grid settings, and using the Set Second Screen command in your grid sets.

Finding your model and serial number

When you contact our Support team, please have your device model and serial number to hand.

This helps us identify the device you are having an issue with. You can find the serial number on the back of your Grid Pad.

The serial number for the device is placed under the stand, lift the stand and the serial number is found here:

