

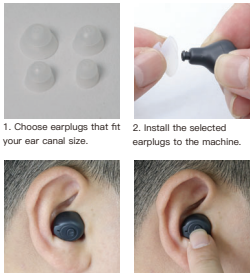








<p>OPERATION GUIDE</p>  <p>01</p>	<p>Before wearing</p>  <p>Please use a cotton swab to clean the ear canal before wearing. If the cerumen is too much, it will block the sound outlet of the product and cause no sound.</p> <p>02</p>	<p>How to wear</p>  <p>03</p>	<p>Note: Be sure to choose earplugs that fit your ear canal and wear them correctly to fit your ears as much as possible. If you don't wear it properly, there will be echoes (whistles/sounds).</p> <p>Problems encountered in wearing</p> <p>1. There is no sound when wearing?</p> <p>① Check whether the volume level is adjusted to suitable for your hearing</p> <p>② Check whether the hearing aid has power</p> <p>③ Whether there is earwax in the ear canal, because the earwax is easy blocked the sound of the hearing aid</p> <p>04</p>
<p>Note: If the hearing aid is silent or the sound becomes low due to the blockage of the sound outlet, please clean the sound outlet.</p> <p>2. Why is there a howling (sharp sound)?</p> <p>Howling (sharp sound) is a manifestation of the normal working of hearing aids. The hearing aid and the ear canal are not completely matched, causing the sound from the receiver to be received by the microphone again, and the sound forms a loop, which will be continuously amplified and produce howling.</p> <p>Solution: The product is equipped with earplugs of different sizes, please choose the earplugs that fit your ear canal</p> <p>05</p>	<p>Turn ON</p> <p>After the hearing aid leaves the charging case for 2 seconds, it will automatically turn on or long press the M button 2 seconds to turn it on, and you will hear a "Power on" sound.</p>  <p>Take out</p> <p>* Please avoid covering the hole when operating the device</p> <p>* The first charge of charging case should be guaranteed to reach 10-12 hours</p> <p>06</p>	<p>Turn OFF</p> <p>After use, directly put the hearing aid back into the charging compartment, it will automatically turn off, or long press the M button for 6 seconds to turn off</p>  <p>Put back</p>  <p>Long press the M key</p> <p>07</p>	<p>Bluetooth mode</p> <p>Bluetooth function:</p> <p>For the first use (never connected with Bluetooth before), pressing the M-key of one of the hearing aids three times in succession after powering on, the beep of "pairing" can be heard.</p> <p>Then turning on the Bluetooth function on the mobile phone, searching for "GM-903B1", and clicking it to connect the Bluetooth</p>  <p>*Tips:</p> <p>If connected with Bluetooth function on the mobile phone before, the hearing aids can be automatically connected with it after powering on.</p> <p>08</p>

<p>Volume control</p>  <p>As shown in the picture: lightly press the M button to adjust the volume, there are 5 levels of volume, and the volume level is 1---5 cyclic adjustment;</p> <p>Touch once, increase the volume by one level, when you hear the "beep" it means the volume has reached the maximum level.</p> <p>* When adjusting the volume, please do not cover the sound pickup hole above</p> <p>09</p>	<p>Charging</p> <p>Charging compartment</p> <p>Display number flashing=charging</p> <p>Display number show 100= charging is completed</p> <p>Hearing aid</p> <p>Green icon flashing= charging</p> <p>Green icon turn off= charging is completed</p> <p>It can work 6 hours if connected with Bluetooth, and work 10 hours if not.</p> <p>The charging compartment can be used to charge one hearing aid 4 times</p>  <p>Charging</p> <p>Charging completed</p> <p>10</p>	<p>Charging</p> <p>1. Charging the charging case: it can be charged through the power adapter or through the USB interface</p> <p>2. Hearing aid charging: Put the hearing aids back to the charging case, it will be automatic charged</p> <p>3. If the hearing aid and charging case are not used for a long time, the battery self-discharge will cause the battery to be seriously insufficient. It should be charged, and the charging time should reach 12 hours. Even if the charging indicator shows that the charging is complete, the charging should not be stopped (guarantee to reach 12 hours)</p> <p>11</p>	<p>4. During use, the hearing aid volume becomes low or when the hearing aid is not in use, the power should be turned off and charged (this will better protect the hearing aid battery). This is very important! Otherwise it will cause excessive discharge</p> <p>5. The hazards of excessive battery discharge will reduce the battery capacity and affect the battery life to varying degrees</p> <p>6. When using AC power to charge, if the charging is over, unplug the power in time</p> <p>12</p>
<p>21-day wear chart</p> <p>THE JOURNEY TO LISTENING</p> <p>Unlike vision correction, which simply wears glasses, hearing correction is a rather detailed process. Once you start your hearing correction journey, you will experience countless "new sounds", and your brain needs to adapt for a period of time. In just a few weeks, you will be shocked by the sound that you have never experienced before.</p> <p>If you have any questions, we will serve you wholeheartedly.</p> <p>13</p>	<p>THE FIRST WEEK</p> <p>Some sounds may be too noisy when wearing a hearing aid for the first time, such as the buzzing of the refrigerator, the wrinkling of paper, or the sound of your footsteps. This is normal, and your brain is relearning how to prioritize these sounds that you don't hear very often. Sometimes you may not like the noise, but please don't give up wearing.</p> <p>Patient</p> <p>At the beginning, wear the hearing aid for 30 minutes a day, and then gradually increase the wearing time until the whole day.</p> <p>14</p>	<p>THE SECOND WEEK</p> <p>The more you wear the hearing aid, the more you can adapt to the sound environment. You even can ask people around you not to speak loudly to you anymore, because you can easily communicate with people.</p> <p>Exercise</p> <p>Listen carefully to a song you are familiar with or watch a TV game show, and guess the answer without looking at the subtitles.</p> <p>15</p>	<p>THE THIRD WEEK</p> <p>After a few weeks, you will feel more energetic because you no longer need to spend a lot of energy to listen what others say. The sound will become clearer and clearer, and you will once again enjoy the sounds that surround you in your life, such as the melodious singing of birds in the backyard, the cheerful laughter of children, and the rustling of leaves in the wind.</p> <p>Participate</p> <p>Have a meal in a restaurant with friends, watch a movie in a theater, join a family gathering, and no longer have to bother to speculate about missed information.</p> <p>16</p>

FCC Warning:

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with minimum distance 0cm between the radiator and your body.