





USER'S MANUAL

Smart Steam Eye Massager

ECP-03

PRODUCT DESCRIPTION

* Products and accessories shall be subject to actual conditions.

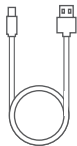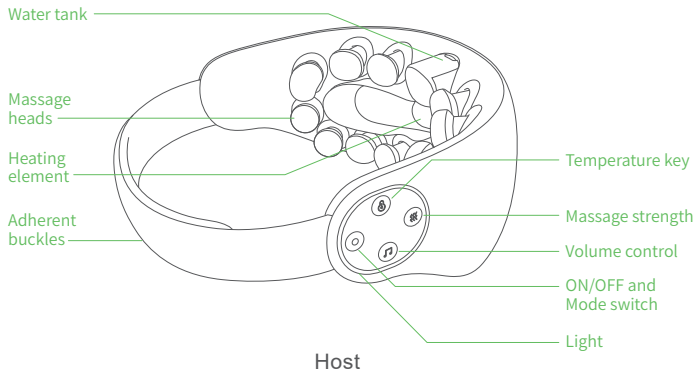

Charging line



water tank



Manual

NOTICE FOR USE

1. In order to make you use this product correctly and safely, please be sure to read the "instruction manual". We will not bear any legal responsibility and related expenses for the consequences caused by not using this product according to the relevant instructions.
2. This product is a personal care massage product and cannot be used for diagnosis, treatment and other commercial purposes.
3. If discomfort occurs around the eyes and eyes during the use of the product, please stop using it immediately and consult a doctor.
4. Please remove glasses, contact lenses and headgear before using the product.
5. Please keep your eyes closed when using the product.
6. This product is not waterproof. Please do not use it in high-temperature and humid environment such as bathroom, and keep away from fire source and corrosive environment.
7. Please use the charging cable equipped with this product. This product cannot use PD charging cable (both sides are Type-C), and the power supply used should meet the national standard.
8. As we are constantly improving our products, please understand that we will not inform you again if there are any changes in the parameters.

Matters need attention:

1. Do not touch the product with sharp objects, and do not hit or beat the product.
2. Do not use it while driving or taking a bath.
3. Do not use this product after makeup to avoid damaging your makeup.
4. Do not massage in places other than the eyes.
5. If you have used eye drops, you should wait at least 30 minutes before using this product.
6. This product has a hot compress function and a hot surface. People who are not sensitive to heat should not use it for a long time to avoid cryogenic burns.
7. Do not use it in the following situations:
 - Feeling unwell and having a fever
 - Skin damage, redness, inflammation, eczema or excessive sunburn
8. The following people should consult a doctor before using this product:
 - Myopia over 800 degrees
 - Just had laser eye surgery.
 - floaters

9. The following people should not use this product:


- Allergic constitution
- Skin diseases or atopic dermatitis
- Pregnant women

(skin discomfort and other problems may be caused by hormonal imbalance)

10. The following people are forbidden to use this product:

- Facial nerve disorder, skin sensory disturbance, skin can't perceive heat
- Eye bleeding, retinal detachment, severe eye inflammation, severe cataract, severe glaucoma

HOW TO USE

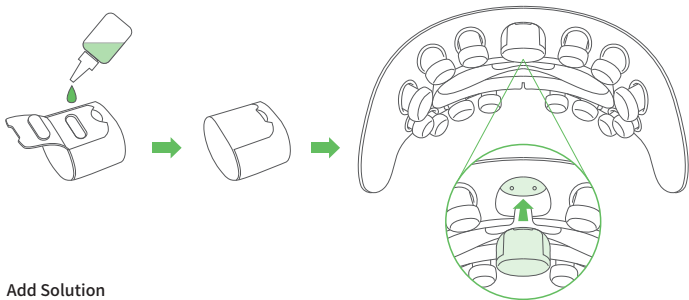

Add Solution


Open the water tank cover, add a proper amount of [Fruit Eye Care Solution] or mineral water.

Close the lid of the water tank and push the water tank into the groove of the host. The feeling of magnetic attraction means that the installation is correct.

Please do not add other solution.

ON/OFF

Press and hold the power-on key for 3 seconds to turn the host on/off.

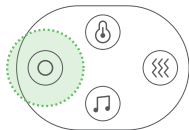


Click switch temperature

No heating	Low heating	Middle heating	High heating
------------	-------------	----------------	--------------

Click switch Mode

Vitality mode
Deeply hydrate
Eye exercises
Sleep mode

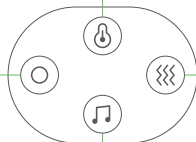

Click switch massage strength

No massage
Gentle knead
Middle knead
hard Knead

Click switch volume control

Bluetooth function

You can use Bluetooth to search for connected devices. After connecting, the phone can adjust the music volume, press the host “🎵” to switch music, click the next song, and double-click the previous song.

Scan code to download App,
more intelligent services



Please keep your eyes closed when using,

The equipment will automatically shut down after working for 10 minutes.

After use, skin should be moisturized, such as applying eye cream and eye care products.



LIGHT

White light: working

Red light: charging

Green light: full

CHARGING

This product uses a rechargeable lithium battery. Please use the equipped charging cable to charge it. This product does not support the use of PD charging cable (double-headed Type-C). The selected power adapter shall conform to the national standard.

Rated voltage: 5V 

Rated power: 10W

Charging time: about 1 hour

(if you use a non-fast charging power adapter, the charging time will be extended)

Tips

1. When using the product for the first time or for long-term reuse, the product should be fully charged before use.
2. Please don't turn it on during charging.
4. Do not place the battery box in the ambient temperature $<+0^{\circ}\text{C}$ or $>+45^{\circ}\text{C}$ for charging.
5. Do not put the product into the fire or heat it.
6. Before discarding this product, the battery must be taken out according to the requirements of laws and regulations, and safely disposed of according to local laws and regulations. Non-professionals please don't disassemble it yourself!

CLEANING

1. Clean products can be wiped with clean soft cloth. The whole product is not waterproof, so it is forbidden to put it in water or wash it with water.
2. Before cleaning, please ensure that the product is turned off and unplugged.
3. Please don't use detergent or detergent to clean.

Matters needing attention

Please change the liquid in the water tank in time and do not store it for a long time. It is forbidden to use microwave ovens, ovens, hair dryers and other drying machines and accessories.

Host clean

Clean soft cloth can be used, which can be soaked in clean water, wrung out and then wiped; Do not use alcohol, paint thinner, gasoline or nail polish remover for cleaning.

Water tank cleaning

1. The water tank can be removed and washed separately.
2. If there is scale, pour lemon juice into the water tank, soak it until the scale disappears, and then wash it repeatedly with water.
3. Please completely dry in the shade before assembling.



STORAGE

1. The cleaned products should be thoroughly dried and stored in an environment that meets the storage conditions. Storage conditions: temperature +5°C ~ +45°C; In an environment with relative humidity $\leq 80\%RH$, which is protected from light, well ventilated and free of corrosive gas, and away from heating devices and open flames.
2. The host should be prevented from falling and bumping during storage.
3. Do not wrap the power cord around the host during storage.
4. If you don't use it for a long time, it is recommended to pour out the remaining water in the water tank and store it to avoid breeding bacteria.
5. Long-term storage should charge the product regularly to avoid shortening the service life of the battery; When using again, charge the main unit first, and clean the water tank before using.



Matters need attention

1. Please store it out of the reach of children.
2. Do not store or use the product in the bathroom or shower room or other high-temperature and humid environment.
3. In case of abnormal phenomenon or failure, please stop using it immediately, disconnect the power supply and contact the after-sales service. Do not modify or repair it without permission.

Or abnormal condition:

Equipment, adapter or charging cable is deformed or abnormally hot.

Burning smell from equipment, adapter or charging cable.

The equipment, adapter or charging cable makes abnormal sound during use or charging.

4. In case of any abnormal phenomenon during use, please follow the "Troubleshooting Table" to solve it first. If it still can't work, please contact the after-sales unit.

RECYCLING INSTRUCTIONS

Please discard it according to local environmental laws and guidelines, and take care of the environment; Be sure to take the battery out of the product before it is discarded.

The lithium battery used in this product is an important resource that can be recycled. When it is discarded, it can't be disassembled, discarded or burned at will. Please give it to the garbage recycling station or a professional organization for treatment.



TROUBLESHOOTING GUIDE

Keep shut down the power off	Low battery: Reminder for charging.
	Others: Please contact after-sales.
No charging	Incorrect charging line or poor contact: Please check the connection of charging cable and power supply.
	Adapter mismatch: check whether the adapter meets the requirements of the manual.
	Others: Please contact after-sales.
No spray	Wrong assembly of water tank: please check the installation direction of water tank.
	Calcium in the water clogs the spray nozzle: 1. Pour lemon juice into the water tank and wait for 3 minutes. 2. Please power on the host, and open mode key, and pour lemon juice out if the host could spray normally. 3. Use the host after cleaning and keep cleaning of the water tank.
	There are bubbles in the water: Pour out and add eye care solution again.
	Others: Please contact after-sales.

PRODUCT PARAMETER

Product Name: Smart Steam Eye Massager Pro

Model: ECP-03

Host Size: 192×127×72mm

Weight: 280g

Rated voltage: 5V 

Rated power: 10W

Applicable Standard: GB 4706.1-2005, GB 4706.10-2008

Communications: BLE4.2 and above

CMIIT ID: 2020DP12319

Temperature usage: Under control +18°C~+40°C for use it

Storage method: Temperature+5°C ~ 45°C

Relative humidity ≤ 80%RH

Working hours: 10Minure/time

Water tank volume: 3.5mL

Charging Time: About 1 hour

SUPPORT

The final interpretation right belongs to Tianjin Fruitech Technology Co., Ltd.

If you have any questions, please contact us

(Email) hipee@fruitech.cn

FCC Statement

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Caution: Any changes or modifications to this device not explicitly approved by manufacturer could void your authority to operate this equipment.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

RF Exposure Information

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

CERTIFICATION

Model: ECP-03

Manufacturer: Tianjin Fruitech Technology Co.,Ltd.

Executive Standard: GB 4706.1-2005, GB 4706.10-2008

FCC ID:

Production Date: Printed on the package

Shelf Life: 5 years

Examiner:



©2022 Tianjin Fruitech Technology Co.,Ltd. All rights reserved.
Tianjin Huayuan Industrial Park,the 6 th Haitai Development Street,
Haitai Green Industrial Base G-1101-15 to G-1101-32.
Tianjin, 300384
Made in China
www.fruitech.cn
(E-mail) hipee@fruitech.cn