


Intelligent speech recognition
beauty eye protection instrument

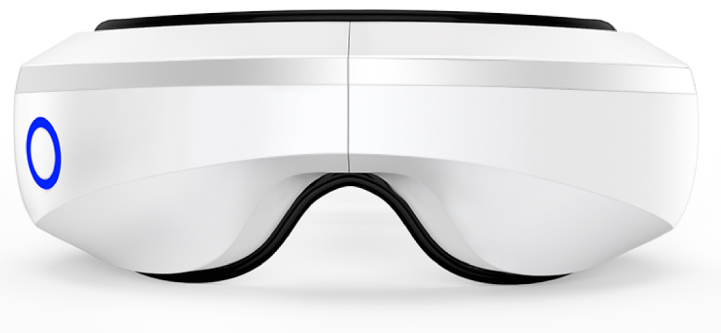

User's manual

Thank you for choosing our product. Please read this manual carefully before use, Especially the precautions to ensure proper use.
Please keep this manual properly for future use.

Catalogue

Product introduction and applicable population ----- 3

Product features and structure----- 4

Preparation and charging before use ----- 5

Method of use -----6-7

Cleaning and storage----- 8

Technical Parameters ----- 9

Thank you for choosing our product. Please read this manual carefully before use, Especially the precautions to ensure proper use.Please keep this manual properly for future use.

Product introduction and applicable population

Introduction

The intelligent wireless eye massager uses low-frequency vibration motor, environmentally-friendly PU airbag and graphemes carbon fiber far-infrared heating wire, and micro-computer chip control technology to help the eye muscles relax and stay away from eye discomfort by performing. The constant temperature hot compress can effectively improve the temperature around the eye. With the micro-motor vibration massage and the airbag simulation like people to acupressure. It can drive eye muscles together to sport give people a full range of comfort.

Scope of application

This product is suitable for people who use computers, mobile phones and other electronic devices for a long time. It can effectively avoid long-term use of the eyes, wearing glasses, insomnia, myopia, hyperopic, dark circles, under-eye bags and other eye problems caused by eye fatigue, blurred vision.


Warning

This product is a personal product of health massage, not a substitute for medicine or medical equipment.

Product features and structure

- Automatic timing
- Far infrared Hot compress
- Intelligent voice
- Portable rechargeable
- Circulating pressure
- 180 degree folding
- Meditation music
- Imitation human massage


Product structure

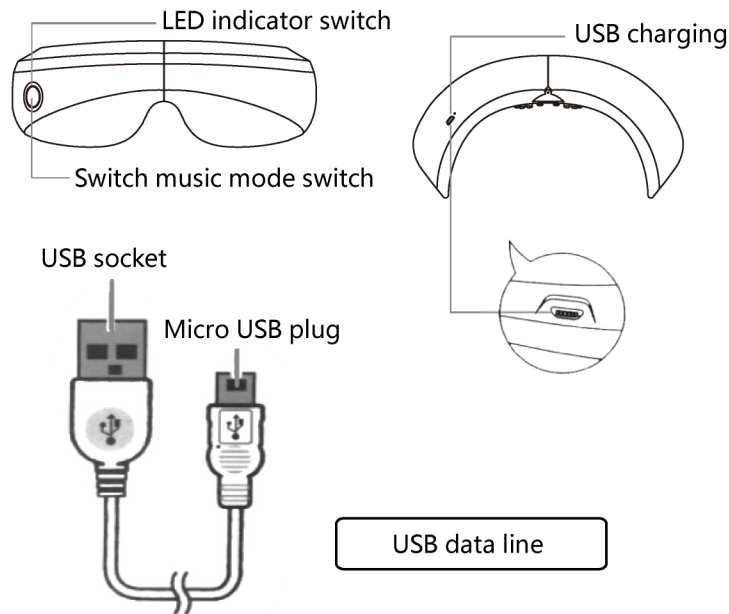

Preparation and charging before use

Battery charging

- 1 Plug the USB plugged into the USB socket of the power adapter.
- 2 Power adapter plugged into a power outlet.
- 3 Insert the Micro USB plug of the USB Data line into the Micro USB socket of the product.

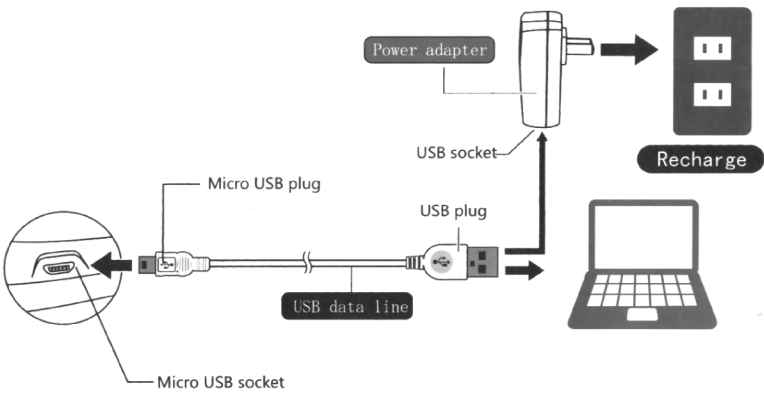

Warning

During the charging process, the indicator light flickers: after the charge is completed,the indicator light is always on. In the process of charging, the product cannot start work.
When the working product is charged, the product will automatically stop.
The charging time is 3 to 5 hours. After the charge is completed , please disconnect the connection of the power adapter.

Method of use(AR206)

Align the eye protector with your eyes, bring the product, adjust the elastic band to the most suitable size, and U,D,L,R to the most comfortable position.



Note:

Please remove the glasses, contact lenses, hairpins, hair accessories, etc. before wearing to avoid discomfort during use.
Don't use this product when your hands are wet. Do not drop and step on the product.

Start up the product

Press" "continuously for 2 seconds,to start the product.
After startup,the machine starts to work and indicator light is always on; the default is automatic mode and the intensity is the weakest.The working time is 15 minutes.

Warning

this product has a timing function, the product will stop working automatically when the time is over.

Method of use(AR206)

Switching mode

Press " " for less than one second to switch mode. Operation once, switch to clear mode; Operate twice, switch to dynamic mode; Operation three times, switch to beauty mode; Operation four times, switch to sleep mode; Operation five times,, switch back to auto mode, and repeat.

Mode introduction

- Operational mode
- Automatic mode = [Change the technique] + [music per 30 seconds]
 - Clear mode = [Vibration massage] + [Hot compress massage] + [music]
 - Vitality mode = [Air pressure] + [vibration massage] + [music]
 - Beauty mode = [Hot compress massage] + [music]
 - Sleep mode = [Air pressure massage]

Music control

Press twice " " button are less than a second to turn off the music. Press" "twice again, less than one second to open music.
Close the product:
In working, press the " " continuously for 2 seconds and shutdown the machine.

Transport and storage

Storage environment: normal temperature, no corrosive gas, well ventilated, dry indoors.
Transportation: This product have easy earthquake protection facilities available, suitable for air, railway, highway and ship transportation. It should avoid the bad conditions such as rain, sun exposure, inversion and collision.
When the storage time exceeds three months, it should be checked out of the box and powered up for half an hour.

Cleaning and storage

Cleaning

Before cleaning, turn off the power and unplug the power adapter from the outlet. Use a soft, dry cloth to clean the surface of the product.

Warning

Be sure to turn off the power when cleaning, . Do not use corrosive cleaning agents such as benzene or thinner to clean the product, as this may cause discoloration, surface aging, structural damage, etc.

Storage method:

If it is not used for a long time, please cover it with a dust cover (self-prepared), and put the product in a dry, ventilated room to avoid sun exposure.

Warning

this product contains battery, never approach a hot place or put it in a fire.

Technical Parameters

Product parameters

Product name : Wireless Music Eye Massager

Product model: AR206

Rated voltage : 5V Rated power. 3W

Executive standard:GB4706.1-2005 GB 4706.10-2008

Packing list

EYE massager*1 USB Data line*1 Instruction manual*1

Reserved

Improper disposal of waste will cause environmental pollution

Product qualification certificate

Product Name: Wireless Music Eye Massager

Product model: AR206

Date of manufacture: _____

Examination clerk: _____

Examination conclusions: qualified

Version number: V1.1

Warning Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.