

REV	CO NO.	DESCRIPTION	REVISED BY
A	CHO-01717	INITIAL RELEASE	B. STEVENSON
B	CHO-02106	REMOVED SYMBOLS: NRTL, IONIZING RADIATION, AND FCC ID. ADDED RX SYMBOL. CHANGED FROM BARCODE TO DATA MATRIX. CORRECTED BLACK COLOR TO CMYK #00000 IN MASTER FILE.	B. STEVENSON
C	CHO-02248	ADDED FCC ID, NON-IONIZING RADIATION SYMBOL, AND PHONE NUMBER. ADDED JAMA REQ. TRACEABILITY FILE.	B. STEVENSON

NOTES, UNLESS OTHERWISE SPECIFIED:

1. ARTWORK SHOWN IS FOR REFERENCE ONLY.
2. DYNAMIC DATA FIELDS ARE REPRESENTED WITH PLACEHOLDER INFORMATION INDICATING FORMAT.

DATE OF MANUFACTURING: YYYY-MM-DD
SERIAL NUMBER: XXXXXX

3. UDI DATA MATRIX MUST BE IN GS1 FORMATTING CONTAINING THE GTIN NUMBER, DATE OF MANUFACTURING, AND SERIAL NUMBER.


GTIN NUMBER: 95200000000033

THIS IS A PLACEHOLDER NUMBER
THE LAST DIGIT IS A CHECK DIGIT

DATE OF MANUFACTURING: YYMMDD
SERIAL NUMBER: XXXXXX

4. UDI DATA MATRIX MUST HAVE HUMAN READABLE NUMBERS ADJACENT TO THE MACHINE READABLE INFORMATION WITH GS1 APPLICATION IDENTIFIERS IN PARENTHESES.

5. REFER TO DC-00219 FOR BARTENDER FILES.

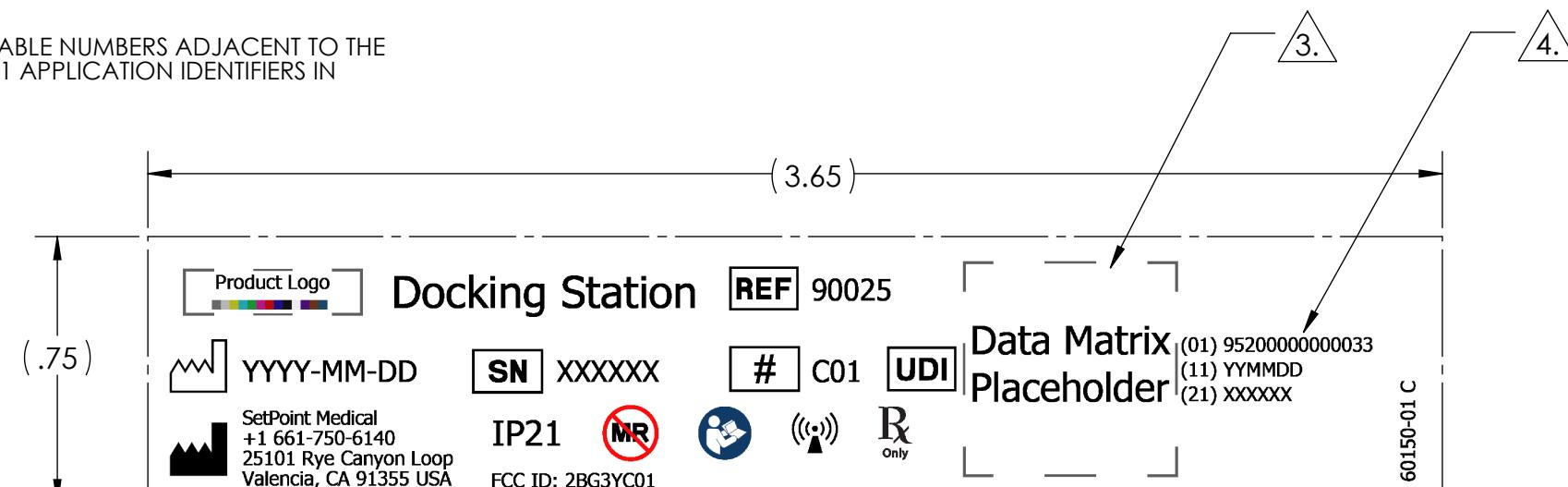

Note: FCC ID is 4-point type size font and is permanently affixed via adhesive

65150-01	C	AW, LABEL, PRODUCT, DOCKING STATION	65150-01 REV C.AI
PART NO.	PART REV.	DESCRIPTION	ARTWORK FILE

This document is **CONFIDENTIAL** and contains designs and/or other information that are the property of SetPoint Medical. This document and the designs and/or other information contained therein, may not, in whole or part, be duplicated or disclosed or used for the manufacture of any part(s) disclosed herein, or for any other purpose, without the prior written permission of SetPoint Medical.

UNLESS OTHERWISE SPECIFIED ALL DIMENSIONS ARE IN INCHES. ALL DIMENSIONS APPLY AFTER FINISH.	SetPoint MEDICAL	25101 RYE CANYON LOOP VALENCIA, CA 91355 (661) 750-6140 (661) 257-3828 FAX
TOLERANCES LINEAR: .X $\pm .1$.XX $\pm .01$.XXX $\pm .005$ ANGULAR: X° $\pm .1$ XX° $\pm .5$		TITLE AW, LABEL, PRODUCT, DOCKING STATION
SURFACE FINISH PER ANSI-B46-1 125/		DRAWN BY B. STEVENSON DATE 21JUN2023 SIZE B DRAWING NUMBER DWG-65150 REV C
DO NOT SCALE SHEET 1 of 2		QAF-002-005 Rev E

NOTES:

-  20. DYNAMIC DATA MUST BE WITHIN THE BOUNDARIES OF THE BOXES SHOWN.
-  21. DATA FIELDS MAY NOT OVERLAP WITH OTHER IMAGE OBJECTS.
-  22. DATA MATRIX X-DIMENSION 0.024IN MINIMUM.

This technical drawing illustrates a label layout for a medical device. The label is divided into several sections with specific dimensions:

- Product Logo:** Located at the top left, with a bounding box of [111, 250, 250, 350].
- Text:** "Docking Station" is centered at [250, 280, 350, 480].
- Reference:** "REF 90025" is located at [250, 510, 350, 630].
- Text:** "YYYY-MM-DD" is at [380, 150, 430, 280].
- Text:** "SN XXXXXX" is at [380, 320, 430, 450].
- Text:** "# C01" is at [380, 510, 430, 600].
- Text:** "UDI" is at [380, 630, 430, 680].
- Text:** "IP21" is at [500, 330, 550, 380].
- Text:** "MR" is at [500, 430, 550, 460] with a red circle around it.
- Icons:** A person icon is at [480, 480, 530, 530], a signal icon is at [480, 550, 530, 600], and a Rx Only icon is at [480, 610, 530, 640].
- Text:** "SetPoint Medical" is at [480, 150, 530, 250].
"+1 661-750-6140" is at [530, 150, 580, 250].
"25101 Rye Canyon Loop" is at [580, 150, 630, 250].
"Valencia, CA 91355 USA" is at [630, 150, 680, 250].
- Text:** "FCC ID: 2BG3YC01" is at [580, 330, 630, 450].
- Text:** "Data Matrix Placeholder" is at [330, 680, 480, 830].
- Text:** "01 95200000000033" is at [330, 840, 380, 950].
"11 YYMMDD" is at [380, 840, 430, 950].
"21 XXXXX" is at [430, 840, 480, 950].
- Text:** "60150-01 C" is at [480, 970, 530, 995].

Dimensions are indicated by dashed lines and arrows:

- 0 to .25
- .05 to .33
- .33 to .43
- .43 to .70
- .70 to .75
- .25 to .85
- .85 to 1.15
- 1.15 to 1.58
- 1.58 to 2.30
- 2.30 to 3.50

Font: REGULAR TAHOMA 6 PT

20. 21. REGULAR
TAHOMA 6 PT

—REGULAR TAHOMA 4 PT

-01 DYNAMIC DATA PLACEMENT

DO NOT SCALE			SIZE	DRAWING NUMBER	REV
SCALE	N/A	SHEET 2 of 2	B	DWG-65150	C