




STERRAD® 100NX® RFID CONFIGURATION



The RFID Transceiver Assembly resides on the interface board within the metal control enclosure assembly.

The RFID Transceiver Assembly interfaces with the RFID antenna board through a coaxial cable. The antenna board is mounted on the delivery system. Peroxide cassettes with an RFID tag are inserted into the delivery system and the system reads and writes to this tag.

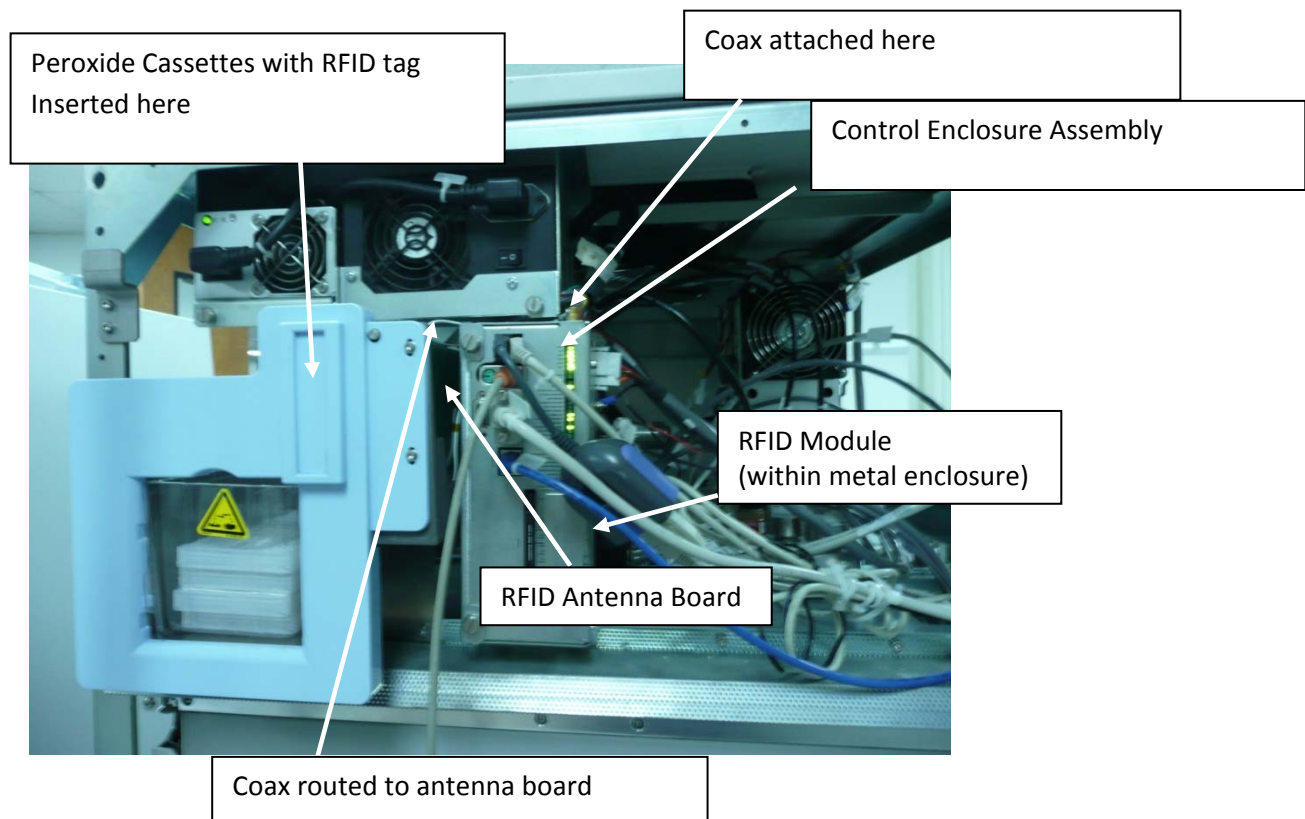

ADVANCED STERILIZATION PRODUCTS®

Division of Ethicon, Inc.
a **Johnson & Johnson** company

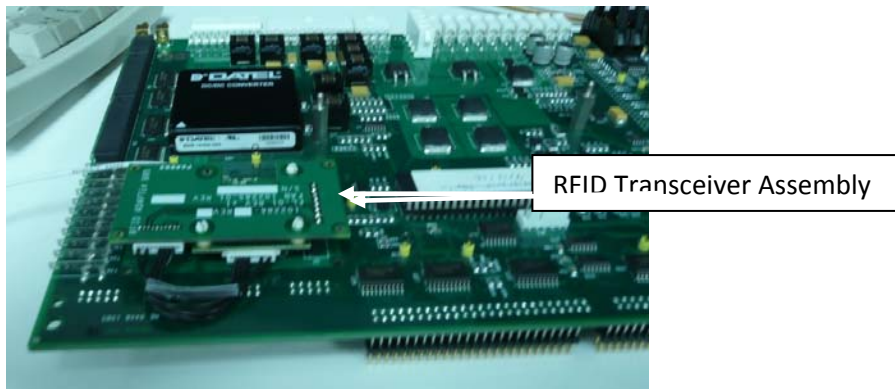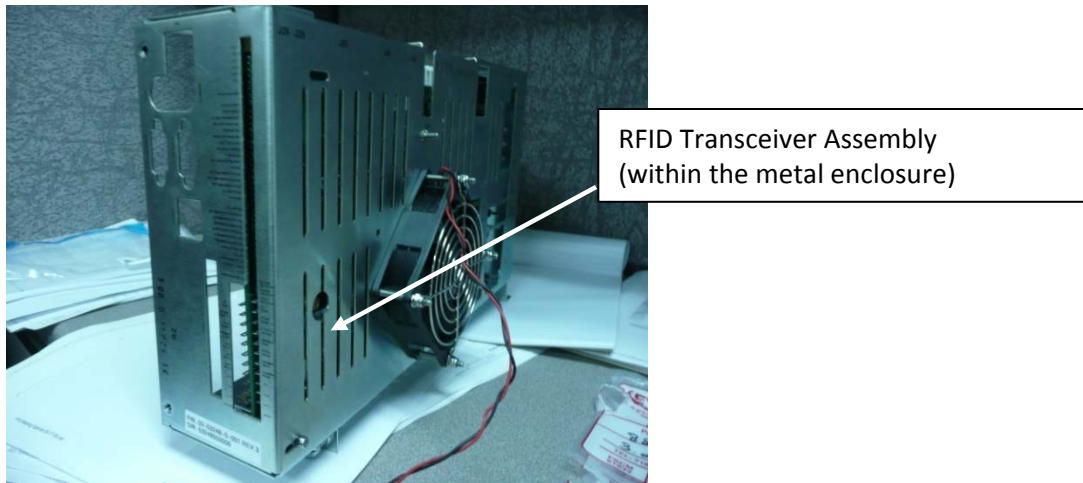

04/19/12

ADVANCED STERILIZATION PRODUCTS®

Division of Ethicon, Inc.

a **Johnson & Johnson** company