## Schmid & Partner **Engineering AG**

Zeughausstrasse 43, 8004 Zurich, Switzerland, Phone +41 1 245 97 00, Fax +41 1 245 97 79

#### **Calibration Certificate**

#### **Dosimetric E-Field Probe**

Туре:	ЕГЗДУ6
Serial Number:	
Place of Calibration:	Zurich
Date of Calibration:	11nty 245; 2002; 114
Calibration Interval:	12 months

Schmid & Partner Engineering AG hereby certifies, that this device has been calibrated on the date indicated above. The calibration was performed in accordance with specifications and procedures of Schmid & Partner Engineering AG.

Wherever applicable, the standards used in the calibration process are traceable to international standards. In all other cases the standards of the Laboratory for EMF and Microwave Electronics at the Swiss Federal Institute of Technology (ETH) in Zurich, Switzerland have been applied.

Calibrated by:

Approved by:

D. Velka)

### Schmid & Partner Engineering AG

Zeughausstrasse 43, 8004 Zurich, Switzerland, Telephone +41 1 245 97 00, Fax +41 1 245 97 79

# Probe ET3DV6

SN:1514

Manufactured:

Last calibration:

Repaired:

Recalibrated:

November 24, 1999

October 25, 2002

June 28, 2002


July 25, 2002


Calibrated for System DASY3


## DASY3 - Parameters of Probe: ET3DV6 SN:1514

Sensitiv	vity in Free	Space		Diode C	ompress	ion	
	NormX	1.70	μV/(V/m) <sup>2</sup>		DCP X	93	mV
	NormY		μV/(V/m) <sup>2</sup>		DCP Y	93	mV
	NormZ		μV/(V/m) <sup>2</sup>		DCP Z	93	mV
Sensiti	vity in Tiss	ue Simu	lating Liquid				
Head	900	MHz	ε <sub>r</sub> = 41.5 ± 5%	σ=	0.97 ± 5% n	nho/m	
Head	835	MHz	e, = 41.5 ± 5%	σ*	0.90 ± 5% n	nho/m	
	ConvF X	6.2	± 9.5% (k=2)		Boundary e	ffect:	
	ConvF Y	6.2	± 9.5% (k=2)		Alpha	0.48	
	ConvF Z	6.2	± 9.5% (k=2)		Depth	2.19	
Head	1800	MHz	e, = 40.0 ± 5%	g=	1.40 ± 5% n	nho/m	
Head	1900		ε, = 40.0 ± 5%		1.40 ± 5% n	nho/m	
	ConvF X	5.2	± 9.5% (k=2)		Boundary e	ffect:	
	ConvF Y	5.2	± 9.5% (k=2)		Alpha	0.57	
	ConvF Z	5.2	± 9.5% (k=2)		Depth	2.30	
Bound	ary Effect						
Head	900	MHz	Typical SAR gradier	nt: 5 % per n	nm		
	Probe Tip to	Boundary			1 mm	2 mm	
			orrection Algorithm		9.9	5.4	
	SAR (%)	With Corre	ection Algorithm		0.2	0.5	
Head	1800	MHz	Typical SAR gradier	nt: 10 % per	mm .		
	Probe Tip to	Boundary			1 mm	2 mm	
			orrection Algorithm		12.7	8.2	
	SAR <sub>be</sub> [%]	With Corre	ection Algorithm		0.2	0.2	
Senso	r Offset						
	Probe Tip to	Sensor Ce	enter	2.7		mm	
	Optical Surf	Optical Surface Detection				mm	


## Receiving Pattern ( $\phi$ ), $\theta = 0^{\circ}$

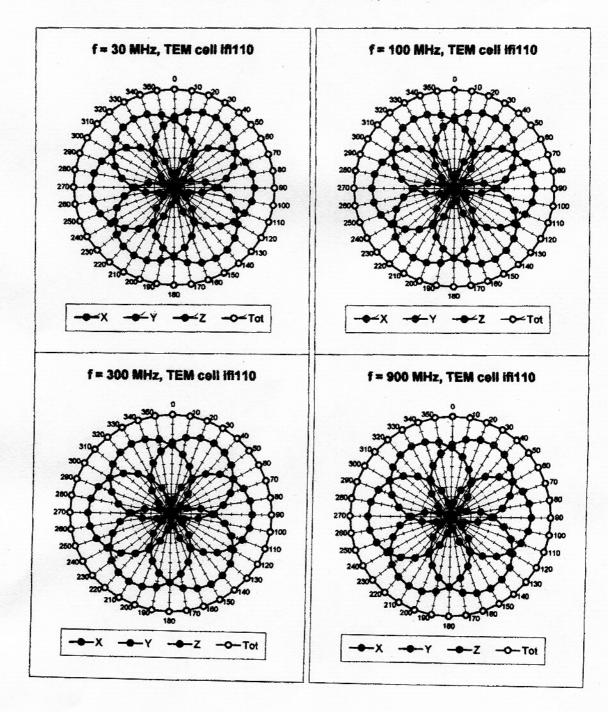

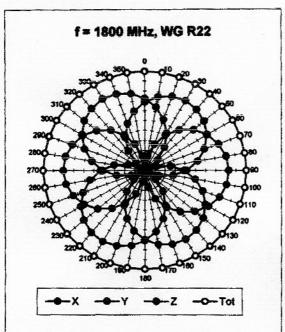

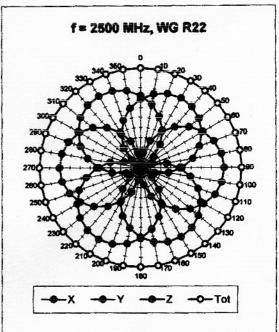


Isotropy Error ( $\phi$ ),  $\theta = 0^{\circ}$ 

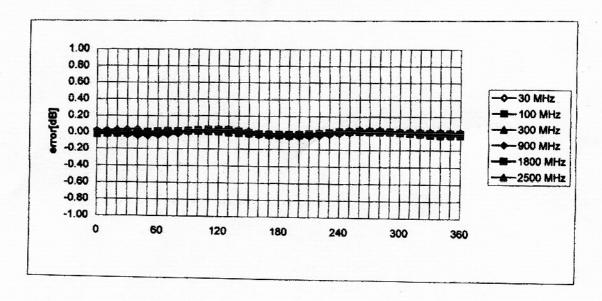


## Frequency Response of E-Field


(TEM-Cell:ifi110, Waveguide R22)

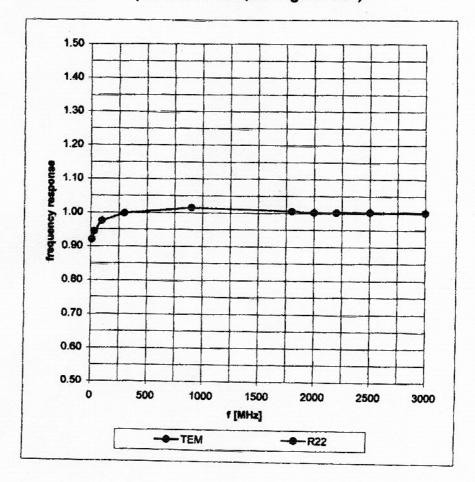


## Dynamic Range f(SAR<sub>brain</sub>)


(Waveguide R22)

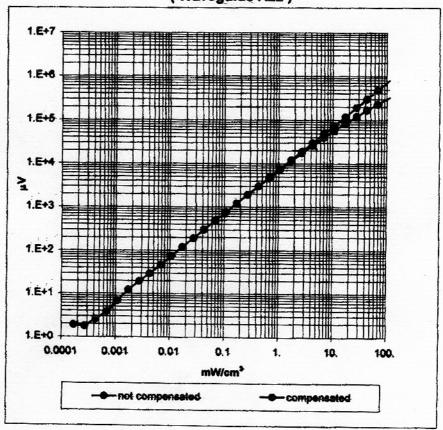

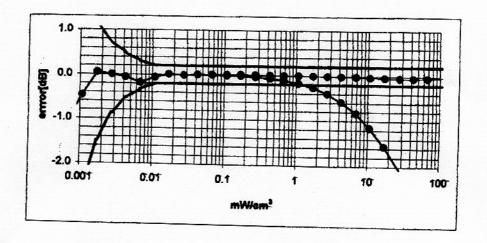

#### **Conversion Factor Assessment**

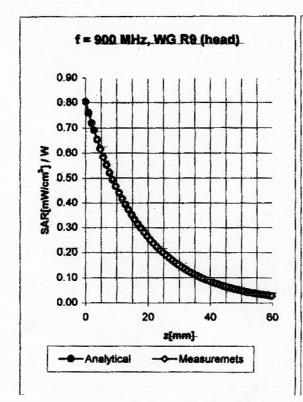

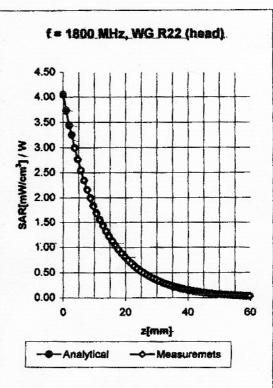


Head	900 MHz	$\varepsilon_r = 41.5 \pm 5\%$	$\sigma = 0.97 \pm 5\% \text{ mho/m}$
Head	835 MHz	e, = 41.5 ± 5%	o = 0.90 ± 5% mho/m
	ConvF X	<b>6.2</b> ± 9.5% (k=2)	Boundary effect:
	ConvF Y	6.2 ± 9.5% (k=2)	Aipha 0.48
	ConvF Z	6.2 ± 9.5% (k=2)	Depth <b>2.19</b>

Head	1800 MHz	$s_r = 40.0 \pm 5\%$	σ = 1.40 ± 5% mho/m
Head	1900 MHz	$\epsilon_{\rm r}$ = 40.0 ± 5%	o = 1.40 ± 5% mho/m
	ConvF X	5.2 ± 9.5% (k=2)	Boundary effect:
	ConvF Y	<b>5.2</b> ± 9.5% (k=2)	Alpha 0.57
	ConvF Z	5.2 ± 9.5% (k=2)	Depth 2.30

#### **Deviation from Isotropy in HSL**

Error  $(\theta,\phi)$ , f = 900 MHz

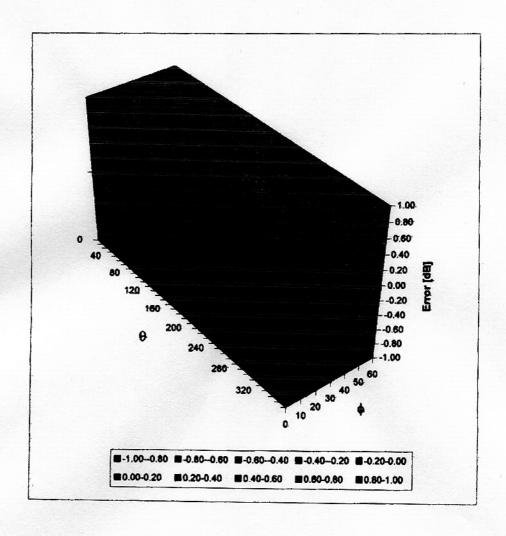

#### Schmid & Partner Engineering AG

Zeughausstrasse 43, 8004 Zurich, Switzerland, Phone +41 1 245 97 00, Fax +41 1 245 97 79

#### **Additional Conversion Factors**

for Dosimetric E-Field Probe

Туре:	ET3DV6
Serial Number:	1514
Place of Assessment:	Zurich
Date of Assessment:	July 26, 2002
Probe Calibration Date:	July 25, 2002

Schmid & Partner Engineering AG hereby certifies that conversion factor(s) of this probe have been evaluated on the date indicated above. The assessment was performed using the FDTD numerical code SEMCAD of Schmid & Partner Engineering AG. Since the evaluation is coupled with measured conversion factors, it has to be recalculated yearly, i.e., following the re-calibration schedule of the probe. The uncertainty of the numerical assessment is based on the extrapolation from measured value at 900 MHz or at 1800 MHz.

Assessed by:

Meno Katz-

#### Dosimetric E-Field Probe ET3DV6 SN:1514

Conversion factor (± standard deviation)

835 MHz	ConvF	6.0 ± 8%	$\varepsilon_r = 55.2 \pm 5\%$ $\sigma = 0.97 \pm 5\% \text{ mho/m}$ (body tissue)
900 MHz	ConvF	5.9 ± 8%	$\varepsilon_r = 55.0 \pm 5\%$ $\sigma = 1.05 \pm 5\% \text{ mho/m}$ (body tissue)
1800 MHz	ConvF	4.8 ± 8%	$\varepsilon_r = 53.3 \pm 5\%$ $\sigma = 1.52 \pm 5\% \text{ mho/m}$ (body tissue)
1900 MHz	ConvF	4.6 ± 8%	$\varepsilon_r = 53.3 \pm 5\%$ $\sigma = 1.52 \pm 5\% \text{ mho/m}$ (body tissue)
1950 MHz	ConvF	4.5 ± 8%	$\varepsilon_r = 53.3 \pm 5\%$ $\sigma = 1.52 \pm 5\% \text{ mho/m}$ (body tissue)
1950 MHz	ConvF	4.9 ± 8%	$g_r = 40.0 \pm 5\%$ $\sigma = 1.40 \pm 5\% \text{ mho/m}$ (head tissue)