### **NCL CALIBRATION LABORATORIES**

Calibration File No: DC-916 Project Number: APREL-ALSAS 10U

# CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the **NCL CALIBRATION LABORATORIES** by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

**APREL Validation Dipole** 

Manufacturer: APREL Laboratories
Part number: ALS-D-1900-S-2
Frequency: 1900 MHz
Serial No: 1900-210-00704

Customer: APREL

Calibrated: 1<sup>st</sup> March 2009 Released on: 1<sup>st</sup> march 2009

Palescad By	
neleased by.	



51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4162

· -	$\boldsymbol{\smallfrown}$	n	~	111	$\sim$	ns
u	u		u	ILI	u	113

D. Brooks

Dipole 1900-210-00704 was a re-calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} \, +/- \, 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} \, +/- \, 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this device has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

This page has been reviewed for content and attested to by signature within this document.

## **Calibration Results Summary**

The following results relate the Calibrated Dipole and should be used as a quick reference for the user.

#### **Mechanical Dimensions**


**Length:** 68.0 mm **Height:** 39.5 mm

## **Electrical Specification**

SWR: 1.10 U Return Loss: -26.3 dB Impedance:  $48.4 \Omega$ 

### **System Validation Results**

Frequency	1 Gram	10 Gram	Peak
1900 MHz	37.96	19.81	70.56

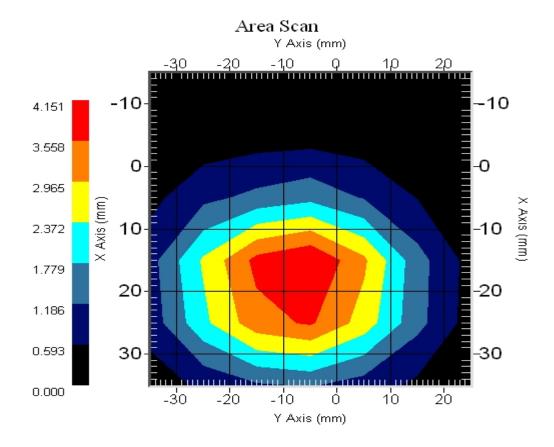

#### Introduction

This Calibration Report has been produced in line with the SSI Dipole Calibration Procedure SSI-TP-018-ALSAS. The results contained within this report are for Validation Dipole 1900-210-00704. The calibration routine consisted of a three-step process. Step 1 was a mechanical verification of the dipole to ensure that it meets the mechanical specifications. Step 2 was an Electrical Calibration for the Validation Dipole, where the SWR, Impedance, and the Return loss were assessed. Step 3 involved a System Validation using the ALSAS-10U, along with APREL E-020 130 MHz to 26 GHz E-Field Probe Serial Number 212.

#### References

SSI-TP-018-ALSAS Dipole Calibration Procedure
SSI-TP-016 Tissue Calibration Procedure
IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average
Specific Absorption Rate (SAR) in the Human Body Due to Wireless
Communications Devices: Experimental Techniques"

### **Conditions**

Dipole 1900-210-00704 was a re-calibration.

**Ambient Temperature of the Laboratory:** 22 °C +/- 0.5 °C **Temperature of the Tissue:** 20 °C +/- 0.5 °C

# **Dipole Calibration Results**

### **Mechanical Verification**

APREL	APREL	Measured	Measured
Length	Height	Length	Height
68.0 mm	39.5 mm	69.1mm	40.3 mm

### **Tissue Validation**


Head Tissue 1900 MHz	Measured
Dielectric constant, ε <sub>r</sub>	40.0
Conductivity, σ [S/m]	1.40

### **Electrical Calibration**


Test	Result
S11 R/L	-26.3 dB
SWR	1.10 U
Impedance	48.4 Ω

The Following Graphs are the results as displayed on the Vector Network Analyzer.


#### **S11 Parameter Return Loss**

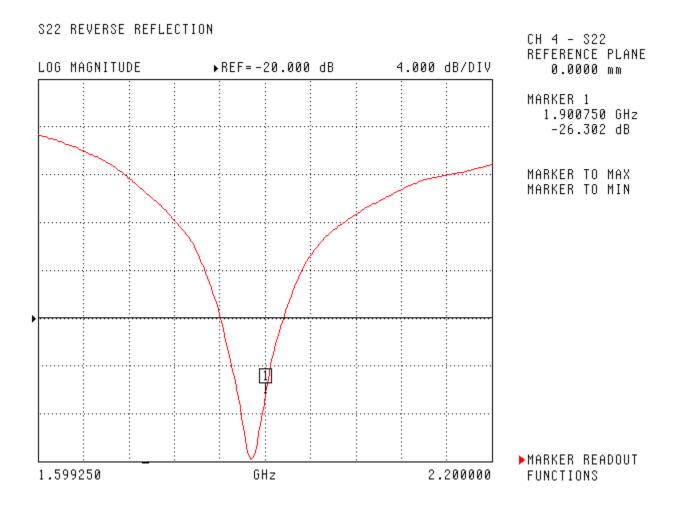

## SWR

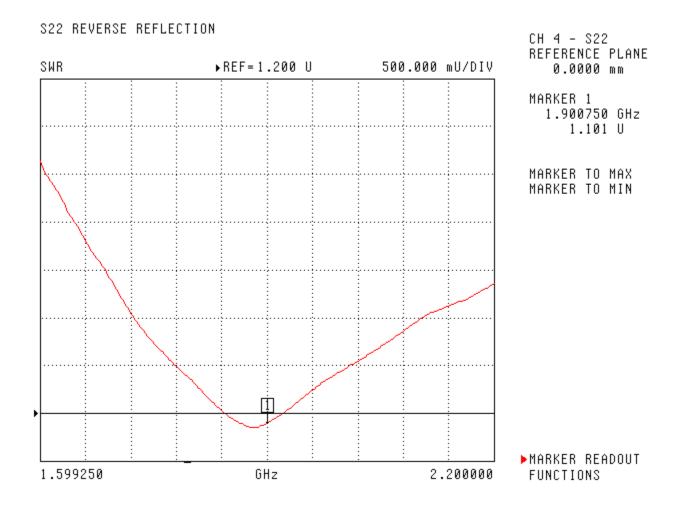


# **Smith Chart Dipole Impedance**

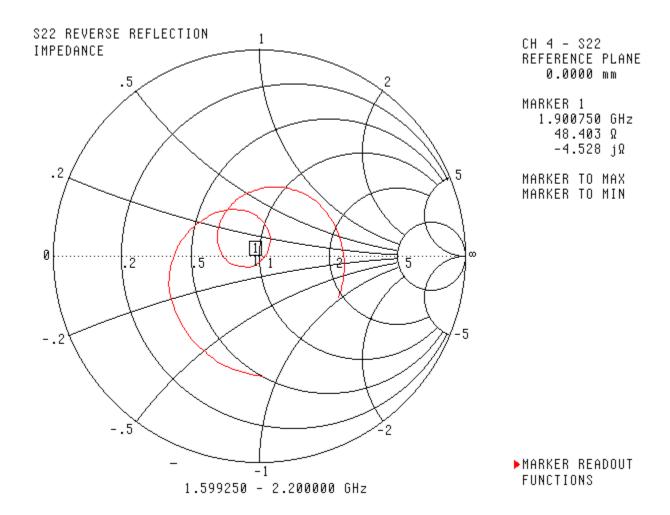

## **System Validation Results Using the Electrically Calibrated Dipole**

Head Tissue Frequency	1 Gram	10 Gram	Peak Above Feed Point
1900 MHz	37.96	19.81	70.56

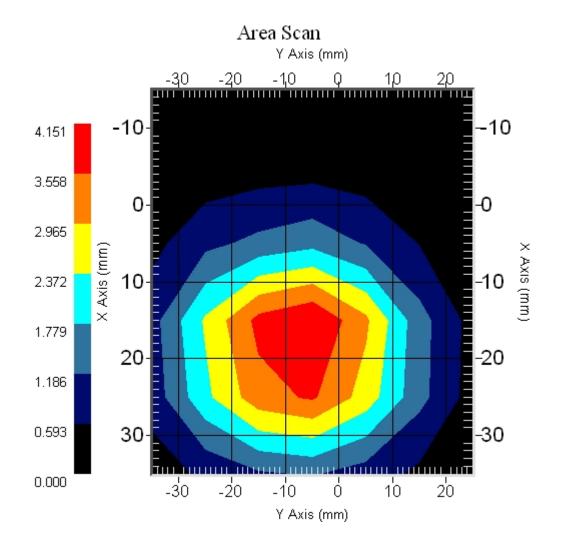

# **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2007.

## **NCL CALIBRATION LABORATORIES**

Calibration File No.: CP-913

Client.: APREL

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the **NCL CALIBRATION LABORATORIES** by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

Equipment: Miniature Isotropic RF Probe 835 MHz

Manufacturer: APREL Laboratories

Model No.: E-020

Serial No.: 225

Calibration in Body Tissue

Calibration Procedure: SSI/DRB-TP-D01-032-E020-V2

Project No: Internal APREL

Calibrated: 3<sup>rd</sup> May 2009 Released on: 3<sup>rd</sup> May 2009

This Calibration Certificate is Incomplete Unless Accompanied with the Calibration Results Summary

Released By:

NCL CALIBRATION LABORATORIES

51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4161

#### Introduction

This Calibration Report reproduces the results of the calibration performed in line with the SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure. The results contained within this report are for APREL E-Field Probe E-020 226.

### References

SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average Specific Absorption Rate (SAR) in the Human Body Due to Wireless Communications Devices: Experimental Techniques" SSI-TP-011 Tissue Calibration Procedure

#### **Conditions**

Probe 225 was a re-calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this probe has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Jesse Hones

## **Calibration Results Summary**

**Probe Type**: E-Field Probe E-020

Serial Number: 226

Frequency: 835 MHz

Sensor Offset: 1.56 mm

Sensor Length: 2.5 mm

**Tip Enclosure:** Ertalyte\*

**Tip Diameter:** <2.9 mm

**Tip Length:** 60 mm

Total Length: 290 mm

## Sensitivity in Air

 Channel X:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Y:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Z:
  $1.2 \, \mu V/(V/m)^2$ 

**Diode Compression Point:** 95 mV

<sup>\*</sup>Resistive to recommended tissue recipes per IEEE-1528

## **Sensitivity in Body Tissue**

Frequency: 835 MHz

**Epsilon:** 56.1 (+/-5%) **Sigma:** 0.95 S/m (+/-5%)

ConvF

Channel X: 6.6

Channel Y: 6.6

Channel Z: 6.6

Tissue sensitivity values were calculated using the load impedance of the APREL Laboratories Daq-Paq.


## **Boundary Effect:**

Uncertainty resulting from the boundary effect is less than 2% for the distance between the tip of the probe and the tissue boundary, when less than 1.44mm.


## **Spatial Resolution:**


The measured probe tip diameter is 2.9 mm (+/- 0.01 mm) and therefore meets the requirements of SSI/DRB-TP-D01-032 for spatial resolution.

# Receiving Pattern 835 MHz (Air)

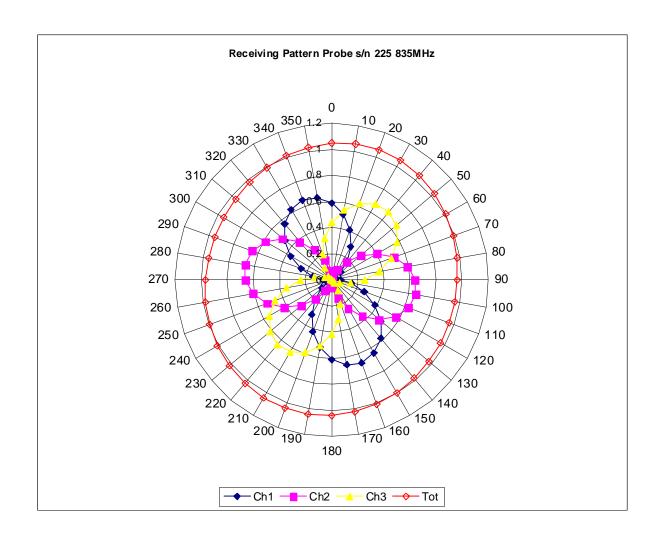

# **Isotropy Error 835 MHz (Air)**

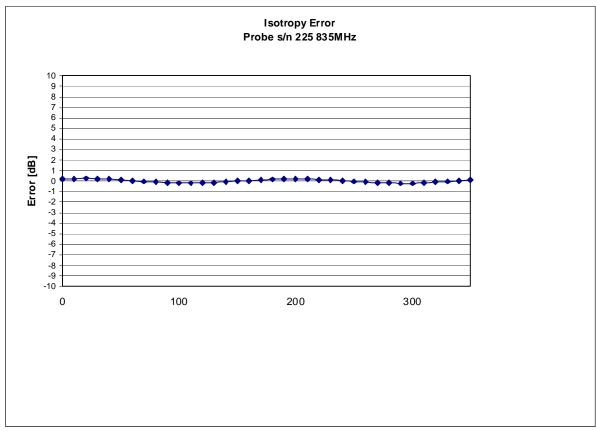

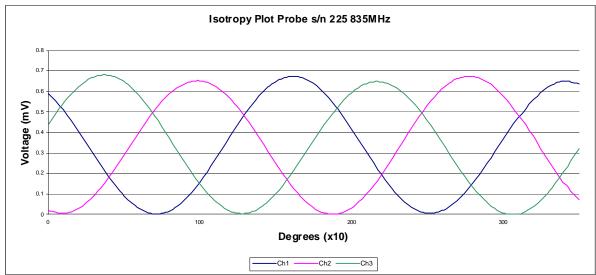


**Isotropicity in Tissue:** 

0.10 dB


Page 6 of 10

# **Dynamic Range**

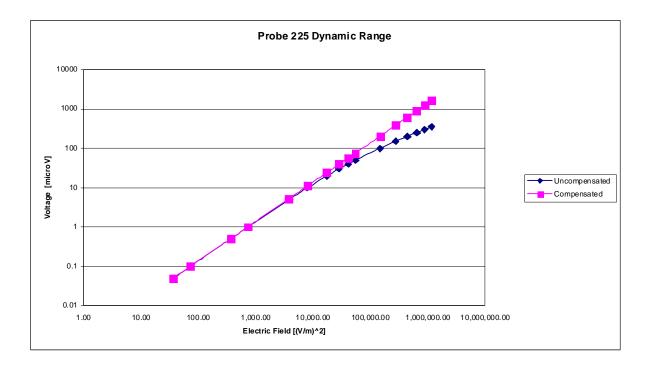

## **Video Bandwidth**

### **Probe Frequency Characteristics**

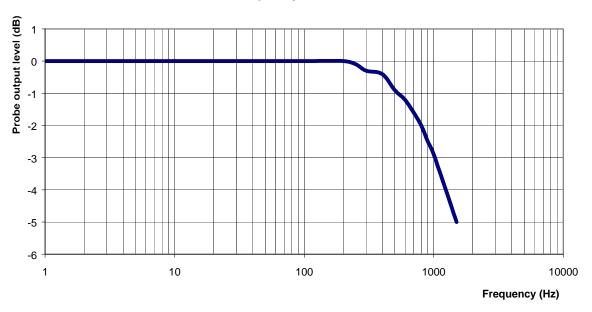

Video Bandwidth at 500 Hz 1 dB Video Bandwidth at 1.02 KHz: 3 dB

## **Conversion Factor Uncertainty Assessment**

Frequency: 835MHz

**Epsilon:** 56.1 (+/-5%) **Sigma:** 0.95 S/m (+/-5%)

ConvF

**Channel X:** 6.6 7%(K=2)

**Channel Y:** 6.6 7%(K=2)

**Channel Z:** 6.6 7%(K=2)

To minimize the uncertainty calculation all tissue sensitivity values were calculated using a load impedance of 5 M $\Omega$ .

### **Boundary Effect:**

For a distance of 1.4mm the evaluated uncertainty (increase in the probe sensitivity) is less than 2%.

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2007.

### **NCL CALIBRATION LABORATORIES**

Calibration File No.: CP-925

Client.: APREL

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the NCL CALIBRATION LABORATORIES by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

Equipment: Miniature Isotropic RF Probe 1800 MHz

Manufacturer: APREL Laboratories

Model No.: E-020

Serial No.: 206

**BODY Calibration** 

Calibration Procedure: SSI/DRB-TP-D01-032-E020 Project No: RFEL-Probe-206-Loan-5366

> Calibrated: 18<sup>th</sup> September 2008 Released on: 19<sup>th</sup> September 2008

This Calibration Certificate is Incomplete Unless Accompanied with the Calibration Results Summary

This Calibration Certificate is Incomplete Unless Accomplanied with the Calibration Results Summary

Released By:

NCL CALIBRATION LABORATORIES

17 Bentley Avenue NEPEAN, ONTARIO CANADA K2E 6T7 CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4161 FAX: (613) 820-4161

#### Introduction

This Calibration Report reproduces the results of the calibration performed in line with the SSI/DRB-TP-D01-032-E020 E-Field Probe Calibration Procedure. The results contained within this report are for APREL E-Field Probe E-020 206.

#### References

SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure

IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average Specific Absorption Rate (SAR) in the Human Head Due to Wireless Communications Devices: Experimental Techniques"

SSI-TP-011 Tissue Calibration Procedure

IEC 62209 "Human exposure to radio frequency fields from hand-held and Head-mounted wireless communication devices – Human models, instrumentation, and procedures –Part 1 & 2: Procedure to determine the Specific Absorption Rate (SAR) for hand-held devices used in close proximity of the ear (frequency range of 300 MHz to 3 GHz)"

IEEE 1309 Draft Standard for Calibration of Electromagnetic Field Sensors and Probes, Excluding Antennas, from 9kHz to 40GHz

### **Conditions**

Probe 206 is a re-calibration.

Ambient Temperature of the Laboratory: 22 °C +/- 0.5 °C Temperature of the Tissue: 21 °C +/- 0.5 °C

We the undersigned attest that to the best of our knowledge the calibration of

this probe has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Jesse Hones

## **Calibration Results Summary**

**Probe Type**: E-Field Probe E-020

Serial Number: 206

Frequency: 1800 MHz

Sensor Offset: 1.56 mm

Sensor Length: 2.5 mm

Tip Enclosure: Ertalyte\*

**Tip Diameter:** 5 mm

**Tip Length:** 60 mm

Total Length: 290 mm

## **Sensitivity in Air**

 Channel X:
  $0.79 \, \mu V/(V/m)^2$  

 Channel Y:
  $0.79 \, \mu V/(V/m)^2$  

 Channel Z:
  $0.79 \, \mu V/(V/m)^2$ 

**Diode Compression Point:** 101 mV

NOTE:

Amplifier impedance for this probe is not normalized in software.

<sup>\*</sup>Resistive to recommended tissue recipes per IEEE-1528

## **Sensitivity in Body Tissue**

Frequency: 1800 MHz

**Epsilon:** 54.0 (+/-5%) **Sigma:** 1.40 S/m (+/-10%)

ConvF

Channel X: 4.29

Channel Y: 4.29

**Channel Z:** 4.29

Tissue sensitivity values were calculated using the load impedance of the APREL Laboratories Daq-Paq.


## **Boundary Effect:**

Uncertainty resulting from the boundary effect is less than 2% for the distance between the tip of the probe and the tissue boundary, when less than 2.65mm.


## **Spatial Resolution:**


The measured probe tip diameter is 5 mm (+/- 0.01 mm) and therefore meets the requirements of SSI/DRB-TP-D01-032 for spatial resolution.

# Receiving Pattern 1800 MHz (Air)

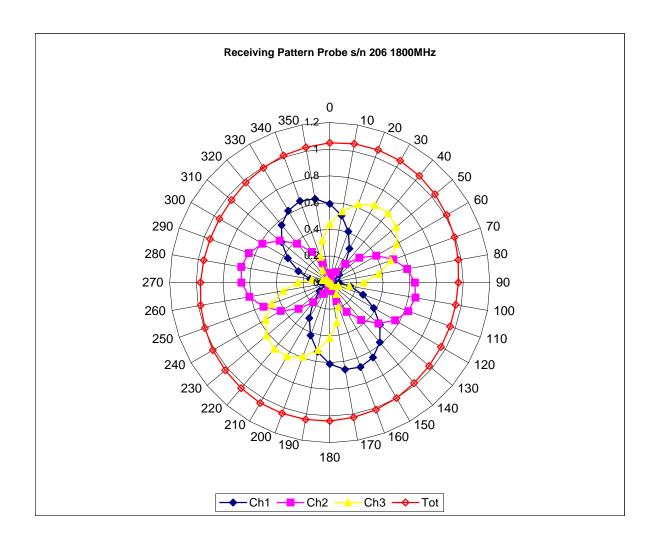


# Isotropy Error 1800 MHz (Air)

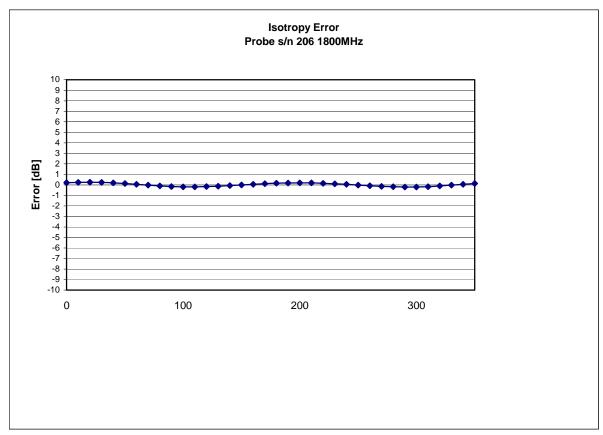

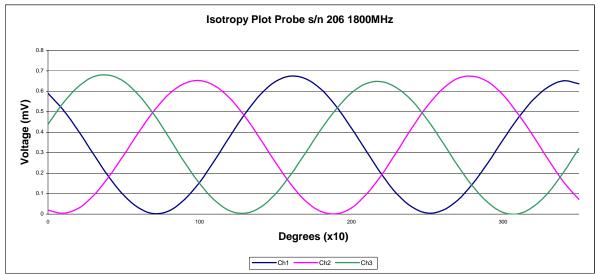


**Isotropicity:** 0.10 dB

# **Dynamic Range**

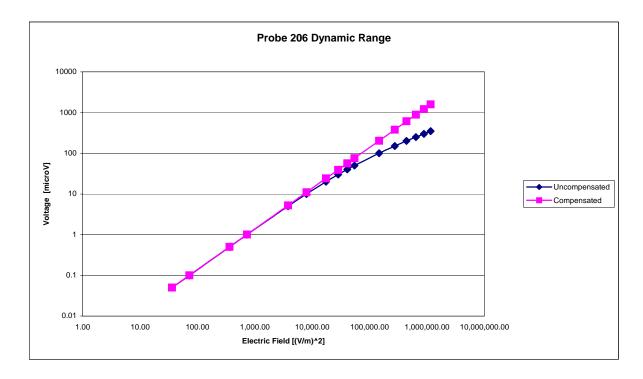

## **Video Bandwidth**

## **Probe Frequency Characteristics**

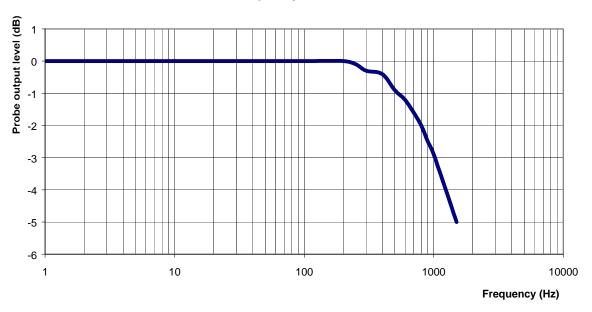

Video Bandwidth at 500 Hz 1 dB Video Bandwidth at 1.02 KHz: 3 dB

## **Conversion Factor Uncertainty Assessment**

Frequency: 1800MHz

**Epsilon:** 54.0 (+/-5%) **Sigma:** 1.40 S/m (+/-10%)

ConvF

**Channel X:** 4.29 7%(K=2)

**Channel Y:** 4.29 7%(K=2)

**Channel Z:** 4.29 7%(K=2)

To minimize the uncertainty calculation all tissue sensitivity values were calculated using a load impedance of 5 M $\Omega$ .

## **Boundary Effect:**

For a distance of 2.6mm the evaluated uncertainty (increase in the probe sensitivity) is less than 2%.

## **NCL Calibration Laboratories**

Division of APREL Laboratories.

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2008.

### **NCL CALIBRATION LABORATORIES**

Calibration File No.: CP-914

Client.: APREL

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the **NCL CALIBRATION LABORATORIES** by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

Equipment: Miniature Isotropic RF Probe 1900 MHz

Manufacturer: APREL Laboratories

Model No.: E-020 Serial No.: 225

Calibration in Body Tissue

Calibration Procedure: SSI/DRB-TP-D01-032-E020-V2

Project No: Internal APREL

Calibrated: 3<sup>rd</sup> May 2009 Released on: 3<sup>rd</sup> May 2009

This Calibration Certificate is Incomplete Unless Accompanied with the Calibration Results Summary

Released By:

NCL CALIBRATION LABORATORIES

51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4161

#### Introduction

This Calibration Report reproduces the results of the calibration performed in line with the SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure. The results contained within this report are for APREL E-Field Probe E-020 226.

### References

SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average Specific Absorption Rate (SAR) in the Human Body Due to Wireless Communications Devices: Experimental Techniques" SSI-TP-011 Tissue Calibration Procedure

#### **Conditions**

Probe 225 was a re-calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this probe has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Jesse Hones

## **Calibration Results Summary**

**Probe Type**: E-Field Probe E-020

Serial Number: 226

Frequency: 1900 MHz

Sensor Offset: 1.56 mm

Sensor Length: 2.5 mm

**Tip Enclosure:** Ertalyte\*

**Tip Diameter:** <2.9 mm

**Tip Length:** 60 mm

Total Length: 290 mm

## Sensitivity in Air

 Channel X:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Y:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Z:
  $1.2 \, \mu V/(V/m)^2$ 

**Diode Compression Point:** 95 mV

<sup>\*</sup>Resistive to recommended tissue recipes per IEEE-1528

## **Sensitivity in Body Tissue**

Frequency: 1900 MHz

**Epsilon:** 54.0 (+/-5%) **Sigma:** 1.45 S/m (+/-5%)

ConvF

Channel X: 4.9

Channel Y: 4.9

Channel Z: 4.9

Tissue sensitivity values were calculated using the load impedance of the APREL Laboratories Daq-Paq.


## **Boundary Effect:**

Uncertainty resulting from the boundary effect is less than 2% for the distance between the tip of the probe and the tissue boundary, when less than 1.44mm.


## **Spatial Resolution:**


The measured probe tip diameter is 2.9 mm (+/- 0.01 mm) and therefore meets the requirements of SSI/DRB-TP-D01-032 for spatial resolution.

# Receiving Pattern 1900 MHz (Air)

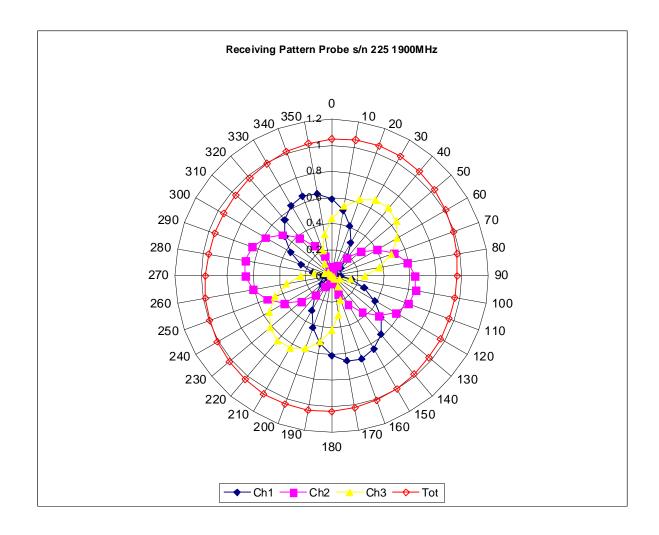

# Isotropy Error 1900 MHz (Air)

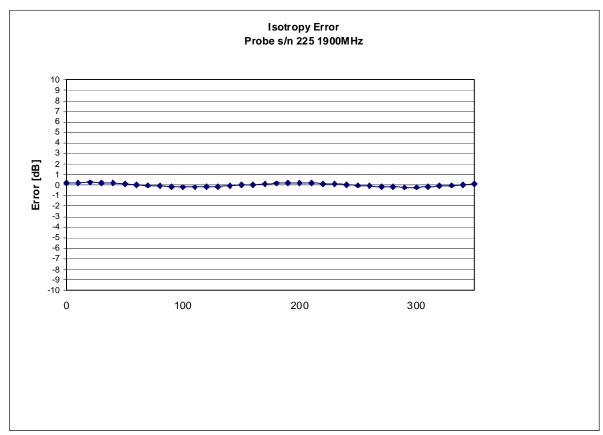

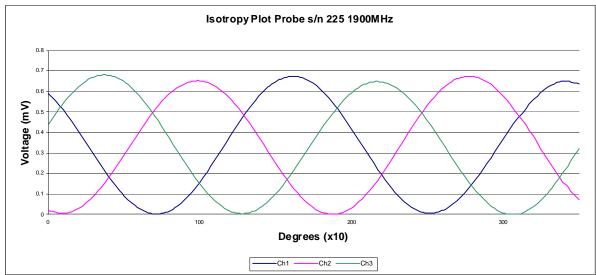


**Isotropicity in Tissue:** 

0.10 dB


Page 6 of 10

# **Dynamic Range**

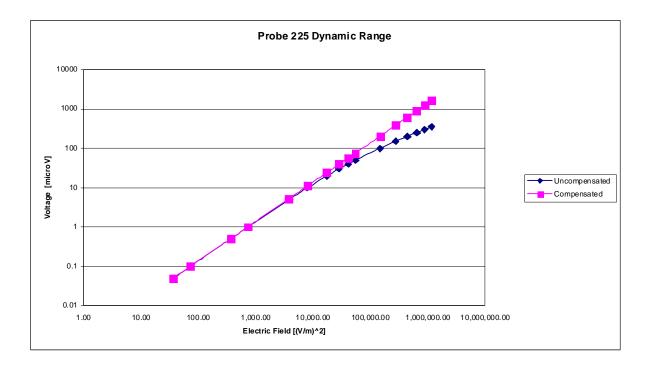

## **Video Bandwidth**

### **Probe Frequency Characteristics**

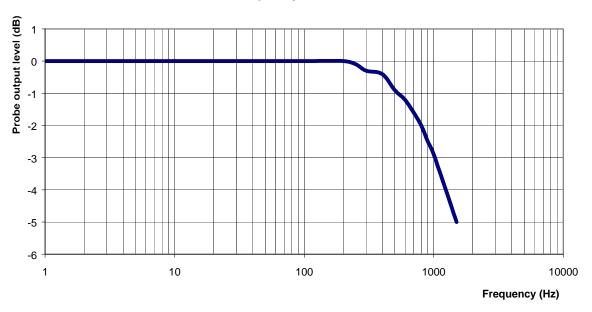

Video Bandwidth at 500 Hz 1 dB Video Bandwidth at 1.02 KHz: 3 dB

## **Conversion Factor Uncertainty Assessment**

Frequency: 1900MHz

**Epsilon:** 54.0 (+/-5%) **Sigma:** 1.45 S/m (+/-5%)

ConvF

**Channel X:** 4.9 7%(K=2)

**Channel Y:** 4.9 7%(K=2)

**Channel Z:** 4.9 7%(K=2)

To minimize the uncertainty calculation all tissue sensitivity values were calculated using a load impedance of 5 M $\Omega$ .

### **Boundary Effect:**

For a distance of 1.4mm the evaluated uncertainty (increase in the probe sensitivity) is less than 2%.

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2007.

### **NCL CALIBRATION LABORATORIES**

Calibration File No.: CP-885

Client.: APREL

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the NCL CALIBRATION LABORATORIES by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

Equipment: Miniature Isotropic RF Probe 2450 MHz

Manufacturer: APREL Laboratories

Model No.: E-030 Serial No.: 222

Calibration in Body Tissue

Calibration Procedure: SSI/DRB-TP-D01-032-E020-V2

Project No: Internal APREL

Calibrated: 3<sup>rd</sup> May 2009 Released on: 3<sup>rd</sup> May 2009

This Calibration Certificate is Incomplete Unless Accompanied with the Calibration Results Summary

Released By:

NCL CALIBRATION LABORATORIES

51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4161

#### Introduction

This Calibration Report reproduces the results of the calibration performed in line with the SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure. The results contained within this report are for APREL E-Field Probe E-030 222.

#### References

SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average Specific Absorption Rate (SAR) in the Human Body Due to Wireless Communications Devices: Experimental Techniques" SSI-TP-011 Tissue Calibration Procedure

#### **Conditions**

Probe 222 was a new probe taken from stock prior to calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this probe has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Jesse Hones

## **Calibration Results Summary**

**Probe Type**: E-Field Probe E-030

Serial Number: 222

Frequency: 2450 MHz

Sensor Offset: 0.44 mm

Sensor Length: 2.5 mm

**Tip Enclosure:** Ertalyte\*

**Tip Diameter:** <2.9 mm

**Tip Length:** 60 mm

Total Length: 290 mm

# Sensitivity in Air

 Channel X:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Y:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Z:
  $1.2 \, \mu V/(V/m)^2$ 

**Diode Compression Point:** 95 mV

<sup>\*</sup>Resistive to recommended tissue recipes per IEEE-1528

### **Sensitivity in Body Tissue**

Frequency: 2450 MHz

**Epsilon:** 52.7 (+/-5%) **Sigma:** 1.95 S/m (+/-5%)

ConvF

Channel X: 4.01

Channel Y: 4.01

**Channel Z:** 4.01

Tissue sensitivity values were calculated using the load impedance of the APREL Laboratories Daq-Paq and corrected for broadband calibration factor.

### **Boundary Effect:**

Uncertainty resulting from the boundary effect is less than 2% for the distance between the tip of the probe and the tissue boundary, when less than 1.44mm.

## **Spatial Resolution:**


The measured probe tip diameter is 2.9 mm (+/- 0.01 mm) and therefore meets the requirements of SSI/DRB-TP-D01-032 for spatial resolution.

#### **Broad Band Calibration:**


The probe was assessed for sensitivity and conversion factor using a +/- 40MHz deviation from the centre frequency.


Deviation at -40MHz: -1.56% Deviation at +40MHz: +1.3%

# Receiving Pattern 2450 MHz (Air)

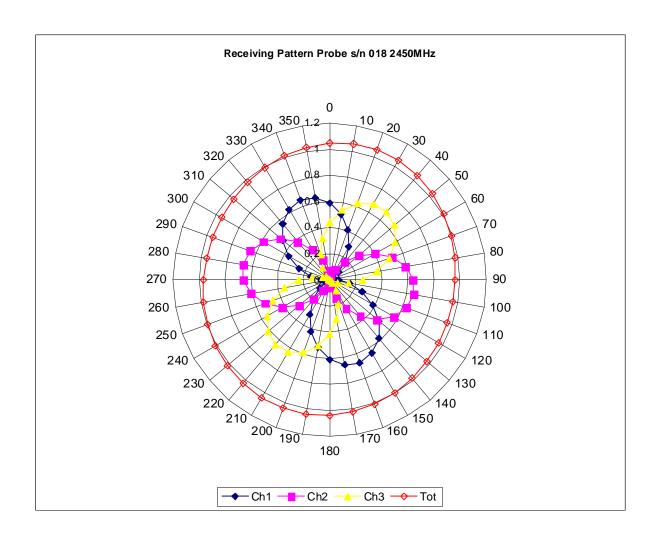

# Isotropy Error 2450 MHz (Air)

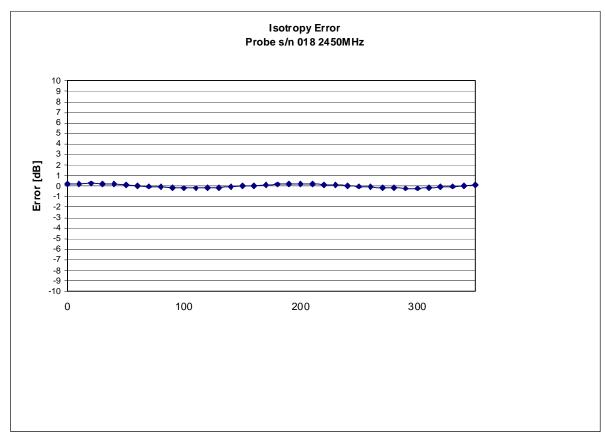

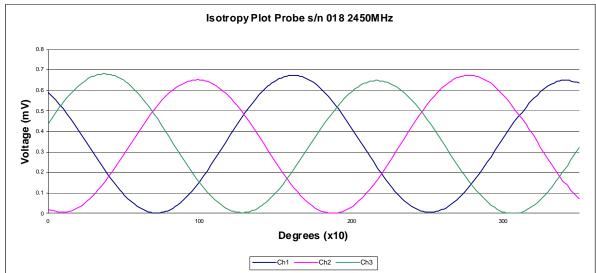


**Isotropicity in Tissue:** 

0.10 dB


Page 6 of 10

# **Dynamic Range**

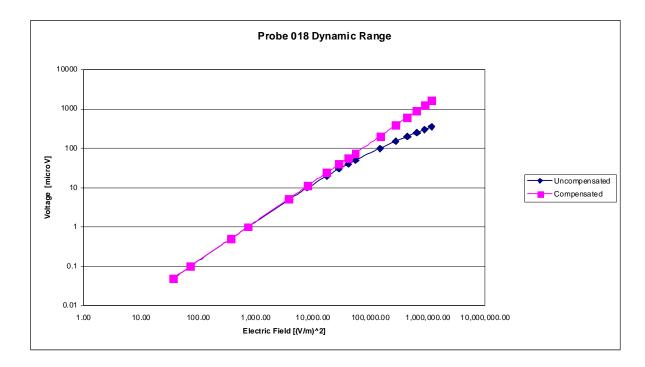

## **Video Bandwidth**

### **Probe Frequency Characteristics**

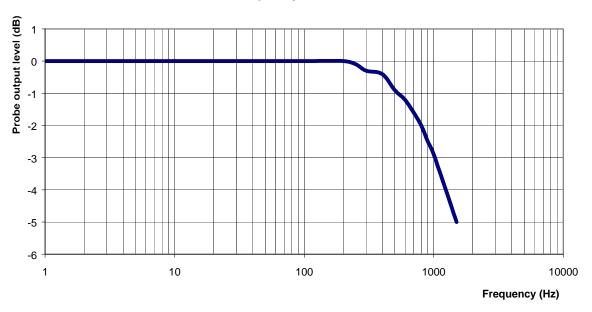

Video Bandwidth at 500 Hz 1 dB Video Bandwidth at 1.02 KHz: 3 dB

## **Conversion Factor Uncertainty Assessment**

Frequency: 2450MHz

**Epsilon:** 52.7 (+/-5%) **Sigma:** 1.95 S/m (+/-5%)

ConvF

**Channel X:** 4.01 7%(K=2)

**Channel Y:** 4.01 7%(K=2)

**Channel Z:** 4.01 7%(K=2)

To minimize the uncertainty calculation all tissue sensitivity values were calculated using a load impedance of 5 M $\Omega$ .

### **Boundary Effect:**

For a distance of 0.4mm the evaluated uncertainty (increase in the probe sensitivity) is less than 2%.

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2008.

### **NCL CALIBRATION LABORATORIES**

Calibration File No.: CP-987

Client.: APREL

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the NCL CALIBRATION LABORATORIES by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

Equipment: Miniature Isotropic RF Probe 2600 MHz

Manufacturer: APREL Laboratories

Model No.: E-030

Serial No.: 222

Calibration in Body Tissue

Calibration Procedure: SSI/DRB-TP-D01-032-E020-V2

Project No: Internal APREL

Calibrated: 9<sup>th</sup> January 2009 Released on: 10<sup>th</sup> January 2009

This Calibration Certificate is Incomplete Unless Accompanied with the Calibration Results Summary

Released By:

NCL CALIBRATION LABORATORIES

51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4161

#### Introduction

This Calibration Report reproduces the results of the calibration performed in line with the SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure. The results contained within this report are for APREL E-Field Probe E-030 222.

#### References

SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average Specific Absorption Rate (SAR) in the Human Body Due to Wireless Communications Devices: Experimental Techniques" SSI-TP-011 Tissue Calibration Procedure

#### **Conditions**

Probe 222 is a reassigned lab serial number for a new probe.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this probe has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Jesse Hones

## **Calibration Results Summary**

**Probe Type**: E-Field Probe E-030

Serial Number: 222

Frequency: 2600 MHz

Sensor Offset: 0.44 mm

Sensor Length: 2.5 mm

**Tip Enclosure:** Ertalyte\*

**Tip Diameter:** <2.9 mm

**Tip Length:** 60 mm

Total Length: 290 mm

# Sensitivity in Air

 Channel X:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Y:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Z:
  $1.2 \, \mu V/(V/m)^2$ 

**Diode Compression Point:** 95 mV

<sup>\*</sup>Resistive to recommended tissue recipes per IEEE-1528

### **Sensitivity in Body Tissue Measured**

Frequency: 2600 MHz

**Epsilon:** 51.30 (+/-5%) **Sigma:** 2.19 S/m (+/-5%)

ConvF

Channel X: 3.8

Channel Y: 3.8

Channel Z: 3.8

Tissue sensitivity values were calculated using the load impedance of the APREL Laboratories Daq-Paq and corrected for broadband calibration factor.

## **Boundary Effect:**

Uncertainty resulting from the boundary effect is less than 2% for the distance between the tip of the probe and the tissue boundary, when less than 1.44mm.

## **Spatial Resolution:**

The measured probe tip diameter is 2.9 mm (+/- 0.01 mm) and therefore meets the requirements of SSI/DRB-TP-D01-032 for spatial resolution.

### **Broad Band Calibration Measured:**

Frequency: 2580MHz

**Epsilon:** 51.37 (+/-5%) **Sigma:** 2.16 S/m (+/-5%)

Frequency: 2620 MHz


**Epsilon:** 51.22 (+/-5%) **Sigma:** 2.22 S/m (+/-5%)

The probe was assessed for sensitivity and conversion factor using a +/- 20MHz deviation from the centre frequency.

Deviation at -20MHz: -1.56% Deviation at +20MHz: +1.3%


#### **Modulation Result**

The probe is assessed with a square wave signal modulated at 20% and 40% duty cycles using a 10MHz and 20MHz span to assess the probe linearity and modulation response.

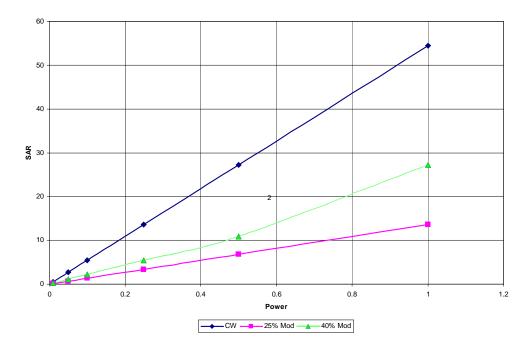



Duty Cycle (factor) 2 @ 20% Modulation: 3.2% Duty Cycle (factor) 4 @ 40% Modulation: 4.3%

# Receiving Pattern 2600 MHz (Air)

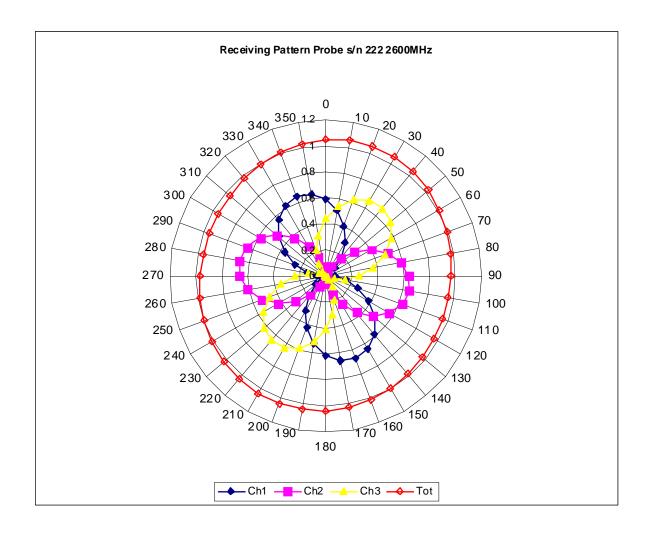

# Isotropy Error 2600 MHz (Air)

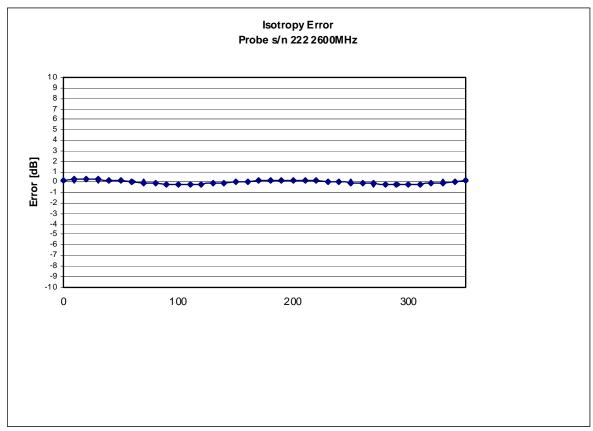

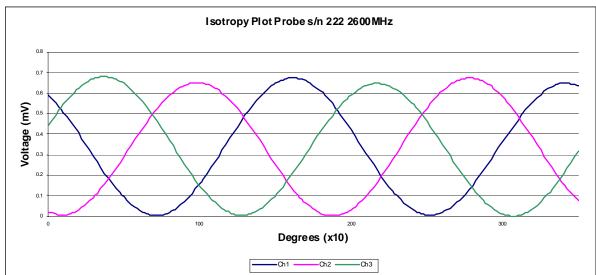


**Isotropicity in Tissue:** 

0.10 dB


Page 7 of 10

# **Dynamic Range**

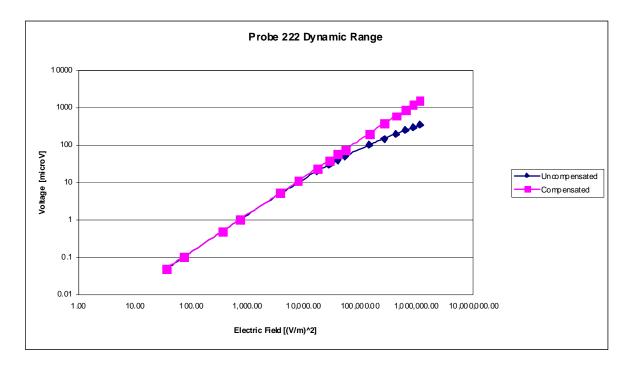

## **Video Bandwidth**

### **Probe Frequency Characteristics**

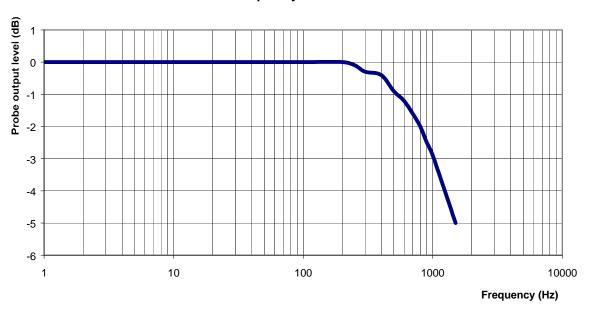

Video Bandwidth at 500 Hz 1 dB Video Bandwidth at 1.02 KHz: 3 dB

## **Conversion Factor Uncertainty Assessment**

Frequency: 2600 MHz

**Epsilon:** 51.30 (+/-5%) **Sigma:** 2.19 S/m (+/-5%)

ConvF

**Channel X:** 3.8 7%(K=2)

**Channel Y:** 3.8 7%(K=2)

**Channel Z:** 3.8 7%(K=2)

To minimize the uncertainty calculation all tissue sensitivity values were calculated using a load impedance of 5 M $\Omega$ .

### **Boundary Effect:**

For a distance of 0.4mm the evaluated uncertainty (increase in the probe sensitivity) is less than 2%.

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2009.

### **NCL CALIBRATION LABORATORIES**

Calibration File No.: CP-886

Client.: APREL

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the NCL CALIBRATION LABORATORIES by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

Equipment: Miniature Isotropic RF Probe 5200 MHz

Manufacturer: APREL Laboratories

Model No.: E-030 Serial No.: 222

Calibration in Body Tissue

Calibration Procedure: SSI/DRB-TP-D01-032-E020-V2

Project No: Internal APREL

Calibrated: 3<sup>rd</sup> May 2009 Released on: 3<sup>rd</sup> May 2009

This Calibration Certificate is Incomplete Unless Accompanied with the Calibration Results Summary

Released By:

NCL CALIBRATION LABORATORIES

51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4161

#### Introduction

This Calibration Report reproduces the results of the calibration performed in line with the SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure. The results contained within this report are for APREL E-Field Probe E-030 222.

#### References

SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average Specific Absorption Rate (SAR) in the Human Body Due to Wireless Communications Devices: Experimental Techniques" SSI-TP-011 Tissue Calibration Procedure

#### **Conditions**

Probe 222 was a new probe taken from stock prior to calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this probe has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Jesse Hones

## **Calibration Results Summary**

**Probe Type**: E-Field Probe E-030

Serial Number: 222

Frequency: 5200 MHz

Sensor Offset: 0.44 mm

Sensor Length: 2.5 mm

Tip Enclosure: Ertalyte\*

**Tip Diameter:** <2.9 mm

**Tip Length:** 60 mm

Total Length: 290 mm

# Sensitivity in Air

 Channel X:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Y:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Z:
  $1.2 \, \mu V/(V/m)^2$ 

**Diode Compression Point:** 95 mV

<sup>\*</sup>Resistive to recommended tissue recipes per IEEE-1528

### **Sensitivity in Body Tissue**

Frequency: 5200 MHz

**Epsilon:** 43.0 (+/-10%) **Sigma:** 5.75 S/m (+/-10%)

ConvF

Channel X: 3.3

Channel Y: 3.3

Channel Z: 3.3

Tissue sensitivity values were calculated using the load impedance of the APREL Laboratories Daq-Paq.

### **Boundary Effect:**

Uncertainty resulting from the boundary effect is less than 2% for the distance between the tip of the probe and the tissue boundary, when less than 2.44mm.

## **Spatial Resolution:**


The measured probe tip diameter is 2.9 mm (+/- 0.01 mm) and therefore meets the requirements of SSI/DRB-TP-D01-032 for spatial resolution.

#### **Broad Band Calibration:**


The probe was assessed for sensitivity and conversion factor using a +/- 40MHz deviation from the centre frequency.


Deviation at -40MHz: -4.16% Deviation at +40MHz: +2.78%

# Receiving Pattern 5200 MHz (Air)

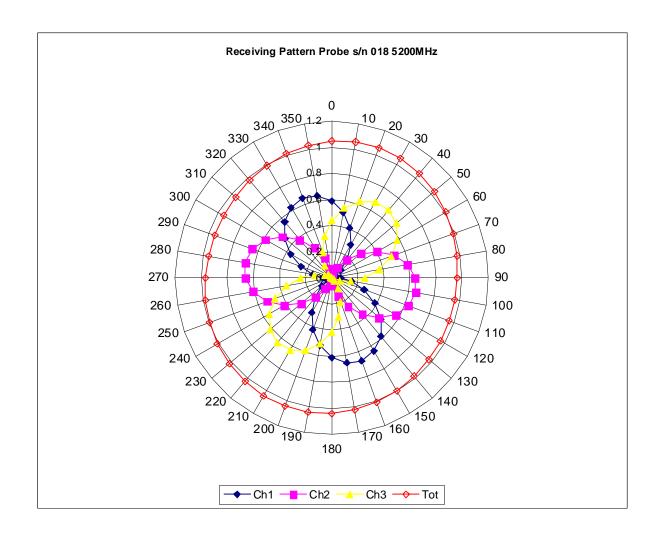

# Isotropy Error 5200 MHz (Air)

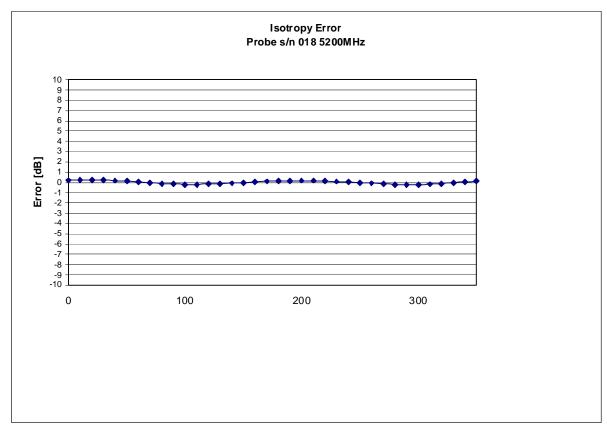


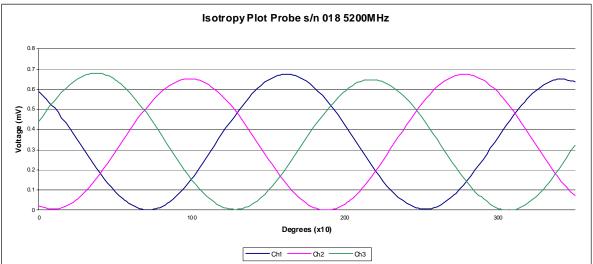

**Isotropicity in Tissue:** 


0.10 dB

# **Dynamic Range**

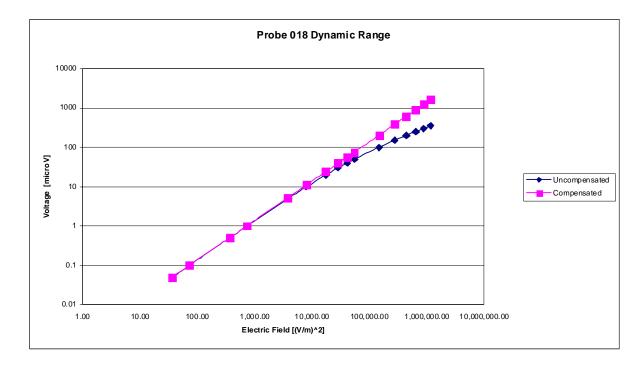

## **Video Bandwidth**

### **Probe Frequency Characteristics**

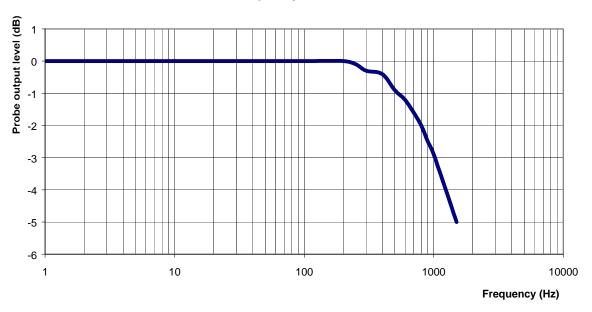

Video Bandwidth at 500 Hz 1 dB Video Bandwidth at 1.02 KHz: 3 dB

## **Conversion Factor Uncertainty Assessment**

Frequency: 5200MHz

**Epsilon:** 43.0 (+/-10%) **Sigma:** 5.75 S/m (+/-10%)

ConvF

**Channel X:** 3.3 7%(K=2)

**Channel Y:** 3.3 7%(K=2)

**Channel Z:** 3.3 7%(K=2)

To minimize the uncertainty calculation all tissue sensitivity values were calculated using a load impedance of 5 M $\Omega$ .

### **Boundary Effect:**

For a distance of 0.4mm the evaluated uncertainty (increase in the probe sensitivity) is less than 2%.

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2008.

### **NCL CALIBRATION LABORATORIES**

Calibration File No.: CP-887

Client.: APREL

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the NCL CALIBRATION LABORATORIES by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

Equipment: Miniature Isotropic RF Probe 5600 MHz

Manufacturer: APREL Laboratories

Model No.: E-030 Serial No.: 222

Calibration in Body Tissue

Calibration Procedure: SSI/DRB-TP-D01-032-E020-V2

Project No: Internal APREL

Calibrated: 3<sup>rd</sup> May 2009 Released on: 3<sup>rd</sup> May 2009

This Calibration Certificate is Incomplete Unless Accompanied with the Calibration Results Summary

Released By:

NCL CALIBRATION LABORATORIES

51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4161

#### Introduction

This Calibration Report reproduces the results of the calibration performed in line with the SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure. The results contained within this report are for APREL E-Field Probe E-030 222.

#### References

SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average Specific Absorption Rate (SAR) in the Human Body Due to Wireless Communications Devices: Experimental Techniques" SSI-TP-011 Tissue Calibration Procedure

#### **Conditions**

Probe 222 was a new probe taken from stock prior to calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this probe has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Jesse Hones

## **Calibration Results Summary**

**Probe Type**: E-Field Probe E-030

Serial Number: 222

Frequency: 5600 MHz

Sensor Offset: 0.44 mm

Sensor Length: 2.5 mm

Tip Enclosure: Ertalyte\*

**Tip Diameter:** <2.9 mm

Tip Length: 60 mm

Total Length: 290 mm

## Sensitivity in Air

 Channel X:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Y:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Z:
  $1.2 \, \mu V/(V/m)^2$ 

**Diode Compression Point:** 95 mV

<sup>\*</sup>Resistive to recommended tissue recipes per IEEE-1528

### **Sensitivity in Body Tissue Measured**

Frequency: 5600 MHz

**Epsilon:** 46.0 (+/-10%) **Sigma:** 5.85 S/m (+/-10%)

ConvF

Channel X: 3.0

Channel Y: 3.0

Channel Z: 3.0

Tissue sensitivity values were calculated using the load impedance of the APREL Laboratories Daq-Paq.

### **Boundary Effect:**

Uncertainty resulting from the boundary effect is less than 2% for the distance between the tip of the probe and the tissue boundary, when less than 2.44mm.

## **Spatial Resolution:**


The measured probe tip diameter is 2.9 mm (+/- 0.01 mm) and therefore meets the requirements of SSI/DRB-TP-D01-032 for spatial resolution.

#### **Broad Band Calibration:**


The probe was assessed for sensitivity and conversion factor using a +/- 40MHz deviation from the centre frequency.


Deviation at -40MHz: -3.77% Deviation at +40MHz: +4.28%

## Receiving Pattern 5600 MHz (Air)

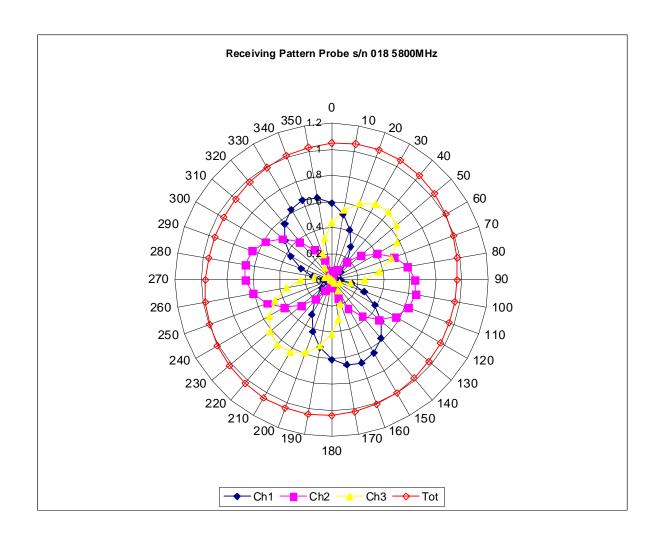

# Isotropy Error 5600 MHz (Air)

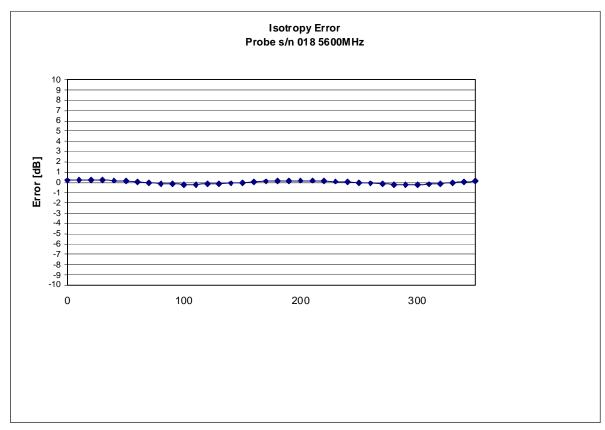


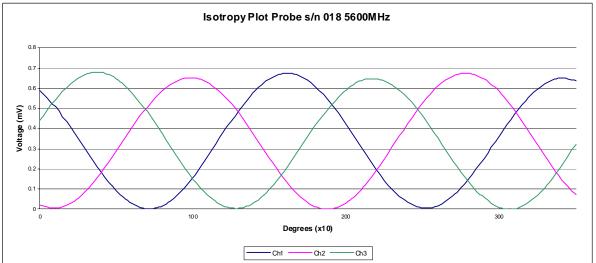

**Isotropicity in Tissue:** 


0.10 dB

# **Dynamic Range**

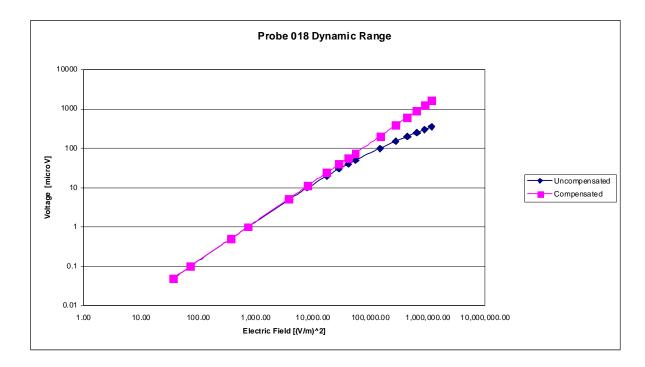

## **Video Bandwidth**

#### **Probe Frequency Characteristics**

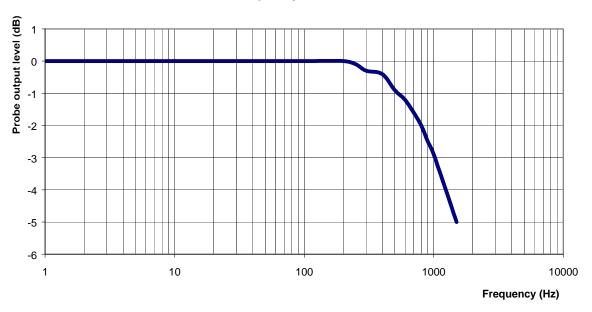

Video Bandwidth at 500 Hz 1 dB Video Bandwidth at 1.02 KHz: 3 dB

## **Conversion Factor Uncertainty Assessment**

## **Sensitivity in Body Tissue Measured**

Frequency: 5600 MHz

**Epsilon:** 46.0 (+/-10%) **Sigma:** 5.85 S/m (+/-10%)

ConvF

**Channel X:** 3.0 7%(K=2)

**Channel Y:** 3.0 7%(K=2)

**Channel Z:** 3.0 7%(K=2)

To minimize the uncertainty calculation all tissue sensitivity values were calculated using a load impedance of 5 M $\Omega$ .

#### **Boundary Effect:**

For a distance of 0.4mm the evaluated uncertainty (increase in the probe sensitivity) is less than 2%.

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2008.

#### **NCL CALIBRATION LABORATORIES**

Calibration File No.: CP-888

Client.: APREL

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the NCL CALIBRATION LABORATORIES by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

Equipment: Miniature Isotropic RF Probe 5800 MHz

Manufacturer: APREL Laboratories

Model No.: E-030 Serial No.: 222

Calibration in Body Tissue

Calibration Procedure: SSI/DRB-TP-D01-032-E020-V2

Project No: Internal APREL

Calibrated: 3<sup>rd</sup> May 2009 Released on: 3<sup>rd</sup> May 2009

This Calibration Certificate is Incomplete Unless Accompanied with the Calibration Results Summary

Released By:

NCL CALIBRATION LABORATORIES

51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4161

#### Introduction

This Calibration Report reproduces the results of the calibration performed in line with the SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure. The results contained within this report are for APREL E-Field Probe E-030 222.

#### References

SSI/DRB-TP-D01-032-E020-V2 E-Field Probe Calibration Procedure IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average Specific Absorption Rate (SAR) in the Human Body Due to Wireless Communications Devices: Experimental Techniques" SSI-TP-011 Tissue Calibration Procedure

#### **Conditions**

Probe 222 was a new probe taken from stock prior to calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this probe has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Jesse Hones

## **Calibration Results Summary**

**Probe Type**: E-Field Probe E-030

Serial Number: 222

Frequency: 5800 MHz

Sensor Offset: 0.44 mm

Sensor Length: 2.5 mm

Tip Enclosure: Ertalyte\*

**Tip Diameter:** <2.9 mm

**Tip Length:** 60 mm

Total Length: 290 mm

## Sensitivity in Air

 Channel X:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Y:
  $1.2 \, \mu V/(V/m)^2$  

 Channel Z:
  $1.2 \, \mu V/(V/m)^2$ 

**Diode Compression Point:** 95 mV

<sup>\*</sup>Resistive to recommended tissue recipes per IEEE-1528

### **Sensitivity in Body Tissue**

Frequency: 5800 MHz

**Epsilon:** 48.2 (+/-10%) **Sigma:** 6.0 S/m (+/-10%)

ConvF

Channel X: 3.2

Channel Y: 3.2

**Channel Z:** 3.2

Tissue sensitivity values were calculated using the load impedance of the APREL Laboratories Daq-Paq and corrected for broadband calibration factor.

### **Boundary Effect:**

Uncertainty resulting from the boundary effect is less than 2% for the distance between the tip of the probe and the tissue boundary, when less than 2.44mm.

## **Spatial Resolution:**


The measured probe tip diameter is 2.9 mm (+/- 0.01 mm) and therefore meets the requirements of SSI/DRB-TP-D01-032 for spatial resolution.

#### **Broad Band Calibration:**


The probe was assessed for sensitivity and conversion factor using a +/- 40MHz deviation from the centre frequency.


Deviation at -40MHz: -3.07% Deviation at +40MHz: +3.22%

## Receiving Pattern 5800 MHz (Air)

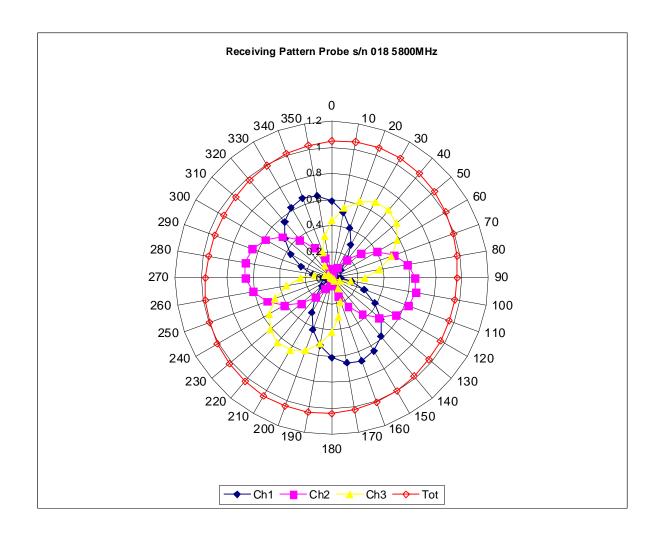

# Isotropy Error 5800 MHz (Air)

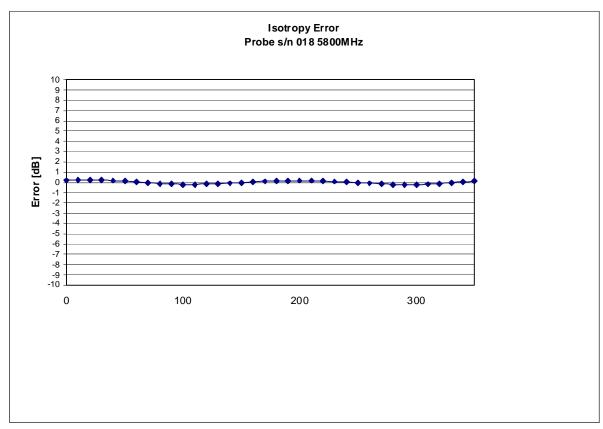


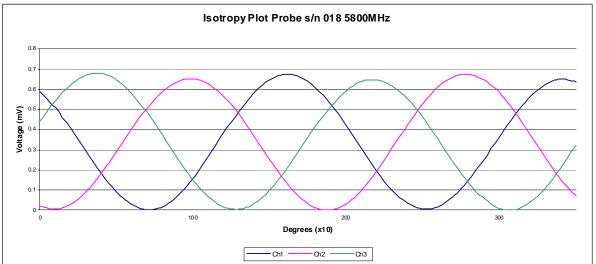

**Isotropicity in Tissue:** 


0.10 dB

# **Dynamic Range**

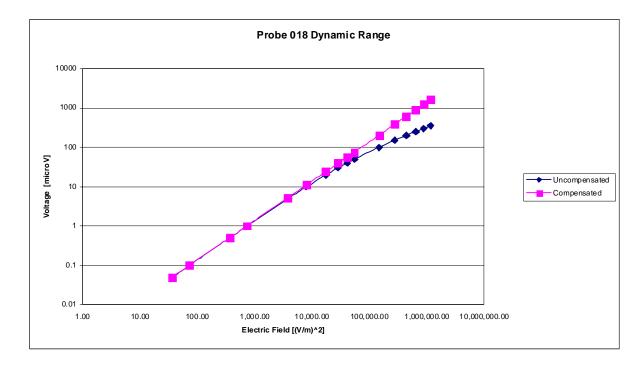

## **Video Bandwidth**

#### **Probe Frequency Characteristics**

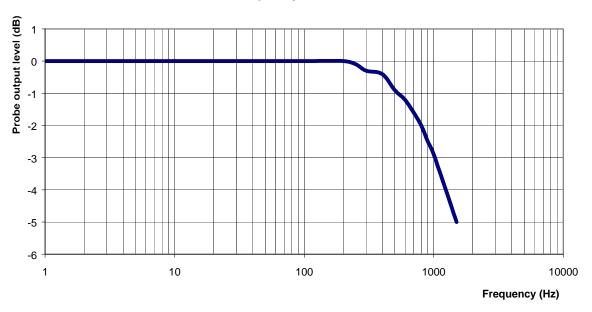

Video Bandwidth at 500 Hz 1 dB Video Bandwidth at 1.02 KHz: 3 dB

## **Conversion Factor Uncertainty Assessment**

Frequency: 5800MHz

**Epsilon:** 48.2 (+/-10%) **Sigma:** 6.0 S/m (+/-10%)

ConvF

**Channel X:** 3.2 7%(K=2)

**Channel Y:** 3.2 7%(K=2)

**Channel Z:** 3.2 7%(K=2)

To minimize the uncertainty calculation all tissue sensitivity values were calculated using a load impedance of 5 M $\Omega$ .

#### **Boundary Effect:**

For a distance of 0.4mm the evaluated uncertainty (increase in the probe sensitivity) is less than 2%.

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2008.

#### **NCL CALIBRATION LABORATORIES**

Calibration File No: DC-889
Project Number: APREL-ALSAS10U

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the **NCL CALIBRATION LABORATORIES** by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

**APREL Validation Dipole** 

Manufacturer: APREL Laboratories
Part number: ALS-D-2450-S-2
Frequency: 2450 MHz
Serial No: 301581

Customer: APREL

Calibrated: 4<sup>th</sup> May 2009 Released on: 4<sup>th</sup> May 2009

Released By:	



51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4162

· -	$\boldsymbol{\smallfrown}$	n	~	111	$\sim$	ns
u	u		u	ILI	u	113

Dipole 301581 was new and taken from stock prior to calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} \, +/- \, 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} \, +/- \, 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this device has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Constantine Teodorian

## **Calibration Results Summary**

The following results relate the Calibrated Dipole and should be used as a quick reference for the user.

#### **Mechanical Dimensions**


**Length:** 52.4 mm **Height:** 30.3 mm

### **Electrical Specification**

SWR: 1.056 U Return Loss: -32.0 dB Impedance:  $50.2 \Omega$ 

### **System Validation Results**

Frequency	1 Gram	10 Gram	Peak
2450 MHz	53.1	24.4	101.8

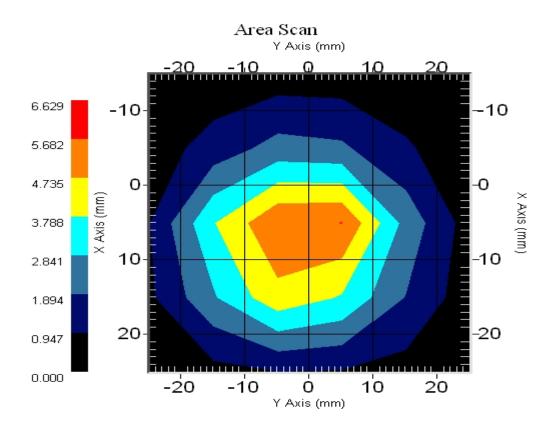

#### Introduction

This Calibration Report has been produced in line with the SSI Dipole Calibration Procedure SSI-TP-018-ALSAS. The results contained within this report are for Validation Dipole 301581. The calibration routine consisted of a three-step process. Step 1 was a mechanical verification of the dipole to ensure that it meets the mechanical specifications. Step 2 was an Electrical Calibration for the Validation Dipole, where the SWR, Impedance, and the Return loss were assessed. Step 3 involved a System Validation using the ALSAS-10U, along with APREL E-020 130 MHz to 26 GHz E-Field Probe Serial Number 212.

#### References

SSI-TP-018-ALSAS Dipole Calibration Procedure

SSI-TP-016 Tissue Calibration Procedure

IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average Specific Absorption Rate (SAR) in the Human Body Due to Wireless Communications Devices: Experimental Techniques"

IEC-62209 "Human exposure to radio frequency fields from hand-held and body-mounted wireless communication devices – Human models, instrumentation, and procedures"

Part 1: "Procedure to determine the Specific Absorption Rate (SAR) for hand-held devices used in close proximity of the ear (frequency range of 300 MHz to 3 GHz)" IEC-62209 "Human exposure to radio frequency fields from hand-held and body-mounted wireless communication devices – Human models, instrumentation, and procedures"

Part 2 *Draft*: "Procedure to determine the Specific Absorption Rate (SAR) for handheld devices used in close proximity of the ear (frequency range of 30 MHz to 6 GHz)"

#### **Conditions**

Dipole 301581 was new taken from stock.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $20 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

# **Dipole Calibration Results**

### **Mechanical Verification**

APREL	APREL APREL		Measured
Length	Height	Length	Height
51.5 mm	30.4 mm	52.4 mm	30.3 mm

#### **Tissue Validation**


Head Tissue 2450 MHz	Measured
Dielectric constant, ε <sub>r</sub>	39.2
Conductivity, σ [S/m]	1.80

#### **Electrical Calibration**


Test	Result
S11 R/L	-32.0 dB
SWR	1.05 U
Impedance	50.2 Ω

The Following Graphs are the results as displayed on the Vector Network Analyzer.


#### **S11 Parameter Return Loss**

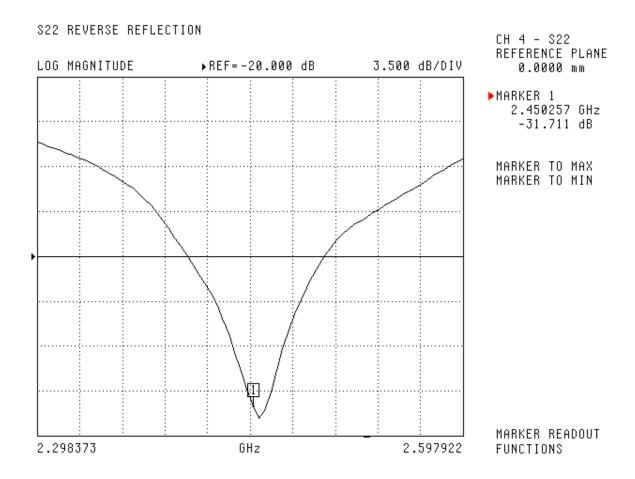

#### SWR

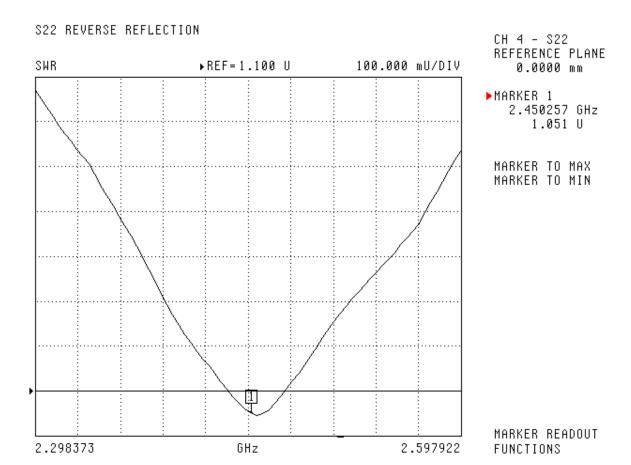


## **Smith Chart Dipole Impedance**

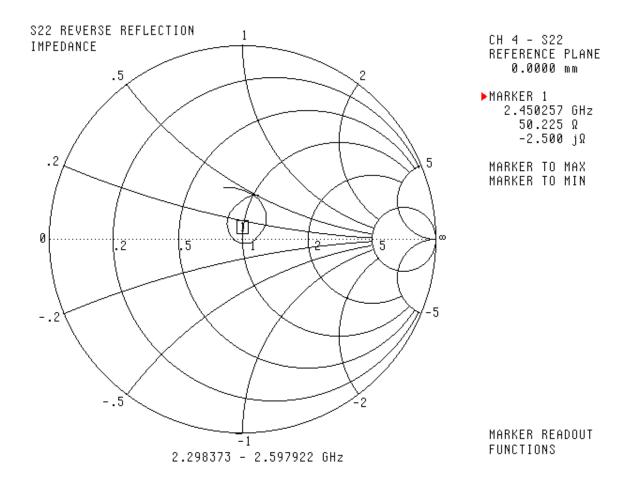

## **System Validation Results Using the Electrically Calibrated Dipole**

Frequency	1 Gram	10 Gram	Peak
2450 MHz	53.1	24.4	101.8

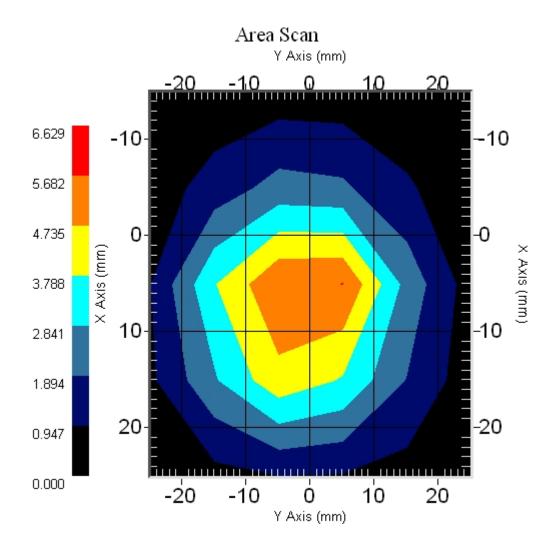

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2008.

#### **NCL CALIBRATION LABORATORIES**

Calibration File No: DC-890 Project Number: APREL-ALSAS 10U

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the **NCL CALIBRATION LABORATORIES** by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

**APREL Validation Dipole** 

Manufacturer: APREL Laboratories
Part number: ALS-D-5258-S-2
Frequency: 5.2GHz to 5.8GHz
Serial No: 5258-235-00802

Customer: APREL

Serial Number: ALS-BB-001

Calibrated: 24<sup>th</sup> May 2008 Released on: 24<sup>th</sup> May 2008

Released By:	



51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4162

· -	$\boldsymbol{\smallfrown}$	n	~	111	$\sim$	ns
u	u		u	ILI	u	113

Dipole 5258-235-00802 was new and taken from stock prior to calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} \, +/- \, 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} \, +/- \, 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this device has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

Constantine Teodorian

This page has been reviewed for content and attested to by signature within this document.

## **Calibration Results Summary**

The following results relate the Calibrated Dipole and should be used as a quick reference for the user.

#### **Mechanical Dimensions**


Length: 23.3 mm Height: 20.3 mm

### **Electrical Specification**


SWR: 1.22 U Return Loss: -20.0 dB Impedance:  $50.0 \Omega$ 

#### **System Validation Results**


Frequency	1 Gram	10 Gram	Peak
5200 MHz	51.9	17.9	223.1

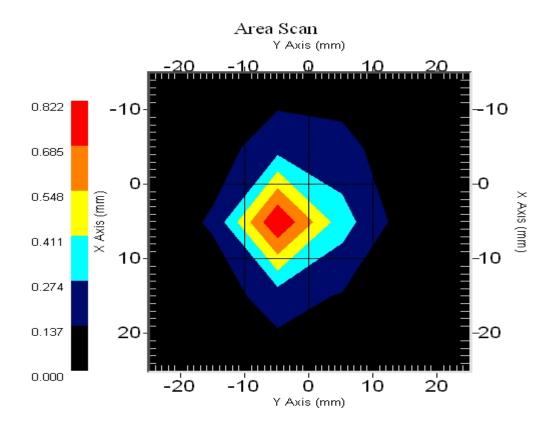

Frequency	1 Gram	10 Gram	Peak
5600 MHz	52.97	18.2	243.1

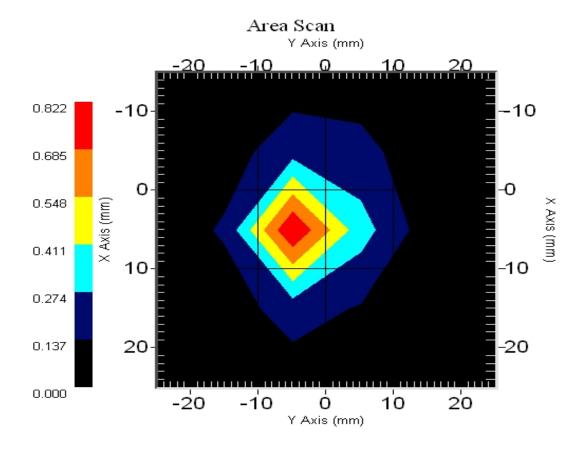

Frequency	1 Gram	10 Gram	Peak
5800 MHz	48.97	17.2	207.1

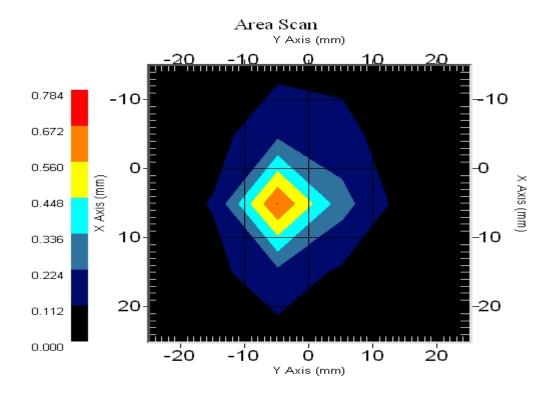

#### Introduction

This Calibration Report has been produced in line with the SSI Dipole Calibration Procedure SSI-TP-018-ALSAS. The results contained within this report are for Validation Dipole 5258-235-00802. The calibration routine consisted of a three-step process. Step 1 was a mechanical verification of the dipole to ensure that it meets the mechanical specifications. Step 2 was an Electrical Calibration for the Validation Dipole, where the SWR, Impedance, and the Return loss were assessed. Step 3 involved a System Validation using the ALSAS-10U, along with APREL E-030 018 E-Field Probe.

#### References

SSI-TP-018-ALSAS Dipole Calibration Procedure
SSI-TP-016 Tissue Calibration Procedure
IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average
Specific Absorption Rate (SAR) in the Human Body Due to Wireless
Communications Devices: Experimental Techniques"

#### **Conditions**

Dipole 5258-235-00802 was new taken from stock.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $20 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

# **Dipole Calibration Results**

#### **Tissue Validation**

Head Tissue 5200 MHz	Measured
Dielectric constant, ε <sub>r</sub>	47.0
Conductivity, σ [S/m]	5.30

Head Tissue 5600 MHz	Measured
Dielectric constant, ε <sub>r</sub>	46.1
Conductivity, σ [S/m]	5.78

Head Tissue 5800 MHz	Measured
Dielectric constant, ε <sub>r</sub>	46.7
Conductivity, σ [S/m]	6.22

### **Mechanical Verification**


APREL Length	APREL Height	Measured Length	Measured Height
23.1 mm	20.7 mm	23.3 mm	20.3 mm

### **Electrical Calibration**


S11	5200MHz	5800MHz
RL (dB)	-21.16	-22.34
SWR	1.2	1.17
Impedance (ohms)	51.38	43.92

The Following Graphs are the results as displayed on the Vector Network Analyzer.


#### **S11 Parameter Return Loss**

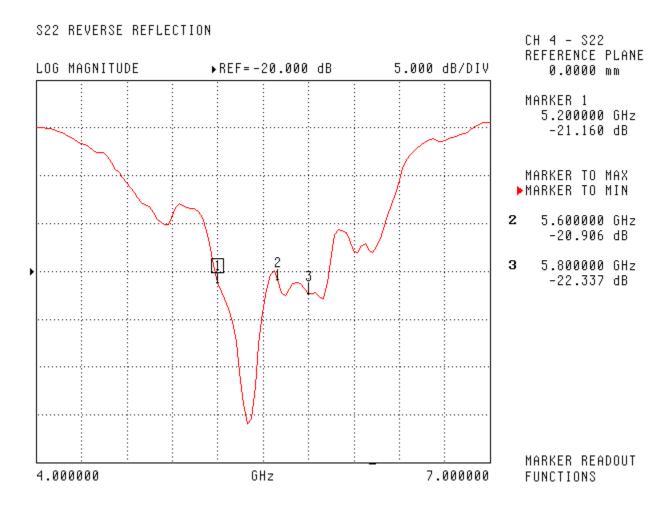

### SWR

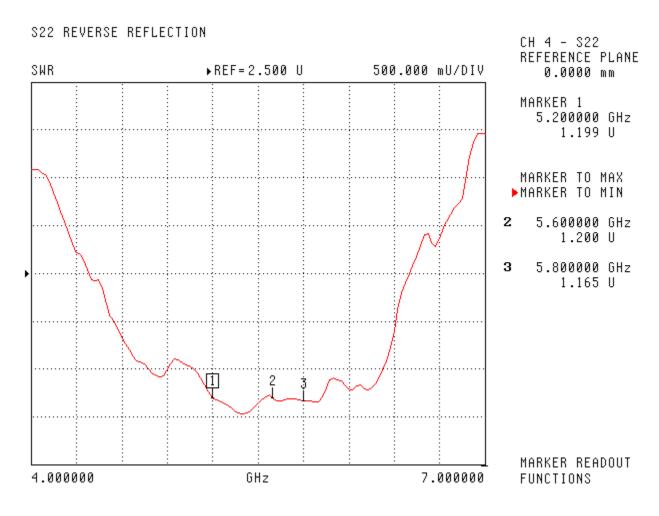

## **Smith Chart Dipole Impedance**

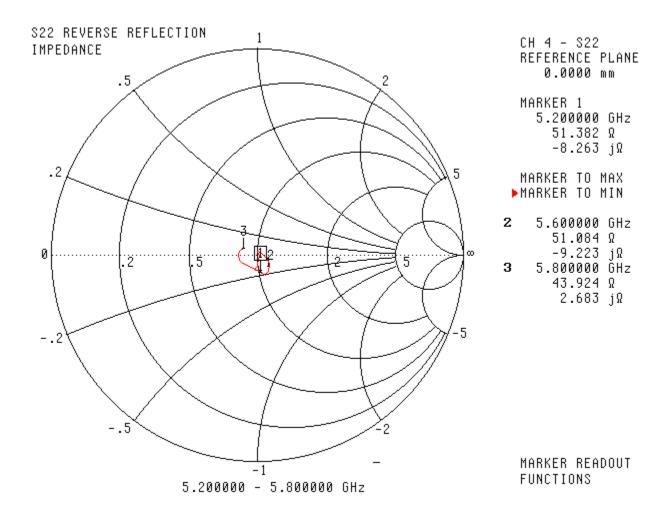

### **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2008.

#### **NCL CALIBRATION LABORATORIES**

Calibration File No: DC-915 Project Number: APREL-ALSAS-10U

## CERTIFICATE OF CALIBRATION

It is certified that the equipment identified below has been calibrated in the **NCL CALIBRATION LABORATORIES** by qualified personnel following recognized procedures and using transfer standards traceable to NRC/NIST.

**APREL Validation Dipole** 

Manufacturer: APREL Laboratories
Part number: ALS-D-835-S-2
Frequency: 835 MHz
Serial No: 835-180-00554

Customer: APREL

Calibrated: 1<sup>st</sup> March 2009 Released on: 1<sup>st</sup> March 2009

Released Bv:		
Released by		



51 SPECTRUM WAY NEPEAN, ONTARIO CANADA K2R 1E6 Division of APREL Lab. TEL: (613) 820-4988 FAX: (613) 820-4162

_				• 4		
•	$\smallfrown$	n	М	ПŤ	$\mathbf{I}$	ns
	LJ					

D. Brooks

Dipole 835-180-00554 was a re-calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} \, +/- \, 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $21 \,^{\circ}\text{C} \, +/- \, 0.5 \,^{\circ}\text{C}$ 

We the undersigned attest that to the best of our knowledge the calibration of this device has been accurately conducted and that all information contained within this report has been reviewed for accuracy.

Stuart Nicol

This page has been reviewed for content and attested to by signature within this document.

### **Calibration Results Summary**

The following results relate the Calibrated Dipole and should be used as a quick reference for the user.

#### **Mechanical Dimensions**


**Length:** 161.0 mm **Height:** 89.8 mm

### **Electrical Specification**

SWR: 1.01 U Return Loss: -48.6 dB Impedance:  $49.6 \Omega$ 

#### **System Validation Results**

Frequency	1 Gram	10 Gram	Peak
835 MHz	9.49	6.1	14.21

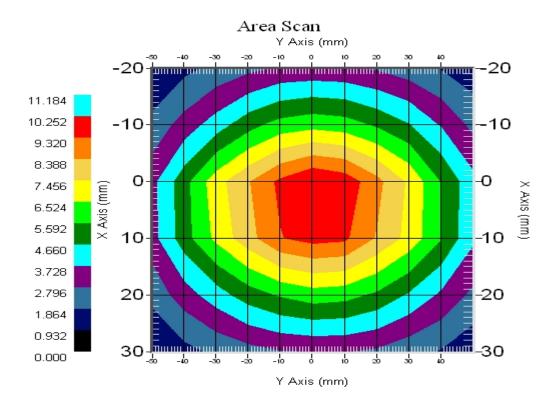

#### Introduction

This Calibration Report has been produced in line with the SSI Dipole Calibration Procedure SSI-TP-018-ALSAS. The results contained within this report are for Validation Dipole 835-180-00554. The calibration routine consisted of a three-step process. Step 1 was a mechanical verification of the dipole to ensure that it meets the mechanical specifications. Step 2 was an Electrical Calibration for the Validation Dipole, where the SWR, Impedance, and the Return loss were assessed. Step 3 involved a System Validation using the ALSAS-10U, along with APREL E-020 130 MHz to 26 GHz E-Field Probe Serial Number 212.

#### References

SSI-TP-018-ALSAS Dipole Calibration Procedure
SSI-TP-016 Tissue Calibration Procedure
IEEE 1528 "Recommended Practice for Determining the Peak Spatial-Average
Specific Absorption Rate (SAR) in the Human Body Due to Wireless
Communications Devices: Experimental Techniques"

#### **Conditions**

Dipole 835-180-00554 was a re-calibration.

Ambient Temperature of the Laboratory:  $22 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ Temperature of the Tissue:  $20 \,^{\circ}\text{C} + /- 0.5 \,^{\circ}\text{C}$ 

# **Dipole Calibration Results**

#### **Mechanical Verification**

APREL	APREL	Measured	Measured
Length	Height	Length	Height
161.0 mm	89.8 mm	162.2 mm	90.3 mm

#### **Tissue Validation**


Head Tissue 835MHz	Measured
Dielectric constant, ε <sub>r</sub>	41.5
Conductivity, σ [S/m]	0.90

#### **Electrical Calibration**

Test	Result	
S11 RL	-48.6 dB	
SWR	1.01 U	
Impedance	49.6 Ω	


The Following Graphs are the results as displayed on the Vector Network Analyzer.

#### **S11 Parameter Return Loss**

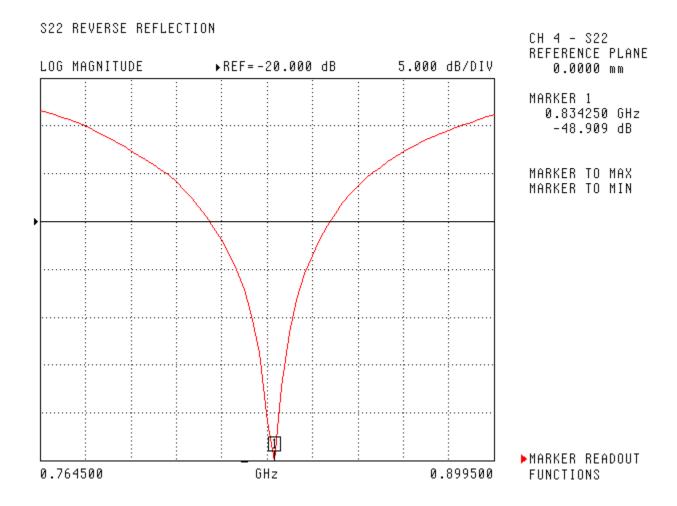


### **SWR**



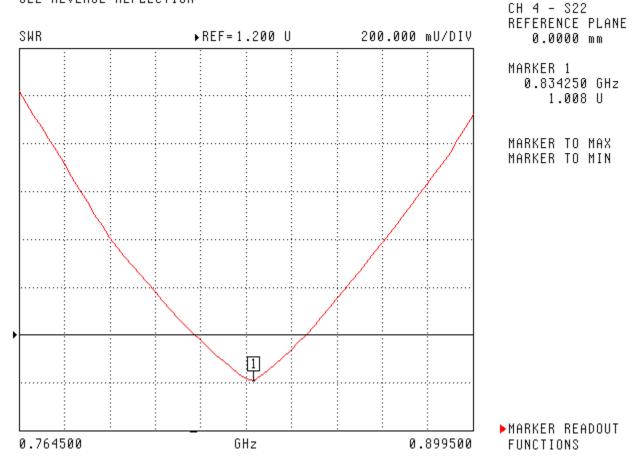


## **Smith Chart Dipole Impedance**

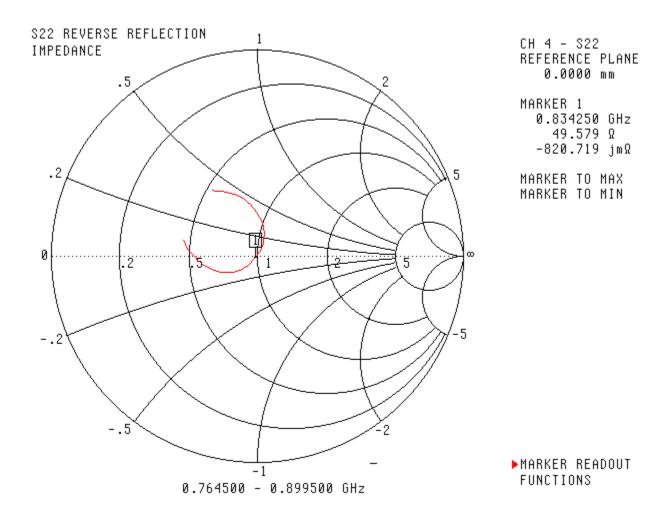

### System Validation Results Using the Electrically Calibrated Dipole

Head Tissue Frequency	1 Gram	10 Gram	Peak Above Feed Point
835 MHz	9.49	6.1	14.21

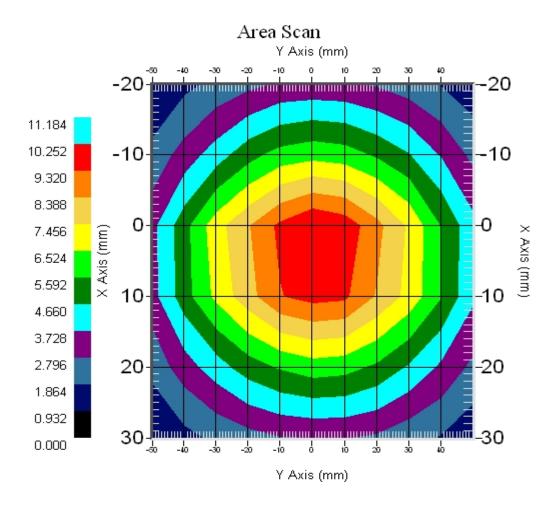

## **Test Equipment**

The test equipment used during Probe Calibration, manufacturer, model number and, current calibration status are listed and located on the main APREL server R:\NCL\Calibration Equipment\Instrument List May 2007.