
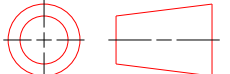


NOTES:

- 1. JOIN (3) TO (2) BY SONIC WELDING.
- 2. ASSEMBLE (1) TO (2) & (3) AFTER STEP 1 IS COMPLETE.
- 3. REFER TO PROE MODELS: 3983-000.ASM, 3983-MST.PRT.

REVISION HISTORY				
REV	EN	DESCRIPTION	DATE	EN ORIGINATOR
01	--	INITIAL ENG RELEASE	4-21-03	--
02	--	ADDED LABELS, ITEMS 5 - 9	6-12-03	--
03	--	REMOVED: 4194-000 & 4355-000. RENUMBERED ITEMS	9-10-03	--

UNLESS OTHERWISE SPECIFIED DIMENSIONS ARE IN INCHES TOLERANCES ON: .X ± 0.1 .XX ± 0.01 .XXX ± 0.001 ANGLES ± 0.5	DRAWN BY ADB	DATE 4-21-03	<div> NONIN MEDICAL, INC. PLYMOUTH, MN</div> <div>CONFIDENTIAL The information disclosed herein is confidential and is proprietary to Nonin Medical, Inc. Nonin reserves all patent, proprietary and reproduction rights therein. This information is not to be used or disclosed to any persons without the express written authorization of a representative of Nonin Medical, Inc.</div>		
	CHECKED BY	DATE			
	ELECTRONIC FILE NAME 3983-000.DRW		TITLE FINAL ASSY, PATIENT OXIMETER, BLUETOOTH		
	DRAWING CONFORMS TO ANSI Y 14.5M-1994		SIZE B	DWG NO 3983-000	REV 003
THIRD ANGLE PROJECTION <div></div>	DO NOT SCALE DRAWING		SCALE N/A	SHEET 1	OF 1