

PAM™3000

Patient Assessment Monitor

Label Position

1 Introduction

This document is a reference for label placement on the PAM 3000 Bed Sensor Panel. The Label appearance and its placement must be consistent for engineering and production units.

2 Label appearance

The compliancy PAM 3000 label (Wireless 2000 P/N 621.001.10.00) is a matt silver aluminium foil adhesive type, made from 3M-7940 material.

It must contain the follow information:

- a) Wireless 2000 logo
- b) Model name
- c) CSA regulatory notice of using unique accessories (Power supply W2K P/N 375.001.10.00)
- d) CSA regulatory notice of proper grounding
- e) Manufacturer name
- f) CSA Master Contract number
- g) Type BF applied part logo as per IEC 60417-5333
- h) FCC ID number



Fig. 1 PAM 3000 Label

PAM™3000

Patient Assessment Monitor

Label Position

3 Label Placement

The label should be placed on the lower part of the PAM 3000 Bed Sensor Panel (power connector facing down) from the left side of a power jack. The lower edge of the label should be not closer than 2cm from the lower edge of the Bed Sensor. The left edge of the label should be as far as 7 cm from the center line of the Bed Sensor, referring to the power jack.

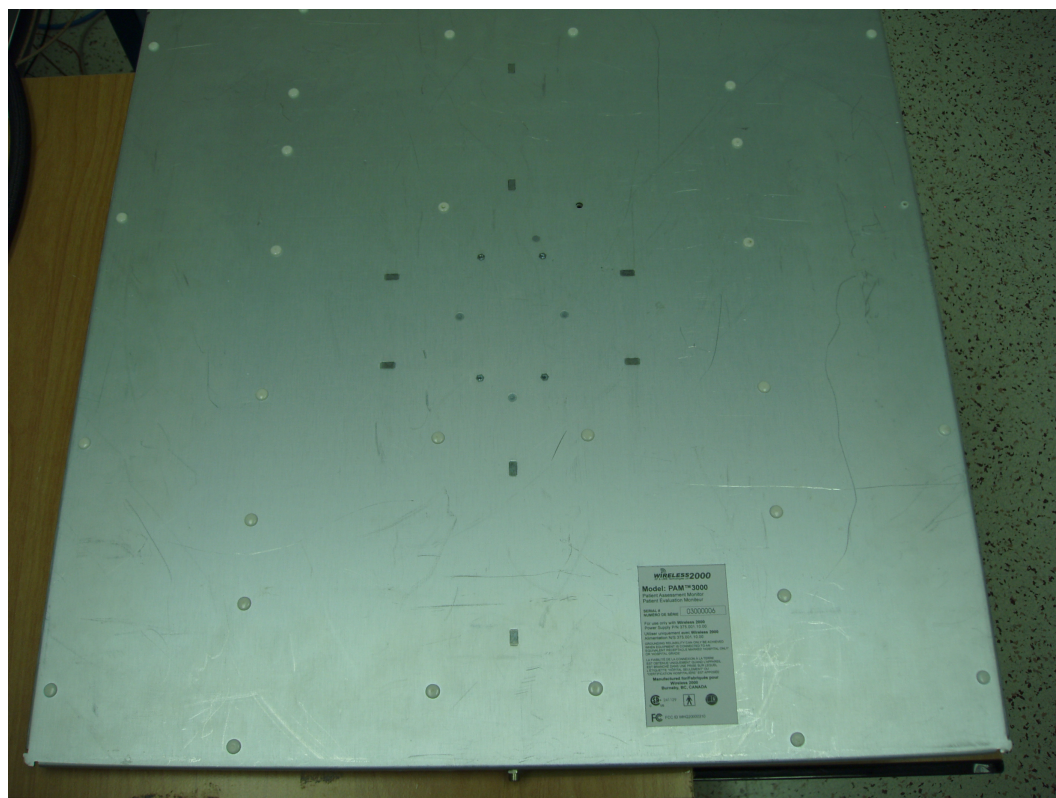


Fig.2 Label Placement